



TWO-YEAR WARRANTY

This Dazzlepro product is warranted to be free of manufacturer's defects in materials or workmanship for one year from the date of purchase. Damage or wear resulting from an accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair are not covered under this warranty.

Should this product require warranty service (or replacement at our discretion) please contact client service to obtain a Return Merchandise Authorization number (RMA) and return instructions, proof of purchase is required. Products returned without a Dazzlepro generated RMA number will not be accepted and the sender will not receive a refund, replacement, or repaired product.



info@dazzlepro.com
1.888.358.0330

There are no express warranties except as listed above. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

DO NOT RETURN TO RETAILER, PLEASE CONTACT US DIRECTLY FOR WARRANTY SERVICE OR REPAIR. THANK YOU.

Thank you for choosing
Dazzlepro
in the spirit of welcoming you
we would like to offer ...



Your next purchase!
Just enter this exclusive code



at checkout when you visit
dazzlepro.com

© 2015, Dazzlepro. All Rights Reserved

DazzleproTM
CLEAN CASE
UV Sanitizing Chamber



Welcome to Dazzlepro™

Please take a moment to read these instructions thoroughly so that you can enjoy all the features of your Dazzlepro Clean Case UV Sanitizing Chamber.

The Dazzlepro Clean Case harnesses the power of germicidal UV light to kill microorganisms, mold, fungi and bacteria in just 10 minutes with no water, soaps, powders or chemicals needed.

The uniquely designed tray is perfect for cleaning dentures, retainers, mouth guards, night guards, whitening trays and other dental devices, leaving them sanitized and ready for use whenever you need them. Now you can enjoy clean devices with confidence, and that's something to smile about!

Say Hello to a Beautiful Smile.




Inserting the Batteries

- 1 Turn the unit over and remove the battery cover.
- 2 Insert 2 AA alkaline batteries, being sure to position them correctly according to the polarity indicator. Replace the battery cover.
- 3 For best results, always use fresh batteries.



Cleaning the Unit

Use a soft, slightly damp cloth to wipe the outside and inside chamber of unit clean. Do not use harsh or abrasive chemicals or materials which may damage the unit or cause it to malfunction.

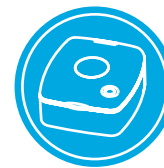

Using the Clean Case

Rinse your dental device in water and pat dry to remove any residual build-up before sanitizing.

Open the lid and insert your device.

Close the lid completely then push the Power button on top of the lid to activate the sanitization process. The Power button will light up to indicate the sanitization process has begun.

Sanitization will automatically stop if the lid is opened at any time during the process. It is important to complete the sanitization process to insure that all microorganisms are destroyed.

The unit will shut off automatically after approximately 10 minutes when the sanitization process is complete. Open the lid and remove your device for immediate use, or store your device in the chamber to keep it clean for future use.

