



Good Brushing Technique:

While proper brushing should take at least two minutes, most adults don't come close to brushing that long. To make the process easier, the Elite Sonic Toothbrush is equipped with a Two-Minute Smart Auto Timer and Area Change Alert which reminds you to brush each section of your mouth for 30 seconds. Once you have achieved the desired brushing time, it will automatically shut off.

To experience the superior cleaning power of the Elite Sonic Toothbrush, place the bristles gently against your teeth and let the sonic movement go to work. Don't press too hard or scrub, simply move the brush around your mouth for the allotted amount of time. Pay extra attention to the gum-line, hard-to-reach back teeth, and areas around crowns, fillings, or other dental work.


After brushing, switch to the Interdental Brush Head for a deeper clean between your teeth where the standard head can't reach. Similar to flossing, The Interdental Brush Head effectively removes plaque and food debris, promotes gum health, and freshens breath. It's also ideal for cleaning between braces, crowns, and implants.

How Important is the Toothpaste I Use?

It's important that you choose a toothpaste that's right for your specific oral needs. A wide selection of toothpastes are available and designed for nearly every condition, so ask your dental professional which toothpaste is right for you.

Your Elite Sonic Toothbrush

- 1 Brush Head
- 2 Power Button
- 3 Brush Mode Button
- 4 Mode Indicators
- 5 Brush Handle
- 6 Charge Indicator
- 7 Charging Base
- 8 UV Chamber
- 9 UV Chamber Cover
- 10 Removable Drip Cup

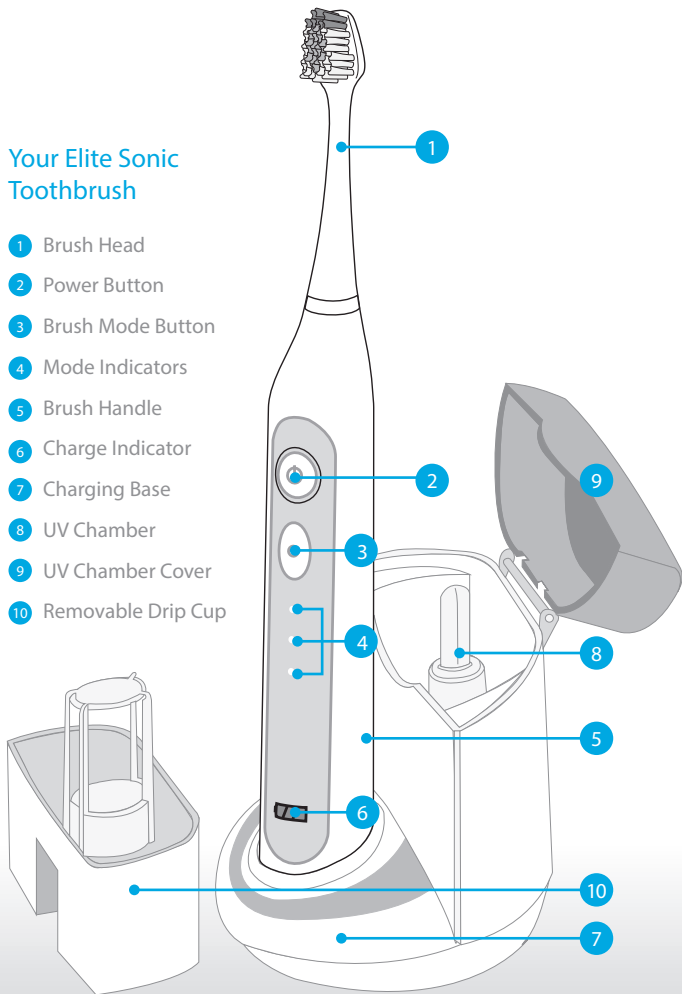



TWO - YEAR WARRANTY

This Dazzlepro product is warranted to be free of manufacturer's defects in materials or workmanship for two years from the date of purchase. Damage or wear resulting from an accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair are not covered under this warranty.

Should this product require warranty service (or replacement at our discretion) please contact client service to obtain a Return Merchandise Authorization number (RMA) and return instructions, proof of purchase is required. Products returned without a Dazzlepro generated RMA number will not be accepted and the sender will not receive a refund, replacement, or repaired product.

Dazzlepro

info@dazzlepro.com

1.888.358.0330

There are no express warranties except as listed above. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

DO NOT RETURN TO RETAILER. PLEASE CONTACT US DIRECTLY
FOR WARRANTY SERVICE OR REPAIR. THANK YOU.



Safety & Precautions

Always follow basic safety precautions when using electronic devices.

The Dazzlepro Elite Sonic Toothbrush has a waterproof handle, is electrically safe, and is specifically designed for use in the bathroom. If the device fails to function, do not attempt to dismantle, exchange, or manipulate any part of it, in order to avoid risk of electric shock. This product is not intended for use by children under the age of three.



A white electric toothbrush with a grey and white bristle head is shown diagonally. The handle has a grey section with a power button. The background features a dynamic splash of blue water.

Dazzlepro™

ELITE SONIC TOOTHBRUSH

Platinum Edition with UV Sanitizing Charging Base

20150225

Welcome to

Dazzlepro™

Please take a moment to read these instructions thoroughly, so that you can enjoy all the features that make the Dazzlepro Elite Sonic Toothbrush our most powerful toothbrush yet.

At 40,000 brush strokes per minute, the Dazzlepro Elite Sonic Toothbrush removes more stains for an overall whiter smile. The Two Minute Smart Auto Timer and Area Change Alert help you achieve the optimal brush time, and the UV Sanitizing Charging Base effectively removes 99.9%* of tested bacteria and germs from your toothbrush heads while recharging it for your next use. Choose between 3 Operating Modes to fit your oral care needs.

Say Hello to a Beautiful Smile.

* Sanitizes and kills 99.9% of Staphylococcus aureus, Enterococcus hirae, Pseudomonas aeruginosa, Escherichia coli.



3 Brushing Modes

Choose the Brushing Mode that matches your oral care needs:

- Clean - gentle cleaning mode that's great for sensitive areas
- Dazzle - powerful and thorough mode
- Pulse - massage and light stimulation mode for healthy gums

Selecting Your Brushing Mode

- Turn the brush on by pressing the "Power" button (2).
- Press the "Brush Mode" button (3) to cycle through, and select the desired mode: Clean, Dazzle, or Pulse.
- The "Mode Indicators" (4) will light up to show which mode you are selecting.
- When powered on, your Elite Sonic will start in the mode you last used.
- The brush will automatically turn off after 2 minutes. To turn your brush off before the 2 minutes, press the "Power" button.

It is recommended to store your brush heads in the UV Sanitizing Chamber in order to keep them clean and germ free* before each use.

* Sanitizes and kills 99.9% of *Staphylococcus aureus*, *Enterococcus hirae*, *Pseudomonas aeruginosa*, and *Escherichia coli* after 8 minute timer is complete.



Sanitizing Your Brush Heads


- Open the UV Chamber Cover (9) and fit the brush heads onto the small pegs at the bottom of the Drip Cup (10). Make sure the bristles of brush heads are facing the UV tube to achieve effective results.
- Close the UV Chamber Cover to automatically activate the UV light. After 7 minutes, when the sanitization process is complete, the UV light will automatically shut-off. If the cover is opened before the end of the cycle, the light will also turn off.

Subscribe & Save Program

To keep your Dazzlepro toothbrush in top condition, dentists recommend replacing your brush heads every two to three months.

Our Subscribe & Save Program delivers new brush heads right to your door, exactly when you need them.

New brush heads are available at dazzlepro.com

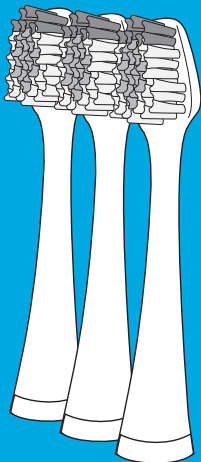



Two-Minute Smart Auto Timer

The Elite Sonic Toothbrush is programmed to run through a complete brushing cycle in two minutes, so you can easily follow ADA recommendations. The cycle is broken down into four, 30-second intervals for each area of your mouth. When the 30 seconds is up, your brush will alert you and briefly pause, allowing you to switch. Once two minutes has passed, and all four quadrants are sparkly clean, the brush will automatically shut off. It's that simple.



Connecting & Charging

It's recommended to fully charge your toothbrush prior to first use. A full charge takes approx. 24 hours and can provide enough power for up to two weeks of regular brushing.

- Plug in the adaptor to an electrical outlet.
- Place the Elite Sonic Brush Handle (5) onto the Charging Base (7), so that the Charge Indicator (6) lights up.
- For everyday use, store the Elite Sonic Brush Handle on the Charging Base to maintain its full power and charge.



Cleaning & Recommendations

- After brushing, rinse your Brush Head thoroughly under running water.
- To clean the Handle, remove the Brush Head, then wipe the the Handle clean with a soft cloth. Do not submerge the Handle in water.
- Wipe the Handle dry before placing it onto the Charging Base, and wipe the Brush Head dry before placing it into the UV Chamber.
- To clean the Removable Drip Cup, open the UV Chamber Cover, lift out the Drip Cup, and rinse with water and a mild soap.
- Before cleaning your Charging Base, make sure to unplug the power adaptor. Only use a dry cloth to wipe the surface clean.
- Place the charging station in a safe, dry place.

IMPORTANT:

Never place the UV Sanitizing Charging Base under running water or submerge it into water.



This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of this product in household waste at the end of its useful life. Disposal can take place at appropriate collection points provided in your area.