



Care@Home™

Cellular and Ethernet
Version

C7000 Series Getting Started







The Best Location for the C7000



Good cellular reception

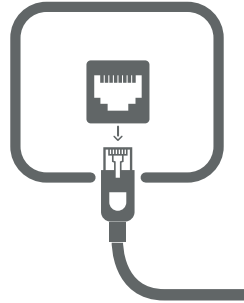

Near network port

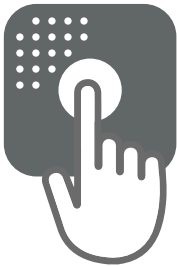

Convenient for resident to access



Where speakerphone can be heard throughout premises




Setting Up the C7000




1 Place thumb in oval hole and lift off the cover



3 Insert the backup battery

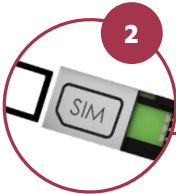

2 Insert the SIM

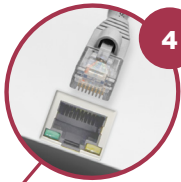

4 Plug in Ethernet cable and attach to network port




6 Turn the power switch on



5 Connect the power cable and plug into a power source. Then close the back cover



Using the C7000

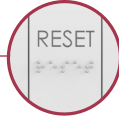

Press to cancel an emergency call



Press to start an emergency call. Red light:

- Blinking: alarm in progress
- On: message acknowledged



Shows power status

- On: AC power connected
- Blinking: using battery
- Off: no power



Shows communication status

- **Green:** if on, ready for call. If blinking, message being sent
- **Red:** if on, communication failure
- **Orange:** active conversation



Shows problem status

- On: malfunction in the C7000 or a peripheral
- Blinking: supervision lost, or low battery in a peripheral

For more information, see: <http://www.essence-grp.com/resources/smart-care-download.html>

Legal notice: http://www.essence-grp.com/html/concerning_documents.html

www.essence-grp.com • info@essence-grp.com

