





Version

C7000 and Emergency Pendant Plus Kit

Getting Started





The Best Location for the C7000



Near a phone jack



Near a network port



Convenient for resident to access







Where speakerphone can be heard throughout premises



Setting Up the C7000

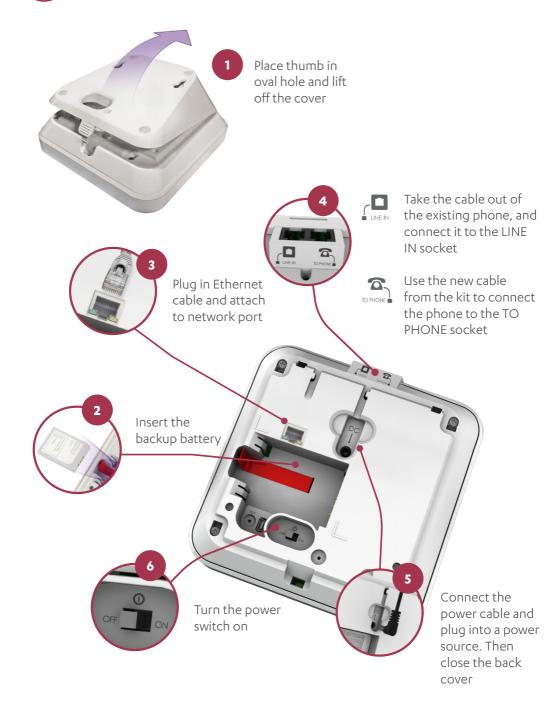

Using the C7000

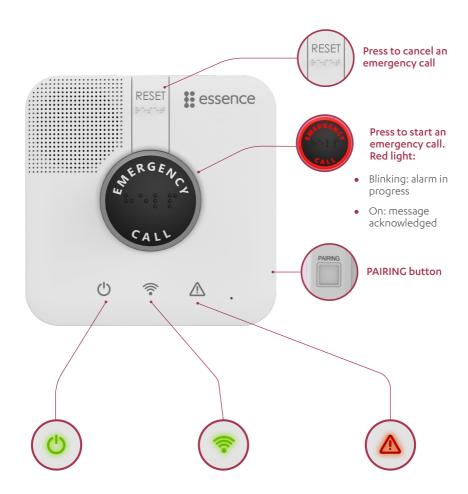

Shows power status

- On: AC power connected
- Blinking: using battery
- Off: no power

Shows communication status

- Green: if on, ready for call.
 If blinking, message being sent
- Red: if on, communication failure
- Orange: active conversation

Shows problem status

- On: malfunction in the C7000 or a peripheral
- Blinking: supervision lost, or low battery in a peripheral



Setting Up the Pendant




Press the PAIRING button for 5 seconds.





Press the button on the pendant.

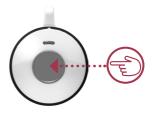

The panel beeps and the EMERGENCY button blinks blue.





The Pendant in Action





The EPA blinks red when the panic button is pressed



The EPA lights green when a fall is detected or the battery is replaced

For more information, see: http://www.essence-grp.com/resources/smart-care-download.html

Legal notice: http://www.essence-grp.com/html/concerning_documents.html

www.essence-grp.com • info@essence-grp.com

