




# Care@Home™

## Control Panel Getting Started





# The Best Location for the Control Panel

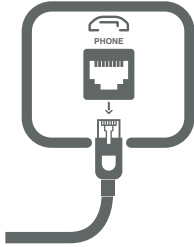

Near a phone jack  
(PSTN version only)



Good cellular reception  
(Cellular version only)



Convenient for  
resident to access



Where speakerphone can be  
heard throughout premises



Control Panel

At least 2 m (6 ft.)  
from other devices



Motion detector



# Setting Up the Control Panel

1 Push here

2 Lift the cover off

4 Insert the backup battery

3 Insert the SIM

Cellular version only

7 Turn the power switch on

PSTN version only



5 Take the cable out of the existing phone, and connect it to the LINE IN socket




Use the new cable from the kit to connect the phone to the TO PHONE socket

6 Connect the power cable and plug into a power source. Then close the back cover



# Using the Control Panel

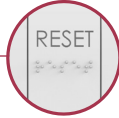

Press to cancel an emergency call



Press to start an emergency call. Red light:

- Blinking: alarm in progress
- On: message acknowledged



## Shows power status

- On: AC power connected
- Blinking: using battery
- Off: no power



## Shows communication status

- **Green:** if on, ready for call. If blinking, message being sent
- **Red:** if on, communication failure
- **Orange:** active conversation



## Shows problem status

- On: malfunction in the control panel or a peripheral
- Blinking: supervision lost, or low battery in a peripheral

For more information, see: <http://www.essence-grp.com/resources/smart-care-download.html>

Legal notice: [http://www.essence-grp.com/html/concerning\\_documents.html](http://www.essence-grp.com/html/concerning_documents.html)

[www.essence-grp.com](http://www.essence-grp.com) • [info@essence-grp.com](mailto:info@essence-grp.com)

