

Multiple Well Plates - 6, 12, 24, 48 & 96 Well



- Available in tissue culture treated or non-treated versions
- 1** Plates are labeled with alpha-numeric markings for easy identification of individual wells
- 2** Each well has a raised rim which nests into the rings on the single direction lid to reduce evaporation and limit the risk of cross contamination
- 3** Lids have a stacking feature for added stability when handling tall stacks of plates

MULTIPLE WELL PLATES - TISSUE CULTURE TREATED S, R/D F, NP, USP, PS

Part No.	Product Description	Bottom Shape	Cell Growth Area	Pack Description	Qty/ Pack	Qty/ Case
229105	6 Well Tissue Culture Plate	Flat	9.60cm ²	Individual		50
229106	6 Well Tissue Culture Plate	Flat	9.60cm ²	Individual	50	100
229111	12 Well Tissue Culture Plate	Flat	3.85cm ²	Individual		50
229112	12 Well Tissue Culture Plate	Flat	3.85cm ²	Individual	50	100
229123	24 Well Tissue Culture Plate	Flat	1.93cm ²	Individual		50
229124	24 Well Tissue Culture Plate	Flat	1.93cm ²	Individual	50	100
229147	48 Well Tissue Culture Plate	Flat	0.84cm ²	Individual		50
229148	48 Well Tissue Culture Plate	Flat	0.84cm ²	Individual	50	100
229195	96 Well Tissue Culture Plate	Flat	0.33cm ²	Individual		50
229196	96 Well Tissue Culture Plate	Flat	0.33cm ²	Individual	50	100
229197	96 Well Tissue Culture Plate	Flat	0.33cm ²	Bulk	5	100
229190	96 Well Tissue Culture Plate	Round	0.33cm ²	Individual	50	100

MULTIPLE WELL PLATES - NON-TREATED S, R/D F, NP, USP, PS

Part No.	Product Description	Bottom Shape	Working Volume	Pack Description	Qty/ Pack	Qty/ Case
229506	6 Well Non-treated Plate	Flat	2 - 8mL	Individual	50	100
229512	12 Well Non-treated Plate	Flat	2 - 4mL	Individual	50	100
229524	24 Well Non-treated Plate	Flat	0.5 - 1.5mL	Individual	50	100
229548	48 Well Non-treated Plate	Flat	0.3 - 1.0mL	Individual	50	100
229590	96 Well Non-treated Plate	Round	50 - 300µL	Individual	50	100
229591	96 Well Non-treated Plate without Lid	Round	50 - 300µL	Individual	50	100
229592	96 Well Non-treated Plate without Lid	Flat	50 - 300µL	Individual	50	100
229596	96 Well Non-treated Plate	Flat	50 - 300µL	Individual	50	100

NEW Multiple Well Plates - 1 & 4 Well


- Available in tissue culture treated or non-treated configurations
- 1 and 4 well plates have all the same standard features as 6 - 96 Multiple Well Plates (see page 24)

1 Well Plates

- Useful for high volume and robotic screening of DNA libraries as well as manual plating
- Same footprint as standard multiple well plates
- Suitable for dot-blotting bacterial clones on 11.4 x 7.3cm (4.5 x 2.875in) nylon membranes

4 Well Plates

- Expanded spacing of wells helps reduce cross contamination
- Well diameter similar to 24 well plate
- Smaller footprint than standard multiple well plates

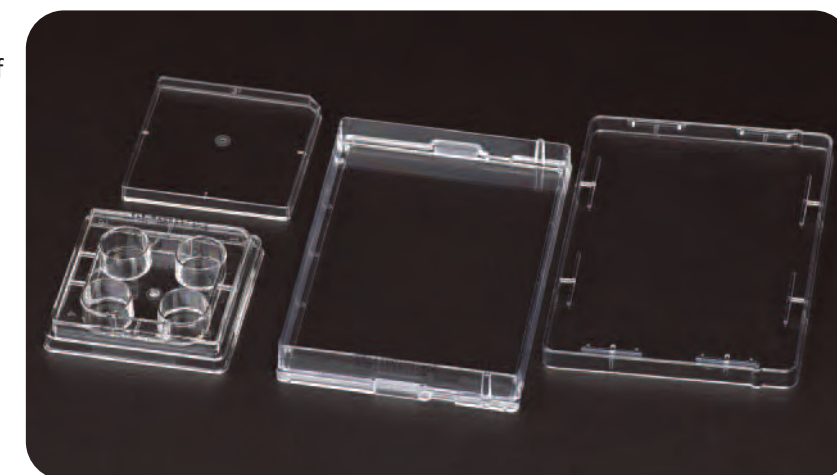


MULTIPLE WELL PLATES - TISSUE CULTURE TREATED S, R/D F, NP, USP, PS

Part No.	Product Description	Cell Growth Area	Pack Description	Qty/Case
229101	1 Well Tissue Culture Plate with Lid	97.4cm ²	Individual	50
229103	4 Well Tissue Culture Plate with Lid	1.96cm ²	Individual	50

MULTIPLE WELL PLATES - NON-TREATED S, R/D F, NP, USP, PS

Part No.	Product Description	Working Volume	Pack Description	Qty/Case
229501	1 Well Non-treated Plate with Lid	19.5 - 48.7mL	Individual	50
229503	4 Well Non-treated Plate with Lid	500 - 800µL	Individual	50

96 Deep Well Storage Plates

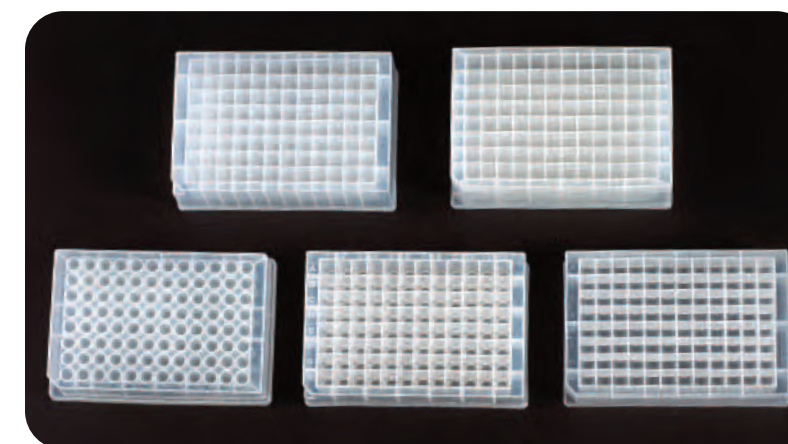

- 96 Deep Well Storage Plates can be used for the growth of bacterial cultures or storage of compounds
- Square to round transition of V-Bottom maximizes working volume
- Pyramid-Bottom eases sample retrieval and improves mixing by creating a sloshing effect rather than a vortex when agitated
- Polypropylene offers a low binding surface to prevent samples from sticking to the side walls during elution, and is chemically inert for combinatorial chemistry applications

96 DEEP WELL STORAGE PLATES N, R/D F, NP, USP, A, PP

	Part No.	Product Description	Well Volume	Well Shape	Bottom Shape	Pack Description	Qty/ Case	Sterile (Y/N)
NEW	229570	96 Deep Well Storage Plate	0.5mL	Round	V-Bottom	5/Bag	25	N
	229571	96 Deep Well Storage Plate	1.1mL	Square	V-Bottom	5/Bag	25	N
NEW	229573	96 Deep Well Storage Plate	1.5mL	Square	Round Bottom	5/Bag	25	N
	229572	96 Deep Well Storage Plate	2.0mL	Square	Pyramid Bottom	5/Bag	25	N
NEW	229574	96 Deep Well Storage Plate	2.0mL	Square	V-Bottom	5/Bag	25	N