

Guideline for Disinfection and Sterilization in Healthcare Facilities, 2008

William A. Rutala, Ph.D., M.P.H.^{1,2}, David J. Weber, M.D., M.P.H.^{1,2}, and the Healthcare Infection Control Practices Advisory Committee (HICPAC)³

¹Hospital Epidemiology University of North Carolina Health Care System Chapel Hill, NC 27514

²Division of Infectious Diseases University of North Carolina School of Medicine Chapel Hill, NC 27599-7030

³HICPAC Members

Robert A. Weinstein, MD (Chair) Cook County Hospital Chicago, IL

Jane D. Siegel, MD (Co-Chair) University of Texas Southwestern Medical Center Dallas, TX

Michele L. Pearson, MD (Executive Secretary) Centers for Disease Control and Prevention Atlanta, GA

Raymond Y.W. Chinn, MD Sharp Memorial Hospital San Diego, CA

Alfred DeMaria, Jr, MD Massachusetts Department of Public Health Jamaica Plain, MA

James T. Lee, MD, PhD University of Minnesota Minneapolis, MN

William A. Rutala, PhD, MPH University of North Carolina Health Care System Chapel Hill, NC

William E. Scheckler, MD University of Wisconsin Madison, WI

Beth H. Stover, RN Kosair Children's Hospital Louisville, KY

Marjorie A. Underwood, RN, BSN CIC Mt. Diablo Medical Center Concord, CA

This guideline discusses use of products by healthcare personnel in healthcare settings such as hospitals, ambulatory care and home care; the recommendations are not intended for consumer use of the products discussed.

Table	of	Contents
Iabic	UI.	Contents

Executive Summary	7
Introduction	8
Methods	8
Definition of Terms	8
A Rational Approach to Disinfection and Sterilization	. 10
Critical Items	. 10
Semicritical Items	. 10
Noncritical Items	. 11
Changes in Disinfection and Sterilization Since 1981	. 12
Disinfection of Healthcare Equipment	. 13
Concerns about Implementing the Spaulding Scheme	. 13
Reprocessing of Endoscopes	. 13
Laparoscopes and Arthroscopes	. 17
Tonometers, Cervical Diaphragm Fitting Rings, Cryosurgical Instruments, and Endocavitary Probes	18
Dental Instruments	. 20
Disinfection of HBV-, HCV-, HIV- or TB-Contaminated Devices	. 21
Disinfection in the Hemodialysis Unit	
Inactivation of Clostridium difficile	
OSHA Bloodborne Pathogen Standard	. 22
Emerging Pathogens (Cryptosporidium, Helicobacter pylori, Escherichia coli O157:H7,	
Rotavirus, Human Papilloma Virus, Norovirus, Severe Acute Respiratory Syndrome [SARS] Coronavirus)	. 23
Inactivation of Bioterrorist Agents	
Toxicological, Environmental and Occupational Concerns	
Disinfection in Ambulatory Care, Home Care, and the Home	
Susceptibility of Antibiotic-Resistant Bacteria to Disinfectants	
Surface Disinfection	
Air Disinfection	
Microbial Contamination of Disinfectants	
Factors Affecting the Efficacy of Disinfection and Sterilization	
Number and Location of Microorganisms	
Innate Resistance of Microorganisms	
Concentration and Potency of Disinfectants	
Physical and Chemical Factors	
Organic and Inorganic Matter	
Duration of Exposure	
Biofilms	
Cleaning	
Disinfection	
Chemical Disinfectants	
Alcohol	
Overview	
Mode of Action	
Microbicidal Activity.	
Uses	
Chlorine and Chlorine Compounds	
Overview	

Mode of Action	
Microbicidal Activity.	
Uses	
Formaldehyde	
Overview	
Mode of Action	
Microbicidal Activity.	
Uses	
Glutaraldehyde	
Overview	
Mode of Action	
Microbicidal Activity	
Uses	
Hydrogen Peroxide	
Overview	
Mode of Action	
Microbicidal Activity.	
Uses	
lodophors	
Overview	
Mode of Action	
Microbicidal Activity	
Uses	
Ortho-phthalaldehyde (OPA)	
Overview	
Mode of Action	
Microbicidal Activity.	
Uses	
Peracetic Acid	
Overview	
Mode of Action	
Microbicidal Activity.	
•	
Peracetic Acid and Hydrogen Peroxide	
Overview	
Microbicidal Activity	
Uses	
Phenolics	51
Overview	
Mode of Action	
Microbicidal Activity.	
Uses	
Quaternary Ammonium Compounds	
Overview.	
Mode of Action	
Microbicidal Activity.	
Uses	
scellaneous Inactivating Agents	
Other Germicides	

Metals as Microbicides	54
Ultraviolet Radiation (UV)	54
Pasteurization	55
Flushing- and Washer-Disinfectors	55
The Regulatory Framework for Disinfectants and Sterilants	56
EPA and FDA	56
CDC	56
Test Methods	57
Neutralization of Germicides	57
Sterilization	58
Steam Sterilization	58
Overview	58
Microbicidal Activity	60
Mode of Action.	60
Uses	60
Flash Sterilization	60
Overview	60
Uses	61
Low-Temperature Sterilization Technologies	61
Ethylene Oxide "Gas" Sterilization	62
Overview	62
Mode of Action.	63
Microbicidal Activity	63
Uses	63
Hydrogen Peroxide Gas Plasma	
Overview.	63
Mode of Action	
Microbicidal Activity	
Uses	
Peracetic Acid Sterilization	
Overview.	
Mode of Action.	
Microbicidal Activity	
Uses	
Microbicidal Activity of Low-Temperature Sterilization Technologies	
Bioburden of Surgical Devices	
Effect of Cleaning on Sterilization Efficacy	
Other Sterilization Methods	
Ionizing Radiation	
Dry-Heat Sterilizers	
Liquid Chemicals	
Performic Acid	
Filtration.	
Microwave	
Glass Bead "Sterilizer"	
Vaporized Hydrogen Peroxide (VHP®).	
Ozone	
Formaldehyde Steam	
Gaseous chlorine dioxide	72

Vaµ	porized Peracetic Acid	72	
Infr	Infrared radiation		
Sterilizi	ng Practices	72	
Ove	erview	72	
Ste	rilization Cycle Verification	73	
Phy	vsical Facilities	73	
Cle	aning	73	
Pac	skaging	74	
Loa	nding	75	
Sto	rage	75	
	nitoring		
	ngle-Use Medical Devices		
Conclusion.		81	
Web-Based	Disinfection and Sterilization Resources	82	
	dations for Disinfection and Sterilization in Healthcare Facilities		
A. Ratio	nale	83	
	sings		
C. Reco	ommendations	83	
1.	Occupational Health and Exposure		
2.	Cleaning of Patient-Care Devices	83	
3.	Indications for Sterilization, High-Level Disinfection, and Low-Level Disinfection		
4.	Selection and Use of Low-Level Disinfectants for Noncritical Patient-Care Devices		
5.	Cleaning and Disinfecting Environmental Surfaces in Healthcare Facilities	85	
6.	Disinfectant Fogging	86	
7.	High-Level Disinfection of Endoscopes	86	
8.	Management of Equipment and Surfaces in Dentistry	89	
9.	Processing Patient-Care Equipment Contaminated with Bloodborne Pathogens (HBV, Hepatitis C Virus, HIV), Antibiotic-Resistant Bacteria (e.g., Vancomycin-		
	Resistant Enterococci, Methicillin-Resistant Staphylococcus aureus, Multidrug		
	Resistant Tuberculosis), or Emerging Pathogens (e.g., Cryptosporidium, Helicobacter		
	pylori, Escherichia coli O157:H7, Clostridium difficile, Mycobacterium tuberculosis,	00	
10	Severe Acute Respiratory Syndrome Coronavirus), or Bioterrorist Agents		
	Disinfection Strategies for Other Semicritical Devices		
	Disinfection by Healthcare Personnel in Ambulatory Care and Home Care		
	Microbial Contamination of Disinfectants		
	Flash Sterilization		
	Methods of Sterilization		
	Packaging		
	Monitoring of Sterilizers		
	Load Configuration.		
	Storage of Sterile Items		
	Quality Control		
	Reuse of Single-Use Medical Devices		
	e Indicators		
	gements		
-	Figure		
	of Financial Interests and Relationships (2000- July 2004)		
17010101062		. 120	

Executive Summary

The Guideline for Disinfection and Sterilization in Healthcare Facilities, 2008, presents evidencebased recommendations on the preferred methods for cleaning, disinfection and sterilization of patientcare medical devices and for cleaning and disinfecting the healthcare environment. This document supercedes the relevant sections contained in the 1985 Centers for Disease Control (CDC) Guideline for Handwashing and Environmental Control.¹ Because maximum effectiveness from disinfection and sterilization results from first cleaning and removing organic and inorganic materials, this document also reviews cleaning methods. The chemical disinfectants discussed for patient-care equipment include alcohols, glutaraldehyde, formaldehyde, hydrogen peroxide, iodophors, ortho-phthalaldehyde, peracetic acid, phenolics, quaternary ammonium compounds, and chlorine. The choice of disinfectant, concentration, and exposure time is based on the risk for infection associated with use of the equipment and other factors discussed in this guideline. The sterilization methods discussed include steam sterilization, ethylene oxide (ETO), hydrogen peroxide gas plasma, and liquid peracetic acid. When properly used, these cleaning, disinfection, and sterilization processes can reduce the risk for infection associated with use of invasive and noninvasive medical and surgical devices. However, for these processes to be effective, health-care workers should adhere strictly to the cleaning, disinfection, and sterilization recommendations in this document and to instructions on product labels.

In addition to updated recommendations, new topics addressed in this guideline include

- 1. inactivation of antibiotic-resistant bacteria, bioterrorist agents, emerging pathogens, and bloodborne pathogens;
- toxicologic, environmental, and occupational concerns associated with disinfection and sterilization practices;
- 3. disinfection of patient-care equipment used in ambulatory settings and home care;
- 4. new sterilization processes, such as hydrogen peroxide gas plasma and liquid peracetic acid; and
- 5. disinfection of complex medical instruments (e.g., endoscopes).

Introduction

In the United States, approximately 46.5 million surgical procedures and even more invasive medical procedures—including approximately 5 million gastrointestinal endoscopies—are performed each year.² Each procedure involves contact by a medical device or surgical instrument with a patient's sterile tissue or mucous membranes. A major risk of all such procedures is the introduction of pathogens that can lead to infection. Failure to properly disinfect or sterilize equipment carries not only risk associated with breach of host barriers but also risk for person-to-person transmission (e.g., hepatitis B virus) and transmission of environmental pathogens (e.g., *Pseudomonas aeruginosa*).

Disinfection and sterilization are essential for ensuring that medical and surgical instruments do not transmit infectious pathogens to patients. Because sterilization of all patient-care items is not necessary, health-care policies must identify, primarily on the basis of the items' intended use, whether cleaning, disinfection, or sterilization is indicated.

Multiple studies in many countries have documented lack of compliance with established guidelines for disinfection and sterilization. ³⁻⁶ Failure to comply with scientifically-based guidelines has led to numerous outbreaks. ⁶⁻¹² This guideline presents a pragmatic approach to the judicious selection and proper use of disinfection and sterilization processes; the approach is based on well-designed studies assessing the efficacy (through laboratory investigations) and effectiveness (through clinical studies) of disinfection and sterilization procedures.

Methods

This guideline resulted from a review of all MEDLINE articles in English listed under the MeSH headings of *disinfection* or *sterilization* (focusing on health-care equipment and supplies) from January 1980 through August 2006. References listed in these articles also were reviewed. Selected articles published before 1980 were reviewed and, if still relevant, included in the guideline. The three major peer-reviewed journals in infection control—*American Journal of Infection Control, Infection Control and Hospital Epidemiology,* and *Journal of Hospital Infection*—were searched for relevant articles published from January 1990 through August 2006. Abstracts presented at the annual meetings of the Society for Healthcare Epidemiology of America and Association for professionals in Infection Control and Epidemiology, Inc. during 1997–2006 also were reviewed; however, abstracts were not used to support the recommendations.

Definition of Terms

Sterilization describes a process that destroys or eliminates all forms of microbial life and is carried out in health-care facilities by physical or chemical methods. Steam under pressure, dry heat, EtO gas, hydrogen peroxide gas plasma, and liquid chemicals are the principal sterilizing agents used in health-care facilities. Sterilization is intended to convey an absolute meaning; unfortunately, however, some health professionals and the technical and commercial literature refer to "disinfection" as "sterilization" and items as "partially sterile." When chemicals are used to destroy all forms of microbiologic life, they can be called chemical sterilants. These same germicides used for shorter exposure periods also can be part of the disinfection process (i.e., high-level disinfection).

Disinfection describes a process that eliminates many or all pathogenic microorganisms, except bacterial spores, on inanimate objects (Tables 1 and 2). In health-care settings, objects usually are disinfected by liquid chemicals or wet pasteurization. Each of the various factors that affect the efficacy of

disinfection can nullify or limit the efficacy of the process.

Factors that affect the efficacy of both disinfection and sterilization include prior cleaning of the object; organic and inorganic load present; type and level of microbial contamination; concentration of and exposure time to the germicide; physical nature of the object (e.g., crevices, hinges, and lumens); presence of biofilms; temperature and pH of the disinfection process; and in some cases, relative humidity of the sterilization process (e.g., ethylene oxide).

Unlike sterilization, disinfection is not sporicidal. A few disinfectants will kill spores with prolonged exposure times (3–12 hours); these are called *chemical sterilants*. At similar concentrations but with shorter exposure periods (e.g., 20 minutes for 2% glutaraldehyde), these same disinfectants will kill all microorganisms except large numbers of bacterial spores; they are called *high-level disinfectants*. *Low-level disinfectants* can kill most vegetative bacteria, some fungi, and some viruses in a practical period of time (≤10 minutes). *Intermediate-level disinfectants* might be cidal for mycobacteria, vegetative bacteria, most viruses, and most fungi but do not necessarily kill bacterial spores. Germicides differ markedly, primarily in their antimicrobial spectrum and rapidity of action.

Cleaning is the removal of visible soil (e.g., organic and inorganic material) from objects and surfaces and normally is accomplished manually or mechanically using water with detergents or enzymatic products. Thorough cleaning is essential before high-level disinfection and sterilization because inorganic and organic materials that remain on the surfaces of instruments interfere with the effectiveness of these processes. *Decontamination* removes pathogenic microorganisms from objects so they are safe to handle, use, or discard.

Terms with the suffix *cide* or *cidal* for killing action also are commonly used. For example, a germicide is an agent that can kill microorganisms, particularly pathogenic organisms ("germs"). The term *germicide* includes both antiseptics and disinfectants. *Antiseptics* are germicides applied to living tissue and skin; *disinfectants* are antimicrobials applied only to inanimate objects. In general, antiseptics are used only on the skin and not for surface disinfection, and disinfectants are not used for skin antisepsis because they can injure skin and other tissues. Virucide, fungicide, bactericide, sporicide, and tuberculocide can kill the type of microorganism identified by the prefix. For example, a bactericide is an agent that kills bacteria. ¹³⁻¹⁸

A Rational Approach to Disinfection and Sterilization

More than 30 years ago, Earle H. Spaulding devised a rational approach to disinfection and sterilization of patient-care items and equipment.¹⁴ This classification scheme is so clear and logical that it has been retained, refined, and successfully used by infection control professionals and others when planning methods for disinfection or sterilization. ^{1, 13, 15, 17, 19, 20} Spaulding believed the nature of disinfection could be understood readily if instruments and items for patient care were categorized as critical, semicritical, and noncritical according to the degree of risk for infection involved in use of the items. The CDC *Guideline for Handwashing and Hospital Environmental Control* ²¹, *Guidelines for the Prevention of Transmission of Human Immunodeficiency Virus (HIV) and Hepatitis B Virus (HBV) to Health-Care and Public-Safety Workers*²², and *Guideline for Environmental Infection Control in Health-Care Facilities*²³ employ this terminology.

Critical Items

Critical items confer a high risk for infection if they are contaminated with any microorganism. Thus, objects that enter sterile tissue or the vascular system must be sterile because any microbial contamination could transmit disease. This category includes surgical instruments, cardiac and urinary catheters, implants, and ultrasound probes used in sterile body cavities. Most of the items in this category should be purchased as sterile or be sterilized with steam if possible. Heat-sensitive objects can be treated with EtO, hydrogen peroxide gas plasma; or if other methods are unsuitable, by liquid chemical sterilants. Germicides categorized as chemical sterilants include ≥2.4% glutaraldehyde-based formulations, 0.95% glutaraldehyde with 1.64% phenol/phenate, 7.5% stabilized hydrogen peroxide, 7.35% hydrogen peroxide with 0.23% peracetic acid, 0.2% peracetic acid, and 0.08% peracetic acid with 1.0% hydrogen peroxide. Liquid chemical sterilants reliably produce sterility only if cleaning precedes treatment and if proper guidelines are followed regarding concentration, contact time, temperature, and pH.

Semicritical Items

Semicritical items contact mucous membranes or nonintact skin. This category includes respiratory therapy and anesthesia equipment, some endoscopes, laryngoscope blades ²⁴, esophageal manometry probes, cystoscopes ²⁵, anorectal manometry catheters, and diaphragm fitting rings. These medical devices should be free from all microorganisms; however, small numbers of bacterial spores are permissible. Intact mucous membranes, such as those of the lungs and the gastrointestinal tract, generally are resistant to infection by common bacterial spores but susceptible to other organisms, such as bacteria, mycobacteria, and viruses. Semicritical items minimally require high-level disinfection using chemical disinfectants. Glutaraldehyde, hydrogen peroxide, *ortho*-phthalaldehyde, and peracetic acid with hydrogen peroxide are cleared by the Food and Drug Administration (FDA) and are dependable high-level disinfectants provided the factors influencing germicidal procedures are met (Table 1). When a disinfectant is selected for use with certain patient-care items, the chemical compatibility after extended use with the items to be disinfected also must be considered.

High-level disinfection traditionally is defined as complete elimination of all microorganisms in or on an instrument, except for small numbers of bacterial spores. The FDA definition of high-level disinfection is a sterilant used for a shorter contact time to achieve a 6-log₁₀ kill of an appropriate *Mycobacterium* species. Cleaning followed by high-level disinfection should eliminate enough pathogens to prevent transmission of infection. ^{26, 27}

Laparoscopes and arthroscopes entering sterile tissue ideally should be sterilized between patients. However, in the United States, this equipment sometimes undergoes only high-level disinfection between patients. ²⁸⁻³⁰ As with flexible endoscopes, these devices can be difficult to clean and high-level disinfect or sterilize because of intricate device design (e.g., long narrow lumens, hinges). Meticulous cleaning must precede any high-level disinfection or sterilization process. Although sterilization is preferred, no reports have been published of outbreaks resulting from high-level disinfection of these scopes when they are properly cleaned and high-level disinfected. Newer models of these instruments can withstand steam sterilization that for critical items would be preferable to high-level disinfection.

Rinsing endoscopes and flushing channels with sterile water, filtered water, or tap water will prevent adverse effects associated with disinfectant retained in the endoscope (e.g., disinfectant-induced colitis). Items can be rinsed and flushed using sterile water after high-level disinfection to prevent contamination with organisms in tap water, such as nontuberculous mycobacteria, ^{10, 31, 32} *Legionella*, ³³⁻³⁵ or gramnegative bacilli such as *Pseudomonas*. ^{1, 17, 36-38} Alternatively, a tapwater or filtered water (0.2µ filter) rinse should be followed by an alcohol rinse and forced air drying. ^{28, 38-40} Forced-air drying markedly reduces bacterial contamination of stored endoscopes, most likely by removing the wet environment favorable for bacterial growth. ³⁹ After rinsing, items should be dried and stored (e.g., packaged) in a manner that protects them from recontamination.

Some items that may come in contact with nonintact skin for a brief period of time (i.e., hydrotherapy tanks, bed side rails) are usually considered noncritical surfaces and are disinfected with intermediate-level disinfectants (i.e., phenolic, iodophor, alcohol, chlorine)²³. Since hydrotherapy tanks have been associated with spread of infection, some facilities have chosen to disinfect them with recommended levels of chlorine ^{23, 41}.

In the past, high-level disinfection was recommended for mouthpieces and spirometry tubing (e.g., glutaraldehyde) but cleaning the interior surfaces of the spirometers was considered unnecessary. ⁴² This was based on a study that showed that mouthpieces and spirometry tubing become contaminated with microorganisms but there was no bacterial contamination of the surfaces inside the spirometers. Filters have been used to prevent contamination of this equipment distal to the filter; such filters and the proximal mouthpiece are changed between patients.

Noncritical Items

Noncritical items are those that come in contact with intact skin but not mucous membranes. Intact skin acts as an effective barrier to most microorganisms; therefore, the sterility of items coming in contact with intact skin is "not critical." In this guideline, noncritical items are divided into noncritical patient care items and noncritical environmental surfaces ^{43, 44}. Examples of noncritical patient-care items are bedpans, blood pressure cuffs, crutches and computers ⁴⁵. In contrast to critical and some semicritical items, most noncritical reusable items may be decontaminated where they are used and do not need to be transported to a central processing area. Virtually no risk has been documented for transmission of infectious agents to patients through noncritical items ³⁷ when they are used as noncritical items and do not contact non-intact skin and/or mucous membranes. Table 1 lists several low-level disinfectants that may be used for noncritical items. Most Environmental Protection Agency (EPA)-registered disinfectants have a 10-minute label claim. However, multiple investigators have demonstrated the effectiveness of these disinfectants against vegetative bacteria (e.g., Listeria, Escherichia coli, Salmonella, vancomycin-resistant Enterococci, methicillin-resistant Staphylococcus aureus), yeasts (e.g., Candida), mycobacteria (e.g., Mycobacterium tuberculosis), and viruses (e.g. poliovirus) at exposure times of 30-60 seconds⁴⁶⁻⁶⁴ Federal law requires all applicable label instructions on EPA-registered products to be followed (e.g., use-dilution, shelf life, storage, material compatibility, safe use, and disposal). If the user selects exposure conditions (e.g., exposure time) that differ from those on the EPA-registered products label, the user assumes liability for any injuries resulting from offlabel use and is potentially subject to enforcement action under Federal Insecticide, Fungicide, and Rodenticide Act (FIFRA) 65.

Noncritcal environmental surfaces include bed rails, some food utensils, bedside tables, patient furniture and floors. Noncritical environmental surfaces frequently touched by hand (e.g., bedside tables,

bed rails) potentially could contribute to secondary transmission by contaminating hands of health-care workers or by contacting medical equipment that subsequently contacts patients ^{13, 46-48, 51, 66, 67}. Mops and reusable cleaning cloths are regularly used to achieve low-level disinfection on environmental surfaces. However, they often are not adequately cleaned and disinfected, and if the water-disinfectant mixture is not changed regularly (e.g., after every three to four rooms, at no longer than 60-minute intervals), the mopping procedure actually can spread heavy microbial contamination throughout the health-care facility ⁶⁸. In one study, standard laundering provided acceptable decontamination of heavily contaminated mopheads but chemical disinfection with a phenolic was less effective. ⁶⁸ Frequent laundering of mops (e.g., daily), therefore, is recommended. Single-use disposable towels impregnated with a disinfectant also can be used for low-level disinfection when spot-cleaning of noncritical surfaces is needed⁴⁵.

Changes in Disinfection and Sterilization Since 1981

The Table in the CDC Guideline for Environmental Control prepared in 1981 as a guide to the appropriate selection and use of disinfectants has undergone several important changes (Table 1).¹⁵ First, formaldehyde-alcohol has been deleted as a recommended chemical sterilant or high-level disinfectant because it is irritating and toxic and not commonly used. Second, several new chemical sterilants have been added, including hydrogen peroxide, peracetic acid ^{58, 69, 70}, and peracetic acid and hydrogen peroxide in combination. Third, 3% phenolics and iodophors have been deleted as high-level disinfectants because of their unproven efficacy against bacterial spores, M. tuberculosis, and/or some fungi. ^{55, 71} Fourth, isopropyl alcohol and ethyl alcohol have been excluded as high-level disinfectants ¹⁵ because of their inability to inactivate bacterial spores and because of the inability of isopropyl alcohol to inactivate hydrophilic viruses (i.e., poliovirus, coxsackie virus).72 Fifth, a 1:16 dilution of 2.0% glutaraldehyde-7.05% phenol-1.20% sodium phenate (which contained 0.125% glutaraldehyde, 0.440% phenol, and 0.075% sodium phenate when diluted) has been deleted as a high-level disinfectant because this product was removed from the marketplace in December 1991 because of a lack of bactericidal activity in the presence of organic matter; a lack of fungicidal, tuberculocidal and sporicidal activity; and reduced virucidal activity.^{49, 55, 56, 71, 73-79} Sixth, the exposure time required to achieve high-level disinfection has been changed from 10-30 minutes to 12 minutes or more depending on the FDA-cleared label claim and the scientific literature. 27, 55, 69, 76, 80-84 A glutaraldehyde and an ortho-phthalaldehyde have an FDA-cleared label claim of 5 minutes when used at 35°C and 25°C, respectively, in an automated endoscope reprocessor with FDA-cleared capability to maintain the solution at the appropriate temperature. 85

In addition, many new subjects have been added to the guideline. These include inactivation of emerging pathogens, bioterrorist agents, and bloodborne pathogens; toxicologic, environmental, and occupational concerns associated with disinfection and sterilization practices; disinfection of patient-care equipment used in ambulatory and home care; inactivation of antibiotic-resistant bacteria; new sterilization processes, such as hydrogen peroxide gas plasma and liquid peracetic acid; and disinfection of complex medical instruments (e.g., endoscopes).

Disinfection of Healthcare Equipment

Concerns about Implementing the Spaulding Scheme

One problem with implementing the aforementioned scheme is oversimplification. For example, the scheme does not consider problems with reprocessing of complicated medical equipment that often is heat-sensitive or problems of inactivating certain types of infectious agents (e.g., prions, such as Creutzfeldt-Jakob disease [CJD] agent). Thus, in some situations, choosing a method of disinfection remains difficult, even after consideration of the categories of risk to patients. This is true particularly for a few medical devices (e.g., arthroscopes, laparoscopes) in the critical category because of controversy about whether they should be sterilized or high-level disinfected. ^{28, 86} Heat-stable scopes (e.g., many rigid scopes) should be steam sterilized. Some of these items cannot be steam sterilized because they are heat-sensitive; additionally, sterilization using ethylene oxide (EtO) can be too time-consuming for routine use between patients (new technologies, such as hydrogen peroxide gas plasma and peracetic acid reprocessor, provide faster cycle times). However, evidence that sterilization of these items improves patient care by reducing the infection risk is lacking^{29, 87-91}. Many newer models of these instruments can withstand steam sterilization, which for critical items is the preferred method.

Another problem with implementing the Spaulding scheme is processing of an instrument in the semicritical category (e.g., endoscope) that would be used in conjunction with a critical instrument that contacts sterile body tissues. For example, is an endoscope used for upper gastrointestinal tract investigation still a semicritical item when used with sterile biopsy forceps or in a patient who is bleeding heavily from esophageal varices? Provided that high-level disinfection is achieved, and all microorganisms except bacterial spores have been removed from the endoscope, the device should not represent an infection risk and should remain in the semicritical category ⁹²⁻⁹⁴. Infection with spore-forming bacteria has not been reported from appropriately high-level disinfected endoscopes.

An additional problem with implementation of the Spaulding system is that the optimal contact time for high-level disinfection has not been defined or varies among professional organizations, resulting in different strategies for disinfecting different types of semicritical items (e.g., endoscopes, applanation tonometers, endocavitary transducers, cryosurgical instruments, and diaphragm fitting rings). Until simpler and effective alternatives are identified for device disinfection in clinical settings, following this guideline, other CDC guidelines ^{1, 22, 95, 96} and FDA-cleared instructions for the liquid chemical sterilants/high-level disinfectants would be prudent.

Reprocessing of Endoscopes

Physicians use endoscopes to diagnose and treat numerous medical disorders. Even though endoscopes represent a valuable diagnostic and therapeutic tool in modern medicine and the incidence of infection associated with their use reportedly is very low (about 1 in 1.8 million procedures) ⁹⁷, more healthcare–associated outbreaks have been linked to contaminated endoscopes than to any other medical device ^{6-8, 12, 98}. To prevent the spread of health-care–associated infections, all heat-sensitive endoscopes (e.g., gastrointestinal endoscopes, bronchoscopes, nasopharygoscopes) must be properly cleaned and, at a minimum, subjected to high-level disinfection after each use. High-level disinfection can be expected to destroy all microorganisms, although when high numbers of bacterial spores are present, a few spores might survive.

Because of the types of body cavities they enter, flexible endoscopes acquire high levels of microbial contamination (bioburden) during each use ⁹⁹. For example, the bioburden found on flexible gastrointestinal endoscopes after use has ranged from 10⁵ colony forming units (CFU)/mL to 10¹⁰ CFU/mL, with the highest levels found in the suction channels ⁹⁹⁻¹⁰². The average load on bronchoscopes before cleaning was 6.4x10⁴ CFU/mL. Cleaning reduces the level of microbial contamination by 4–6 log₁₀ ^{83, 103}. Using human immunovirus (HIV)-contaminated endoscopes, several investigators have shown that cleaning completely eliminates the microbial contamination on the scopes ^{104, 105}. Similarly, other investigators found that EtO sterilization or soaking in 2% glutaraldehyde for 20 minutes was effective only when the device first was properly cleaned ¹⁰⁶.

FDA maintains a list of cleared liquid chemical sterilants and high-level disinfectants that can be used to reprocess heat-sensitive medical devices, such as flexible endoscopes [This link is no longer active: http://www.fda.gov/cdrh/ode/germlab.html. The current version of this document may differ from original version: FDA-Cleared Sterilants and High Level Disinfectants with General Claims for Processing Reusable Medical and Dental Devices - March 2015]

(http://www.fda.gov/MedicalDevices/DeviceRegulationandGuidance/ReprocessingofReusableMedicalDevi ces/ucm437347.htm). At this time, the FDA-cleared and marketed formulations include: ≥2.4% glutaraldehyde, 0.55% ortho-phthalaldehyde (OPA), 0.95% glutaraldehyde with 1.64% phenol/phenate, 7.35% hydrogen peroxide with 0.23% peracetic acid, 1.0% hydrogen peroxide with 0.08% peracetic acid, and 7.5% hydrogen peroxide ⁸⁵. These products have excellent antimicrobial activity; however, some oxidizing chemicals (e.g., 7.5% hydrogen peroxide, and 1.0% hydrogen peroxide with 0.08% peracetic acid [latter product is no longer marketed]) reportedly have caused cosmetic and functional damage to endoscopes 69. Users should check with device manufacturers for information about germicide compatibility with their device. If the germicide is FDA-cleared, then it is safe when used according to label directions; however, professionals should review the scientific literature for newly available data regarding human safety or materials compatibility. EtO sterilization of flexible endoscopes is infrequent because it requires a lengthy processing and aeration time (e.g., 12 hours) and is a potential hazard to staff and patients. The two products most commonly used for reprocessing endoscopes in the United States are glutaraldehyde and an automated, liquid chemical sterilization process that uses peracetic acid ¹⁰⁷. The American Society for Gastrointestinal Endoscopy (ASGE) recommends glutaraldehyde solutions that do not contain surfactants because the soapy residues of surfactants are difficult to remove during rinsing ¹⁰⁸. ortho-phthalaldehyde has begun to replace glutaraldehyde in many health-care facilities because it has several potential advantages over glutaraldehyde: is not known to irritate the eyes and nasal passages, does not require activation or exposure monitoring, and has a 12-minute high-level disinfection claim in the United States ⁶⁹. Disinfectants that are not FDA-cleared and should not be used for reprocessing endoscopes include iodophors, chlorine solutions, alcohols, guaternary ammonium compounds, and phenolics. These solutions might still be in use outside the United States, but their use should be strongly discouraged because of lack of proven efficacy against all microorganisms or materials incompatibility.

FDA clearance of the contact conditions listed on germicide labeling is based on the manufacturer's test results [This link is no longer active: http://www.fda.gov/cdrh/ode/germlab.html. The current version of this document may differ from original version: FDA-Cleared Sterilants and High Level Disinfectants with General Claims for Processing Reusable Medical and Dental Devices - March 2015 (http://www.fda.gov/MedicalDevices/DeviceRegulationandGuidance/ReprocessingofReusableMedicalDevi ces/ucm437347.htm)]. Manufacturers test the product under worst-case conditions for germicide formulation (i.e., minimum recommended concentration of the active ingredient), and include organic soil. Typically manufacturers use 5% serum as the organic soil and hard water as examples of organic and inorganic challenges. The soil represents the organic loading to which the device is exposed during actual use and that would remain on the device in the absence of cleaning. This method ensures that the contact conditions completely eliminate the test mycobacteria (e.g., 10⁵ to 10⁶ Mycobacteria tuberculosis in organic soil and dried on a scope) if inoculated in the most difficult areas for the disinfectant to penetrate and contact in the absence of cleaning and thus provides a margin of safety ¹⁰⁹. For 2.4% glutaraldehyde that requires a 45-minute immersion at 25°C to achieve high-level disinfection (i.e., 100% kill of M. tuberculosis). FDA itself does not conduct testing but relies solely on the disinfectant manufacturer's data. Data suggest that *M. tuberculosis* levels can be reduced by at least 8 log₁₀ with cleaning (4 log₁₀) ^{83, 101, 102,} ¹¹⁰, followed by chemical disinfection for 20 minutes at 20°C (4 to 6 log₁₀) ^{83, 93, 111, 112}. On the basis of these data, APIC ¹¹³, the Society of Gastroenterology Nurses and Associates (SGNA) ^{38, 114, 115}, the ASGE ¹⁰⁸. American College of Chest Physicians ¹², and a multi-society guideline ¹¹⁶ recommend alternative contact conditions with 2% glutaraldehyde to achieve high-level disinfection (e.g., that equipment be immersed in 2% glutaraldehyde at 20°C for at least 20 minutes for high-level disinfection). Federal regulations are to follow the FDA-cleared label claim for high-level disinfectants. The FDA-cleared labels for high-level disinfection with >2% glutaraldehvde at 25°C range from 20-90 minutes, depending upon the product based on three tier testing which includes AOAC sporicidal tests, simulated use testing with mycobacterial and in-use testing. The studies supporting the efficacy of >2% glutaraldehyde for 20 minutes at 20°C assume adequate cleaning prior to disinfection, whereas the FDA-cleared label claim incorporates an

added margin of safety to accommodate possible lapses in cleaning practices. Facilities that have chosen to apply the 20 minute duration at 20°C have done so based on the IA recommendation in the July 2003 SHEA position paper, "Multi-society Guideline for Reprocessing Flexible Gastrointestinal Endoscopes"^{19, 57, 83, 94, 108, 111, 116-121}.

Flexible GI Endoscope Reprocessing Update [June 2011]: <u>Multisociety guideline on reprocessing</u> <u>flexible gastrointestinal endoscopes: 2011</u> (<u>http://www.asge.org/uploadedFiles/Publications_and_Products/Practice_Guidelines/Multisociety%20</u> guideline%20on%20reprocessing%20flexible%20gastrointestinal.pdf).

Flexible endoscopes are particularly difficult to disinfect ¹²² and easy to damage because of their intricate design and delicate materials. ¹²³ Meticulous cleaning must precede any sterilization or high-level disinfection of these instruments. Failure to perform good cleaning can result in sterilization or disinfection failure, and outbreaks of infection can occur. Several studies have demonstrated the importance of cleaning in experimental studies with the duck hepatitis B virus (HBV) ^{106, 124}, HIV ¹²⁵and *Helicobacter pylori*. ¹²⁶

An examination of health-care-associated infections related only to endoscopes through July 1992 found 281 infections transmitted by gastrointestinal endoscopy and 96 transmitted by bronchoscopy. The clinical spectrum ranged from asymptomatic colonization to death. Salmonella species and Pseudomonas aeruginosa repeatedly were identified as causative agents of infections transmitted by gastrointestinal endoscopy, and *M. tuberculosis*, atypical mycobacteria, and *P. aeruginosa* were the most common causes of infections transmitted by bronchoscopy ¹². Major reasons for transmission were inadequate cleaning, improper selection of a disinfecting agent, and failure to follow recommended cleaning and disinfection procedures ^{6, 8, 37, 98}, and flaws in endoscope design ^{127, 128} or automated endoscope reprocessors. ^{7, 98} Failure to follow established guidelines has continued to result in infections associated with gastrointestinal endoscopes ⁸ and bronchoscopes ^{7, 12}. Potential device-associated problems should be reported to the FDA Center for Devices and Radiologic Health. One multistate investigation found that 23.9% of the bacterial cultures from the internal channels of 71 gastrointestinal endoscopes grew ≥100,000 colonies of bacteria after completion of all disinfection and sterilization procedures (nine of 25 facilities were using a product that has been removed from the marketplace [six facilities using 1:16 glutaraldehyde phenate], is not FDA-cleared as a high-level disinfectant [an iodophor] or no disinfecting agent) and before use on the next patient¹²⁹. The incidence of postendoscopic procedure infections from an improperly processed endoscope has not been rigorously assessed.

Automated endoscope reprocessors (AER) offer several advantages over manual reprocessing: they automate and standardize several important reprocessing steps¹³⁰⁻¹³², reduce the likelihood that an essential reprocessing step will be skipped, and reduce personnel exposure to high-level disinfectants or chemical sterilants. Failure of AERs has been linked to outbreaks of infections ¹³³ or colonization ^{7, 134}, and the AER water filtration system might not be able to reliably provide "sterile" or bacteria-free rinse water^{135, 136}. Establishment of correct connectors between the AER and the device is critical to ensure complete flow of disinfectants and rinse water ^{7, 137}. In addition, some endoscopes such as the duodenoscopes (e.g., endoscopic retrograde cholangiopancreatography [ERCP]) contain features (e.g., elevator-wire channel) that require a flushing pressure that is not achieved by most AERs and must be reprocessed manually using a 2- to 5-mL syringe, until new duodenoscopes equipped with a wider elevator-channel that AERs can reliably reprocess become available ¹³². Outbreaks involving removable endoscope parts ^{138, 139} such as suction valves and endoscopic accessories designed to be inserted through flexible endoscopes such as biopsy forceps emphasize the importance of cleaning to remove all foreign matter before high-level disinfection or sterilization. ¹⁴⁰ Some types of valves are now available as single-use, disposable products (e.g., bronchoscope valves) or steam sterilizable products (e.g., gastrointestinal endoscope valves).

AERs need further development and redesign ^{7, 141}, as do endoscopes ^{123, 142}, so that they do not represent a potential source of infectious agents. Endoscopes employing disposable components (e.g., protective barrier devices or sheaths) might provide an alternative to conventional liquid chemical high-level disinfection/sterilization^{143, 144}. Another new technology is a swallowable camera-in-a-capsule that travels through the digestive tract and transmits color pictures of the small intestine to a receiver worn outside the body. This capsule currently does not replace colonoscopies.

Published recommendations for cleaning and disinfecting endoscopic equipment should be strictly followed ^{12, 38, 108, 113-116, 145-148}. Unfortunately, audits have shown that personnel do not consistently adhere to guidelines on reprocessing ¹⁴⁹⁻¹⁵¹ and outbreaks of infection continue to occur. ¹⁵²⁻¹⁵⁴ To ensure reprocessing personnel are properly trained, each person who reprocesses endoscopic instruments should receive initial and annual competency testing ^{38, 155}.

In general, endoscope disinfection or sterilization with a liquid chemical sterilant involves five steps after leak testing:

- 1. Clean: mechanically clean internal and external surfaces, including brushing internal channels and flushing each internal channel with water and a detergent or enzymatic cleaners (leak testing is recommended for endoscopes before immersion).
- Disinfect: immerse endoscope in high-level disinfectant (or chemical sterilant) and perfuse (eliminates air pockets and ensures contact of the germicide with the internal channels) disinfectant into all accessible channels, such as the suction/biopsy channel and air/water channel and expose for a time recommended for specific products.
- 3. Rinse: rinse the endoscope and all channels with sterile water, filtered water (commonly used with AERs) or tap water (i.e., high-quality potable water that meets federal clean water standards at the point of use).
- 4. Dry: rinse the insertion tube and inner channels with alcohol, and dry with forced air after disinfection and before storage.
- 5. Store: store the endoscope in a way that prevents recontamination and promotes drying (e.g., hung vertically).

Drving the endoscope (steps 3 and 4) is essential to greatly reduce the chance of recontamination of the endoscope by microorganisms that can be present in the rinse water ^{116, 156}. One study demonstrated that reprocessed endoscopes (i.e., air/water channel, suction/biopsy channel) generally were negative (100% after 24 hours; 90% after 7 days [1 CFU of coagulase-negative Staphylococcus in one channel]) for bacterial growth when stored by hanging vertically in a ventilated cabinet¹⁵⁷. Other investigators found all endoscopes were bacteria-free immediately after high-level disinfection, and only four of 135 scopes were positive during the subsequent 5-day assessment (skin bacteria cultured from endoscope surfaces). All flush-through samples remained sterile ¹⁵⁸. Because tapwater can contain low levels of microorganisms¹⁵⁹, some researchers have suggested that only sterile water (which can be prohibitively expensive) ¹⁶⁰ or AER filtered water be used. The suggestion to use only sterile water or filtered water is not consistent with published guidelines that allow tapwater with an alcohol rinse and forced air-drying ^{38, 108, 113} or the scientific literature. ^{39, 93} In addition, no evidence of disease transmission has been found when a tap water rinse is followed by an alcohol rinse and forced-air drving. AERs produce filtered water by passage through a bacterial filter (e.g., 0.2 µ). Filtered rinse water was identified as a source of bacterial contamination in a study that cultured the accessory and suction channels of endoscopes and the internal chambers of AERs during 1996-2001 and reported 8.7% of samples collected during 1996-1998 had bacterial growth, with 54% being Pseudomonas species. After a system of hot water flushing of the piping (60°C for 60 minutes daily) was introduced, the frequency of positive cultures fell to approximately 2% with only rare isolation of >10 CFU/mL ¹⁶¹. In addition to the endoscope reprocessing steps, a protocol should be developed that ensures the user knows whether an endoscope has been appropriately cleaned and disinfected (e.g., using a room or cabinet for processed endoscopes only) or has not been reprocessed. When users leave endoscopes on movable carts, confusion can result about whether the endoscope has been processed. Although one guideline recommended endoscopes (e.g., duodenoscopes) be reprocessed immediately before use ¹⁴⁷, other guidelines do not require this activity ^{38, 108, 115} and except for the Association of periOperative Registered Nurses (AORN), professional organizations do not recommended that reprocessing be repeated as long as the original processing is done correctly. As part of a quality assurance program, healthcare facility personnel can consider random bacterial surveillance cultures of processed endoscopes to ensure high-level disinfection or sterilization^{7, 162-164}. Reprocessed endoscopes should be free of microbial pathogens except for small numbers of relatively avirulent microbes that represent exogenous environmental contamination (e.g., coagulase-negative Staphylococcus, Bacillus species, diphtheroids). Although recommendations exist for the final rinse water used during endoscope

Guideline for Disinfection and Sterilization in Healthcare Facilities (2008)

reprocessing to be microbiologically cultured at least monthly ¹⁶⁵, a microbiologic standard has not been set, and the value of routine endoscope cultures has not been shown ¹⁶⁶. In addition, neither the routine culture of reprocessed endoscopes nor the final rinse water has been validated by correlating viable counts on an endoscope to infection after an endoscopic procedure. If reprocessed endoscopes were cultured, sampling the endoscope would assess water quality and other important steps (e.g., disinfectant effectiveness, exposure time, cleaning) in the reprocessing procedure. A number of methods for sampling endoscopes and water have been described ^{23, 157, 161, 163, 167, 168}. Novel approaches (e.g., detection of adenosine triphosphate [ATP]) to evaluate the effectiveness of endoscope cleaning ^{169, 170} or endoscope reprocessing ¹⁷¹ also have been evaluated, but no method has been established as a standard for assessing the outcome of endoscope reprocessing.

The carrying case used to transport clean and reprocessed endoscopes outside the health-care environment should not be used to store an endoscope or to transport the instrument within the health-care facility. A contaminated endoscope should never be placed in the carrying case because the case can also become contaminated. When the endoscope is removed from the case, properly reprocessed, and put back in the case, the case could recontaminate the endoscope. A contaminated carrying case should be discarded (Olympus America, June 2002, written communication).

Infection-control professionals should ensure that institutional policies are consistent with national guidelines and conduct infection-control rounds periodically (e.g., at least annually) in areas where endoscopes are reprocessed to ensure policy compliance. Breaches in policy should be documented and corrective action instituted. In incidents in which endoscopes were not exposed to a high-level disinfection process, patients exposed to potentially contaminated endoscopes have been assessed for possible acquisition of HIV, HBV, and hepatitis C virus (HCV). A 14-step method for managing a failure incident associated with high-level disinfection or sterilization has been described [Rutala WA, 2006 #12512]. The possible transmission of bloodborne and other infectious agents highlights the importance of rigorous infection control^{172, 173}.

Laparoscopes and Arthroscopes

Although high-level disinfection appears to be the minimum standard for processing laparoscopes and arthroscopes between patients 28, 86, 174, 175, this practice continues to be debated 89, 90, 176. However, neither side in the high-level disinfection versus sterilization debate has sufficient data on which to base its conclusions. Proponents of high-level disinfection refer to membership surveys ²⁹ or institutional experiences ⁸⁷ involving more than 117,000 and 10,000 laparoscopic procedures, respectively, that cite a low risk for infection (<0.3%) when high-level disinfection is used for gynecologic laparoscopic equipment. Only one infection in the membership survey was linked to spores. In addition, growth of common skin microorganisms (e.g., Staphylococcus epidermidis, diphtheroids) has been documented from the umbilical area even after skin preparation with povidone-iodine and ethyl alcohol. Similar organisms were recovered in some instances from the pelvic serosal surfaces or from the laparoscopic telescopes, suggesting that the microorganisms probably were carried from the skin into the peritoneal cavity ^{177, 178}. Proponents of sterilization focus on the possibility of transmitting infection by spore-forming organisms. Researchers have proposed several reasons why sterility was not necessary for all laparoscopic equipment: only a limited number of organisms (usually <10) are introduced into the peritoneal cavity during laparoscopy; minimal damage is done to inner abdominal structures with little devitalized tissue; the peritoneal cavity tolerates small numbers of spore-forming bacteria; equipment is simple to clean and disinfect; surgical sterility is relative; the natural bioburden on rigid lumened devices is low¹⁷⁹; and no evidence exists that high-level disinfection instead of sterilization increases the risk for infection ^{87, 89, 90}. With the advent of laparoscopic cholecystectomy, concern about high-level disinfection is justifiable because the degree of tissue damage and bacterial contamination is greater than with laparoscopic procedures in gynecology. Failure to completely dissemble, clean, and high-level disinfect laparoscope parts has led to infections in patients¹⁸⁰. Data from one study suggested that disassembly, cleaning, and proper reassembly of laparoscopic equipment used in gynecologic procedures before steam sterilization presents no risk for infection¹⁸¹.

As with laparoscopes and other equipment that enter sterile body sites, arthroscopes ideally should be sterilized before used. Older studies demonstrated that these instruments were commonly (57%) only high-level disinfected in the United States ^{28, 86}. A later survey (with a response rate of only 5%) reported that high-level disinfection was used by 31% and a sterilization process in the remainder of the health-care facilities³⁰ High-level disinfection rather than sterilization presumably has been used because the incidence of infection is low and the few infections identified probably are unrelated to the use of high-level disinfection rather than sterilization. A retrospective study of 12,505 arthroscopic procedures found an infection rate of 0.04% (five infections) when arthroscopes were soaked in 2% glutaraldehyde for 15–20 minutes. Four infections were caused by *S. aureus*; the fifth was an anaerobic streptococcal infection ⁸⁸. Because these organisms are very susceptible to high-level disinfectants, such as 2% glutaraldehyde, the infections most likely originated from the patient's skin. Two cases of *Clostridium perfringens* arthritis have been reported when the arthroscope was disinfected with glutaraldehyde for an exposure time that is not effective against spores ^{182, 183}.

Although only limited data are available, the evidence does not demonstrate that high-level disinfection of arthroscopes and laparoscopes poses an infection risk to the patient. For example, a prospective study that compared the reprocessing of arthroscopes and laparoscopes (per 1,000 procedures) with EtO sterilization to high-level disinfection with glutaraldehyde found no statistically significant difference in infection risk between the two methods (i.e., EtO, 7.5/1,000 procedures; glutaraldehyde, 2.5/1,000 procedures)⁸⁹. Although the debate for high-level disinfection versus sterilization of laparoscopes and arthroscopes will go unsettled until well-designed, randomized clinical trials are published, this guideline should be followed ^{1, 17}. That is, laparoscopes, arthroscopes, and other scopes that enter normally sterile tissue should be sterilized before each use; if this is not feasible, they should receive at least high-level disinfection.

Tonometers, Cervical Diaphragm Fitting Rings, Cryosurgical Instruments, and Endocavitary Probes

Disinfection strategies vary widely for other semicritical items (e.g., applanation tonometers, rectal/vaginal probes, cryosurgical instruments, and diaphragm fitting rings). FDA requests that device manufacturers include at least one validated cleaning and disinfection/sterilization protocol in the labeling for their devices. As with all medications and devices, users should be familiar with the label instructions. One study revealed that no uniform technique was in use for disinfection of applanation tonometers, with disinfectant contact times varying from <15 sec to 20 minutes ²⁸. In view of the potential for transmission of viruses (e.g., herpes simplex virus [HSV], adenovirus 8, or HIV) ¹⁸⁴ by tonometer tips, CDC recommended that the tonometer tips be wiped clean and disinfected for 5-10 minutes with either 3% hydrogen peroxide, 5000 ppm chlorine, 70% ethyl alcohol, or 70% isopropyl alcohol ⁹⁵. However, more recent data suggest that 3% hydrogen peroxide and 70% isopropyl alcohol are not effective against adenovirus capable of causing epidemic keratoconjunctivitis and similar viruses and should not be used for disinfecting applanation tonometers ^{49, 185, 186}. Structural damage to Schiotz tonometers has been observed with a 1:10 sodium hypochlorite (5,000 ppm chlorine) and 3% hydrogen peroxide¹⁸⁷. After disinfection, the tonometer should be thoroughly rinsed in tapwater and air dried before use. Although these disinfectants and exposure times should kill pathogens that can infect the eyes, no studies directly support this ^{188, 189}. The guidelines of the American Academy of Ophthalmology for preventing infections in ophthalmology focus on only one potential pathogen: HIV. ¹⁹⁰ Because a short and simple decontamination procedure is desirable in the clinical setting, swabbing the tonometer tip with a 70% isopropyl alcohol wipe sometimes is practiced. ¹⁸⁹ Preliminary reports suggest that wiping the tonometer tip with an alcohol swab and then allowing the alcohol to evaporate might be effective in eliminating HSV, HIV, and adenovirus^{189, 191, 192}. However, because these studies involved only a few replicates and were conducted in a controlled laboratory setting, further studies are needed before this technique can be recommended. In addition, two reports have found that disinfection of pneumotonometer tips between uses with a 70% isopropyl alcohol wipe contributed to outbreaks of epidemic keratoconjunctivitis caused

by adenovirus type 8^{193, 194}.

Limited studies have evaluated disinfection techniques for other items that contact mucous membranes, such as diaphragm fitting rings, cryosurgical probes, transesophageal echocardiography probes ¹⁹⁵, flexible cystoscopes ¹⁹⁶ or vaginal/rectal probes used in sonographic scanning. Lettau, Bond, and McDougal of CDC supported the recommendation of a diaphragm fitting ring manufacturer that involved using a soap-and-water wash followed by a 15-minute immersion in 70% alcohol⁹⁶. This disinfection method should be adequate to inactivate HIV, HBV, and HSV even though alcohols are not classified as high-level disinfectants because their activity against picornaviruses is somewhat limited⁷². No data are available regarding inactivation of human papillomavirus (HPV) by alcohol or other disinfectants because in vitro replication of complete virions has not been achieved. Thus, even though alcohol for 15 minutes should kill pathogens of relevance in gynecology, no clinical studies directly support this practice.

Vaginal probes are used in sonographic scanning. A vaginal probe and all endocavitary probes without a probe cover are semicritical devices because they have direct contact with mucous membranes (e.g., vagina, rectum, pharynx). While use of the probe cover could be considered as changing the category, this guideline proposes use of a new condom/probe cover for the probe for each patient, and because condoms/probe covers can fail ^{195, 197-199}, the probe also should be high-level disinfected. The relevance of this recommendation is reinforced with the findings that sterile transvaginal ultrasound probe covers have a very high rate of perforations even before use (0%, 25%, and 65% perforations from three suppliers). ¹⁹⁹ One study found, after oocyte retrieval use, a very high rate of perforations in used endovaginal probe covers from two suppliers (75% and 81%)¹⁹⁹, other studies demonstrated a lower rate of perforations after use of condoms (2.0% and 0.9%) ^{197 200}. Condoms have been found superior to commercially available probe covers for covering the ultrasound probe (1.7% for condoms versus 8.3% leakage for probe covers)²⁰¹. These studies underscore the need for routine probe disinfection between examinations. Although most ultrasound manufacturers recommend use of 2% glutaraldehyde for highlevel disinfection of contaminated transvaginal transducers, the this agent has been questioned ²⁰² because it might shorten the life of the transducer and might have toxic effects on the gametes and embryos ²⁰³. An alternative procedure for disinfecting the vaginal transducer involves the mechanical removal of the gel from the transducer, cleaning the transducer in soap and water, wiping the transducer with 70% alcohol or soaking it for 2 minutes in 500 ppm chlorine, and rinsing with tap water and air drying²⁰⁴. The effectiveness of this and other methods ²⁰⁰ has not been validated in either rigorous laboratory experiments or in clinical use. High-level disinfection with a product (e.g., hydrogen peroxide) that is not toxic to staff, patients, probes, and retrieved cells should be used until the effectiveness of alternative procedures against microbes of importance at the cavitary site is demonstrated by welldesigned experimental scientific studies. Other probes such as rectal, cryosurgical, and transesophageal probes or devices also should be high-level disinfected between patients.

Ultrasound probes used during surgical procedures also can contact sterile body sites. These probes can be covered with a sterile sheath to reduce the level of contamination on the probe and reduce the risk for infection. However, because the sheath does not completely protect the probe, the probes should be sterilized between each patient use as with other critical items. If this is not possible, at a minimum the probe should be high-level disinfected and covered with a sterile probe cover.

Some cryosurgical probes are not fully immersible. During reprocessing, the tip of the probe should be immersed in a high-level disinfectant for the appropriate time; any other portion of the probe that could have mucous membrane contact can be disinfected by immersion or by wrapping with a cloth soaked in a high-level disinfectant to allow the recommended contact time. After disinfection, the probe should be rinsed with tap water and dried before use. Health-care facilities that use nonimmersible probes should replace them as soon as possible with fully immersible probes.

As with other high-level disinfection procedures, proper cleaning of probes is necessary to ensure the success of the subsequent disinfection ²⁰⁵. One study demonstrated that vegetative bacteria

inoculated on vaginal ultrasound probes decreased when the probes were cleaned with a towel ²⁰⁶. No information is available about either the level of contamination of such probes by potential viral pathogens such as HBV and HPV or their removal by cleaning (such as with a towel). Because these pathogens might be present in vaginal and rectal secretions and contaminate probes during use, high-level disinfection of the probes after such use is recommended.

Dental Instruments

Scientific articles and increased publicity about the potential for transmitting infectious agents in dentistry have focused attention on dental instruments as possible agents for pathogen transmission^{207, 208}. The American Dental Association recommends that surgical and other instruments that normally penetrate soft tissue or bone (e.g., extraction forceps, scalpel blades, bone chisels, periodontal scalers, and surgical burs) be classified as critical devices that should be sterilized after each use or discarded. Instruments not intended to penetrate oral soft tissues or bone (e.g., amalgam condensers, and air/water syringes) but that could contact oral tissues are classified as semicritical, but sterilization after each use is recommended if the instruments are heat-tolerant ^{43, 209}. If a semicritical item is heat–sensitive, it should, at a minimum, be processed with high-level disinfection ^{43, 210}. Handpieces can be contaminated internally with patient material and should be heat sterilized after each patient. Handpieces that cannot be heat sterilized should not be used. ²¹¹ Methods of sterilization that can be used for critical or semicritical dental instruments and materials that are heat-stable include steam under pressure (autoclave), chemical (formaldehyde) vapor, and dry heat (e.g., 320°F for 2 hours). Dental professionals most commonly use the steam sterilized hand pieces ²¹³. Heat-tolerant alternatives are available for most clinical dental applications and are preferred⁴³.

CDC has divided noncritical surfaces in dental offices into clinical contact and housekeeping surfaces⁴³. Clinical contact surfaces are surfaces that might be touched frequently with gloved hands during patient care or that might become contaminated with blood or other potentially infectious material and subsequently contact instruments, hands, gloves, or devices (e.g., light handles, switches, dental X-ray equipment, chair-side computers). Barrier protective coverings (e.g., clear plastic wraps) can be used for these surfaces, particularly those that are difficult to clean (e.g., light handles, chair switches). The coverings should be changed when visibly soiled or damaged and routinely (e.g., between patients). Protected surfaces should be disinfected at the end of each day or if contamination is evident. If not barrier-protected, these surfaces should be disinfected between patients with an intermediate-disinfectant (i.e., EPA-registered hospital disinfectant with tuberculocidal claim) or low-level disinfectant (i.e., EPA-registered hospital disinfectant with an HBV and HIV label claim) ^{43, 214, 215}.

Most housekeeping surfaces need to be cleaned only with a detergent and water or an EPAregistered hospital disinfectant, depending of the nature of the surface and the type and degree of contamination. When housekeeping surfaces are visibly contaminated by blood or body substances, however, prompt removal and surface disinfection is a sound infection control practice and required by the Occupational Safety and Health Administration (OSHA) ^{43, 214}.

Several studies have demonstrated variability among dental practices while trying to meet these recommendations^{216, 217}. For example, 68% of respondents believed they were sterilizing their instruments but did not use appropriate chemical sterilants or exposure times and 49% of respondents did not challenge autoclaves with biological indicators²¹⁶. Other investigators using biologic indicators have found a high proportion (15%–65%) of positive spore tests after assessing the efficacy of sterilizers used in dental offices. In one study of Minnesota dental offices, operator error, rather than mechanical malfunction²¹⁸, caused 87% of sterilization failures. Common factors in the improper use of sterilizers include chamber overload, low temperature setting, inadequate exposure time, failure to preheat the sterilizer, and interruption of the cycle.

Mail-return sterilization monitoring services use spore strips to test sterilizers in dental clinics, but

delay caused by mailing to the test laboratory could potentially cause false-negatives results. Studies revealed, however, that the post-sterilization time and temperature after a 7-day delay had no influence on the test results²¹⁹. Delays (7 days at 27°C and 37°C, 3-day mail delay) did not cause any predictable pattern of inaccurate spore tests ²²⁰.

Disinfection of HBV-, HCV-, HIV- or TB-Contaminated Devices

The CDC recommendation for high-level disinfection of HBV-, HCV-, HIV- or TB-contaminated devices is appropriate because experiments have demonstrated the effectiveness of high-level disinfectants to inactivate these and other pathogens that might contaminate semicritical devices ^{61, 62, 73, 81, 105, 121, 125, 221-238}. Nonetheless, some healthcare facilities have modified their disinfection procedures when endoscopes are used with a patient known or suspected to be infected with HBV, HIV, or *M. tuberculosis* ^{28, 239}. This is inconsistent with the concept of Standard Precautions that presumes all patients are potentially infected with bloodborne pathogens²²⁸. Several studies have highlighted the inability to distinguish HBV- or HIV-infected patients from noninfected patients on clinical grounds²⁴⁰⁻²⁴². In addition, mycobacterial infection is unlikely to be clinically apparent in many patients. In most instances, hospitals that altered their disinfection procedure used EtO sterilization on the endoscopic instruments because they believed this practice reduced the risk for infection ^{28, 239}. EtO is not routinely used for endoscope sterilization because of the lengthy processing time. Endoscopes and other semicritical devices should be managed the same way regardless of whether the patient is known to be infected with HBV, HCV, HIV or *M. tuberculosis*.

An evaluation of a manual disinfection procedure to eliminate HCV from experimentally contaminated endoscopes provided some evidence that cleaning and 2% glutaraldehyde for 20 minutes should prevent transmission ²³⁶. A study that used experimentally contaminated hysteroscopes detected HCV by polymerase chain reaction (PCR) in one (3%) of 34 samples after cleaning with a detergent, but no samples were positive after treatment with a 2% glutaraldehyde solution for 20 minutes ¹²⁰. Another study demonstrated complete elimination of HCV (as detected by PCR) from endoscopes used on chronically infected patients after cleaning and disinfection for 3–5 minutes in glutaraldehyde ¹¹⁸. Similarly, PCR was used to demonstrate complete elimination of HCV after standard disinfection of experimentally contaminated endoscopes ²³⁶ and endoscopes used on HCV-antibody– positive patients had no detectable HCV RNA after high-level disinfection ²⁴³. The inhibitory activity of a phenolic and a chlorine compound on HCV showed that the phenolic inhibited the binding and replication of HCV, but the chlorine was ineffective, probably because of its low concentration and its neutralization in the presence of organic matter ²⁴⁴.

Disinfection in the Hemodialysis Unit

Hemodialysis systems include hemodialysis machines, water supply, water-treatment systems, and distribution systems. During hemodialysis, patients have acquired bloodborne viruses and pathogenic bacteria ²⁴⁵⁻²⁴⁷. Cleaning and disinfection are important components of infection control in a hemodialysis center. EPA and FDA regulate disinfectants used to reprocess hemodialyzers, hemodialysis machines, and water-treatment systems.

Noncritical surfaces (e.g., dialysis bed or chair, countertops, external surfaces of dialysis machines, and equipment [scissors, hemostats, clamps, blood pressure cuffs, stethoscopes]) should be disinfected with an EPA-registered disinfectant unless the item is visibly contaminated with blood; in that case a tuberculocidal agent (or a disinfectant with specific label claims for HBV and HIV) or a 1:100 dilution of a hypochlorite solution (500–600 ppm free chlorine) should be used ^{246, 248}. This procedure accomplishes two goals: it removes soil on a regular basis and maintains an environment that is consistent with good patient care. Hemodialyzers are disinfected with peracetic acid, formaldehyde, glutaraldehyde, heat pasteurization with citric acid, and chlorine-containing compounds ²⁴⁹. Hemodialysis systems usually are disinfected by chlorine-based disinfectants (e.g., sodium hypochlorite), aqueous

formaldehyde, heat pasteurization, ozone, or peracetic acid ^{250, 251}. All products must be used according to the manufacturers' recommendations. Some dialysis systems use hot-water disinfection to control microbial contamination.

At its high point, 82% of U.S. chronic hemodialysis centers were reprocessing (i.e., reusing) dialyzers for the same patient using high-level disinfection ²⁴⁹. However, one of the large dialysis organizations has decided to phase out reuse and, by 2002 the percentage of dialysis facilities reprocessing hemodialyzers had decreased to 63% ²⁵². The two commonly used disinfectants to reprocess dialyzers were peracetic acid and formaldehyde; 72% used peracetic acid and 20% used formaldehyde to disinfect hemodialyzers. Another 4% of the facilities used either glutaraldehyde or heat pasteurization in combination with citric acid ²⁵². Infection-control recommendations, including disinfection and sterilization and the use of dedicated machines for hepatitis B surface antigen (HBsAg)-positive patients, in the hemodialysis setting were detailed in two reviews ^{245, 246}. The Association for the Advancement of Medical Instrumentation(AAMI) has published recommendations for the reuse of hemodialyzers²⁵³.

Inactivation of Clostridium difficile

The source of health-care-associated acquisition of Clostridium difficile in nonepidemic settings has not been determined. The environment and carriage on the hands of health-care personnel have been considered possible sources of infection ^{66, 254}. Carpeted rooms occupied by a patient with C. difficile were more heavily contaminated with C. difficile than were noncarpeted rooms ²⁵⁵. Because C. difficile spore-production can increase when exposed to nonchlorine-based cleaning agents and the spores are more resistant than vegetative cells to commonly used surface disinfectants²⁵⁶, some investigators have recommended use of dilute solutions of hypochlorite (1,600 ppm available chlorine) for routine environmental disinfection of rooms of patients with C. difficile-associated diarrhea or colitis 257, to reduce the incidence of C. difficile diarrhea ²⁵⁸, or in units with high C. difficile rates. ²⁵⁹ Stool samples of patients with symptomatic C. difficile colitis contain spores of the organism, as demonstrated by ethanol treatment of the stool to reduce the overgrowth of fecal flora when isolating C. difficile in the laboratory^{260, 261}, C. difficile-associated diarrhea rates were shown to have decreased markedly in a bone-marrow transplant unit (from 8.6 to 3.3 cases per 1,000 patient-days) during a period of bleach disinfection (1:10 dilution) of environmental surfaces compared with cleaning with a guaternary ammonium compound. Because no EPA-registered products exist that are specific for inactivating C. difficile spores, use of diluted hypochlorite should be considered in units with high C. difficile rates. Acidified bleach and regular bleach (5000 ppm chlorine) can inactivate 10⁶ C. difficile spores in ≤10 minutes ²⁶². However, studies have shown that asymptomatic patients constitute an important reservoir within the health-care facility and that person-to-person transmission is the principal means of transmission between patients. Thus, combined use of hand washing, barrier precautions, and meticulous environmental cleaning with an EPA-registered disinfectant (e.g., germicidal detergent) should effectively prevent spread of the organism ²⁶³.

Contaminated medical devices, such as colonoscopes and thermometers, can be vehicles for transmission of *C. difficile* spores ²⁶⁴. For this reason, investigators have studied commonly used disinfectants and exposure times to assess whether current practices can place patients at risk. Data demonstrate that 2% glutaraldehyde ^{79, 265-267} and peracetic acid ^{267, 268} reliably kill *C. difficile* spores using exposure times of 5–20 minutes. *ortho*-Phthalaldehyde and $\ge 0.2\%$ peracetic acid (WA Rutala, personal communication, April 2006) also can inactivate $\ge 10^4$ *C. difficile* spores in 10–12 minutes at 20°C ²⁶⁸. Sodium dichloroisocyanurate at a concentration of 1000 ppm available chlorine achieved lower log₁₀ reduction factors against *C. difficile* spores at 10 min, ranging from 0.7 to 1.5, than 0.26% peracetic acid with log₁₀ reduction factors ranging from 2.7 to 6.0²⁶⁸.

OSHA Bloodborne Pathogen Standard

In December 1991, OSHA promulgated a standard entitled "Occupational Exposure to Bloodborne

Pathogens" to eliminate or minimize occupational exposure to bloodborne pathogens ²¹⁴. One component of this requirement is that all equipment and environmental and working surfaces be cleaned and decontaminated with an appropriate disinfectant after contact with blood or other potentially infectious materials. Even though the OSHA standard does not specify the type of disinfectant or procedure, the OSHA original compliance document ²⁶⁹ suggested that a germicide must be tuberculocidal to kill the HBV. To follow the OSHA compliance document a tuberculocidal disinfectant (e.g., phenolic, and chlorine) would be needed to clean a blood spill. However, in February 1997, OSHA amended its policy and stated that EPA-registered disinfectants labeled as effective against HIV and HBV would be considered as appropriate disinfectants ". . . provided such surfaces have not become contaminated with agent(s) or volumes of or concentrations of agent(s) for which higher level disinfection is recommended." When bloodborne pathogens other than HBV or HIV are of concern, OSHA continues to require use of EPA-registered tuberculocidal disinfectants or hypochlorite solution (diluted 1:10 or 1:100 with water) ^{215, 228}. Studies demonstrate that, in the presence of large blood spills, a 1:10 final dilution of EPA-registered hypochlorite solution initially should be used to inactivate bloodborne viruses ^{63, 235} to minimize risk for infection to health-care personnel from percutaneous injury during cleanup.

Emerging Pathogens (*Cryptosporidium, Helicobacter pylori, Escherichia coli* O157:H7, Rotavirus, Human Papilloma Virus, Norovirus, Severe Acute Respiratory Syndrome [SARS] Coronavirus)

Emerging pathogens are of growing concern to the general public and infection-control professionals. Relevant pathogens include *Cryptosporidium parvum*, *Helicobacter pylori*, *E. coli* O157:H7, HIV, HCV, rotavirus, norovirus, severe acute respiratory syndrome (SARS) coronavirus, multidrug-resistant *M. tuberculosis*, and nontuberculous mycobacteria (e.g., *M. chelonae*). The susceptibility of each of these pathogens to chemical disinfectants and sterilants has been studied. With the exceptions discussed below, all of these emerging pathogens are susceptible to currently available chemical disinfectants and sterilants ²⁷⁰.

Cryptosporidium is resistant to chlorine at concentrations used in potable water. *C. parvum* is not completely inactivated by most disinfectants used in healthcare including ethyl alcohol ²⁷¹, glutaraldehyde ^{271, 272}, 5.25% hypochlorite ²⁷¹, peracetic acid ²⁷¹, ortho-phthalaldehyde ²⁷¹, phenol ^{271, 272}, povidone-iodine ^{271, 272}, and quaternary ammonium compounds²⁷¹. The only chemical disinfectants and sterilants able to inactivate greater than 3 log₁₀ of *C. parvum* were 6% and 7.5% hydrogen peroxide ²⁷¹. Sterilization methods will fully inactivate *C. parvum*, including steam ²⁷¹, EtO ^{271, 273}, and hydrogen peroxide gas plasma²⁷¹. Although most disinfectants are ineffective against *C. parvum*, current cleaning and disinfection practices appear satisfactory to prevent healthcare-associated transmission. For example, endoscopes are unlikely to be an important vehicle for transmitting *C. parvum* because the results of bacterial studies indicate mechanical cleaning will remove approximately 10⁴ organisms, and drying results in rapid loss of *C. parvum* viability (e.g., 30 minutes, 2.9 log₁₀ decrease; and 60 minutes, 3.8 log₁₀ decrease) ²⁷¹.

Chlorine at ~1 ppm has been found capable of eliminating approximately 4 log₁₀ of *E. coli* O157:H7 within 1 minute in a suspension test⁶⁴. Electrolyzed oxidizing water at 23°C was effective in 10 minutes in producing a 5-log₁₀ decrease in *E. coli* O157:H7 inoculated onto kitchen cutting boards²⁷⁴. The following disinfectants eliminated >5 log₁₀ of *E. coli* O157:H7 within 30 seconds: a quaternary ammonium compound, a phenolic, a hypochlorite (1:10 dilution of 5.25% bleach), and ethanol⁵³. Disinfectants including chlorine compounds can reduce *E. coli* O157:H7 experimentally inoculated onto alfalfa seeds or sprouts ^{275, 276} or beef carcass surfaces²⁷⁷.

Data are limited on the susceptibility of *H. pylori* to disinfectants. Using a suspension test, one study assessed the effectiveness of a variety of disinfectants against nine strains of *H. pylori*⁶⁰. Ethanol (80%) and glutaraldehyde (0.5%) killed all strains within 15 seconds; chlorhexidine gluconate (0.05%, 1.0%), benzalkonium chloride (0.025%, 0.1%), alkyldiaminoethylglycine hydrochloride (0.1%), povidone-iodine (0.1%), and sodium hypochlorite (150 ppm) killed all strains within 30 seconds. Both ethanol (80%) and

glutaraldehyde (0.5%) retained similar bactericidal activity in the presence of organic matter; the other disinfectants showed reduced bactericidal activity. In particular, the bactericidal activity of povidone-iodine (0.1%) and sodium hypochlorite (150 ppm) markedly decreased in the presence of dried yeast solution with killing times increased to 5 - 10 minutes and 5 - 30 minutes, respectively.

Immersing biopsy forceps in formalin before obtaining a specimen does not affect the ability to culture *H. pylori* from the biopsy specimen ²⁷⁸. The following methods are ineffective for eliminating *H.* pylori from endoscopes: cleaning with soap and water ^{119, 279}, immersion in 70% ethanol for 3 minutes²⁸⁰, instillation of 70% ethanol¹²⁶, instillation of 30 ml of 83% methanol²⁷⁹, and instillation of 0.2% Hyamine solution²⁸¹. The differing results with regard to the efficacy of ethyl alcohol against Helicobacter are unexplained. Cleaning followed by use of 2% alkaline glutaraldehyde (or automated peracetic acid) has been demonstrated by culture to be effective in eliminating H. pylori ^{119, 279, 282}. Epidemiologic investigations of patients who had undergone endoscopy with endoscopes mechanically washed and disinfected with 2.0%-2.3% glutaraldehyde have revealed no evidence of person-toperson transmission of H. pylori ^{126, 283}. Disinfection of experimentally contaminated endoscopes using 2% glutaraldehyde (10-minute, 20-minute, 45-minute exposure times) or the peracetic acid system (with and without active peracetic acid) has been demonstrated to be effective in eliminating H. pylori¹¹⁹. H. pylori DNA has been detected by PCR in fluid flushed from endoscope channels after cleaning and disinfection with 2% glutaraldehyde 284. The clinical significance of this finding is unclear. In vitro experiments have demonstrated a >3.5-log₁₀ reduction in *H. pylori* after exposure to 0.5 mg/L of free chlorine for 80 seconds²⁸⁵.

An outbreak of healthcare-associated rotavirus gastroenteritis on a pediatric unit has been reported ²⁸⁶. Person to person through the hands of health-care workers was proposed as the mechanism of transmission. Prolonged survival of rotavirus on environmental surfaces (90 minutes to >10 days at room temperature) and hands (>4 hours) has been demonstrated. Rotavirus suspended in feces can survive longer ^{287, 288}. Vectors have included hands, fomites, air, water, and food ^{288, 289}. Products with demonstrated efficacy (>3 log₁₀ reduction in virus) against rotavirus within 1 minute include: 95% ethanol, 70% isopropanol, some phenolics, 2% glutaraldehyde, 0.35% peracetic acid, and some quaternary ammonium compounds ^{59, 290-293}. In a human challenge study, a disinfectant spray (0.1% orthophenylphenol and 79% ethanol), sodium hypochlorite (800 ppm free chlorine), and a phenol-based product (14.7% phenol diluted 1:256 in tapwater) when sprayed onto contaminated stainless steel disks, were effective in interrupting transfer of a human rotavirus from stainless steel disk to fingerpads of volunteers after an exposure time of 3- 10 minutes. A quaternary ammonium product (7.05% quaternary ammonium compound diluted 1:128 in tapwater) and tapwater allowed transfer of virus ⁵².

No data exist on the inactivation of HPV by alcohol or other disinfectants because in vitro replication of complete virions has not been achieved. Similarly, little is known about inactivation of noroviruses (members of the family Caliciviridae and important causes of gastroenteritis in humans) because they cannot be grown in tissue culture. Improper disinfection of environmental surfaces contaminated by feces or vomitus of infected patients is believed to play a role in the spread of noroviruses in some settings ²⁹⁴⁻²⁹⁶. Prolonged survival of a norovirus surrogate (i.e., feline calicivirus virus [FCV], a closely related cultivable virus) has been demonstrated (e.g., at room temperature, FCV in a dried state survived for 21–18 days) ²⁹⁷. Inactivation studies with FCV have shown the effectiveness of chlorine, glutaraldehyde, and iodine-based products whereas the quaternary ammonium compound, detergent, and ethanol failed to inactivate the virus completely. ²⁹⁷ An evaluation of the effectiveness of several disinfectants against the feline calicivirus found that bleach diluted to 1000 ppm of available chlorine reduced infectivity of FCV by 4.5 logs in 1 minute. Other effective (\log_{10} reduction factor of >4 in virus) disinfectants included accelerated hydrogen peroxide, 5,000 ppm (3 min); chlorine dioxide, 1,000 ppm chlorine (1 min); a mixture of four quaternary ammonium compounds, 2,470 ppm (10 min); 79% ethanol with 0.1% quaternary ammonium compound (3 min); and 75% ethanol (10 min) ²⁹⁸. A quaternary ammonium compound exhibited activity against feline calicivirus supensions dried on hard surface carriers in 10 minutes ²⁹⁹. Seventy percent ethanol and 70% 1-propanol reduced FCV by a 3–4-log₁₀ reduction in 30 seconds ³⁰⁰.

CDC announced that a previously unrecognized human virus from the coronavirus family is the leading hypothesis for the cause of a described syndrome of SARS ³⁰¹. Two coronaviruses that are known to infect humans cause one third of common colds and can cause gastroenteritis. The virucidal efficacy of chemical germicides against coronavirus has been investigated. A study of disinfectants against coronavirus 229E found several that were effective after a 1-minute contact time; these included sodium hypochlorite (at a free chlorine concentration of 1,000 ppm and 5,000 ppm), 70% ethyl alcohol, and povidone-iodine (1% iodine) ¹⁸⁶. In another study, 70% ethanol, 50% isopropanol, 0.05% benzalkonium chloride, 50 ppm iodine in iodophor, 0.23% sodium chlorite, 1% cresol soap and 0.7% formaldehyde inactivated >3 logs of two animal coronaviruses (mouse hepatitis virus, canine coronavirus) after a 10-minute exposure time ³⁰². The activity of povidone-iodine has been demonstrated against human coronaviruses 229E and OC43 303. A study also showed complete inactivation of the SARS coronavirus by 70% ethanol and povidone-iodine with an exposure times of 1 minute and 2.5% glutaraldehyde with an exposure time of 5 minute ³⁰⁴. Because the SARS coronavirus is stable in feces and urine at room temperature for at least 1-2 days [The current version of this document may differ from original: First data on stability and resistance of SARS coronavirus compiled by members of WHO laboratory network (http://www.who.int/csr/sars/survival 2003 05 04/en/)], surfaces might be a possible source of contamination and lead to infection with the SARS coronavirus and should be disinfected. Until more precise information is available, environments in which SARS patients are housed should be considered heavily contaminated, and rooms and equipment should be thoroughly disinfected daily and after the patient is discharged. EPA-registered disinfectants or 1:100 dilution of household bleach and water should be used for surface disinfection and disinfection on noncritical patient-care equipment. Highlevel disinfection and sterilization of semicritical and critical medical devices, respectively, does not need to be altered for patients with known or suspected SARS.

Free-living amoeba can be pathogenic and can harbor agents of pneumonia such as *Legionella pneumophila*. Limited studies have shown that 2% glutaraldehyde and peracetic acid do not completely inactivate *Acanthamoeba polyphaga* in a 20-minute exposure time for high-level disinfection. If amoeba are found to contaminate instruments and facilitate infection, longer immersion times or other disinfectants may need to be considered ³⁰⁵.

Inactivation of Bioterrorist Agents

Publications have highlighted concerns about the potential for biological terrorism^{306, 307}. CDC has categorized several agents as "high priority" because they can be easily disseminated or transmitted from person to person, cause high mortality, and are likely to cause public panic and social disruption ³⁰⁸. These agents include *Bacillus anthracis* (the cause of anthrax), *Yersinia pestis* (plague), variola major (smallpox), *Clostridium botulinum* toxin (botulism), *Francisella tularensis* (tularemia), filoviruses (Ebola hemorrhagic fever, Marburg hemorrhagic fever); and arenaviruses (Lassa [Lassa fever], Junin [Argentine hemorrhagic fever]), and related viruses³⁰⁸.

A few comments can be made regarding the role of sterilization and disinfection of potential agents of bioterrorism³⁰⁹. First, the susceptibility of these agents to germicides *in vitro* is similar to that of other related pathogens. For example, variola is similar to vaccinia ^{72, 310, 311} and *B. anthracis* is similar to *B. atrophaeus* (formerly *B. subtilis*)^{312, 313}. *B. subtilis* spores, for instance, proved as resistant as, if not more resistant than, *B. anthracis* spores (>6 log₁₀ reduction of *B. anthracis* spores in 5 minutes with acidified bleach [5,250 ppm chlorine])³¹³. Thus, one can extrapolate from the larger database available on the susceptibility of genetically similar organisms³¹⁴. Second, many of the potential bioterrorist agents are stable enough in the environment that contaminated environmental surfaces or fomites could lead to transmission of agents such as *B. anthracis*, *F. tularensis*, variola major, *C. botulinum* toxin, and *C. burnetti* ³¹⁵. Third, data suggest that current disinfection and sterilization practices are appropriate for managing patient-care equipment and environmental surfaces when potentially contaminated patients are evaluated and/or admitted in a health-care facility after exposure to a bioterrorist agent. For example, sodium hypochlorite can be used for surface disinfection (see [This link is no longer active:

http://www.epa.gov/pesticides/factsheets/chemicals/bleachfactsheet.htm.]). In instances where the health-care facility is the site of a bioterrorist attack, environmental decontamination might require special decontamination procedures (e.g., chlorine dioxide gas for *B. anthracis* spores). Because no antimicrobial products are registered for decontamination of biologic agents after a bioterrorist attack, EPA has granted a crises exemption for each product (see [This link is no longer active: http://www.epa.gov/pesticides/factsheets/chemicals/bleachfactsheet.htm.]). Of only theoretical concern is the possibility that a bioterrorist agent could be engineered to be less susceptible to disinfection and sterilization processes ³⁰⁹.

Toxicological, Environmental and Occupational Concerns

Health hazards associated with the use of germicides in healthcare vary from mucous membrane irritation to death, with the latter involving accidental injection by mentally disturbed patients³¹⁶. Although their degrees of toxicity vary ³¹⁷⁻³²⁰, all disinfectants should be used with the proper safety precautions ³²¹ and only for the intended purpose.

Key factors associated with assessing the health risk of a chemical exposure include the duration, intensity (i.e., how much chemical is involved), and route (e.g., skin, mucous membranes, and inhalation) of exposure. Toxicity can be acute or chronic. Acute toxicity usually results from an accidental spill of a chemical substance. Exposure is sudden and often produces an emergency situation. Chronic toxicity results from repeated exposure to low levels of the chemical over a prolonged period. Employers are responsible for informing workers about the chemical hazards in the workplace and implementing control measures. The OSHA Hazard Communication Standard (29 CFR 1910.1200, 1915.99, 1917.28, 1918.90, 1926.59, and 1928.21) requires manufacturers and importers of hazardous chemicals to develop Material Safety Data Sheets (MSDS) for each chemical or mixture of chemicals. Employers must have these data sheets readily available to employees who work with the products to which they could be exposed.

Exposure limits have been published for many chemicals used in health care to help provide a safe environment and, as relevant, are discussed in each section of this guideline. Only the exposure limits published by OSHA carry the legal force of regulations. OSHA publishes a limit as a time-weighted average (TWA), that is, the average concentration for a normal 8-hour work day and a 40-hour work week to which nearly all workers can be repeatedly exposed to a chemical without adverse health effects. For example, the permissible exposure limit (PEL) for EtO is 1.0 ppm, 8 hour TWA. The CDC National Institute for Occupational Safety and Health (NIOSH) develops recommended exposure limits (RELs). RELs are occupational exposure limits recommended by NIOSH as being protective of worker health and safety over a working lifetime. This limit is frequently expressed as a 40-hour TWA exposure for up to 10 hours per day during a 40-hour work week. These exposure limits are designed for inhalation exposures. Irritant and allergic effects can occur below the exposure limits, and skin contact can result in dermal effects or systemic absorption without inhalation. The American Conference on Governmental Industrial Hygienists (ACGIN) also provides guidelines on exposure limits ³²². Information about workplace exposures and methods to reduce them (e.g., work practices, engineering controls, PPE) is available on the OSHA (https://www.osha.gov/) and NIOSH (https://www.cdc.gov/niosh/) websites.

Some states have excluded or limited concentrations of certain chemical germicides (e.g., glutaraldehyde, formaldehyde, and some phenols) from disposal through the sewer system. These rules are intended to minimize environmental harm. If health-care facilities exceed the maximum allowable concentration of a chemical (e.g., \geq 5.0 mg/L), they have three options. First, they can switch to alternative products; for example, they can change from glutaraldehyde to another disinfectant for high-level disinfection or from phenolics to quaternary ammonium compounds for low-level disinfection. Second, the health-care facility can collect the disinfectant and dispose of it as a hazardous chemical. Third, the facility can use a commercially available small-scale treatment method (e.g., neutralize

glutaraldehyde with glycine).

Safe disposal of regulated chemicals is important throughout the medical community. For disposal of large volumes of spent solutions, users might decide to neutralize the microbicidal activity before disposal (e.g., glutaraldehyde). Solutions can be neutralized by reaction with chemicals such as sodium bisulfite ^{323, 324} or glycine ³²⁵.

European authors have suggested that instruments and ventilation therapy equipment should be disinfected by heat rather than by chemicals. The concerns for chemical disinfection include toxic side effects for the patient caused by chemical residues on the instrument or object, occupational exposure to toxic chemicals, and recontamination by rinsing the disinfectant with microbially contaminated tap water ³²⁶.

Disinfection in Ambulatory Care, Home Care, and the Home

With the advent of managed healthcare, increasing numbers of patients are now being cared for in ambulatory-care and home settings. Many patients in these settings might have communicable diseases, immunocompromising conditions, or invasive devices. Therefore, adequate disinfection in these settings is necessary to provide a safe patient environment. Because the ambulatory-care setting (i.e., outpatient facility) provides the same risk for infection as the hospital, the Spaulding classification scheme described in this guideline should be followed (Table 1).¹⁷

The home environment should be much safer than hospitals or ambulatory care. Epidemics should not be a problem, and cross-infection should be rare. The healthcare provider is responsible for providing the responsible family member information about infection-control procedures to follow in the home, including hand hygiene, proper cleaning and disinfection of equipment, and safe storage of cleaned and disinfected devices. Among the products recommended for home disinfection of reusable objects are bleach, alcohol, and hydrogen peroxide. APIC recommends that reusable objects (e.g., tracheostomy tubes) that touch mucous membranes be disinfected by immersion in 70% isopropyl alcohol for 5 minutes or in 3% hydrogen peroxide for 30 minutes. Additionally, a 1:50 dilution of 5.25%–6.15% sodium hypochlorite (household bleach) for 5 minutes should be effective ³²⁷⁻³²⁹. Noncritical items (e.g., blood pressure cuffs, crutches) can be cleaned with a detergent. Blood spills should be handled according to OSHA regulations as previously described (see section on OSHA Bloodborne Pathogen Standard). In general, sterilization of critical items is not practical in homes but theoretically could be accomplished by chemical sterilants or boiling. Single-use disposable items can be used or reusable items sterilized in a hospital ^{330, 331}.

Some environmental groups advocate "environmentally safe" products as alternatives to commercial germicides in the home-care setting. These alternatives (e.g., ammonia, baking soda, vinegar, Borax, liquid detergent) are not registered with EPA and should not be used for disinfecting because they are ineffective against *S. aureus*. Borax, baking soda, and detergents also are ineffective against *Salmonella* Typhi and *E.coli*; however, undiluted vinegar and ammonia are effective against *S*. Typhi and *E.coli*^{53, 332, 333}. Common commercial disinfectants designed for home use also are effective against selected antibiotic-resistant bacteria ⁵³.

Public concerns have been raised that the use of antimicrobials in the home can promote development of antibiotic-resistant bacteria ^{334, 335}. This issue is unresolved and needs to be considered further through scientific and clinical investigations. The public health benefits of using disinfectants in the home are unknown. However, some facts are known: many sites in the home kitchen and bathroom are microbially contaminated ³³⁶, use of hypochlorites markedly reduces bacteria ³³⁷, and good standards of hygiene (e.g., food hygiene, hand hygiene) can help reduce infections in the home ^{338, 339}. In addition, laboratory studies indicate that many commercially prepared household disinfectants are effective against common pathogens ⁵³ and can interrupt surface-to-human transmission of pathogens

⁴⁸. The "targeted hygiene concept"—which means identifying situations and areas (e.g., foodpreparation surfaces and bathroom) where risk exists for transmission of pathogens—may be a reasonable way to identify when disinfection might be appropriate ³⁴⁰.

Susceptibility of Antibiotic-Resistant Bacteria to Disinfectants

As with antibiotics, reduced susceptibility (or acquired "resistance") of bacteria to disinfectants can arise by either chromosomal gene mutation or acquisition of genetic material in the form of plasmids or transposons ^{338, 341-343, 344, 345, 346}. When changes occur in bacterial susceptibility that renders an antibiotic ineffective against an infection previously treatable by that antibiotic, the bacteria are referred to as "resistant." In contrast, reduced susceptibility to disinfectants does not correlate with failure of the disinfectant because concentrations used in disinfection still greatly exceed the cidal level. Thus, the word "resistance" when applied to these changes is incorrect, and the preferred term is "reduced susceptibility" or "increased tolerance"^{344, 347}. No data are available that show that antibiotic-resistant bacteria are less sensitive to the liquid chemical germicides than antibiotic-sensitive bacteria at currently used germicide contact conditions and concentrations.

MRSA and vancomycin-resistant *Enterococcus* (VRE) are important health-care–associated agents. Some antiseptics and disinfectants have been known for years to be, because of MICs, somewhat less inhibitory to *S. aureus* strains that contain a plasmid-carrying gene encoding resistance to the antibiotic gentamicin³⁴⁴. For example, gentamicin resistance has been shown to also encode reduced susceptibility to propamidine, quaternary ammonium compounds, and ethidium bromide ³⁴⁸, and MRSA strains have been found to be less susceptible than methicillin-sensitive *S. aureus* (MSSA) strains to chlorhexidine, propamidine, and the quaternary ammonium compound cetrimide ³⁴⁹. In other studies, MRSA and MSSA strains have been equally sensitive to phenols and chlorhexidine, but MRSA strains were slightly more tolerant to quaternary ammonium compounds ³⁵⁰. Two gene families (*qac*CD [now referred to as *smr*] and *qac*AB) are involved in providing protection against agents that are components of disinfectant formulations such as quaternary ammonium compounds. Staphylococci have been proposed to evade destruction because the protein specified by the *qac*A determinant is a cytoplasmic-membrane–associated protein involved in an efflux system that actively reduces intracellular accumulation of toxicants, such as quaternary ammonium compounds, to intracellular targets. ³⁵¹

Other studies demonstrated that plasmid-mediated formaldehyde tolerance is transferable from *Serratia marcescens* to *E. coli*³⁵² and plasmid-mediated quaternary ammonium tolerance is transferable from *S. aureus* to *E. coli*.³⁵³. Tolerance to mercury and silver also is plasmid borne. ^{341, 343-346}

Because the concentrations of disinfectants used in practice are much higher than the MICs observed, even for the more tolerant strains, the clinical relevance of these observations is questionable. Several studies have found antibiotic-resistant hospital strains of common healthcare-associated pathogens (i.e., *Enterococcus, P. aeruginosa, Klebsiella pneumoniae, E. coli, S. aureus,* and *S. epidermidis*) to be equally susceptible to disinfectants as antibiotic-sensitive strains ^{53, 354-356}. The susceptibility of glycopeptide-intermediate *S. aureus* was similar to vancomycin-susceptible, MRSA ³⁵⁷. On the basis of these data, routine disinfection and housekeeping protocols do not need to be altered because of antibiotic resistance provided the disinfection method is effective ^{358, 359}. A study that evaluated the efficacy of selected cleaning methods (e.g., QUAT-sprayed cloth, and QUAT-immersed cloth) for eliminating VRE found that currently used disinfection must involve contact with all contaminated surfaces ³⁵⁸. A new method using an invisible flurorescent marker to objectively evaluate the thoroughness of cleaning activities in patient rooms might lead to improvement in cleaning of all objects and surfaces but needs further evaluation. ³⁶⁰

Lastly, does the use of antiseptics or disinfectants facilitate the development of disinfectant-tolerant organisms? Evidence and reviews indicate enhanced tolerance to disinfectants can be developed in

response to disinfectant exposure ^{334, 335, 346, 347, 361}. However, the level of tolerance is not important in clinical terms because it is low and unlikely to compromise the effectiveness of disinfectants of which much higher concentrations are used ^{347, 362}.

The issue of whether low-level tolerance to germicides selects for antibiotic-resistant strains is unsettled but might depend on the mechanism by which tolerance is attained. For example, changes in the permeability barrier or efflux mechanisms might affect susceptibility to both antibiotics and germicides, but specific changes to a target site might not. Some researchers have suggested that use of disinfectants or antiseptics (e.g., triclosan) could facilitate development of antibiotic-resistant microorganisms ^{334, 335, 363}. Although evidence in laboratory studies indicates low-level resistance to triclosan, the concentrations of triclosan in these studies were low (generally <1 µg/mL) and dissimilar from the higher levels used in antimicrobial products (2,000–20,000 µg/mL) ^{364, 365}. Thus, researchers can create laboratory-derived mutants that demonstrate reduced susceptibility to antiseptics or disinfectants. In some experiments, such bacteria have demonstrated reduced susceptibility to certain antibiotics 335. There is no evidence that using antiseptics or disinfectants selects for antibiotic-resistant organisms in nature or that such mutants survive in nature³⁶⁶.). In addition, the action of antibiotics and the action of disinfectants differ fundamentally. Antibiotics are selectively toxic and generally have a single target site in bacteria, thereby inhibiting a specific biosynthetic process. Germicides generally are considered nonspecific antimicrobials because of a multiplicity of toxic-effect mechanisms or target sites and are broader spectrum in the types of microorganisms against which they are effective ^{344, 347}.

The rotational use of disinfectants in some environments (e.g., pharmacy production units) has been recommended and practiced in an attempt to prevent development of resistant microbes ^{367, 368}. There have been only rare case reports that appropriately used disinfectants have resulted in a clinical problem arising from the selection or development of nonsusceptible microorganisms ³⁶⁹.

Surface Disinfection

The effective use of disinfectants is part of a multibarrier strategy to prevent health-care–associated infections. Surfaces are considered noncritical items because they contact intact skin. Use of noncritical items or contact with noncritical surfaces carries little risk of causing an infection in patients or staff. Thus, the routine use of germicidal chemicals to disinfect hospital floors and other noncritical items is controversial ³⁷⁰⁻³⁷⁵. A 1991 study expanded the Spaulding scheme by dividing the noncritical environmental surfaces into housekeeping surfaces and medical equipment surfaces ³⁷⁶. The classes of disinfectants used on housekeeping and medical equipment surfaces can be similar. However, the frequency of decontaminating can vary (see Recommendations). Medical equipment surfaces (e.g., blood pressure cuffs, stethoscopes, hemodialysis machines, and X-ray machines) can become contaminated with infectious agents and contribute to the spread of health-care–associated infections ^{248, 375}. For this reason, noncritical medical equipment surfaces should be disinfected with an EPA-registered low- or intermediate-level disinfectant. Use of a disinfectant will provide antimicrobial activity that is likely to be achieved with minimal additional cost or work.

Environmental surfaces (e.g., bedside table) also could potentially contribute to cross-transmission by contamination of health-care personnel from hand contact with contaminated surfaces, medical equipment, or patients ^{50, 375, 377}. A paper reviews the epidemiologic and microbiologic data (Table 3) regarding the use of disinfectants on noncritical surfaces ³⁷⁸.

Of the seven reasons to usie a disinfectant on noncritical surfaces, five are particularly noteworthy and support the use of a germicidal detergent. First, hospital floors become contaminated with microorganisms from settling airborne bacteria: by contact with shoes, wheels, and other objects; and occasionally by spills. The removal of microbes is a component in controling health-care–associated infections. In an investigation of the cleaning of hospital floors, the use of soap and water (80% reduction) was less effective in reducing the numbers of bacteria than was a phenolic disinfectant (94%–

99.9% reduction) ³⁷⁹. However, a few hours after floor disinfection, the bacterial count was nearly back to the pretreatment level. Second, detergents become contaminated and result in seeding the patient's environment with bacteria. Investigators have shown that mop water becomes increasingly dirty during cleaning and becomes contaminated if soap and water is used rather than a disinfectant. For example, in one study, bacterial contamination in soap and water without a disinfectant increased from 10 CFU/mL to 34,000 CFU/mL after cleaning a ward, whereas contamination in a disinfectant solution did not change (20 CFU/mL) ³⁸⁰. Contamination of surfaces close to the patient that are frequently touched by the patient or staff (e.g., bed rails) could result in patient exposures0 ³⁸¹. In a study, using of detergents on floors and patient room furniture, increased bacterial contamination of the patients' environmental surfaces was found after cleaning (average increase = 103.6 CFU/24cm²) ³⁸². In addition, a P. aeruginosa outbreak was reported in a hematology-oncology unit associated with contamination of the surface cleaning equipment when nongermicidal cleaning solutions instead of disinfectants were used to decontaminate the patients' environment ³⁸³ and another study demonstrated the role of environmental cleaning in controlling an outbreak of Acinetobacter baumannii 384. Studies also have shown that, in situations where the cleaning procedure failed to eliminate contamination from the surface and the cloth is used to wipe another surface, the contamination is transferred to that surface and the hands of the person holding the cloth^{381, 385}. Third, the CDC Isolation Guideline recommends that noncritical equipment contaminated with blood, body fluids, secretions, or excretions be cleaned and disinfected after use. The same guideline recommends that, in addition to cleaning, disinfection of the bedside equipment and environmental surfaces (e.g., bedrails, bedside tables, carts, commodes, door-knobs, and faucet handles) is indicated for certain pathogens, e.g., enterococci, which can survive in the inanimate environment for prolonged periods ³⁸⁶. Fourth, OSHA requires that surfaces contaminated with blood and other potentially infectious materials (e.g., amniotic, pleural fluid) be disinfected. Fifth, using a single product throughout the facility can simplify both training and appropriate practice.

Reasons also exist for using a detergent alone on floors because noncritical surfaces contribute minimally to endemic health-care-associated infections ³⁸⁷, and no differences have been found in healthcare-associated infections rates when floors are cleaned with detergent rather than disinfectant ^{382, 388, 389}. However, these studies have been small and of short duration and suffer from low statistical power because the outcome-healthcare-associated infections-is of low frequency. The low rate of infections makes the efficacy of an intervention statistically difficult to demonstrate. Because housekeeping surfaces are associated with the lowest risk for disease transmission, some researchers have suggested that either detergents or a disinfectant/detergent could be used ³⁷⁶. No data exist that show reduced healthcare-associated infection rates with use of surface disinfection of floors, but some data demonstrate reduced microbial load associated with the use of disinfectants. Given this information; other information showing that environmental surfaces (e.g., bedside table, bed rails) close to the patient and in outpatient settings ³⁹⁰ can be contaminated with epidemiologically important microbes (such as VRE and MRSA)^{47, 390-394}; and data showing these organisms survive on various hospital surfaces ^{395, 396}: some researchers have suggested that such surfaces should be disinfected on a regular schedule ³⁷⁸. Spot decontamination on fabrics that remain in hospitals or clinic rooms while patients move in and out (e.g., privacy curtains) also should be considered. One study demonstrated the effectiveness of spraying the fabric with 3% hydrogen peroxide ³⁹⁷. Future studies should evaluate the level of contamination on noncritical environmental surfaces as a function of high and low hand contact and whether some surfaces (e.g., bed rails) near the patient with high contact frequencies require more frequent disinfection. Regardless of whether a detergent or disinfectant is used on surfaces in a health-care facility, surfaces should be cleaned routinely and when dirty or soiled to provide an aesthetically pleasing environment and to prevent potentially contaminated objects from serving as a source for health-careassociated infections. ³⁹⁸ The value of designing surfaces (e.g. hexyl-polyvinylpyridine) that kill bacteria on contact ³⁹⁹ or have sustained antimicrobial activity ⁴⁰⁰ should be further evaluated.

Several investigators have recognized heavy microbial contamination of wet mops and cleaning cloths and the potential for spread of such contamination ^{68, 401}. They have shown that wiping hard surfaces with contaminated cloths can contaminate hands, equipment, and other surfaces ^{68, 402}. Data have been published that can be used to formulate effective policies for decontamination and

Guideline for Disinfection and Sterilization in Healthcare Facilities (2008)

maintenance of reusable cleaning cloths. For example, heat was the most reliable treatment of cleaning cloths as a detergent washing followed by drying at 80°C for 2 hours produced elimination of contamination. However, the dry heating process might be a fire hazard if the mop head contains petroleum-based products or lint builds up within the equipment or vent hose (American Health Care Association, personal communication, March 2003). Alternatively, immersing the cloth in hypochlorite (4,000 ppm) for 2 minutes produced no detectable surviving organisms in 10 of 13 cloths ⁴⁰³. If reusable cleaning cloths or mops are used, they should be decontaminated regularly to prevent surface contamination during cleaning with subsequent transfer of organisms from these surfaces to patients or equipment by the hands of health-care workers. Some hospitals have begun using a new mopping technique involving microfiber materials to clean floors. Microfibers are densely constructed, polyester and polyamide (nylon) fibers, that are approximately 1/16 the thickness of a human hair. The positively charged microfibers attract dust (which has a negative charge) and are more absorbent than a conventional, cotton-loop mop. Microfiber materials also can be wet with disinfectants, such as quaternary ammonium compounds. In one study, the microfiber system tested demonstrated superior microbial removal compared with conventional string mops when used with a detergent cleaner (94% vs 68%). The use of a disinfectant did not improve the microbial elimination demonstrated by the microfiber system (95% vs 94%). However, use of disinfectant significantly improved microbial removal when a conventional string mop was used (95% vs 68%) (WA Rutala, unpublished data, August 2006). The microfiber system also prevents the possibility of transferring microbes from room to room because a new microfiber pad is used in each room.

An important issue concerning use of disinfectants for noncritical surfaces in health-care settings is that the contact time specified on the label of the product is often too long to be practically followed. The labels of most products registered by EPA for use against HBV, HIV, or *M. tuberculosis* specify a contact time of 10 minutes. Such a long contact time is not practical for disinfection of environmental surfaces in a health-care setting because most health-care facilities apply a disinfectant and allow it to dry (~1 minute). Multiple scientific papers have demonstrated significant microbial reduction with contact times of 30 to 60 seconds^{46-56, 58-64}. In addition, EPA will approve a shortened contact time for any product for which the manufacturers will submit confirmatory efficacy data.

Currently, some EPA-registered disinfectants have contact times of one to three minutes. By law, users must follow all applicable label instructions for EPA-registered products. Ideally, product users should consider and use products that have the shortened contact time. However, disinfectant manufacturers also need to obtain EPA approval for shortened contact times so these products will be used correctly and effectively in the health-care environment.

Air Disinfection

Disinfectant spray-fog techniques for antimicrobial control in hospital rooms has been used. This technique of spraying of disinfectants is an unsatisfactory method of decontaminating air and surfaces and is not recommended for general infection control in routine patient-care areas³⁸⁶. Disinfectant fogging is rarely, if ever, used in U.S. healthcare facilities for air and surface disinfection in patient-care areas. Methods (e.g., filtration, ultraviolet germicidal irradiation, chlorine dioxide) to reduce air contamination in the healthcare setting are discussed in another guideline ²³.

Microbial Contamination of Disinfectants

Contaminated disinfectants and antiseptics have been occasional vehicles of health-care infections and pseudoepidemics for more than 50 years. Published reports describing contaminated disinfectants and antiseptic solutions leading to health-care-associated infections have been summarized ⁴⁰⁴. Since this summary additional reports have been published ⁴⁰⁵⁻⁴⁰⁸. An examination of reports of disinfectants contaminated with microorganisms revealed noteworthy observations. Perhaps most importantly, high-

Guideline for Disinfection and Sterilization in Healthcare Facilities (2008)

level disinfectants/liquid chemical sterilants have not been associated with outbreaks due to intrinsic or extrinsic contamination.Members of the genus *Pseudomonas* (e.g., *P. aeruginosa*) are the most frequent isolates from contaminated disinfectants—recovered from 80% of contaminated products. Their ability to remain viable or grow in use-dilutions of disinfectants is unparalleled. This survival advantage for *Pseudomonas* results presumably from their nutritional versatility, their unique outer membrane that constitutes an effective barrier to the passage of germicides, and/or efflux systems ⁴⁰⁹. Although the concentrated solutions of the disinfectants have not been demonstrated to be contaminated at the point of manufacture, an undiluted phenolic can be contaminated by a *Pseudomonas* sp. during use ⁴¹⁰. In most of the reports that describe illness associated with contaminated disinfectants, the product was used to disinfect patient-care equipment, such as cystoscopes, cardiac catheters, and thermometers. Germicides used as disinfectants that were reported to have been contaminated include chlorhexidine, quaternary ammonium compounds, phenolics, and pine oil.

The following control measures should be instituted to reduce the frequency of bacterial growth in disinfectants and the threat of serious healthcare–associated infections from the use of such contaminated products ⁴⁰⁴. First, some disinfectants should not be diluted; those that are diluted must be prepared correctly to achieve the manufacturers' recommended use-dilution. Second, infection-control professionals must learn from the literature what inappropriate activities result in extrinsic contamination (i.e., at the point of use) of germicides and train users to prevent recurrence. Common sources of extrinsic contamination of germicides in the reviewed literature are the water to make working dilutions, contaminated containers, and general contamination of the hospital areas where the germicides are prepared and/or used. Third, stock solutions of germicides must be stored as indicated on the product label. EPA verifies manufacturers' efficacy claims against microorganisms. These measures should provide assurance that products meeting the EPA registration requirements can achieve a certain level of antimicrobial activity when used as directed.

Factors Affecting the Efficacy of Disinfection and Sterilization

The activity of germicides against microorganisms depends on a number of factors, some of which are intrinsic qualities of the organism, others of which are the chemical and external physical environment. Awareness of these factors should lead to better use of disinfection and sterilization processes and will be briefly reviewed. More extensive consideration of these and other factors is available elsewhere ^{13, 14, 16, 411-413}.

Number and Location of Microorganisms

All other conditions remaining constant, the larger the number of microbes, the more time a germicide needs to destroy all of them. Spaulding illustrated this relation when he employed identical test conditions and demonstrated that it took 30 minutes to kill 10 *B. atrophaeus* (formerly *Bacillus subtilis*) spores but 3 hours to kill 100,000 *Bacillus atrophaeus* spores. This reinforces the need for scrupulous cleaning of medical instruments before disinfection and sterilization. Reducing the number of microorganisms that must be inactivated through meticulous cleaning, increases the margin of safety when the germicide is used according to the labeling and shortens the exposure time required to kill the entire microbial load. Researchers also have shown that aggregated or clumped cells are more difficult to inactivate than monodispersed cells⁴¹⁴.

The location of microorganisms also must be considered when factors affecting the efficacy of germicides are assessed. Medical instruments with multiple pieces must be disassembled and equipment such as endoscopes that have crevices, joints, and channels are more difficult to disinfect than are flatsurface equipment because penetration of the disinfectant of all parts of the equipment is more difficult. Only surfaces that directly contact the germicide will be disinfected, so there must be no air pockets and the equipment must be completely immersed for the entire exposure period. Manufacturers should be encouraged to produce equipment engineered for ease of cleaning and disinfection.

Innate Resistance of Microorganisms

Microorganisms vary greatly in their resistance to chemical germicides and sterilization processes (Figure 1) ³⁴² Intrinsic resistance mechanisms in microorganisms to disinfectants vary. For example, spores are resistant to disinfectants because the spore coat and cortex act as a barrier, mycobacteria have a waxy cell wall that prevents disinfectant entry, and gram-negative bacteria possess an outer membrane that acts as a barrier to the uptake of disinfectants ^{341, 343-345}. Implicit in all disinfection strategies is the consideration that the most resistant microbial subpopulation controls the sterilization or disinfection time. That is, to destroy the most resistant types of microorganisms (i.e., bacterial spores), the user needs to employ exposure times and a concentration of germicide needed to achieve complete destruction. Except for prions, bacterial spores possess the highest innate resistance to chemical germicides, followed by coccidia (e.g., Cryptosporidium), mycobacteria (e.g., M. tuberculosis), nonlipid or small viruses (e.g., poliovirus, and coxsackievirus), fungi (e.g., Aspergillus, and Candida), vegetative bacteria (e.g., Staphylococcus, and Pseudomonas) and lipid or medium-size viruses (e.g., herpes, and HIV). The germicidal resistance exhibited by the gram-positive and gram-negative bacteria is similar with some exceptions (e.g., P. aeruginosa which shows greater resistance to some disinfectants) ^{369, 415, 416}. P. aeruginosa also is significantly more resistant to a variety of disinfectants in its "naturally occurring" state than are cells subcultured on laboratory media ^{415, 417}. *Rickettsiae*, *Chlamydiae*, and mycoplasma cannot be placed in this scale of relative resistance because information about the efficacy of germicides against these agents is limited ⁴¹⁸. Because these microorganisms contain lipid and are similar in structure and composition to other bacteria, they can be predicted to be inactivated by the same germicides that destroy lipid viruses and vegetative bacteria. A known exception to this supposition is Coxiella burnetti, which has demonstrated resistance to disinfectants ⁴¹⁹.

Concentration and Potency of Disinfectants

With other variables constant, and with one exception (iodophors), the more concentrated the disinfectant, the greater its efficacy and the shorter the time necessary to achieve microbial kill. Generally not recognized, however, is that all disinfectants are not similarly affected by concentration adjustments. For example, quaternary ammonium compounds and phenol have a concentration exponent of 1 and 6, respectively; thus, halving the concentration of a quaternary ammonium compound requires doubling its disinfecting time, but halving the concentration of a phenol solution requires a 64-fold (i.e., 2⁶) increase in its disinfecting time ^{365, 413, 420}.

Considering the length of the disinfection time, which depends on the potency of the germicide, also is important. This was illustrated by Spaulding who demonstrated using the mucin-loop test that 70% isopropyl alcohol destroyed 10⁴ *M. tuberculosis* in 5 minutes, whereas a simultaneous test with 3% phenolic required 2–3 hours to achieve the same level of microbial kill ¹⁴.

Physical and Chemical Factors

Several physical and chemical factors also influence disinfectant procedures: temperature, pH, relative humidity, and water hardness. For example, the activity of most disinfectants increases as the temperature increases, but some exceptions exist. Furthermore, too great an increase in temperature causes the disinfectant to degrade and weakens its germicidal activity and thus might produce a potential health hazard.

An increase in pH improves the antimicrobial activity of some disinfectants (e.g., glutaraldehyde, quaternary ammonium compounds) but decreases the antimicrobial activity of others (e.g., phenols, hypochlorites, and iodine). The pH influences the antimicrobial activity by altering the disinfectant molecule or the cell surface ⁴¹³.

Relative humidity is the single most important factor influencing the activity of gaseous disinfectants/sterilants, such as EtO, chlorine dioxide, and formaldehyde.

Water hardness (i.e., high concentration of divalent cations) reduces the rate of kill of certain disinfectants because divalent cations (e.g., magnesium, calcium) in the hard water interact with the disinfectant to form insoluble precipitates ^{13, 421}.

Organic and Inorganic Matter

Organic matter in the form of serum, blood, pus, or fecal or lubricant material can interfere with the antimicrobial activity of disinfectants in at least two ways. Most commonly, interference occurs by a chemical reaction between the germicide and the organic matter resulting in a complex that is less germicidal or nongermicidal, leaving less of the active germicide available for attacking microorganisms. Chlorine and iodine disinfectants, in particular, are prone to such interaction. Alternatively, organic material can protect microorganisms from attack by acting as a physical barrier ^{422, 423}.

The effects of inorganic contaminants on the sterilization process were studied during the 1950s and 1960s ^{424, 425}. These and other studies show the protection by inorganic contaminants of microorganisms to all sterilization processes results from occlusion in salt crystals ^{426, 427}. This further emphasizes the importance of meticulous cleaning of medical devices before any sterilization or disinfection procedure because both organic and inorganic soils are easily removed by washing ⁴²⁶.

Duration of Exposure

Items must be exposed to the germicide for the appropriate minimum contact time. Multiple investigators have demonstrated the effectiveness of low-level disinfectants against vegetative bacteria (e.g., *Listeria, E. coli, Salmonella*, VRE, MRSA), yeasts (e.g., *Candida*), mycobacteria (e.g., *M. tuberculosis*), and viruses (e.g., poliovirus) at exposure times of 30–60 seconds ⁴⁶⁻⁶⁴. By law, all

applicable label instructions on EPA-registered products must be followed. If the user selects exposure conditions that differ from those on the EPA-registered product label, the user assumes liability for any injuries resulting from off-label use and is potentially subject to enforcement action under the Federal Insecticide, Fungicide, and Rodenticide Act (FIFRA).

All lumens and channels of endoscopic instruments must contact the disinfectant. Air pockets interfere with the disinfection process, and items that float on the disinfectant will not be disinfected. The disinfectant must be introduced reliably into the internal channels of the device. The exact times for disinfecting medical items are somewhat elusive because of the effect of the aforementioned factors on disinfection efficacy. Certain contact times have proved reliable (Table 1), but, in general, longer contact times are more effective than shorter contact times.

Biofilms

Microorganisms may be protected from disinfectants by production of thick masses of cells ⁴²⁸ and extracellular materials, or biofilms ⁴²⁹⁻⁴³⁵. Biofilms are microbial communities that are tightly attached to surfaces and cannot be easly removed. Once these masses form, microbes within them can be resistant to disinfectants by multiple mechanisms, including physical characteristics of older biofilms, genotypic variation of the bacteria, microbial production of neutralizing enzymes, and physiologic gradients within the biofilm (e.g., pH). Bacteria within biofilms are up to 1,000 times more resistant to antimicrobials than are the same bacteria in suspension ⁴³⁶. Although new decontamination methods ⁴³⁷ are being investigated for removing biofilms, chlorine and monochloramines can effectively inactivate biofilm bacteria ⁴³¹ ⁴³⁸. Investigators have hypothesized that the glycocalyx-like cellular masses on the interior walls of polyvinyl chloride pipe would protect embedded organisms from some disinfectants and be a reservoir for continuous contamination ^{429, 430, 439}. Biofilms have been found in whirlpools ⁴⁴⁰, dental unit waterlines⁴⁴¹, and numerous medical devices (e.g., contact lenses, pacemakers, hemodialysis systems, urinary catheters, central venous catheters, endoscopes) 434, 436, 438, 442. Their presence can have serious implications for immunocompromised patients and patients who have indwelling medical devices. Some enzymes ^{436, 443, 444} and detergents ⁴³⁶ can degrade biofilms or reduce numbers of viable bacteria within a biofilm, but no products are EPA-registered or FDA-cleared for this purpose.

Cleaning

Cleaning is the removal of foreign material (e.g., soil, and organic material) from objects and is normally accomplished using water with detergents or enzymatic products. Thorough cleaning is required before high-level disinfection and sterilization because inorganic and organic materials that remain on the surfaces of instruments interfere with the effectiveness of these processes. Also, if soiled materials dry or bake onto the instruments, the removal process becomes more difficult and the disinfection or sterilization process less effective or ineffective. Surgical instruments should be presoaked or rinsed to prevent drying of blood and to soften or remove blood from the instruments.

Cleaning is done manually in use areas without mechanical units (e.g., ultrasonic cleaners or washerdisinfectors) or for fragile or difficult-to-clean instruments. With manual cleaning, the two essential components are friction and fluidics. Friction (e.g., rubbing/scrubbing the soiled area with a brush) is an old and dependable method. Fluidics (i.e., fluids under pressure) is used to remove soil and debris from internal channels after brushing and when the design does not allow passage of a brush through a channel ⁴⁴⁵. When a washer-disinfector is used, care should be taken in loading instruments: hinged instruments should be opened fully to allow adequate contact with the detergent solution; stacking of instruments in washers should be avoided; and instruments should be disassembled as much as possible.

The most common types of mechanical or automatic cleaners are ultrasonic cleaners, washerdecontaminators, washer-disinfectors, and washer-sterilizers. Ultrasonic cleaning removes soil by cavitation and implosion in which waves of acoustic energy are propagated in aqueous solutions to disrupt the bonds that hold particulate matter to surfaces. Bacterial contamination can be present in used ultrasonic cleaning solutions (and other used detergent solutions) because these solutions generally do not make antibacterial label claims ⁴⁴⁶. Even though ultrasound alone does not significantly inactivate bacteria, sonication can act synergistically to increase the cidal efficacy of a disinfectant 447. Users of ultrasonic cleaners should be aware that the cleaning fluid could result in endotoxin contamination of surgical instruments, which could cause severe inflammatory reactions ⁴⁴⁸. Washer-sterilizers are modified steam sterilizers that clean by filling the chamber with water and detergent through which steam passes to provide agitation. Instruments are subsequently rinsed and subjected to a short steamsterilization cycle. Another washer-sterilizer employs rotating spray arms for a wash cycle followed by a steam sterilization cycle at 285°F ^{449, 450}. Washer-decontaminators/disinfectors act like a dishwasher that uses a combination of water circulation and detergents to remove soil. These units sometimes have a cycle that subjects the instruments to a heat process (e.g., 93°C for 10 minutes) ⁴⁵¹. Washer-disinfectors are generally computer-controlled units for cleaning, disinfecting, and drying solid and hollow surgical and medical equipment. In one study, cleaning (measured as 5-6 log₁₀ reduction) was achieved on surfaces that had adequate contact with the water flow in the machine ⁴⁵². Detailed information about cleaning and preparing supplies for terminal sterilization is provided by professional organizations ^{453, 454} and books ⁴⁵⁵. Studies have shown that manual and mechanical cleaning of endoscopes achieves approximately a 4log₁₀ reduction of contaminating organisms ^{83, 104, 456, 457}. Thus, cleaning alone effectively reduces the number of microorganisms on contaminated equipment. In a quantitative analysis of residual protein contamination of reprocessed surgical instruments, median levels of residual protein contamination per instrument for five trays were 267, 260, 163, 456, and 756 µg ⁴⁵⁸. In another study, the median amount of protein from reprocessed surgical instruments from different hospitals ranged from 8 µg to 91 µg ⁴⁵⁹. When manual methods were compared with automated methods for cleaning reusable accessory devices used for minimally invasive surgical procedures, the automated method was more efficient for cleaning biopsy forceps and ported and nonported laparoscopic devices and achieved a >99% reduction in soil parameters (i.e., protein, carbohydrate, hemoglobin) in the ported and nonported laparoscopic devices ^{460,} . 461

For instrument cleaning, a neutral or near-neutral pH detergent solution commonly is used because such solutions generally provide the best material compatibility profile and good soil removal. Enzymes,

usually proteases, sometimes are added to neutral pH solutions to assist in removing organic material. Enzymes in these formulations attack proteins that make up a large portion of common soil (e.g., blood, pus). Cleaning solutions also can contain lipases (enzymes active on fats) and amylases (enzymes active on starches). Enzymatic cleaners are not disinfectants, and proteinaceous enzymes can be inactivated by germicides. As with all chemicals, enzymes must be rinsed from the equipment or adverse reactions (e.g., fever, residual amounts of high-level disinfectants, proteinaceous residue) could result ^{462, 463}. Enzyme solutions should be used in accordance with manufacturer's instructions, which include proper dilution of the enzymatic detergent and contact with equipment for the amount of time specified on the label 463. Detergent enzymes can result in asthma or other allergic effects in users. Neutral pH detergent solutions that contain enzymes are compatible with metals and other materials used in medical instruments and are the best choice for cleaning delicate medical instruments, especially flexible endoscopes ⁴⁵⁷. Alkalinebased cleaning agents are used for processing medical devices because they efficiently dissolve protein and fat residues ⁴⁶⁴; however, they can be corrosive ⁴⁵⁷. Some data demonstrate that enzymatic cleaners are more effective than neutral detergents ^{465, 466} in removing microorganisms from surfaces but two more recent studies found no difference in cleaning efficiency between enzymatic and alkaline-based cleaners ^{443, 464}. Another study found no significant difference between enzymatic and non-enzymatic cleaners in terms of microbial cleaning efficacy ⁴⁶⁷. A new non-enzyme, hydrogen peroxide-based formulation (not FDA-cleared) was as effective as enzymatic cleaners in removing protein, blood, carbohydrate, and endotoxin from surface test carriers⁴⁶⁸ In addition, this product effected a 5-log₁₀ reduction in microbial loads with a 3-minute exposure at room temperature ⁴⁶⁸.

Although the effectiveness of high-level disinfection and sterilization mandates effective cleaning, no "real-time" tests exist that can be employed in a clinical setting to verify cleaning. If such tests were commercially available they could be used to ensure an adequate level of cleaning.⁴⁶⁹⁻⁴⁷² The only way to ensure adequate cleaning is to conduct a reprocessing verification test (e.g., microbiologic sampling), but this is not routinely recommended ⁴⁷³. Validation of the cleaning processes in a laboratory-testing program is possible by microorganism detection, chemical detection for organic contaminants, radionuclide tagging, and chemical detection for specific ions ^{426, 471}. During the past few years, data have been published describing use of an artificial soil, protein, endotoxin, X-ray contrast medium, or blood to verify the manual or automated cleaning process ^{169, 452, 474-478} and adenosine triphosphate bioluminescence and microbiologic sampling to evaluate the effectiveness of environmental surface cleaning^{170, 479}. At a minimum, all instruments should be individually inspected and be visibly clean.

Disinfection

Many disinfectants are used alone or in combinations (e.g., hydrogen peroxide and peracetic acid) in the health-care setting. These include alcohols, chlorine and chlorine compounds, formaldehyde, glutaraldehyde, *ortho*-phthalaldehyde, hydrogen peroxide, iodophors, peracetic acid, phenolics, and quaternary ammonium compounds. Commercial formulations based on these chemicals are considered unique products and must be registered with EPA or cleared by FDA. In most instances, a given product is designed for a specific purpose and is to be used in a certain manner. Therefore, users should read labels carefully to ensure the correct product is selected for the intended use and applied efficiently.

Disinfectants are not interchangeable, and incorrect concentrations and inappropriate disinfectants can result in excessive costs. Because occupational diseases among cleaning personnel have been associated with use of several disinfectants (e.g., formaldehyde, glutaraldehyde, and chlorine), precautions (e.g., gloves and proper ventilation) should be used to minimize exposure ^{318, 480, 481}. Asthma and reactive airway disease can occur in sensitized persons exposed to any airborne chemical, including germicides. Clinically important asthma can occur at levels below ceiling levels regulated by OSHA or recommended by NIOSH. The preferred method of control is elimination of the chemical (through engineering controls or substitution) or relocation of the worker.

The following overview of the performance characteristics of each provides users with sufficient information to select an appropriate disinfectant for any item and use it in the most efficient way.

Chemical Disinfectants

Alcohol

Overview. In the healthcare setting, "alcohol" refers to two water-soluble chemical compounds—ethyl alcohol and isopropyl alcohol—that have generally underrated germicidal characteristics ⁴⁸². FDA has not cleared any liquid chemical sterilant or high-level disinfectant with alcohol as the main active ingredient. These alcohols are rapidly bactericidal rather than bacteriostatic against vegetative forms of bacteria; they also are tuberculocidal, fungicidal, and virucidal but do not destroy bacterial spores. Their cidal activity drops sharply when diluted below 50% concentration, and the optimum bactericidal concentration is 60%–90% solutions in water (volume/volume) ^{483, 484}.

Mode of Action. The most feasible explanation for the antimicrobial action of alcohol is denaturation of proteins. This mechanism is supported by the observation that absolute ethyl alcohol, a dehydrating agent, is less bactericidal than mixtures of alcohol and water because proteins are denatured more quickly in the presence of water ^{484, 485}. Protein denaturation also is consistent with observations that alcohol destroys the dehydrogenases of *Escherichia coli* ⁴⁸⁶, and that ethyl alcohol increases the lag phase of *Enterobacter aerogenes* ⁴⁸⁷ and that the lag phase effect could be reversed by adding certain amino acids. The bacteriostatic action was believed caused by inhibition of the production of metabolites essential for rapid cell division.

Microbicidal Activity. Methyl alcohol (methanol) has the weakest bactericidal action of the alcohols and thus seldom is used in healthcare ⁴⁸⁸. The bactericidal activity of various concentrations of ethyl alcohol (ethanol) was examined against a variety of microorganisms in exposure periods ranging from 10 seconds to 1 hour ⁴⁸³. *Pseudomonas aeruginosa* was killed in 10 seconds by all concentrations of ethanol from 30% to 100% (v/v), and *Serratia marcescens, E, coli* and *Salmonella typhosa* were killed in 10 seconds by all concentrations of ethanol from 40% to 100%. The gram-positive organisms *Staphylococcus aureus* and *Streptococcus pyogenes* were slightly more resistant, being killed in 10 seconds by ethyl alcohol concentrations of 60%–95%. Isopropyl alcohol (isopropanol) was slightly more bactericidal than ethyl alcohol for *E. coli* and *S. aureus* ⁴⁸⁹.

Ethyl alcohol, at concentrations of 60%–80%, is a potent virucidal agent inactivating all of the lipophilic viruses (e.g., herpes, vaccinia, and influenza virus) and many hydrophilic viruses (e.g.,

adenovirus, enterovirus, rhinovirus, and rotaviruses but not hepatitis A virus (HAV) ⁵⁸ or poliovirus) ⁴⁹. Isopropyl alcohol is not active against the nonlipid enteroviruses but is fully active against the lipid viruses ⁷². Studies also have demonstrated the ability of ethyl and isopropyl alcohol to inactivate the hepatitis B virus(HBV) ^{224, 225} and the herpes virus, ⁴⁹⁰ and ethyl alcohol to inactivate human immunodeficiency virus (HIV) ²²⁷, rotavirus, echovirus, and astrovirus ⁴⁹¹.

In tests of the effect of ethyl alcohol against *M. tuberculosis*, 95% ethanol killed the tubercle bacilli in sputum or water suspension within 15 seconds ⁴⁹². In 1964, Spaulding stated that alcohols were the germicide of choice for tuberculocidal activity, and they should be the standard by which all other tuberculocides are compared. For example, he compared the tuberculocidal activity of iodophor (450 ppm), a substituted phenol (3%), and isopropanol (70%/volume) using the mucin-loop test (10⁶ *M. tuberculosis* per loop) and determined the contact times needed for complete destruction were 120–180 minutes, 45–60 minutes, and 5 minutes, respectively. The mucin-loop test is a severe test developed to produce long survival times. Thus, these figures should not be extrapolated to the exposure times needed when these germicides are used on medical or surgical material ⁴⁸².

Ethyl alcohol (70%) was the most effective concentration for killing the tissue phase of *Cryptococcus neoformans, Blastomyces dermatitidis, Coccidioides immitis,* and *Histoplasma capsulatum* and the culture phases of the latter three organisms aerosolized onto various surfaces. The culture phase was more resistant to the action of ethyl alcohol and required about 20 minutes to disinfect the contaminated surface, compared with <1 minute for the tissue phase ^{493, 494}.

Isopropyl alcohol (20%) is effective in killing the cysts of *Acanthamoeba culbertsoni* (560) as are chlorhexidine, hydrogen peroxide, and thimerosal ⁴⁹⁶.

Uses. Alcohols are not recommended for sterilizing medical and surgical materials principally because they lack sporicidal action and they cannot penetrate protein-rich materials. Fatal postoperative wound infections with *Clostridium* have occurred when alcohols were used to sterilize surgical instruments contaminated with bacterial spores ⁴⁹⁷. Alcohols have been used effectively to disinfect oral and rectal thermometers^{498, 499}, hospital pagers ⁵⁰⁰, scissors ⁵⁰¹, and stethoscopes ⁵⁰². Alcohols have been used to disinfect fiberoptic endoscopes ^{503, 504} but failure of this disinfectant have lead to infection ^{280, 505}. Alcohol towelettes have been used for years to disinfect small surfaces such as rubber stoppers of multiple-dose medication vials or vaccine bottles. Furthermore, alcohol occasionally is used to disinfect external surfaces of equipment (e.g., stethoscopes, ventilators, manual ventilation bags) ⁵⁰⁶, CPR manikins ⁵⁰⁷, ultrasound instruments ⁵⁰⁸ or medication preparation areas. Two studies demonstrated the effectiveness of 70% isopropyl alcohol to disinfect reusable transducer heads in a controlled environment ^{509, 510}. In contrast, three bloodstream infection outbreaks have been described when alcohol was used to disinfect transducer heads in an intensive-care setting ⁵¹¹.

The documented shortcomings of alcohols on equipment are that they damage the shellac mountings of lensed instruments, tend to swell and harden rubber and certain plastic tubing after prolonged and repeated use, bleach rubber and plastic tiles ⁴⁸² and damage tonometer tips (by deterioration of the glue) after the equivalent of 1 working year of routine use ⁵¹². Tonometer biprisms soaked in alcohol for 4 days developed rough front surfaces that potentially could cause corneal damage; this appeared to be caused by weakening of the cementing substances used to fabricate the biprisms ⁵¹³. Corneal opacification has been reported when tonometer tips were swabbed with alcohol immediately before measurement of intraocular pressure ⁵¹⁴. Alcohols are flammable and consequently must be stored in a cool, well-ventilated area. They also evaporate rapidly, making extended exposure time difficult to achieve unless the items are immersed.

Chlorine and Chlorine Compounds

Overview. Hypochlorites, the most widely used of the chlorine disinfectants, are available as liquid (e.g., sodium hypochlorite) or solid (e.g., calcium hypochlorite). The most prevalent chlorine products in

the United States are aqueous solutions of 5.25%-6.15% sodium hypochlorite (see glossary), usually called household bleach. They have a broad spectrum of antimicrobial activity, do not leave toxic residues, are unaffected by water hardness, are inexpensive and fast acting ³²⁸, remove dried or fixed organisms and biofilms from surfaces⁴⁶⁵, and have a low incidence of serious toxicity ⁵¹⁵⁻⁵¹⁷. Sodium hypochlorite at the concentration used in household bleach (5.25-6.15%) can produce ocular irritation or oropharyngeal, esophageal, and gastric burns ^{318, 518-522}. Other disadvantages of hypochlorites include corrosiveness to metals in high concentrations (>500 ppm), inactivation by organic matter, discoloring or "bleaching" of fabrics, release of toxic chlorine gas when mixed with ammonia or acid (e.g., household cleaning agents) ⁵²³⁻⁵²⁵, and relative stability ³²⁷. The microbicidal activity of chlorine is attributed largely to undissociated hypochlorous acid (HOCI). The dissociation of HOCI to the less microbicidal form (hypochlorite ion OCI) depends on pH. The disinfecting efficacy of chlorine decreases with an increase in pH that parallels the conversion of undissociated HOCI to OCI-329, 526. A potential hazard is production of the carcinogen bis(chloromethyl) ether when hypochlorite solutions contact formaldehyde 527 and the production of the animal carcinogen trihalomethane when hot water is hyperchlorinated ⁵²⁸. After reviewing environmental fate and ecologic data, EPA has determined the currently registered uses of hypochlorites will not result in unreasonable adverse effects to the environment ⁵²⁹.

Alternative compounds that release chlorine and are used in the health-care setting include demandrelease chlorine dioxide, sodium dichloroisocyanurate, and chloramine-T. The advantage of these compounds over the hypochlorites is that they retain chlorine longer and so exert a more prolonged bactericidal effect. Sodium dichloroisocyanurate tablets are stable, and for two reasons, the microbicidal activity of solutions prepared from sodium dichloroisocyanurate tablets might be greater than that of sodium hypochlorite solutions containing the same total available chlorine. First, with sodium dichloroisocyanurate, only 50% of the total available chlorine is free (HOCI and OCI-), whereas the remainder is combined (monochloroisocyanurate or dichloroisocyanurate), and as free available chlorine is used up, the latter is released to restore the equilibrium. Second, solutions of sodium dichloroisocvanurate are acidic, whereas sodium hypochlorite solutions are alkaline, and the more microbicidal type of chlorine (HOCI) is believed to predominate ⁵³⁰⁻⁵³³. Chlorine dioxide-based disinfectants are prepared fresh as required by mixing the two components (base solution [citric acid with preservatives and corrosion inhibitors] and the activator solution [sodium chlorite]). In vitro suspension tests showed that solutions containing about 140 ppm chlorine dioxide achieved a reduction factor exceeding 10⁶ of S. aureus in 1 minute and of Bacillus atrophaeus spores in 2.5 minutes in the presence of 3 g/L bovine albumin. The potential for damaging equipment requires consideration because long-term use can damage the outer plastic coat of the insertion tube ⁵³⁴. In another study, chlorine dioxide solutions at either 600 ppm or 30 ppm killed Mycobacterium avium-intracellulare within 60 seconds after contact but contamination by organic material significantly affected the microbicidal properties⁵³⁵.

The microbicidal activity of a new disinfectant, "superoxidized water," has been examined The concept of electrolyzing saline to create a disinfectant or antiseptics is appealing because the basic materials of saline and electricity are inexpensive and the end product (i.e., water) does not damage the environment. The main products of this water are hypochlorous acid (e.g., at a concentration of about 144 mg/L) and chlorine. As with any germicide, the antimicrobial activity of superoxidized water is strongly affected by the concentration of the active ingredient (available free chlorine) ⁵³⁶. One manufacturer generates the disinfectant at the point of use by passing a saline solution over coated titanium electrodes at 9 amps. The product generated has a pH of 5.0-6.5 and an oxidation-reduction potential (redox) of >950 mV. Although superoxidized water is intended to be generated fresh at the point of use, when tested under clean conditions the disinfectant was effective within 5 minutes when 48 hours old ⁵³⁷. Unfortunately, the equipment required to produce the product can be expensive because parameters such as pH, current, and redox potential must be closely monitored. The solution is nontoxic to biologic tissues. Although the United Kingdom manufacturer claims the solution is noncorrosive and nondamaging to endoscopes and processing equipment, one flexible endoscope manufacturer (Olympus Key-Med, United Kingdom) has voided the warranty on the endoscopes if superoxidized water is used to disinfect them ⁵³⁸. As with any germicide formulation, the user should check with the device

manufacturer for compatibility with the germicide. Additional studies are needed to determine whether this solution could be used as an alternative to other disinfectants or antiseptics for hand washing, skin antisepsis, room cleaning, or equipment disinfection (e.g., endoscopes, dialyzers) ^{400, 539, 540}. In October 2002, the FDA cleared superoxidized water as a high-level disinfectant (FDA, personal communication, September 18, 2002).

Mode of Action. The exact mechanism by which free chlorine destroys microorganisms has not been elucidated. Inactivation by chlorine can result from a number of factors: oxidation of sulfhydryl enzymes and amino acids; ring chlorination of amino acids; loss of intracellular contents; decreased uptake of nutrients; inhibition of protein synthesis; decreased oxygen uptake; oxidation of respiratory components; decreased adenosine triphosphate production; breaks in DNA; and depressed DNA synthesis ^{329, 347}. The actual microbicidal mechanism of chlorine might involve a combination of these factors or the effect of chlorine on critical sites ³⁴⁷.

Microbicidal Activity. Low concentrations of free available chlorine (e.g., HOCI, OCI⁻, and elemental chlorine-Cl₂) have a biocidal effect on mycoplasma (25 ppm) and vegetative bacteria (<5 ppm) in seconds in the absence of an organic load ^{329, 418}. Higher concentrations (1,000 ppm) of chlorine are required to kill *M. tuberculosis* using the Association of Official Analytical Chemists (AOAC) tuberculocidal test ⁷³. A concentration of 100 ppm will kill ≥99.9% of *B. atrophaeus* spores within 5 minutes ^{541, 542} and destroy mycotic agents in <1 hour ³²⁹. Acidified bleach and regular bleach (5,000 ppm chlorine) can inactivate 10⁶ *Clostridium difficile* spores in ≤10 minutes ²⁶². One study reported that 25 different viruses were inactivated in 10 minutes with 200 ppm available chlorine ⁷². Several studies have demonstrated the effectiveness of diluted sodium hypochlorite and other disinfectants to inactivate HIV ⁶¹. Chlorine (500 ppm) showed inhibition of *Candida* after 30 seconds of exposure ⁵⁴. In experiments using the AOAC Use-Dilution Method, 100 ppm of free chlorine killed 10⁶–10⁷ *S. aureus*, *Salmonella choleraesuis*, and *P. aeruginosa* in <10 minutes ³²⁷. Because household bleach contains 5.25%–6.15% sodium hypochlorite, or 52,500–61,500 ppm available chlorine, a 1:1,000 dilution provides about 53–62 ppm available chlorine, and a 1:10 dilution of household bleach provides about 5250–6150 ppm.

Data are available for chlorine dioxide that support manufacturers' bactericidal, fungicidal, sporicidal, tuberculocidal, and virucidal label claims ⁵⁴³⁻⁵⁴⁶. A chlorine dioxide generator has been shown effective for decontaminating flexible endoscopes ⁵³⁴ but it is not currently FDA-cleared for use as a high-level disinfectant ⁸⁵. Chlorine dioxide can be produced by mixing solutions, such as a solution of chlorine with a solution of sodium chlorite ³²⁹. In 1986, a chlorine dioxide product was voluntarily removed from the market when its use caused leakage of cellulose-based dialyzer membranes, which allowed bacteria to migrate from the dialysis fluid side of the dialyzer to the blood side ⁵⁴⁷.

Sodium dichloroisocyanurate at 2,500 ppm available chlorine is effective against bacteria in the presence of up to 20% plasma, compared with 10% plasma for sodium hypochlorite at 2,500 ppm ⁵⁴⁸.

"Superoxidized water" has been tested against bacteria, mycobacteria, viruses, fungi, and spores ^{537, 539, 549}. Freshly generated superoxidized water is rapidly effective (<2 minutes) in achieving a 5-log₁₀ reduction of pathogenic microorganisms (i.e., *M. tuberculosis, M. chelonae*, poliovirus, HIV, multidrug-resistant *S. aureus, E. coli, Candida albicans, Enterococcus faecalis, P. aeruginosa*) in the absence of organic loading. However, the biocidal activity of this disinfectant decreased substantially in the presence of organic material (e.g., 5% horse serum) ^{537, 549, 550}. No bacteria or viruses were detected on artificially contaminated endoscopes after a 5-minute exposure to superoxidized water ⁵⁵¹ and HBV-DNA was not detected from any endoscope experimentally contaminated with HBV-positive mixed sera after a disinfectant exposure time of 7 minutes⁵⁵².

Uses. Hypochlorites are widely used in healthcare facilities in a variety of settings. ³²⁸ Inorganic chlorine solution is used for disinfecting tonometer heads ¹⁸⁸ and for spot-disinfection of countertops and floors. A 1:10–1:100 dilution of 5.25%–6.15% sodium hypochlorite (i.e., household bleach) ^{22, 228, 553, 554} or

an EPA-registered tuberculocidal disinfectant ¹⁷has been recommended for decontaminating blood spills. For small spills of blood (i.e., drops of blood) on noncritical surfaces, the area can be disinfected with a 1:100 dilution of 5.25%-6.15% sodium hypochlorite or an EPA-registered tuberculocidal disinfectant. Because hypochlorites and other germicides are substantially inactivated in the presence of blood ^{63, 548,} ^{555, 556}, large spills of blood require that the surface be cleaned before an EPA-registered disinfectant or a 1:10 (final concentration) solution of household bleach is applied ⁵⁵⁷. If a sharps injury is possible, the surface initially should be decontaminated ^{69, 318}, then cleaned and disinfected (1:10 final concentration) ⁶³. Extreme care always should be taken to prevent percutaneous injury. At least 500 ppm available chlorine for 10 minutes is recommended for decontaminating CPR training manikins ⁵⁵⁸. Full-strength bleach has been recommended for self-disinfection of needles and syringes used for illicit-drug injection when needle-exchange programs are not available. The difference in the recommended concentrations of bleach reflects the difficulty of cleaning the interior of needles and syringes and the use of needles and syringes for parenteral injection ⁵⁵⁹. Clinicians should not alter their use of chlorine on environmental surfaces on the basis of testing methodologies that do not simulate actual disinfection practices ^{560, 561}. Other uses in healthcare include as an irrigating agent in endodontic treatment ⁵⁶² and as a disinfectant for manikins, laundry, dental appliances, hydrotherapy tanks ^{23, 41}, regulated medical waste before disposal ³²⁸, and the water distribution system in hemodialysis centers and hemodialysis machines ⁵⁶³.

Chlorine long has been used as the disinfectant in water treatment. Hyperchlorination of a *Legionella*contaminated hospital water system ²³ resulted in a dramatic decrease (from 30% to 1.5%) in the isolation of *L. pneumophila* from water outlets and a cessation of healthcare-associated Legionnaires' disease in an affected unit ^{528, 564}. Water disinfection with monochloramine by municipal water-treatment plants substantially reduced the risk for healthcare–associated Legionnaires disease ^{565, 566}. Chlorine dioxide also has been used to control *Legionella* in a hospital water supply. ⁵⁶⁷ Chloramine T ⁵⁶⁸ and hypochlorites ⁴¹ have been used to disinfect hydrotherapy equipment.

Hypochlorite solutions in tap water at a pH >8 stored at room temperature (23°C) in closed, opaque plastic containers can lose up to 40%–50% of their free available chlorine level over 1 month. Thus, if a user wished to have a solution containing 500 ppm of available chlorine at day 30, he or she should prepare a solution containing 1,000 ppm of chlorine at time 0. Sodium hypochlorite solution does not decompose after 30 days when stored in a closed brown bottle ³²⁷.

The use of powders, composed of a mixture of a chlorine-releasing agent with highly absorbent resin, for disinfecting spills of body fluids has been evaluated by laboratory tests and hospital ward trials. The inclusion of acrylic resin particles in formulations markedly increases the volume of fluid that can be soaked up because the resin can absorb 200–300 times its own weight of fluid, depending on the fluid consistency. When experimental formulations containing 1%, 5%, and 10% available chlorine were evaluated by a standardized surface test, those containing 10% demonstrated bactericidal activity. One problem with chlorine-releasing granules is that they can generate chlorine fumes when applied to urine ⁵⁶⁹.

Formaldehyde

Overview. Formaldehyde is used as a disinfectant and sterilant in both its liquid and gaseous states. Liquid formaldehyde will be considered briefly in this section, and the gaseous form is reviewed elsewhere ⁵⁷⁰. Formaldehyde is sold and used principally as a water-based solution called formalin, which is 37% formaldehyde by weight. The aqueous solution is a bactericide, tuberculocide, fungicide, virucide and sporicide ^{72, 82, 571-573}. OSHA indicated that formaldehyde should be handled in the workplace as a potential carcinogen and set an employee exposure standard for formaldehyde that limits an 8-hour time-weighted average exposure concentration of 0.75 ppm ^{574, 575}. The standard includes a second permissible exposure limit in the form of a short-term exposure limit (STEL) of 2 ppm that is the maximum exposure allowed during a 15-minute period ⁵⁷⁶. Ingestion of formaldehyde can be fatal, and long-term exposure to low levels in the air or on the skin can cause asthma-like respiratory problems and skin irritation, such as dermatitis and itching. For these reasons, employees should have limited direct contact

with formaldehyde, and these considerations limit its role in sterilization and disinfection processes. Key provisions of the OSHA standard that protects workers from exposure to formaldehyde appear in Title 29 of the Code of Federal Regulations (CFR) Part 1910.1048 (and equivalent regulations in states with OSHA-approved state plans) ⁵⁷⁷.

Mode of Action. Formaldehyde inactivates microorganisms by alkylating the amino and sulfhydral groups of proteins and ring nitrogen atoms of purine bases ³⁷⁶.

Microbicidal Activity. Varying concentrations of aqueous formaldehyde solutions destroy a wide range of microorganisms. Inactivation of poliovirus in 10 minutes required an 8% concentration of formalin, but all other viruses tested were inactivated with 2% formalin ⁷². Four percent formaldehyde is a tuberculocidal agent, inactivating 10⁴ *M. tuberculosis* in 2 minutes ⁸², and 2.5% formaldehyde inactivated about 10⁷ *Salmonella* Typhi in 10 minutes in the presence of organic matter ⁵⁷². The sporicidal action of formaldehyde was slower than that of glutaraldehyde in comparative tests with 4% aqueous formaldehyde and 2% glutaraldehyde against the spores of *B. anthracis* ⁸². The formaldehyde solution required 2 hours of contact to achieve an inactivation factor of 10⁴, whereas glutaraldehyde required only 15 minutes.

Uses. Although formaldehyde-alcohol is a chemical sterilant and formaldehyde is a high-level disinfectant, the health-care uses of formaldehyde are limited by its irritating fumes and its pungent odor even at very low levels (<1 ppm). For these reasons and others—such as its role as a suspected human carcinogen linked to nasal cancer and lung cancer ⁵⁷⁸, this germicide is excluded from Table 1. When it is used, , direct exposure to employees generally is limited; however, excessive exposures to formaldehyde have been documented for employees of renal transplant units ^{574, 579}, and students in a gross anatomy laboratory ⁵⁸⁰. Formaldehyde is used in the health-care setting to prepare viral vaccines (e.g., poliovirus and influenza); as an embalming agent; and to preserve anatomic specimens; and historically has been used to sterilize surgical instruments, especially when mixed with ethanol. A 1997 survey found that formaldehyde was used for reprocessing hemodialyzers by 34% of U.S. hemodialysis centers—a 60% decrease from 1983 ^{249, 581}. If used at room temperature, a concentration of 4% with a minimum exposure of 24 hours is required to disinfect disposable hemodialyzers reused on the same patient ^{582, 583}. Aqueous formaldehyde solutions (1%–2%) also have been used to disinfect the internal fluid pathways of dialysis machines ⁵⁸³. To minimize a potential health hazard to dialysis patients, the dialysis equipment must be thoroughly rinsed and tested for residual formaldehyde before use.

Paraformaldehyde, a solid polymer of formaldehyde, can be vaporized by heat for the gaseous decontamination of laminar flow biologic safety cabinets when maintenance work or filter changes require access to the sealed portion of the cabinet.

Glutaraldehyde

Overview. Glutaraldehyde is a saturated dialdehyde that has gained wide acceptance as a high-level disinfectant and chemical sterilant ¹⁰⁷. Aqueous solutions of glutaraldehyde are acidic and generally in this state are not sporicidal. Only when the solution is "activated" (made alkaline) by use of alkalinating agents to pH 7.5–8.5 does the solution become sporicidal. Once activated, these solutions have a shelf-life of minimally 14 days because of the polymerization of the glutaraldehyde molecules at alkaline pH levels. This polymerization blocks the active sites (aldehyde groups) of the glutaraldehyde molecules that are responsible for its biocidal activity.

Novel glutaraldehyde formulations (e.g., glutaraldehyde-phenol-sodium phenate, potentiated acid glutaraldehyde, stabilized alkaline glutaraldehyde) produced in the past 30 years have overcome the problem of rapid loss of activity (e.g., use-life 28–30 days) while generally maintaining excellent microbicidal activity ⁵⁸⁴⁻⁵⁸⁸. However, antimicrobial activity depends not only on age but also on use conditions, such as dilution and organic stress. Manufacturers' literature for these preparations suggests the neutral or alkaline glutaraldehydes possess microbicidal and anticorrosion properties superior to

those of acid glutaraldehydes, and a few published reports substantiate these claims ^{542, 589, 590}. However, two studies found no difference in the microbicidal activity of alkaline and acid glutaraldehydes ^{73, 591}. The use of glutaraldehyde-based solutions in health-care facilities is widespread because of their advantages, including excellent biocidal properties; activity in the presence of organic matter (20% bovine serum); and noncorrosive action to endoscopic equipment, thermometers, rubber, or plastic equipment (Tables 4 and 5).

Mode of Action. The biocidal activity of glutaraldehyde results from its alkylation of sulfhydryl, hydroxyl, carboxyl, and amino groups of microorganisms, which alters RNA, DNA, and protein synthesis. The mechanism of action of glutaraldehydes are reviewed extensively elsewhere ^{592, 593}.

Microbicidal Activity. The in vitro inactivation of microorganisms by glutaraldehydes has been extensively investigated and reviewed ^{592, 593}. Several investigators showed that $\geq 2\%$ aqueous solutions of glutaraldehyde, buffered to pH 7.5–8.5 with sodium bicarbonate effectively killed vegetative bacteria in <2 minutes; *M. tuberculosis*, fungi, and viruses in <10 minutes; and spores of *Bacillus* and *Clostridium* species in 3 hours ^{542, 592-597}. Spores of *C. difficile* are more rapidly killed by 2% glutaraldehyde than are spores of other species of *Clostridium* and *Bacillus* ^{79, 265, 266}. Microorganisms with substantial resistance to glutaraldehyde have been reported, including some mycobacteria (*M. chelonae, Mycobacterium avium-intracellulare, M. xenopi*) ⁵⁹⁸⁻⁶⁰¹, *Methylobacterium mesophilicum* ⁶⁰², *Trichosporon*, fungal ascospores (e.g., *Microascus cinereus*, *Cheatomium globosum*), and *Cryptosporidium*^{271, 603}. *M. chelonae* persisted in a 0.2% glutaraldehyde solution used to store porcine prosthetic heart valves ⁶⁰⁴.

Two percent alkaline glutaraldehyde solution inactivated 10⁵ M. tuberculosis cells on the surface of penicylinders within 5 minutes at 18°C 589. However, subsequent studies⁸² questioned the mycobactericidal prowess of glutaraldehydes. Two percent alkaline glutaraldehyde has slow action (20 to >30 minutes) against *M. tuberculosis* and compares unfavorably with alcohols, formaldehydes, iodine, and phenol⁸². Suspensions of *M. avium*, *M. intracellulare*, and *M. gordonae* were more resistant to inactivation by a 2% alkaline glutaraldehyde (estimated time to complete inactivation: ~60 minutes) than were virulent *M. tuberculosis* (estimated time to complete inactivation ~25 minutes) ⁶⁰⁵. The rate of kill was directly proportional to the temperature, and a standardized suspension of *M. tuberculosis* could not be sterilized within 10 minutes ⁸⁴. An FDA-cleared chemical sterilant containing 2.5% glutaraldehyde uses increased temperature (35°C) to reduce the time required to achieve high-level disinfection (5 minutes)^{85,} ⁶⁰⁶, but its use is limited to automatic endoscope reprocessors equipped with a heater. In another study employing membrane filters for measurement of mycobactericidal activity of 2% alkaline glutaraldehyde, complete inactivation was achieved within 20 minutes at 20°C when the test inoculum was 10⁶ M. tuberculosis per membrane ⁸¹. Several investigators ^{55, 57, 73, 76, 80, 81, 84, 605} have demonstrated that glutaraldehyde solutions inactivate 2.4 to >5.0 log₁₀ of *M. tuberculosis* in 10 minutes (including multidrugresistant *M. tuberculosis*) and 4.0-6.4 log₁₀ of *M. tuberculosis* in 20 minutes. On the basis of these data and other studies, 20 minutes at room temperature is considered the minimum exposure time needed to reliably kill Mycobacteria and other vegetative bacteria with ≥2% glutaraldehvde ^{17, 19, 27, 57, 83, 94, 108, 111, 117-} 121, 607

Glutaraldehyde is commonly diluted during use, and studies showed a glutaraldehyde concentration decline after a few days of use in an automatic endoscope washer ^{608, 609}. The decline occurs because instruments are not thoroughly dried and water is carried in with the instrument, which increases the solution's volume and dilutes its effective concentration ⁶¹⁰. This emphasizes the need to ensure that semicritical equipment is disinfected with an acceptable concentration of glutaraldehyde. Data suggest that 1.0%–1.5% glutaraldehyde is the minimum effective concentration for >2% glutaraldehyde solutions when used as a high-level disinfectant ^{76, 589, 590, 609}. Chemical test strips or liquid chemical monitors ^{610, 611} are available for determining whether an effective concentration of glutaraldehyde is present despite repeated use and dilution. The frequency of testing should be based on how frequently the solutions are used (e.g., used daily, test daily; used weekly, test before use; used 30 times per day, test each 10th use), but the strips should not be used to extend the use life beyond the expiration date. Data suggest the

chemicals in the test strip deteriorate with time ⁶¹² and a manufacturer's expiration date should be placed on the bottles. The bottle of test strips should be dated when opened and used for the period of time indicated on the bottle (e.g., 120 days). The results of test strip monitoring should be documented. The glutaraldehyde test kits have been preliminarily evaluated for accuracy and range ⁶¹² but the reliability has been questioned ⁶¹³. To ensure the presence of minimum effective concentration of the high-level disinfectant, manufacturers of some chemical test strips recommend the use of quality-control procedures to ensure the strips perform properly. If the manufacturer of the chemical test strip recommends a qualitycontrol procedure, users should comply with the manufacturer's recommendations. The concentration should be considered unacceptable or unsafe when the test indicates a dilution below the product's minimum effective concentration (MEC) (generally to $\leq 1.0\%$ –1.5% glutaraldehyde) by the indicator not changing color.

A 2.0% glutaraldehyde–7.05% phenol–1.20% sodium phenate product that contained 0.125% glutaraldehyde–0.44% phenol–0.075% sodium phenate when diluted 1:16 is not recommended as a high-level disinfectant because it lacks bactericidal activity in the presence of organic matter and lacks tuberculocidal, fungicidal, virucidal, and sporicidal activity ^{49, 55, 56, 71, 73-79, 614}. In December 1991, EPA issued an order to stop the sale of all batches of this product because of efficacy data showing the product is not effective against spores and possibly other microorganisms or inanimate objects as claimed on the label ⁶¹⁵. FDA has cleared a glutaraldehyde–phenol/phenate concentrate as a high-level disinfectant that contains 1.12% glutaraldehyde with 1.93% phenol/phenate at its use concentration. Other FDA cleared glutaraldehyde sterilants that contain 2.4%–3.4% glutaraldehyde are used undiluted ⁶⁰⁶.

Uses. Glutaraldehyde is used most commonly as a high-level disinfectant for medical equipment such as endoscopes ^{69, 107, 504}, spirometry tubing, dialyzers ⁶¹⁶, transducers, anesthesia and respiratory therapy equipment ⁶¹⁷, hemodialysis proportioning and dialysate delivery systems ^{249, 618}, and reuse of laparoscopic disposable plastic trocars ⁶¹⁹. Glutaraldehyde is noncorrosive to metal and does not damage lensed instruments, rubber. or plastics. Glutaraldehyde should not be used for cleaning noncritical surfaces because it is too toxic and expensive.

Colitis believed caused by glutaraldehyde exposure from residual disinfecting solution in endoscope solution channels has been reported and is preventable by careful endoscope rinsing ^{318, 620-630}. One study found that residual glutaraldehyde levels were higher and more variable after manual disinfection (<0.2 mg/L to 159.5 mg/L) than after automatic disinfection (0.2–6.3 mg/L)⁶³¹. Similarly, keratopathy and corneal decompensation were caused by ophthalmic instruments that were inadequately rinsed after soaking in 2% glutaraldehyde ^{632, 633}.

Healthcare personnel can be exposed to elevated levels of glutaraldehyde vapor when equipment is processed in poorly ventilated rooms, when spills occur, when glutaraldehyde solutions are activated or changed,⁶³⁴, or when open immersion baths are used. Acute or chronic exposure can result in skin irritation or dermatitis, mucous membrane irritation (eye, nose, mouth), or pulmonary symptoms ^{318, 635-639}. Epistaxis, allergic contact dermatitis, asthma, and rhinitis also have been reported in healthcare workers exposed to glutaraldehyde ^{636, 640-647}.

Glutaraldehyde exposure should be monitored to ensure a safe work environment. Testing can be done by four techniques: a silica gel tube/gas chromatography with a flame ionization detector, dinitrophenylhydrazine (DNPH)-impregnated filter cassette/high-performance liquid chromatography (HPLC) with an ultraviolet (UV) detector, a passive badge/HPLC, or a handheld glutaraldehyde air monitor ⁶⁴⁸. The silica gel tube and the DNPH-impregnated cassette are suitable for monitoring the 0.05 ppm ceiling limit. The passive badge, with a 0.02 ppm limit of detection, is considered marginal at the Americal Council of Governmental Industrial Hygienists (ACGIH) ceiling level. The ceiling level is considered too close to the glutaraldehyde meter's 0.03 ppm limit of detection to provide confidence in the readings ⁶⁴⁸. ACGIH does not require a specific monitoring schedule for glutaraldehyde; however, a monitoring schedule is needed to ensure the level is less than the ceiling limit. For example, monitoring

should be done initially to determine glutaraldehyde levels, after procedural or equipment changes, and in response to worker complaints ⁶⁴⁹. In the absence of an OSHA permissible exposure limit, if the glutaraldehyde level is higher than the ACGIH ceiling limit of 0.05 ppm, corrective action and repeat monitoring would be prudent ⁶⁴⁹.

Engineering and work-practice controls that can be used to resolve these problems include ducted exhaust hoods, air systems that provide 7–15 air exchanges per hour, ductless fume hoods with absorbents for the glutaraldehyde vapor, tight-fitting lids on immersion baths, personal protection (e.g., nitrile or butyl rubber gloves but not natural latex gloves, goggles) to minimize skin or mucous membrane contact, and automated endoscope processors ^{7, 650}. If engineering controls fail to maintain levels below the ceiling limit, institutions can consider the use of respirators (e.g., a half-face respirator with organic vapor cartridge ⁶⁴⁰ or a type "C" supplied air respirator with a full facepiece operated in a positive pressure mode) ⁶⁵¹. In general, engineering controls are preferred over work-practice and administrative controls because they do not require active participation by the health-care worker. Even though enforcement of the OSHA ceiling limit was suspended in 1993 by the U.S. Court of Appeals ⁵⁷⁷, limiting employee exposure to 0.05 ppm (according to ACGIH) is prudent because, at this level, glutaraldehyde can irritate the eyes, throat, and nose ^{318, 577, 639, 652}. If glutaraldehyde disposal through the sanitary sewer system is restricted, sodium bisulfate can be used to neutralize the glutaraldehyde and make it safe for disposal.

Hydrogen Peroxide

Overview. The literature contains several accounts of the properties, germicidal effectiveness, and potential uses for stabilized hydrogen peroxide in the health-care setting. Published reports ascribe good germicidal activity to hydrogen peroxide and attest to its bactericidal, virucidal, sporicidal, and fungicidal properties ⁶⁵³⁻⁶⁵⁵. (Tables 4 and 5) The FDA website lists cleared liquid chemical sterilants and high-level disinfectants containing hydrogen peroxide and their cleared contact conditions.

Mode of Action. Hydrogen peroxide works by producing destructive hydroxyl free radicals that can attack membrane lipids, DNA, and other essential cell components. Catalase, produced by aerobic organisms and facultative anaerobes that possess cytochrome systems, can protect cells from metabolically produced hydrogen peroxide by degrading hydrogen peroxide to water and oxygen. This defense is overwhelmed by the concentrations used for disinfection ^{653, 654}.

Microbicidal Activity. Hydrogen peroxide is active against a wide range of microorganisms, including bacteria, yeasts, fungi, viruses, and spores 78, 654. A 0.5% accelerated hydrogen peroxide demonstrated bactericidal and virucidal activity in 1 minute and mycobactericidal and fungicidal activity in 5 minutes 656. Bactericidal effectiveness and stability of hydrogen peroxide in urine has been demonstrated against a variety of health-care-associated pathogens; organisms with high cellular catalase activity (e.g., S. aureus, S. marcescens, and Proteus mirabilis) required 30-60 minutes of exposure to 0.6% hydrogen peroxide for a 10⁸ reduction in cell counts, whereas organisms with lower catalase activity (e.g., E. coli, Streptococcus species, and Pseudomonas species) required only 15 minutes' exposure 657. In an investigation of 3%, 10%, and 15% hydrogen peroxide for reducing spacecraft bacterial populations, a complete kill of 10⁶ spores (i.e., *Bacillus* species) occurred with a 10% concentration and a 60-minute exposure time. A 3% concentration for 150 minutes killed 10⁶ spores in six of seven exposure trials ⁶⁵⁸. A 10% hydrogen peroxide solution resulted in a 10³ decrease in *B*. atrophaeus spores, and a ≥10⁵ decrease when tested against 13 other pathogens in 30 minutes at 20°C ^{659, 660}. A 3.0% hydrogen peroxide solution was ineffective against VRE after 3 and 10 minutes exposure times ⁶⁶¹ and caused only a 2-log₁₀ reduction in the number of Acanthamoeba cysts in approximately 2 hours ⁶⁶². A 7% stabilized hydrogen peroxide proved to be sporicidal (6 hours of exposure), mycobactericidal (20 minutes), fungicidal (5 minutes) at full strength, virucidal (5 minutes) and bactericidal (3 minutes) at a 1:16 dilution when a quantitative carrier test was used 655. The 7% solution of hydrogen peroxide, tested after 14 days of stress (in the form of germ-loaded carriers and respiratory therapy equipment), was sporicidal (>7 log₁₀ reduction in 6 hours), mycobactericidal (>6.5 log₁₀ reduction in 25

minutes), fungicidal (>5 log₁₀ reduction in 20 minutes), bactericidal (>6 log₁₀ reduction in 5 minutes) and virucidal (5 log₁₀ reduction in 5 minutes) ⁶⁶³. Synergistic sporicidal effects were observed when spores were exposed to a combination of hydrogen peroxide (5.9%–23.6%) and peracetic acid ⁶⁶⁴. Other studies demonstrated the antiviral activity of hydrogen peroxide against rhinovirus ⁶⁶⁵. The time required for inactivating three serotypes of rhinovirus using a 3% hydrogen peroxide solution was 6–8 minutes; this time increased with decreasing concentrations (18-20 minutes at 1.5%, 50–60 minutes at 0.75%).

Concentrations of hydrogen peroxide from 6% to 25% show promise as chemical sterilants. The product marketed as a sterilant is a premixed, ready-to-use chemical that contains 7.5% hydrogen peroxide and 0.85% phosphoric acid (to maintain a low pH) ⁶⁹. The mycobactericidal activity of 7.5% hydrogen peroxide has been corroborated in a study showing the inactivation of >10⁵ multidrug-resistant *M. tuberculosis* after a 10-minute exposure ⁶⁶⁶. Thirty minutes were required for >99.9% inactivation of poliovirus and HAV 667. Three percent and 6% hydrogen peroxide were unable to inactivate HAV in 1 minute in a carrier test ⁵⁸. When the effectiveness of 7.5% hydrogen peroxide at 10 minutes was compared with 2% alkaline glutaraldehyde at 20 minutes in manual disinfection of endoscopes, no significant difference in germicidal activity was observed ⁶⁶⁸.). No complaints were received from the nursing or medical staff regarding odor or toxicity. In one study, 6% hydrogen peroxide (unused product was 7.5%) was more effective in the high-level disinfection of flexible endoscopes than was the 2% glutaraldehyde solution ⁴⁵⁶. A new, rapid-acting 13.4% hydrogen peroxide formulation (that is not yet FDA-cleared) has demonstrated sporicidal, mycobactericidal, fungicidal, and virucidal efficacy. Manufacturer data demonstrate that this solution sterilizes in 30 minutes and provides high-level disinfection in 5 minutes⁶⁶⁹. This product has not been used long enough to evaluate material compatibility to endoscopes and other semicritical devices, and further assessment by instrument manufacturers is needed.

Under normal conditions, hydrogen peroxide is extremely stable when properly stored (e.g., in dark containers). The decomposition or loss of potency in small containers is less than 2% per year at ambient temperatures ⁶⁷⁰.

Uses. Commercially available 3% hydrogen peroxide is a stable and effective disinfectant when used on inanimate surfaces. It has been used in concentrations from 3% to 6% for disinfecting soft contact lenses (e.g., 3% for 2–3 hrs) ^{653, 671, 672}, tonometer biprisms ⁵¹³, ventilators ⁶⁷³, fabrics ³⁹⁷, and endoscopes ⁴⁵⁶. Hydrogen peroxide was effective in spot-disinfecting fabrics in patients' rooms ³⁹⁷. Corneal damage from a hydrogen peroxide-soaked tonometer tip that was not properly rinsed has been reported ⁶⁷⁴. Hydrogen peroxide also has been instilled into urinary drainage bags in an attempt to eliminate the bag as a source of bladder bacteriuria and environmental contamination ⁶⁷⁵. Although the instillation of hydrogen peroxide into the bag reduced microbial contamination of the bag, this procedure did not reduce the incidence of catheter-associated bacteriuria ⁶⁷⁵.

A chemical irritation resembling pseudomembranous colitis caused by either 3% hydrogen peroxide or a 2% glutaraldehyde has been reported ⁶²¹. An epidemic of pseudomembrane-like enteritis and colitis in seven patients in a gastrointestinal endoscopy unit also has been associated with inadequate rinsing of 3% hydrogen peroxide from the endoscope ⁶⁷⁶.

As with other chemical sterilants, dilution of the hydrogen peroxide must be monitored by regularly testing the minimum effective concentration (i.e., 7.5%–6.0%). Compatibility testing by Olympus America of the 7.5% hydrogen peroxide found both cosmetic changes (e.g., discoloration of black anodized metal finishes) ⁶⁹ and functional changes with the tested endoscopes (Olympus, written communication, October 15, 1999).

lodophors

Overview. lodine solutions or tinctures long have been used by health professionals primarily as antiseptics on skin or tissue. lodophors, on the other hand, have been used both as antiseptics and

disinfectants. FDA has not cleared any liquid chemical sterilant or high-level disinfectants with iodophors as the main active ingredient. An iodophor is a combination of iodine and a solubilizing agent or carrier; the resulting complex provides a sustained-release reservoir of iodine and releases small amounts of free iodine in aqueous solution. The best-known and most widely used iodophor is povidone-iodine, a compound of polyvinylpyrrolidone with iodine. This product and other iodophors retain the germicidal efficacy of iodine but unlike iodine generally are nonstaining and relatively free of toxicity and irritancy ^{677, 678}.

Several reports that documented intrinsic microbial contamination of antiseptic formulations of povidone-iodine and poloxamer-iodine ⁶⁷⁹⁻⁶⁸¹ caused a reappraisal of the chemistry and use of iodophors⁶⁸². "Free" iodine (I₂) contributes to the bactericidal activity of iodophors and dilutions of iodophors demonstrate more rapid bactericidal action than does a full-strength povidone-iodine solution. The reason for the observation that dilution increases bactericidal activity is unclear, but dilution of povidone-iodine might weaken the iodine linkage to the carrier polymer with an accompanying increase of free iodine in solution ⁶⁸⁰. Therefore, iodophors must be diluted according to the manufacturers' directions to achieve antimicrobial activity.

Mode of Action. Iodine can penetrate the cell wall of microorganisms quickly, and the lethal effects are believed to result from disruption of protein and nucleic acid structure and synthesis.

Microbicidal Activity. Published reports on the in vitro antimicrobial efficacy of iodophors demonstrate that iodophors are bactericidal, mycobactericidal, and virucidal but can require prolonged contact times to kill certain fungi and bacterial spores ^{14, 71-73, 290, 683-686}. Three brands of povidone-iodine solution have demonstrated more rapid kill (seconds to minutes) of *S. aureus* and *M. chelonae* at a 1:100 dilution than did the stock solution ⁶⁸³. The virucidal activity of 75–150 ppm available iodine was demonstrated against seven viruses ⁷². Other investigators have questioned the efficacy of iodophors against poliovirus in the presence of organic matter ⁶⁸⁵and rotavirus SA-11 in distilled or tapwater ²⁹⁰. Manufacturers' data demonstrate that commercial iodophors are not sporicidal, but they are tuberculocidal, fungicidal, virucidal, and bactericidal at their recommended use-dilution.

Uses. Besides their use as an antiseptic, iodophors have been used for disinfecting blood culture bottles and medical equipment, such as hydrotherapy tanks, thermometers, and endoscopes. Antiseptic iodophors are not suitable for use as hard-surface disinfectants because of concentration differences. Iodophors formulated as antiseptics contain less free iodine than do those formulated as disinfectants ³⁷⁶. Iodine or iodine-based antiseptics should not be used on silicone catheters because they can adversely affect the silicone tubing ⁶⁸⁷.

Ortho-phthalaldehyde (OPA)

Overview. Ortho-phthalaldehyde is a high-level disinfectant that received FDA clearance in October 1999. It contains 0.55% 1,2-benzenedicarboxaldehyde (OPA). OPA solution is a clear, pale-blue liquid with a pH of 7.5. (Tables 4 and 5)

Mode of Action. Preliminary studies on the mode of action of OPA suggest that both OPA and glutaraldehyde interact with amino acids, proteins, and microorganisms. However, OPA is a less potent cross-linking agent. This is compensated for by the lipophilic aromatic nature of OPA that is likely to assist its uptake through the outer layers of mycobacteria and gram-negative bacteria ⁶⁸⁸⁻⁶⁹⁰. OPA appears to kill spores by blocking the spore germination process ⁶⁹¹.

Microbicidal Activity. Studies have demonstrated excellent microbicidal activity in vitro ^{69, 100, 271, 400, 692-703}. For example, OPA has superior mycobactericidal activity (5-log₁₀ reduction in 5 minutes) to glutaraldehyde. The mean times required to produce a 6-log₁₀ reduction for *M. bovis* using 0.21% OPA was 6 minutes, compared with 32 minutes using 1.5% glutaraldehyde ⁶⁹³. OPA showed good activity against the mycobacteria tested, including the glutaraldehyde-resistant strains, but 0.5% OPA was not sporicidal with 270 minutes of exposure. Increasing the pH from its unadjusted level (about 6.5) to pH 8

improved the sporicidal activity of OPA ⁶⁹⁴. The level of biocidal activity was directly related to the temperature. A greater than 5-log₁₀ reduction of *B. atrophaeus* spores was observed in 3 hours at 35°C, than in 24 hours at 20°C. Also, with an exposure time \leq 5 minutes, biocidal activity decreased with increasing serum concentration. However, efficacy did not differ when the exposure time was \geq 10 minutes ⁶⁹⁷. In addition, OPA is effective (>5-log₁₀ reduction) against a wide range of microorganisms, including glutaraldehyde-resistant mycobacteria and *B. atrophaeus* spores ⁶⁹⁴.

The influence of laboratory adaptation of test strains, such as *P. aeruginosa*, to 0.55% OPA has been evaluated. Resistant and multiresistant strains increased substantially in susceptibility to OPA after laboratory adaptation (log₁₀ reduction factors increased by 0.54 and 0.91 for resistant and multiresistant strains, respectively) ⁷⁰⁴. Other studies have found naturally occurring cells of *P. aeurginosa* were more resistant to a variety of disinfectants than were subcultured cells ⁷⁰⁵.

Uses. OPA has several potential advantages over glutaraldehyde. It has excellent stability over a wide pH range (pH 3–9), is not a known irritant to the eyes and nasal passages ⁷⁰⁶, does not require exposure monitoring, has a barely perceptible odor, and requires no activation. OPA, like glutaraldehyde, has excellent material compatibility. A potential disadvantage of OPA is that it stains proteins gray (including unprotected skin) and thus must be handled with caution ⁶⁹. However, skin staining would indicate improper handling that requires additional training and/or personal protective equipment (e.g., gloves, eye and mouth protection, and fluid-resistant gowns). OPA residues remaining on inadequately water-rinsed transesophageal echo probes can stain the patient's mouth ⁷⁰⁷. Meticulous cleaning, using the correct OPA exposure time (e.g., 12 minutes) and copious rinsing of the probe with water should eliminate this problem. The results of one study provided a basis for a recommendation that rinsing of instruments disinfected with OPA will require at least 250 mL of water per channel to reduce the chemical residue to a level that will not compromise patient or staff safety (<1 ppm) ⁷⁰⁸. Personal protective equipment should be worn when contaminated instruments, equipment, and chemicals are handled ⁴⁰⁰. In addition, equipment must be thoroughly rinsed to prevent discoloration of a patient's skin or mucous membrane.

In April 2004, the manufacturer of OPA disseminated information to users about patients who reportedly experienced an anaphylaxis-like reaction after cystoscopy where the scope had been reprocessed using OPA. Of approximately 1 million urologic procedures performed using instruments reprocessed using OPA, 24 cases (17 cases in the United States, six in Japan, one in the United Kingdom) of anaphylaxis-like reactions have been reported after repeated cystoscopy (typically after four to nine treatments). Preventive measures include removal of OPA residues by thorough rinsing and not using OPA for reprocessing urologic instrumentation used to treat patients with a history of bladder cancer (Nevine Erian, personal communication, June 4, 2004; Product Notification, Advanced Sterilization Products, April 23, 2004)⁷⁰⁹.

A few OPA clinical studies are available. In a clinical-use study, OPA exposure of 100 endoscopes for 5 minutes resulted in a >5-log₁₀ reduction in bacterial load. Furthermore, OPA was effective over a 14-day use cycle ¹⁰⁰. Manufacturer data show that OPA will last longer in an automatic endoscope reprocessor before reaching its MEC limit (MEC after 82 cycles) than will glutaraldehyde (MEC after 40 cycles) ⁴⁰⁰. High-pressure liquid chromatography confirmed that OPA levels are maintained above 0.3% for at least 50 cycles ^{706, 710}. OPA must be disposed in accordance with local and state regulations. If OPA disposal through the sanitary sewer system is restricted, glycine (25 grams/gallon) can be used to neutralize the OPA and make it safe for disposal.

The high-level disinfectant label claims for OPA solution at 20°C vary worldwide (e.g., 5 minutes in Europe, Asia, and Latin America; 10 minutes in Canada and Australia; and 12 minutes in the United States). These label claims differ worldwide because of differences in the test methodology and requirements for licensure. In an automated endoscope reprocessor with an FDA-cleared capability to maintain solution temperatures at 25°C, the contact time for OPA is 5 minutes.

Peracetic Acid

Overview. Peracetic, or peroxyacetic, acid is characterized by rapid action against all microorganisms. Special advantages of peracetic acid are that it lacks harmful decomposition products (i.e., acetic acid, water, oxygen, hydrogen peroxide), enhances removal of organic material ⁷¹¹, and leaves no residue. It remains effective in the presence of organic matter and is sporicidal even at low temperatures (Tables 4 and 5). Peracetic acid can corrode copper, brass, bronze, plain steel, and galvanized iron but these effects can be reduced by additives and pH modifications. It is considered unstable, particularly when diluted; for example, a 1% solution loses half its strength through hydrolysis in 6 days, whereas 40% peracetic acid loses 1%–2% of its active ingredients per month ⁶⁵⁴.

Mode of Action. Little is known about the mechanism of action of peracetic acid, but it is believed to function similarly to other oxidizing agents—that is, it denatures proteins, disrupts the cell wall permeability, and oxidizes sulfhydryl and sulfur bonds in proteins, enzymes, and other metabolites ⁶⁵⁴.

Microbicidal Activity. Peracetic acid will inactivate gram-positive and gram-negative bacteria, fungi, and yeasts in ≤ 5 minutes at <100 ppm. In the presence of organic matter, 200–500 ppm is required. For viruses, the dosage range is wide (12–2250 ppm), with poliovirus inactivated in yeast extract in 15 minutes with 1,500–2,250 ppm. In one study, 3.5% peracetic acid was ineffective against HAV after 1-minute exposure using a carrier test ⁵⁸. Peracetic acid (0.26%) was effective (log₁₀ reduction factor >5) against all test strains of mycobacteria (*M. tuberculosis, M. avium-intracellulare, M. chelonae, and M. fortuitum*) within 20–30 minutes in the presence or absence of an organic load ^{607, 712}. With bacterial spores, 500–10,000 ppm (0.05%–1%) inactivates spores in 15 seconds to 30 minutes using a spore suspension test ^{654, 659, 713-715}.

Uses. An automated machine using peracetic acid to chemically sterilize medical (e.g., endoscopes, arthroscopes), surgical, and dental instruments is used in the United States⁷¹⁶⁻⁷¹⁸. As previously noted, dental handpieces should be steam sterilized. The sterilant, 35% peracetic acid, is diluted to 0.2% with filtered water at 50°C. Simulated-use trials have demonstrated excellent microbicidal activity ^{111, 718-722}, and three clinical trials have demonstrated both excellent microbial killing and no clinical failures leading to infection^{90, 723, 724}. The high efficacy of the system was demonstrated in a comparison of the efficacies of the system with that of ethylene oxide. Only the peracetic acid system completely killed 6 log₁₀ of *M*. chelonae, E. faecalis, and B. atrophaeus spores with both an organic and inorganic challenge⁷²². An investigation that compared the costs, performance, and maintenance of urologic endoscopic equipment processed by high-level disinfection (with glutaraldehyde) with those of the peracetic acid system reported no clinical differences between the two systems. However, the use of this system led to higher costs than the high-level disinfection, including costs for processing (\$6.11 vs. \$0.45 per cycle), purchasing and training (\$24,845 vs. \$16), installation (\$5,800 vs. \$0), and endoscope repairs (\$6,037 vs. \$445) ⁹⁰. Furthermore, three clusters of infection using the peracetic acid automated endoscope reprocessor were linked to inadequately processed bronchoscopes when inappropriate channel connectors were used with the system ⁷²⁵. These clusters highlight the importance of training, proper model-specific endoscope connector systems, and quality-control procedures to ensure compliance with endoscope manufacturer recommendations and professional organization guidelines. An alternative highlevel disinfectant available in the United Kingdom contains 0.35% peracetic acid. Although this product is rapidly effective against a broad range of microorganisms ^{466, 726, 727}, it tarnishes the metal of endoscopes and is unstable, resulting in only a 24-hour use life 727.

Peracetic Acid and Hydrogen Peroxide

Overview. Two chemical sterilants are available that contain peracetic acid plus hydrogen peroxide (i.e., 0.08% peracetic acid plus 1.0% hydrogen peroxide [no longer marketed]; and 0.23% peracetic acid plus 7.35% hydrogen peroxide (Tables 4 and 5).

Microbicidal Activity. The bactericidal properties of peracetic acid and hydrogen peroxide have been demonstrated ⁷²⁸. Manufacturer data demonstrated this combination of peracetic acid and hydrogen

peroxide inactivated all microorganisms except bacterial spores within 20 minutes. The 0.08% peracetic acid plus 1.0% hydrogen peroxide product effectively inactivated glutaraldehyde-resistant mycobacteria⁷²⁹.

Uses. The combination of peracetic acid and hydrogen peroxide has been used for disinfecting hemodialyzers ⁷³⁰. The percentage of dialysis centers using a peracetic acid-hydrogen peroxide-based disinfectant for reprocessing dialyzers increased from 5% in 1983 to 56% in 1997²⁴⁹. Olympus America does not endorse use of 0.08% peracetic acid plus 1.0% hydrogen peroxide (Olympus America, personal communication, April 15, 1998) on any Olympus endoscope because of cosmetic and functional damage and will not assume liability for chemical damage resulting from use of this product. This product is not currently available. FDA has cleared a newer chemical sterilant with 0.23% peracetic acid and 7.35% hydrogen peroxide (Tables 4 and 5). After testing the 7.35% hydrogen peroxide and 0.23% peracetic acid product, Olympus America concluded it was not compatible with the company's flexible gastrointestinal endoscopes; this conclusion was based on immersion studies where the test insertion tubes had failed because of swelling and loosening of the black polymer layer of the tube (Olympus America, personal communication, September 13, 2000).

Phenolics

Overview. Phenol has occupied a prominent place in the field of hospital disinfection since its initial use as a germicide by Lister in his pioneering work on antiseptic surgery. In the past 30 years, however, work has concentrated on the numerous phenol derivatives or phenolics and their antimicrobial properties. Phenol derivatives originate when a functional group (e.g., alkyl, phenyl, benzyl, halogen) replaces one of the hydrogen atoms on the aromatic ring. Two phenol derivatives commonly found as constituents of hospital disinfectants are *ortho*-phenylphenol and *ortho*-benzyl-*para*-chlorophenol. The antimicrobial properties of these compounds and many other phenol derivatives are much improved over those of the parent chemical. Phenolics are absorbed by porous materials, and the residual disinfectant can irritate tissue. In 1970, depigmentation of the skin was reported to be caused by phenolic germicidal detergents containing *para*-tertiary butylphenol and *para*-tertiary amylphenol⁷³¹.

Mode of Action. In high concentrations, phenol acts as a gross protoplasmic poison, penetrating and disrupting the cell wall and precipitating the cell proteins. Low concentrations of phenol and higher molecular-weight phenol derivatives cause bacterial death by inactivation of essential enzyme systems and leakage of essential metabolites from the cell wall ⁷³².

Microbicidal Activity. Published reports on the antimicrobial efficacy of commonly used phenolics showed they were bactericidal, fungicidal, virucidal, and tuberculocidal ^{14, 61, 71, 73, 227, 416, 573, 732-738}. One study demonstrated little or no virucidal effect of a phenolic against coxsackie B4, echovirus 11, and poliovirus 1 ⁷³⁶. Similarly, 12% *ortho*-phenylphenol failed to inactivate any of the three hydrophilic viruses after a 10-minute exposure time, although 5% phenol was lethal for these viruses ⁷². A 0.5% dilution of a phenolic (2.8% *ortho*-phenylphenol and 2.7% *ortho*-benzyl-*para*-chlorophenol) inactivated HIV ²²⁷ and a 2% solution of a phenolic (15% ortho-phenylphenol and 6.3% para-tertiary-amylphenol) inactivated all but one of 11 fungi tested ⁷¹.

Manufacturers' data using the standardized AOAC methods demonstrate that commercial phenolics are not sporicidal but are tuberculocidal, fungicidal, virucidal, and bactericidal at their recommended usedilution. Attempts to substantiate the bactericidal label claims of phenolics using the AOAC Use-Dilution Method occasionally have failed ^{416, 737}. However, results from these same studies have varied dramatically among laboratories testing identical products.

Uses. Many phenolic germicides are EPA-registered as disinfectants for use on environmental surfaces (e.g., bedside tables, bedrails, and laboratory surfaces) and noncritical medical devices. Phenolics are not FDA-cleared as high-level disinfectants for use with semicritical items but could be used to preclean

or decontaminate critical and semicritical devices before terminal sterilization or high-level disinfection.

The use of phenolics in nurseries has been questioned because of hyperbilirubinemia in infants placed in bassinets where phenolic detergents were used ⁷³⁹. In addition, bilirubin levels were reported to increase in phenolic-exposed infants, compared with nonphenolic-exposed infants, when the phenolic was prepared according to the manufacturers' recommended dilution ⁷⁴⁰. If phenolics are used to clean nursery floors, they must be diluted as recommended on the product label. Phenolics (and other disinfectants) should not be used to clean infant bassinets and incubators while occupied. If phenolics are used to terminally clean infant bassinets and incubators, the surfaces should be rinsed thoroughly with water and dried before reuse of infant bassinets and incubators ¹⁷.

Quaternary Ammonium Compounds

Overview. The quaternary ammonium compounds are widely used as disinfectants. Health-care– associated infections have been reported from contaminated quaternary ammonium compounds used to disinfect patient-care supplies or equipment, such as cystoscopes or cardiac catheters ^{741, 742}. The quaternaries are good cleaning agents, but high water hardness ⁷⁴³ and materials such as cotton and gauze pads can make them less microbicidal because of insoluble precipitates or cotton and gauze pads absorb the active ingredients, respectively. One study showed a significant decline (~40%–50% lower at 1 hour) in the concentration of quaternaries released when cotton rags or cellulose-based wipers were used in the open-bucket system, compared with the nonwoven spunlace wipers in the closed-bucket system ⁷⁴⁴ As with several other disinfectants (e.g., phenolics, iodophors) gram-negative bacteria can survive or grow in them ⁴⁰⁴.

Chemically, the quaternaries are organically substituted ammonium compounds in which the nitrogen atom has a valence of 5, four of the substituent radicals (R1-R4) are alkyl or heterocyclic radicals of a given size or chain length, and the fifth (X⁻) is a halide, sulfate, or similar radical ⁷⁴⁵. Each compound exhibits its own antimicrobial characteristics, hence the search for one compound with outstanding antimicrobial properties. Some of the chemical names of quaternary ammonium compounds used in healthcare are alkyl dimethyl benzyl ammonium chloride, alkyl didecyl dimethyl ammonium chloride, and dialkyl dimethyl ammonium chloride. The newer quaternary ammonium compounds (i.e., fourth generation), referred to as twin-chain or dialkyl quaternaries (e.g. didecyl dimethyl ammonium bromide and dioctyl dimethyl ammonium bromide), purportedly remain active in hard water and are tolerant of anionic residues ⁷⁴⁶.

A few case reports have documented occupational asthma as a result of exposure to benzalkonium chloride ⁷⁴⁷.

Mode of Action. The bactericidal action of the quaternaries has been attributed to the inactivation of energy-producing enzymes, denaturation of essential cell proteins, and disruption of the cell membrane⁷⁴⁶. Evidence exists that supports these and other possibilities ^{745 748}.

Microbicidal Activity. Results from manufacturers' data sheets and from published scientific literature indicate that the quaternaries sold as hospital disinfectants are generally fungicidal, bactericidal, and virucidal against lipophilic (enveloped) viruses; they are not sporicidal and generally not tuberculocidal or virucidal against hydrophilic (nonenveloped) viruses^{14, 54-56, 58, 59, 61, 71, 73, 186, 297, 748, 749}. The poor mycobactericidal activities of quaternary ammonium compounds have been demonstrated ^{55, 73}. Quaternary ammonium compounds (as well as 70% isopropyl alcohol, phenolic, and a chlorine-containing wipe [80 ppm]) effectively (>95%) remove and/or inactivate contaminants (i.e., multidrug-resistant *S. aureus*, vancomycin-resistant *Entercoccus*, *P. aeruginosa*) from computer keyboards with a 5-second application time. No functional damage or cosmetic changes occurred to the computer keyboards after 300 applications of the disinfectants ⁴⁵.

Attempts to reproduce the manufacturers' bactericidal and tuberculocidal claims using the AOAC tests

with a limited number of quaternary ammonium compounds occasionally have failed ^{73, 416, 737}. However, test results have varied extensively among laboratories testing identical products ^{416, 737}.

Uses. The quaternaries commonly are used in ordinary environmental sanitation of noncritical surfaces, such as floors, furniture, and walls. EPA-registered quaternary ammonium compounds are appropriate to use for disinfecting medical equipment that contacts intact skin (e.g., blood pressure cuffs).

Miscellaneous Inactivating Agents

Other Germicides

Several compounds have antimicrobial activity but for various reasons have not been incorporated into the armamentarium of health-care disinfectants. These include mercurials, sodium hydroxide, β -propiolactone, chlorhexidine gluconate, cetrimide-chlorhexidine, glycols (triethylene and propylene), and the Tego disinfectants. Two authoritative references examine these agents in detail ^{16, 412}.

A peroxygen-containing formulation had marked bactericidal action when used as a 1% weight/volume solution and virucidal activity at 3% ⁴⁹, but did not have mycobactericidal activity at concentrations of 2.3% and 4% and exposure times ranging from 30 to 120 minutes ⁷⁵⁰. It also required 20 hours to kill *B. atrophaeus* spores ⁷⁵¹. A powder-based peroxygen compound for disinfecting contaminated spill was strongly and rapidly bactericidal ⁷⁵².

In preliminary studies, nanoemulsions (composed of detergents and lipids in water) showed activity against vegetative bacteria, enveloped viruses and *Candida*. This product represents a potential agent for use as a topical biocidal agent. ⁷⁵³⁻⁷⁵⁵.

New disinfectants that require further evaluation include glucoprotamin⁷⁵⁶, tertiary amines ⁷⁰³. and a light-activated antimicrobial coating ⁷⁵⁷. Several other disinfection technologies might have potential applications in the healthcare setting ⁷⁵⁸.

Metals as Microbicides

Comprehensive reviews of antisepsis ⁷⁵⁹, disinfection⁴²¹, and anti-infective chemotherapy ⁷⁶⁰ barely mention the antimicrobial activity of heavy metals^{761, 762}. Nevertheless, the anti-infective activity of some heavy metals has been known since antiquity. Heavy metals such as silver have been used for prophylaxis of conjunctivitis of the newborn, topical therapy for burn wounds, and bonding to indwelling catheters, and the use of heavy metals as antiseptics or disinfectants is again being explored ⁷⁶³. Inactivation of bacteria on stainless steel surfaces by zeolite ceramic coatings containing silver and zinc ions has also been demonstrated ^{764, 765}.

Metals such as silver, iron, and copper could be used for environmental control, disinfection of water, or reusable medical devices or incorporated into medical devices (e.g., intravascular catheters)^{400, 761-763, 766-770}. A comparative evaluation of six disinfectant formulations for residual antimicrobial activity demonstrated that only the silver disinfectant demonstrated significant residual activity against *S. aureus* and *P. aeruginosa* ⁷⁶³. Preliminary data suggest metals are effective against a wide variety of microorganisms.

Clinical uses of other heavy metals include copper-8-quinolinolate as a fungicide against *Aspergillus*, copper-silver ionization for *Legionella* disinfection ⁷⁷¹⁻⁷⁷⁴, organic mercurials as an antiseptic (e.g., mercurochrome) and preservative/disinfectant (e.g., thimerosal [currently being removed from vaccines]) in pharmaceuticals and cosmetics ⁷⁶².

Ultraviolet Radiation (UV)

The wavelength of UV radiation ranges from 328 nm to 210 nm (3280 A to 2100 A). Its maximum bactericidal effect occurs at 240–280 nm. Mercury vapor lamps emit more than 90% of their radiation at 253.7 nm, which is near the maximum microbicidal activity ⁷⁷⁵. Inactivation of microorganisms results from destruction of nucleic acid through induction of thymine dimers. UV radiation has been employed in the disinfection of drinking water ⁷⁷⁶, air ⁷⁷⁵, titanium implants ⁷⁷⁷, and contact lenses⁷⁷⁸. Bacteria and viruses are more easily killed by UV light than are bacterial spores ⁷⁷⁵. UV radiation has several potential applications, but unfortunately its germicidal effectiveness and use is influenced by organic matter; wavelength; type of suspension; temperature; type of microorganism; and UV intensity, which is affected

by distance and dirty tubes⁷⁷⁹. The application of UV radiation in the health-care environment (i.e., operating rooms, isolation rooms, and biologic safety cabinets) is limited to destruction of airborne organisms or inactivation of microorganisms on surfaces. The effect of UV radiation on postoperative wound infections was investigated in a double-blind, randomized study in five university medical centers. After following 14,854 patients over a 2-year period, the investigators reported the overall wound infection rate was unaffected by UV radiation, although postoperative infection in the "refined clean" surgical procedures decreased significantly (3.8%–2.9%) ⁷⁸⁰. No data support the use of UV lamps in isolation rooms, and this practice has caused at least one epidemic of UV-induced skin erythema and keratoconjunctivitis in hospital patients and visitors ⁷⁸¹.

Pasteurization

Pasteurization is not a sterilization process; its purpose is to destroy all pathogenic microorganisms. However, pasteurization does not destroy bacterial spores. The time-temperature relation for hot-water pasteurization is generally ~70°C (158°F) for 30 minutes. The water temperature and time should be monitored as part of a quality-assurance program ⁷⁸². Pasteurization of respiratory therapy ^{783, 784} and anesthesia equipment ⁷⁸⁵ is a recognized alternative to chemical disinfection. The efficacy of this process has been tested using an inoculum that the authors believed might simulate contamination by an infected patient. Use of a large inoculum (10⁷) of *P. aeruginosa* or *Acinetobacter calcoaceticus* in sets of respiratory tubing before processing demonstrated that machine-assisted chemical processing was more efficient than machine-assisted pasteurization with a disinfection failure rate of 6% and 83%, respectively ⁷⁸³. Other investigators found hot water disinfection to be effective (inactivation factor >5 log₁₀) against multiple bacteria, including multidrug-resistant bacteria, for disinfecting reusable anesthesia or respiratory therapy equipment ⁷⁸⁴⁻⁷⁸⁶.

Flushing- and Washer-Disinfectors

Flushing- and washer-disinfectors are automated and closed equipment that clean and disinfect objects from bedpans and washbowls to surgical instruments and anesthesia tubes. Items such as bedpans and urinals can be cleaned and disinfected in flushing-disinfectors. They have a short cycle of a few minutes. They clean by flushing with warm water, possibly with a detergent, and then disinfect by flushing the items with hot water or with steam. Because this machine empties, cleans, and disinfects, manual cleaning is eliminated, fewer disposable items are needed, and fewer chemical germicides are used. A microbiologic evaluation of one washer/disinfector demonstrated complete inactivation of suspensions of *E. faecalis* or poliovirus ⁷⁸⁷. Other studies have shown that strains of *Enterococcus faecium* can survive the British Standard for heat disinfection of bedpans (80°C for 1 minute). The significance of this finding with reference to the potential for enterococci to survive and disseminate in the health-care environment is debatable ⁷⁸⁸⁻⁷⁹⁰. These machines are available and used in many European countries.

Surgical instruments and anesthesia equipment are more difficult to clean. They are run in washerdisinfectors on a longer cycle of approximately 20–30 minutes with a detergent. These machines also disinfect by hot water at approximately 90°C⁷⁹¹.

The Regulatory Framework for Disinfectants and Sterilants

Before using the guidance provided in this document, health-care workers should be aware of the federal laws and regulations that govern the sale, distribution, and use of disinfectants and sterilants. In particular, health-care workers need to know what requirements pertain to them when they apply these products. Finally, they should understand the relative roles of EPA, FDA, and CDC so the context for the guidance provided in this document is clear.

EPA and FDA

In the United States, chemical germicides formulated as sanitizers, disinfectants, or sterilants are regulated in interstate commerce by the Antimicrobials Division, Office of Pesticides Program, EPA, under the authority of the Federal Insecticide, Fungicide, and Rodenticide Act (FIFRA) of 1947, as amended ⁷⁹². Under FIFRA, any substance or mixture of substances intended to prevent, destroy, repel, or mitigate any pest (including microorganisms but excluding those in or on living humans or animals) must be registered before sale or distribution. To obtain a registration, a manufacturer must submit specific data about the safety and effectiveness of each product. For example, EPA requires manufacturers of sanitizers, disinfectants, or chemical sterilants to test formulations by using accepted methods for microbiocidal activity, stability, and toxicity to animals and humans. The manufacturers submit these data to EPA along with proposed labeling. If EPA concludes the product can be used without causing "unreasonable adverse effects," then the product and its labeling are registered, and the manufacturer can sell and distribute the product in the United States.

FIFRA also requires users of products to follow explicitly the labeling directions on each product. The following standard statement appears on all labels under the "Directions for Use" heading: "It is a violation of federal law to use this product in a manner inconsistent with its labeling." This statement means a health-care worker must follow the safety precautions and use directions on the labeling of each registered product. Failure to follow the specified use-dilution, contact time, method of application, or any other condition of use is considered a misuse of the product and potentially subject to enforcement action under FIFRA.

In general, EPA regulates disinfectants and sterilants used on environmental surfaces, and not those used on critical or semicritical medical devices; the latter are regulated by FDA. In June 1993, FDA and EPA issued a "Memorandum of Understanding" that divided responsibility for review and surveillance of chemical germicides between the two agencies. Under the agreement, FDA regulates liquid chemical sterilants used on critical and semicritical devices, and EPA regulates disinfectants used on noncritical surfaces and gaseous sterilants ⁷⁹³. In 1996, Congress passed the Food Quality Protection Act (FQPA). This act amended FIFRA in regard to several types of products regulated by both EPA and FDA. One provision of FQPA removed regulation of liquid chemical sterilants used on critical and semicritical medical devices from EPA's jurisdiction, and it now rests solely with FDA ^{792, 794}. EPA continues to register nonmedical chemical sterilants. FDA and EPA have considered the impact of FQPA, and in January 2000, FDA published its final guidance document on product submissions and labeling. Antiseptics are considered antimicrobial drugs used on living tissue and thus are regulated by FDA under the Food, Drug and Cosmetic Act. FDA regulates liquid chemical sterilants and high-level disinfectants intended to process critical and semicritical devices. FDA has published recommendations on the types of test methods that manufacturers should submit to FDA for 510[k] clearance for such agents.

CDC

At CDC, the mission of the Coordinating Center for Infections Diseases is to guide the public on how to prevent and respond to infectious diseases in both health-care settings and at home. With respect to disinfectants and sterilants, part of CDC's role is to inform the public (in this case healthcare personnel) of current scientific evidence pertaining to these products, to comment about their safety and efficacy, and to recommend which chemicals might be most appropriate or effective for specific microorganisms and settings.

Test Methods

The methods EPA has used for registration are standardized by the AOAC International; however, a survey of scientific literature reveals a number of problems with these tests that were reported during 1987–1990 ^{58, 76, 80, 428, 736, 737, 795-800} that cause them to be neither accurate nor reproducible ^{416, 737}. As part of their regulatory authority, EPA and FDA support development and validation of methods for assessing disinfection claims ⁸⁰¹⁻⁸⁰³. For example, EPA has supported the work of Dr. Syed Sattar and coworkers who have developed a two-tier quantitative carrier test to assess sporicidal, mycobactericidal, bactericidal, fungicidal, virucidal, and protozoacidal activity of chemical germicides ^{701, 803}. EPA is accepting label claims against hepatitis B virus (HBV) using a surrogate organism, the duck HBV, to quantify disinfectant activity ^{124, 804}. EPA also is accepting labeling claims against hepatitis C virus using the bovine viral diarrhea virus as a surrogate.

For nearly 30 years, EPA also performed intramural preregistration and postregistration efficacy testing of some chemical disinfectants in its own laboratories. In 1982, this was stopped, reportedly for budgetary reasons. At that time, manufacturers did not need to have microbiologic activity claims verified by EPA or an independent testing laboratory when registering a disinfectant or chemical sterilant ⁸⁰⁵. This occurred when the frequency of contaminated germicides and infections secondary to their use had increased ⁴⁰⁴. Investigations demonstrating that interlaboratory reproducibility of test results was poor and manufacturers' label claims were not verifiable 416, 737 and symposia sponsored by the American Society for Microbiology⁸⁰⁰ heightened awareness of these problems and reconfirmed the need to improve the AOAC methods and reinstate a microbiologic activity verification program. A General Accounting Office report entitled Disinfectants: EPA Lacks Assurance They Work⁸⁰⁶ seemed to provide the necessary impetus for EPA to initiate corrective measures, including cooperative agreements to improve the AOAC methods and independent verification testing for all products labeled as sporicidal and disinfectants labeled as tuberculocidal. For example, of 26 sterilant products tested by EPA, 15 were canceled because of product failure. A list of products registered with EPA and labeled for use as sterilants or tuberculocides or against HIV and/or HBV is available through EPA's website at [This link is no longer active: http://www.epa.gov/oppad001/chemregindex.htm. The current version of this document may differ from original version: Selected EPA-registered Disinfectants (https://www.epa.gov/pesticideregistration/selected-epa-registered-disinfectants).]. Organizations (e.g., Organization for Economic Cooperation and Development) are working to standardize requirements for germicide testing and registration.

Neutralization of Germicides

One of the difficulties associated with evaluating the bactericidal activity of disinfectants is prevention of bacteriostasis from disinfectant residues carried over into the subculture media. Likewise, small amounts of disinfectants on environmental surfaces can make an accurate bacterial count difficult to get when sampling of the health-care environment as part of an epidemiologic or research investigation. One way these problems may be overcome is by employing neutralizers that inactivate residual disinfectants ⁸⁰⁷⁻⁸⁰⁹. Two commonly used neutralizing media for chemical disinfectants are Letheen Media and D/E Neutralizing Media. The former contains lecithin to neutralize quaternaries and polysorbate 80 (Tween 80) to neutralize phenolics, hexachlorophene, formalin, and, with lecithin, ethanol. The D/E Neutralizing media will neutralize a broad spectrum of antiseptic and disinfectant chemicals, including quaternary ammonium compounds, phenols, iodine and chlorine compounds, mercurials, formaldehyde, and glutaraldehyde ⁸¹⁰. A review of neutralizers used in germicide testing has been published⁸⁰⁸.

Sterilization

Most medical and surgical devices used in healthcare facilities are made of materials that are heat stable and therefore undergo heat, primarily steam, sterilization. However, since 1950, there has been an increase in medical devices and instruments made of materials (e.g., plastics) that require low-temperature sterilization. Ethylene oxide gas has been used since the 1950s for heat- and moisture-sensitive medical devices. Within the past 15 years, a number of new, low-temperature sterilization systems (e.g., hydrogen peroxide gas plasma, peracetic acid immersion, ozone) have been developed and are being used to sterilize medical devices. This section reviews sterilization technologies used in healthcare and makes recommendations for their optimum performance in the processing of medical devices ^{1, 18, 811-820}.

Sterilization destroys all microorganisms on the surface of an article or in a fluid to prevent disease transmission associated with the use of that item. While the use of inadequately sterilized critical items represents a high risk of transmitting pathogens, documented transmission of pathogens associated with an inadequately sterilized critical item is exceedingly rare ^{821, 822}. This is likely due to the wide margin of safety associated with the sterilization processes used in healthcare facilities. The concept of what constitutes "sterile" is measured as a probability of sterility for each item to be sterilized. This probability is commonly referred to as the sterility assurance level (SAL) of the product and is defined as the probability of a single viable microorganism occurring on a product after sterilization. SAL is normally expressed a 10⁻ⁿ. For example, if the probability of a spore surviving were one in one million, the SAL would be 10^{-6 823, 824}. In short, a SAL is an estimate of lethality of the entire sterilization process and is a conservative calculation. Dual SALs (e.g., 10⁻³ SAL for blood culture tubes, drainage bags; 10⁻⁶ SAL for scalpels, implants) have been used in the United States for many years and the choice of a 10⁻⁶ SAL was strictly arbitrary and not associated with any adverse outcomes (e.g., patient infections) ⁸²³.

Medical devices that have contact with sterile body tissues or fluids are considered critical items. These items should be sterile when used because any microbial contamination could result in disease transmission. Such items include surgical instruments, biopsy forceps, and implanted medical devices. If these items are heat resistant, the recommended sterilization process is steam sterilization, because it has the largest margin of safety due to its reliability, consistency, and lethality. However, reprocessing heat- and moisture-sensitive items requires use of a low-temperature sterilization technology (e.g., ethylene oxide, hydrogen peroxide gas plasma, peracetic acid) ⁸²⁵. A summary of the advantages and disadvantages for commonly used sterilization technologies is presented in Table 6.

Steam Sterilization

Overview. Of all the methods available for sterilization, moist heat in the form of saturated steam under pressure is the most widely used and the most dependable. Steam sterilization is nontoxic, inexpensive ⁸²⁶, rapidly microbicidal, sporicidal, and rapidly heats and penetrates fabrics (Table 6) ⁸²⁷. Like all sterilization processes, steam sterilization has some deleterious effects on some materials, including corrosion and combustion of lubricants associated with dental handpieces²¹²; reduction in ability to transmit light associated with laryngoscopes⁸²⁸; and increased hardening time (5.6 fold) with plaster-cast ⁸²⁹.

The basic principle of steam sterilization, as accomplished in an autoclave, is to expose each item to direct steam contact at the required temperature and pressure for the specified time. Thus, there are four parameters of steam sterilization: steam, pressure, temperature, and time. The ideal steam for sterilization is dry saturated steam and entrained water (dryness fraction ≥97%)^{813, 819}. Pressure serves as a means to obtain the high temperatures necessary to quickly kill microorganisms. Specific temperatures are 121°C (250°F) and 132°C (270°F). These temperatures (and other high temperatures) ⁸³⁰ must be maintained for a minimal time to kill microorganisms. Recognized minimum exposure periods for sterilization of wrapped healthcare supplies are 30 minutes at 121°C (250°F) in a gravity displacement

sterilizer or 4 minutes at 132°C (270°C) in a prevacuum sterilizer (Table 7). At constant temperatures, sterilization times vary depending on the type of item (e.g., metal versus rubber, plastic, items with lumens), whether the item is wrapped or unwrapped, and the sterilizer type.

The two basic types of steam sterilizers (autoclaves) are the gravity displacement autoclave and the high-speed prevacuum sterilizer. In the former, steam is admitted at the top or the sides of the sterilizing chamber and, because the steam is lighter than air, forces air out the bottom of the chamber through the drain vent. The gravity displacement autoclaves are primarily used to process laboratory media, water, pharmaceutical products, regulated medical waste, and nonporous articles whose surfaces have direct steam contact. For gravity displacement sterilizers the penetration time into porous items is prolonged because of incomplete air elimination. This point is illustrated with the decontamination of 10 lbs of microbiological waste, which requires at least 45 minutes at 121°C because the entrapped air remaining in a load of waste greatly retards steam permeation and heating efficiency^{831, 832}. The high-speed prevacuum sterilizers are similar to the gravity displacement sterilizers except they are fitted with a vacuum pump (or ejector) to ensure air removal from the sterilizing chamber and load before the steam is admitted. The advantage of using a vacuum pump is that there is nearly instantaneous steam penetration even into porous loads. The Bowie-Dick test is used to detect air leaks and inadequate air removal and consists of folded 100% cotton surgical towels that are clean and preconditioned. A commercially available Bowie-Dick-type test sheet should be placed in the center of the pack. The test pack should be placed horizontally in the front, bottom section of the sterilizer rack, near the door and over the drain, in an otherwise empty chamber and run at 134°C for 3.5 minutes^{813, 819}. The test is used each day the vacuum-type steam sterilizer is used, before the first processed load. Air that is not removed from the chamber will interfere with steam contact. Smaller disposable test packs (or process challenge devices) have been devised to replace the stack of folded surgical towels for testing the efficacy of the vacuum system in a prevacuum sterilizer. 833 These devices are "designed to simulate product to be sterilized and to constitute a defined challenge to the sterilization process"^{819, 834}. They should be representative of the load and simulate the greatest challenge to the load⁸³⁵. Sterilizer vacuum performance is acceptable if the sheet inside the test pack shows a uniform color change. Entrapped air will cause a spot to appear on the test sheet, due to the inability of the steam to reach the chemical indicator. If the sterilizer fails the Bowie-Dick test, do not use the sterilizer until it is inspected by the sterilizer maintenance personnel and passes the Bowie-Dick test^{813, 819, 836}.

Another design in steam sterilization is a steam flush-pressure pulsing process, which removes air rapidly by repeatedly alternating a steam flush and a pressure pulse above atmospheric pressure. Air is rapidly removed from the load as with the prevacuum sterilizer, but air leaks do not affect this process because the steam in the sterilizing chamber is always above atmospheric pressure. Typical sterilization temperatures and times are 132°C to 135°C with 3 to 4 minutes exposure time for porous loads and instruments^{827, 837}.

Like other sterilization systems, the steam cycle is monitored by mechanical, chemical, and biological monitors. Steam sterilizers usually are monitored using a printout (or graphically) by measuring temperature, the time at the temperature, and pressure. Typically, chemical indicators are affixed to the outside and incorporated into the pack to monitor the temperature or time and temperature. The effectiveness of steam sterilization is monitored with a biological indicator containing spores of *Geobacillus stearothermophilus* (formerly *Bacillus stearothermophilus*). Positive spore test results are a relatively rare event ⁸³⁸ and can be attributed to operator error, inadequate steam delivery⁸³⁹, or equipment malfunction.

Portable (table-top) steam sterilizers are used in outpatient, dental, and rural clinics⁸⁴⁰. These sterilizers are designed for small instruments, such as hypodermic syringes and needles and dental instruments. The ability of the sterilizer to reach physical parameters necessary to achieve sterilization should be monitored by mechanical, chemical, and biological indicators.

Microbicidal Activity. The oldest and most recognized agent for inactivation of microorganisms is heat. D-values (time to reduce the surviving population by 90% or 1 log₁₀) allow a direct comparison of the heat resistance of microorganisms. Because a D-value can be determined at various temperatures, a subscript is used to designate the exposure temperature (i.e., D_{121C}). D_{121C}-values for *Geobacillus stearothermophilus* used to monitor the steam sterilization process range from 1 to 2 minutes. Heat-resistant nonspore-forming bacteria, yeasts, and fungi have such low D_{121C} values that they cannot be experimentally measured⁸⁴¹.

Mode of Action. Moist heat destroys microorganisms by the irreversible coagulation and denaturation of enzymes and structural proteins. In support of this fact, it has been found that the presence of moisture significantly affects the coagulation temperature of proteins and the temperature at which microorganisms are destroyed.

Uses. Steam sterilization should be used whenever possible on all critical and semicritical items that are heat and moisture resistant (e.g., steam sterilizable respiratory therapy and anesthesia equipment), even when not essential to prevent pathogen transmission. Steam sterilizers also are used in healthcare facilities to decontaminate microbiological waste and sharps containers ^{831, 832, 842} but additional exposure time is required in the gravity displacement sterilizer for these items.

Flash Sterilization

Overview. "Flash" steam sterilization was originally defined by Underwood and Perkins as sterilization of an unwrapped object at 132°C for 3 minutes at 27-28 lbs. of pressure in a gravity displacement sterilizer⁸⁴³. Currently, the time required for flash sterilization depends on the type of sterilizer and the type of item (i.e., porous vs non-porous items)(see Table 8). Although the wrapped method of sterilization is preferred for the reasons listed below, correctly performed flash sterilization is an effective process for the sterilization of critical medical devices^{844, 845}. Flash sterilization is a modification of conventional steam sterilization (either gravity, prevacuum, or steam-flush pressure-pulse) in which the flashed item is placed in an open tray or is placed in a specially designed, covered, rigid container to allow for rapid penetration of steam. Historically, it is not recommended as a routine sterilization method because of the lack of timely biological indicators to monitor performance, absence of protective packaging following sterilization, possibility for contamination of processed items during transportation to the operating rooms, and the sterilization cycle parameters (i.e., time, temperature, pressure) are minimal. To address some of these concerns, many healthcare facilities have done the following: placed equipment for flash sterilization in close proximity to operating rooms to facilitate aseptic delivery to the point of use (usually the sterile field in an ongoing surgical procedure); extended the exposure time to ensure lethality comparable to sterilized wrapped items (e.g., 4 minutes at 132°C)^{846, 847}; used biological indicators that provide results in 1 hour for flash-sterilized items^{846, 847}; and used protective packaging that permits steam penetration^{812, 817-819, 845, 848}. Further, some rigid, reusable sterilization container systems have been designed and validated by the container manufacturer for use with flash cycles. When sterile items are open to air, they will eventually become contaminated. Thus, the longer a sterile item is exposed to air, the greater the number of microorganisms that will settle on it. Sterilization cycle parameters for flash sterilization are shown in Table 8.

A few adverse events have been associated with flash sterilization. When evaluating an increased incidence of neurosurgical infections, the investigators noted that surgical instruments were flash sterilized between cases and 2 of 3 craniotomy infections involved plate implants that were flash sterilized⁸⁴⁹. A report of two patients who received burns during surgery from instruments that had been flash sterilized reinforced the need to develop policies and educate staff to prevent the use of instruments hot enough to cause clinical burns⁸⁵⁰. Staff should use precautions to prevent burns with potentially hot instruments (e.g., transport tray using heat-protective gloves). Patient burns may be prevented by either air-cooling the instruments or immersion in sterile liquid (e.g., saline).

Uses. Flash sterilization is considered acceptable for processing cleaned patient-care items that cannot be packaged, sterilized, and stored before use. It also is used when there is insufficient time to sterilize an item by the preferred package method. Flash sterilization should not be used for reasons of convenience, as an alternative to purchasing additional instrument sets, or to save time⁸¹⁷. Because of the potential for serious infections, flash sterilization is not recommended for implantable devices (i.e., devices placed into a surgically or naturally formed cavity of the human body); however, flash sterilization may be unavoidable for some devices (e.g., orthopedic screw, plates). If flash sterilization of an implantable device is unavoidable, recordkeeping (i.e., load identification, patient's name/hospital identifier, and biological indicator result) is essential for epidemiological tracking (e.g., of surgical site infection, tracing results of biological indicators to patients who received the item to document sterility), and for an assessment of the reliability of the sterilization process (e.g., evaluation of biological monitoring records and sterilization maintenance records noting preventive maintenance and repairs with dates).

Low-Temperature Sterilization Technologies

Ethylene oxide (ETO) has been widely used as a low-temperature sterilant since the 1950s. It has been the most commonly used process for sterilizing temperature- and moisture-sensitive medical devices and supplies in healthcare institutions in the United States. Two types of ETO sterilizers are available, mixed gas and 100% ETO. Until 1995, ethylene oxide sterilizers combined ETO with a chloroflourocarbon (CFC) stabilizing agent, most commonly in a ratio of 12% ETO mixed with 88% CFC (referred to as 12/88 ETO).

For several reasons, healthcare personnel have been exploring the use of new low-temperature sterilization technologies^{825, 851}. First, CFCs were phased out in December 1995 under provisions of the Clean Air Act ⁸⁵². CFCs were classified as a Class I substance under the Clean Air Act because of scientific evidence linking them to destruction of the earth's ozone layer. Second, some states (e.g., California, New York, Michigan) require the use of ETO abatement technology to reduce the amount of ETO being released into ambient air from 90 to 99.9% depending on the state. Third, OSHA regulates the acceptable vapor levels of ETO (i.e., 1 ppm averaged over 8 hours) due to concerns that ETO exposure represents an occupational hazard³¹⁸. These constraints have led to the development of alternative technologies for low-temperature sterilization in the healthcare setting.

Alternative technologies to ETO with chlorofluorocarbon that are currently available and cleared by the FDA for medical equipment include 100% ETO; ETO with a different stabilizing gas, such as carbon dioxide or hydrochlorofluorocarbons (HCFC); immersion in peracetic acid; hydrogen peroxide gas plasma; and ozone. Technologies under development for use in healthcare facilities, but not cleared by the FDA, include vaporized hydrogen peroxide, vapor phase peracetic acid, gaseous chlorine dioxide, ionizing radiation, or pulsed light ^{400, 758, 853}. However, there is no guarantee that these new sterilization technologies will receive FDA clearance for use in healthcare facilities.

These new technologies should be compared against the characteristics of an ideal low-temperature (<60°C) sterilant (Table 9). ⁸⁵¹ While it is apparent that all technologies will have limitations (Table 9), understanding the limitations imposed by restrictive device designs (e.g., long, narrow lumens) is critical for proper application of new sterilization technology⁸⁵⁴. For example, the development of increasingly small and complex endoscopes presents a difficult challenge for current sterilization processes. This occurs because microorganisms must be in direct contact with the sterilant for inactivation to occur. Several peer-reviewed scientific publications have data demonstrating concerns about the efficacy of several of the low-temperature sterilization processes (i.e., gas plasma, vaporized hydrogen peroxide, ETO, peracetic acid), particularly when the test organisms are challenged in the presence of serum and salt and a narrow lumen vehicle^{469, 721, 825, 855, 856}. Factors shown to affect the efficacy of sterilization are shown in Table 10.

Ethylene Oxide "Gas" Sterilization

Overview. ETO is a colorless gas that is flammable and explosive. The four essential parameters (operational ranges) are: gas concentration (450 to 1200 mg/l); temperature (37 to 63°C); relative humidity (40 to 80%)(water molecules carry ETO to reactive sites); and exposure time (1 to 6 hours). These influence the effectiveness of ETO sterilization^{814, 857, 858}. Within certain limitations, an increase in gas concentration and temperature may shorten the time necessary for achieving sterilization.

The main disadvantages associated with ETO are the lengthy cycle time, the cost, and its potential hazards to patients and staff; the main advantage is that it can sterilize heat- or moisture-sensitive medical equipment without deleterious effects on the material used in the medical devices (Table 6). Acute exposure to ETO may result in irritation (e.g., to skin, eyes, gastrointestinal or respiratory tracts) and central nervous system depression⁸⁵⁹⁻⁸⁶². Chronic inhalation has been linked to the formation of cataracts, cognitive impairment, neurologic dysfunction, and disabling polyneuropathies^{860, 861, 863-866}. Occupational exposure in healthcare facilities has been linked to hematologic changes ⁸⁶⁷ and an increased risk of spontaneous abortions and various cancers^{318, 868-870}. ETO should be considered a known human carcinogen⁸⁷¹.

The basic ETO sterilization cycle consists of five stages (i.e., preconditioning and humidification, gas introduction, exposure, evacuation, and air washes) and takes approximately 2 1/2 hrs excluding aeration time. Mechanical aeration for 8 to 12 hours at 50 to 60°C allows desorption of the toxic ETO residual contained in exposed absorbent materials. Most modern ETO sterilizers combine sterilization and aeration in the same chamber as a continuous process. These ETO models minimize potential ETO exposure during door opening and load transfer to the aerator. Ambient room aeration also will achieve desorption of the toxic ETO but requires 7 days at 20°C. There are no federal regulations for ETO sterilizer emission; however, many states have promulgated emission-control regulations⁸¹⁴.

The use of ETO evolved when few alternatives existed for sterilizing heat- and moisture-sensitive medical devices; however, favorable properties (Table 6) account for its continued widespread use⁸⁷². Two ETO gas mixtures are available to replace ETO-chlorofluorocarbon (CFC) mixtures for large capacity, tank-supplied sterilizers. The ETO-carbon dioxide (CO₂) mixture consists of 8.5% ETO and 91.5% CO₂. This mixture is less expensive than ETO-hydrochlorofluorocarbons (HCFC), but a disadvantage is the need for pressure vessels rated for steam sterilization, because higher pressures (28-psi gauge) are required. The other mixture, which is a drop-in CFC replacement, is ETO mixed with HCFC. HCFCs are approximately 50-fold less damaging to the earth's ozone layer than are CFCs. The EPA will begin regulation of HCFC in the year 2015 and will terminate production in the year 2030. Two companies provide ETO-HCFC mixtures as drop-in replacement for CFC-12; one mixture consists of 8.6% ETO and 91.4% HCFC, and the other mixture is composed of 10% ETO and 90% HCFC⁸⁷². An alternative to the pressurized mixed gas ETO systems is 100% ETO. The 100% ETO sterilizers using unit-dose cartridges eliminate the need for external tanks.

ETO is absorbed by many materials. For this reason, following sterilization the item must undergo aeration to remove residual ETO. Guidelines have been promulgated regarding allowable ETO limits for devices that depend on how the device is used, how often, and how long in order to pose a minimal risk to patients in normal product use⁸¹⁴.

ETO toxicity has been established in a variety of animals. Exposure to ETO can cause eye pain, sore throat, difficulty breathing and blurred vision. Exposure can also cause dizziness, nausea, headache, convulsions, blisters and vomiting and coughing⁸⁷³. In a variety of *in vitro* and animal studies, ETO has been demonstrated to be carcinogenic. ETO has been linked to spontaneous abortion, genetic damage, nerve damage, peripheral paralysis, muscle weakness, and impaired thinking and memory⁸⁷³. Occupational exposure in healthcare facilities has been linked to an increased risk of spontaneous abortions and various cancers³¹⁸. Injuries (e.g., tissue burns) to patients have been associated with ETO residues in implants used in surgical procedures⁸⁷⁴. Residual ETO in capillary flow dialysis membranes has been shown to be neurotoxic in vitro⁸⁷⁵. OSHA has established a PEL of 1 ppm airborne ETO in the workplace, expressed as a TWA for an 8-hour work shift in a 40-hour work week. The "action level" for

ETO is 0.5 ppm, expressed as an 8-hour TWA, and the short-term excursion limit is 5 ppm, expressed as a 15-minute TWA⁸¹⁴. For details of the requirements in OSHA's ETO standard for occupational exposures, see Title 29 of the Code of Federal Regulations (CFR) Part 1910.1047⁸⁷³. Several personnel monitoring methods (e.g., charcoal tubes and passive sampling devices) are in use⁸¹⁴. OSHA has established a PEL of 5 ppm for ethylene chlorohydrin (a toxic by-product of ETO) in the workplace⁸⁷⁶. Additional information regarding use of ETO in health care facilities is available from NIOSH.

Mode of Action. The microbicidal activity of ETO is considered to be the result of alkylation of protein, DNA, and RNA. Alkylation, or the replacement of a hydrogen atom with an alkyl group, within cells prevents normal cellular metabolism and replication⁸⁷⁷.

Microbicidal Activity. The excellent microbicidal activity of ETO has been demonstrated in several studies ^{469, 721, 722, 856, 878, 879} and summarized in published reports⁸⁷⁷. ETO inactivates all microorganisms although bacterial spores (especially *B. atrophaeus*) are more resistant than other microorganisms. For this reason *B. atrophaeus* is the recommended biological indicator.

Like all sterilization processes, the effectiveness of ETO sterilization can be altered by lumen length, lumen diameter, inorganic salts, and organic materials^{469, 721, 722, 855, 856, 879}. For example, although ETO is not used commonly for reprocessing endoscopes²⁸, several studies have shown failure of ETO in inactivating contaminating spores in endoscope channels ⁸⁵⁵or lumen test units ^{469, 721, 879} and residual ETO levels averaging 66.2 ppm even after the standard degassing time⁴⁵⁶. Failure of ETO also has been observed when dental handpieces were contaminated with *Streptococcus mutans* and exposed to ETO⁸⁸⁰. It is recommended that dental handpieces be steam sterilized.

Uses. ETO is used in healthcare facilities to sterilize critical items (and sometimes semicritical items) that are moisture or heat sensitive and cannot be sterilized by steam sterilization.

Hydrogen Peroxide Gas Plasma

Overview. New sterilization technology based on plasma was patented in 1987 and marketed in the United States in 1993. Gas plasmas have been referred to as the fourth state of matter (i.e., liquids, solids, gases, and gas plasmas). Gas plasmas are generated in an enclosed chamber under deep vacuum using radio frequency or microwave energy to excite the gas molecules and produce charged particles, many of which are in the form of free radicals. A free radical is an atom with an unpaired electron and is a highly reactive species. The proposed mechanism of action of this device is the production of free radicals within a plasma field that are capable of interacting with essential cell components (e.g., enzymes, nucleic acids) and thereby disrupt the metabolism of microorganisms. The type of seed gas used and the depth of the vacuum are two important variables that can determine the effectiveness of this process.

In the late 1980s the first hydrogen peroxide gas plasma system for sterilization of medical and surgical devices was field-tested. According to the manufacturer, the sterilization chamber is evacuated and hydrogen peroxide solution is injected from a cassette and is vaporized in the sterilization chamber to a concentration of 6 mg/l. The hydrogen peroxide vapor diffuses through the chamber (50 minutes), exposes all surfaces of the load to the sterilant, and initiates the inactivation of microorganisms. An electrical field created by a radio frequency is applied to the chamber to create a gas plasma. Microbicidal free radicals (e.g., hydroxyl and hydroperoxyl) are generated in the plasma. The excess gas is removed and in the final stage (i.e., vent) of the process the sterilization chamber is returned to atmospheric pressure by introduction of high-efficiency filtered air. The by-products of the cycle (e.g., water vapor, oxygen) are nontoxic and eliminate the need for aeration. Thus, the sterilized materials can be handled safely, either for immediate use or storage. The process operates in the range of 37-44°C and has a cycle time of 75 minutes. If any moisture is present on the objects the vacuum will not be achieved and the cycle aborts^{856, 881-883}.

A newer version of the unit improves sterilizer efficacy by using two cycles with a hydrogen peroxide

diffusion stage and a plasma stage per sterilization cycle. This revision, which is achieved by a software modification, reduces total processing time from 73 to 52 minutes. The manufacturer believes that the enhanced activity obtained with this system is due in part to the pressure changes that occur during the injection and diffusion phases of the process and to the fact that the process consists of two equal and consecutive half cycles, each with a separate injection of hydrogen peroxide. ^{856, 884, 885} This system and a smaller version ^{400, 882} have received FDA 510[k] clearance with limited application for sterilization of medical devices (Table 6). The biological indicator used with this system is *Bacillus atrophaeus* spores⁸⁵¹. The newest version of the unit, which employs a new vaporization system that removes most of the water from the hydrogen peroxide, has a cycle time from 28-38 minutes (see manufacturer's literature for device dimension restrictions).

Penetration of hydrogen peroxide vapor into long or narrow lumens has been addressed outside the United States by the use of a diffusion enhancer. This is a small, breakable glass ampoule of concentrated hydrogen peroxide (50%) with an elastic connector that is inserted into the device lumen and crushed immediately before sterilization^{470, 885}. The diffusion enhancer has been shown to sterilize bronchoscopes contaminated with *Mycobacteria tuberculosis*⁸⁸⁶. At the present time, the diffusion enhancer is not FDA cleared.

Another gas plasma system, which differs from the above in several important ways, including the use of peracetic acid-acetic acid-hydrogen peroxide vapor, was removed from the marketplace because of reports of corneal destruction to patients when ophthalmic surgery instruments had been processed in the sterilizer^{887, 888}. In this investigation, exposure of potentially wet ophthalmologic surgical instruments with small bores and brass components to the plasma gas led to degradation of the brass to copper and zinc^{888, 889}. The experimenters showed that when rabbit eyes were exposed to the rinsates of the gas plasma-sterilized instruments, corneal decompensation was documented. This toxicity is highly unlikely with the hydrogen peroxide gas plasma process since a toxic, soluble form of copper would not form (LA Feldman, written communication, April 1998).

Mode of Action. This process inactivates microorganisms primarily by the combined use of hydrogen peroxide gas and the generation of free radicals (hydroxyl and hydroproxyl free radicals) during the plasma phase of the cycle.

Microbicidal Activity. This process has the ability to inactivate a broad range of microorganisms, including resistant bacterial spores. Studies have been conducted against vegetative bacteria (including mycobacteria), yeasts, fungi, viruses, and bacterial spores^{469, 721, 856, 881-883, 890-893}. Like all sterilization processes, the effectiveness can be altered by lumen length, lumen diameter, inorganic salts, and organic materials^{469, 721, 855, 856, 890, 891, 893}.

Uses. Materials and devices that cannot tolerate high temperatures and humidity, such as some plastics, electrical devices, and corrosion-susceptible metal alloys, can be sterilized by hydrogen peroxide gas plasma. This method has been compatible with most (>95%) medical devices and materials tested^{884, 895}.

Peracetic Acid Sterilization

Overview. Peracetic acid is a highly biocidal oxidizer that maintains its efficacy in the presence of organic soil. Peracetic acid removes surface contaminants (primarily protein) on endoscopic tubing^{711, 717}. An automated machine using peracetic acid to sterilize medical, surgical, and dental instruments chemically (e.g., endoscopes, arthroscopes) was introduced in 1988. This microprocessor-controlled, low-temperature sterilization method is commonly used in the United States¹⁰⁷. The sterilant, 35% peracetic acid, and an anticorrosive agent are supplied in a single-dose container. The container is punctured at the time of use, immediately prior to closing the lid and initiating the cycle. The concentrated peracetic acid is diluted to 0.2% with filtered water (0.2 μm) at a temperature of approximately 50°C. The diluted peracetic

acid is circulated within the chamber of the machine and pumped through the channels of the endoscope for 12 minutes, decontaminating exterior surfaces, lumens, and accessories. Interchangeable trays are available to permit the processing of up to three rigid endoscopes or one flexible endoscope. Connectors are available for most types of flexible endoscopes for the irrigation of all channels by directed flow. Rigid endoscopes are placed within a lidded container, and the sterilant fills the lumens either by immersion in the circulating sterilant or by use of channel connectors to direct flow into the lumen(s) (see below for the importance of channel connectors). The peracetic acid is discarded via the sewer and the instrument rinsed four times with filtered water. Concern has been raised that filtered water may be inadequate to maintain sterility⁸⁹⁶. Limited data have shown that low-level bacterial contamination may follow the use of filtered water in an AER but no data has been published on AERs using the peracetic acid system¹⁶¹. Clean filtered air is passed through the chamber of the machine and endoscope channels to remove excess water⁷¹⁹. As with any sterilization process, the system can only sterilize surfaces that can be contacted by the sterilant. For example, bronchoscopy-related infections occurred when bronchoscopes were processed using the wrong connector^{155, 725}. Investigation of these incidents revealed that bronchoscopes were inadequately reprocessed when inappropriate channel connectors were used and when there were inconsistencies between the reprocessing instructions provided by the manufacturer of the bronchoscope and the manufacturer of the automatic endoscope reprocessor¹⁵⁵. The importance of channel connectors to achieve sterilization was also shown for rigid lumen devices^{137, 856}.

The manufacturers suggest the use of biological monitors (*G. stearothermophilus* spore strips) both at the time of installation and routinely to ensure effectiveness of the process. The manufacturer's clip must be used to hold the strip in the designated spot in the machine as a broader clamp will not allow the sterilant to reach the spores trapped under it⁸⁹⁷. One investigator reported a 3% failure rate when the appropriate clips were used to hold the spore strip within the machine⁷¹⁸. The use of biological monitors designed to monitor either steam sterilization or ETO for a liquid chemical sterilizer has been questioned for several reasons including spore wash-off from the filter paper strips which may cause less valid monitoring⁸⁹⁸⁻⁹⁰¹. The processor is equipped with a conductivity probe that will automatically abort the cycle if the buffer system is not detected in a fresh container of the peracetic acid solution. A chemical monitoring strip that detects that the active ingredient is >1500 ppm is available for routine use as an additional process control.

Mode of Action. Only limited information is available regarding the mechanism of action of peracetic acid, but it is thought to function as other oxidizing agents, i.e., it denatures proteins, disrupts cell wall permeability, and oxidizes sulfhydral and sulfur bonds in proteins, enzymes, and other metabolites^{654, 726}.

Microbicidal Activity. Peracetic acid will inactivate gram-positive and gram-negative bacteria, fungi, and yeasts in <5 minutes at <100 ppm. In the presence of organic matter, 200-500 ppm is required. For viruses, the dosage range is wide (12-2250 ppm), with poliovirus inactivated in yeast extract in 15 minutes with 1500 to 2250 ppm. Bacterial spores in suspension are inactivated in 15 seconds to 30 minutes with 500 to 10,000 ppm (0.05 to 1%)⁶⁵⁴.

Simulated-use trials have demonstrated microbicidal activity ^{111, 718-722} and three clinical trials have demonstrated both microbial killing and no clinical failures leading to infection^{90, 723, 724}. Alfa and co-workers, who compared the peracetic acid system with ETO, demonstrated the high efficacy of the system. Only the peracetic acid system was able to completely kill 6-log₁₀ of *Mycobacterium chelonae*, *Enterococcus faecalis*, and *B. atrophaeus* spores with both an organic and inorganic challenge⁷²². Like other sterilization processes, the efficacy of the process can be diminished by soil challenges ⁹⁰² and test conditions⁸⁵⁶.

Uses. This automated machine is used to chemically sterilize medical (e.g., GI endoscopes) and surgical (e.g., flexible endoscopes) instruments in the United States. Lumened endoscopes must be connected to an appropriate channel connector to ensure that the sterilant has direct contact with the contaminated lumen. ^{137, 856, 903} Olympus America has not listed this system as a compatible product for

use in reprocessing Olympus bronchoscopes and gastrointestinal endoscopes (Olympus America, January 30, 2002, written communication).

Microbicidal Activity of Low-Temperature Sterilization Technologies

Sterilization processes used in the United States must be cleared by FDA, and they require that sterilizer microbicidal performance be tested under simulated-use conditions⁹⁰⁴. FDA requires that the test article be inoculated with 10⁶ colony-forming units of the most resistant test organism and prepared with organic and inorganic test loads as would occur after actual use. FDA requires manufacturers to use organic soil (e.g., 5% fetal calf serum), dried onto the device with the inoculum, to represent soil remaining on the device following marginal cleaning. However, 5% fetal calf serum as a measure of marginal cleaning has not been validated by measurements of protein load on devices following use and the level of protein removal by various cleaning methods. The inocula must be placed in various locations of the test articles, including those least favorable to penetration and contact with the sterilant (e.a.. lumens). Cleaning before sterilization is not allowed in the demonstration of sterilization efficacy⁹⁰⁴. Several studies have evaluated the relative microbicidal efficacy of these low-temperature sterilization technologies (Table 11). These studies have either tested the activity of a sterilization process against specific microorganisms^{892, 905, 906}, evaluated the microbicidal activity of a singular technology ^{711, 719, 724}, 855, 879, 882-884, 890, 891, 907 or evaluated the comparative effectiveness of several sterilization technologies^{271,} 426, 469, 721, 722, 856, 908, 909. Several test methodologies use stainless steel or porcelain carriers that are inoculated with a test organism. Commonly used test organisms include vegetative bacteria. mycobacteria, and spores of Bacillus species. The available data demonstrate that low-temperature sterilization technologies are able to provide a 6-log₁₀ reduction of microbes when inoculated onto carriers in the absence of salt and serum. However, tests can be constructed such that all of the available sterilization technologies are unable to reliably achieve complete inactivation of a microbial load. 425, 426, 469, 721, 856, 909 For example, almost all of the sterilization processes will fail to reliably inactivate the microbial load in the presence of salt and serum^{469, 721, 909}.

The effect of salts and serums on the sterilization process were studied initially in the 1950s and 1960s^{424, 910}. These studies showed that a high concentration of crystalline-type materials and a low protein content provided greater protection to spores than did serum with a high protein content⁴²⁶. A study by Doyle and Ernst demonstrated resistance of spores by crystalline material applied not only to low-temperature sterilization technology but also to steam and dry heat⁴²⁵. These studies showed that occlusion of *Bacillus atrophaeus* spores in calcium carbonate crystals dramatically increased the time required for inactivation as follows: 10 seconds to 150 minutes for steam (121°C), 3.5 hours to 50 hours for dry heat (121°C), 30 seconds to >2 weeks for ETO (54°C). Investigators have corroborated and extended these findings^{469, 470, 721, 855, 908, 909}. While soils containing both organic and inorganic materials impair microbial killing, soils that contain a high inorganic salt-to-protein ratio favor crystal formation and impair sterilization by occlusion of organisms^{425, 426, 881}.

Alfa and colleagues demonstrated a 6-log₁₀ reduction of the microbial inoculum of porcelain penicylinders using a variety of vegetative and spore-forming organisms (Table 11)⁴⁶⁹. However, if the bacterial inoculum was in tissue-culture medium supplemented with 10% serum, only the ETO 12/88 and ETO-HCFC sterilization mixtures could sterilize 95% to 97% of the penicylinder carriers. The plasma and 100% ETO sterilizer demonstrated significantly reduced activity (Table 11). For all sterilizers evaluated using penicylinder carriers (i.e., ETO 12/88, 100% ETO, hydrogen peroxide gas plasma), there was a 3-to 6-log₁₀ reduction of inoculated bacteria even in the presence of serum and salt. For each sterilizer evaluated, the ability to inactivate microorganisms in the presence of salt and serum was reduced even further when the inoculum was placed in a narrow-lumen test object (3 mm diameter by 125 cm long). Although there was a 2- to 4-log₁₀ reduction in microbial kill, less than 50% of the lumen test objects were sterile when processed using any of the sterilization methods evaluated except the peracetic acid immersion system (Table 11)⁷²¹. Complete killing (or removal) of 6-log₁₀ of *Enterococcus faecalis, Mycobacterium chelonei*, and *Bacillus atrophaeus* spores in the presence of salt and serum and lumen test objects was observed only for the peracetic acid immersion system.

With respect to the results by Alfa and coworkers⁴⁶⁹, Jacobs showed that the use of the tissue culture media created a technique-induced sterilization failure⁴²⁶. Jacobs et al. showed that microorganisms mixed with tissue culture media, used as a surrogate body fluid, formed physical crystals that protected the microorganisms used as a challenge. If the carriers were exposed for 60 sec to nonflowing water, the salts dissolved and the protective effect disappeared. Since any device would be exposed to water for a short period of time during the washing procedure, these protective effects would have little clinical relevance⁴²⁶.

Narrow lumens provide a challenge to some low-temperature sterilization processes. For example, Rutala and colleagues showed that, as lumen size decreased, increased failures occurred with some lowtemperature sterilization technologies. However, some low-temperature processes such as ETO-HCFC and the hydrogen peroxide gas plasma process remained effective even when challenged by a lumen as small as 1 mm in the absence of salt and serum⁸⁵⁶.

The importance of allowing the sterilant to come into contact with the inoculated carrier is demonstrated by comparing the results of two investigators who studied the peracetic acid immersion system. Alfa and coworkers demonstrated excellent activity of the peracetic acid immersion system against three test organisms using a narrow-lumen device. In these experiments, the lumen test object was connected to channel irrigators, which ensured that the sterilant had direct contact with the contaminated carriers⁷²². This effectiveness was achieved through a combination of organism wash-off and peracetic acid sterilant killing the test organisms⁷²². The data reported by Rutala et al. demonstrated failure of the peracetic acid immersion system to eliminate Geobacillus stearothermophilus spores from a carrier placed in a lumen test object. In these experiments, the lumen test unit was not connected to channel irrigators. The authors attributed the failure of the peracetic acid immersion system to eliminate the high levels of spores from the center of the test unit to the inability of the peracetic acid to diffuse into the center of 40-cm long, 3-mm diameter tubes. This may be caused by an air lock or air bubbles formed in the lumen, impeding the flow of the sterilant through the long and narrow lumen and limiting complete access to the Bacillus spores^{137, 856}. Experiments using a channel connector specifically designed for 1-. 2-, and 3-mm lumen test units with the peracetic acid immersion system were completely effective in eliminating an inoculum of 10⁶ Geobacillus stearothermophilus spores⁷. The restricted diffusion environment that exists in the test conditions would not exist with flexible scopes processed in the peracetic acid immersion system, because the scopes are connected to channel irrigators to ensure that the sterilant has direct contact with contaminated surfaces. Alfa and associates attributed the efficacy of the peracetic acid immersion system to the ability of the liquid chemical process to dissolve salts and remove protein and bacteria due to the flushing action of the fluid⁷²².

Bioburden of Surgical Devices

In general, used medical devices are contaminated with a relatively low bioburden of organisms^{179, 911, 912}. Nystrom evaluated medical instruments used in general surgical, gynecological, orthopedic, and earnose-throat operations and found that 62% of the instruments were contaminated with <10¹ organisms after use, 82% with <10², and 91% with <10³. After being washed in an instrument washer, more than 98% of the instruments had <10¹ organisms, and none >10² organisms⁹¹¹. Other investigators have published similar findings^{179, 912}. For example, after a standard cleaning procedure, 72% of 50 surgical instruments contained <10¹ organisms, 86% <10², and only 6% had >3 X 10²⁹¹². In another study of rigid-lumen medical devices, the bioburden on both the inner and outer surface of the lumen ranged from 10¹ to 10⁴ organisms per device. After cleaning, 83% of the devices had a bioburden ≤10² organisms¹⁷⁹. In all of these studies, the contaminating microflora consisted mainly of vegetative bacteria, usually of low pathogenicity (e.g., coagulase-negative *Staphylococcus*)^{179, 911, 912}.

An evaluation of the microbial load on used critical medical devices such as spinal anesthesia needles and angiographic catheters and sheaths demonstrated that mesophilic microorganisms were detected at levels of 10¹ to 10² in only two of five needles. The bioburden on used angiographic catheters

and sheath introducers exceeded 10³ CFUs on 14% (3 of 21) and 21% (6 of 28), respectively⁹⁰⁷.

Effect of Cleaning on Sterilization Efficacy

The effect of salt and serum on the efficacy of low-temperature sterilization technologies has raised concern regarding the margin of safety of these technologies. Experiments have shown that salts have the greatest impact on protecting microorganisms from killing426, 469. However, other studies have suggested that these concerns may not be clinically relevant. One study evaluated the relative rate of removal of inorganic salts, organic soil, and microorganisms from medical devices to better understand the dynamics of the cleaning process426. These tests were conducted by inoculating Alfa soil (tissueculture media and 10% fetal bovine serum) 469 containing 106 G. stearothermophilus spores onto the surface of a stainless-steel scalpel blade. After drying for 30 minutes at 35°C followed by 30 minutes at room temperature, the samples were placed in water at room temperature. The blades were removed at specified times, and the concentration of total protein and chloride ion was measured. The results showed that soaking in deionized water for 60 seconds resulted in a >95% release rate of chloride ion from NaCl solution in 20 seconds, Alfa soil in 30 seconds, and fetal bovine serum in 120 seconds. Thus, contact with water for short periods, even in the presence of protein, rapidly leads to dissolution of salt crystals and complete inactivation of spores by a low-temperature sterilization process (Table 10). Based on these experimental data, cleaning procedures would eliminate the detrimental effect of high salt content on a low-temperature sterilization process.

These articles ^{426, 469, 721} assessing low-temperature sterilization technology reinforce the importance of meticulous cleaning before sterilization. These data support the critical need for healthcare facilities to develop rigid protocols for cleaning contaminated objects before sterilization⁴⁷². Sterilization of instruments and medical devices is compromised if the process is not preceded by meticulous cleaning.

The cleaning of any narrow-lumen medical device used in patient care presents a major challenge to reprocessing areas. While attention has been focused on flexible endoscopes, cleaning issues related to other narrow-lumen medical devices such as sphinctertomes have been investigated⁹¹³. This study compared manual cleaning with that of automated cleaning with a narrow-lumen cleaner and found that only retro-flushing with the narrow lumen cleaner provided adequate cleaning of the three channels. If reprocessing was delayed for more than 24 hours, retro-flush cleaning was no longer effective and ETO sterilization failure was detected when devices were held for 7 days ⁹¹³. In another study involving simulated-use cleaning of laparoscopic devices, Alfa found that minimally the use of retro-flushing should be used during cleaning of non-ported laparoscopic devices⁹¹⁴.

Other Sterilization Methods

Ionizing Radiation. Sterilization by ionizing radiation, primarily by cobalt 60 gamma rays or electron accelerators, is a low-temperature sterilization method that has been used for a number of medical products (e.g., tissue for transplantation, pharmaceuticals, medical devices). There are no FDA-cleared ionizing radiation sterilization processes for use in healthcare facilities. Because of high sterilization costs, this method is an unfavorable alternative to ETO and plasma sterilization in healthcare facilities but is suitable for large-scale sterilization. Some deleterious effects on patient-care equipment associated with gamma radiation include induced oxidation in polyethylene ⁹¹⁵ and delamination and cracking in polyethylene knee bearings⁹¹⁶. Several reviews ^{917, 918} dealing with the sources, effects, and application of ionizing radiation may be referred to for more detail.

Dry-Heat Sterilizers. This method should be used only for materials that might be damaged by moist heat or that are impenetrable to moist heat (e.g., powders, petroleum products, sharp instruments). The advantages for dry heat include the following: it is nontoxic and does not harm the environment; a dry heat cabinet is easy to install and has relatively low operating costs; it penetrates materials; and it is noncorrosive for metal and sharp instruments. The disadvantages for dry heat are the slow rate of heat penetration and microbial killing makes this a time-consuming method. In addition, the high temperatures

are not suitable for most materials⁹¹⁹. The most common time-temperature relationships for sterilization with hot air sterilizers are 170°C (340°F) for 60 minutes, 160°C (320°F) for 120 minutes, and 150°C (300°F) for 150 minutes. *B. atrophaeus* spores should be used to monitor the sterilization process for dry heat because they are more resistant to dry heat than are *G. stearothermophilus* spores. The primary lethal process is considered to be oxidation of cell constituents.

There are two types of dry-heat sterilizers: the static-air type and the forced-air type. The static-air type is referred to as the oven-type sterilizer as heating coils in the bottom of the unit cause the hot air to rise inside the chamber via gravity convection. This type of dry-heat sterilizer is much slower in heating, requires longer time to reach sterilizing temperature, and is less uniform in temperature control throughout the chamber than is the forced-air type. The forced-air or mechanical convection sterilizer is equipped with a motor-driven blower that circulates heated air throughout the chamber at a high velocity, permitting a more rapid transfer of energy from the air to the instruments⁹²⁰.

Liquid Chemicals. Several FDA-cleared liquid chemical sterilants include indications for sterilization of medical devices (Tables 4 and 5)⁶⁹. The indicated contact times range from 3 hours to 12 hours. However, except for a few of the products, the contact time is based only on the conditions to pass the AOAC Sporicidal Test as a sterilant and not on simulated use testing with devices. These solutions are commonly used as high-level disinfectants when a shorter processing time is required. Generally, chemical liquid sterilants cannot be monitored using a biological indicator to verify sterility^{899, 900}.

The survival kinetics for thermal sterilization methods, such as steam and dry heat, have been studied and characterized extensively, whereas the kinetics for sterilization with liquid sterilants are less well understood⁹²¹. The information that is available in the literature suggests that sterilization processes based on liquid chemical sterilants, in general, may not convey the same sterility assurance level as sterilization achieved using thermal or physical methods⁸²³. The data indicate that the survival curves for liquid chemical sterilants may not exhibit log-linear kinetics and the shape of the survivor curve may vary depending of the formulation, chemical nature and stability of the liquid chemical sterilant. In addition, the design of the AOAC Sporicidal Test does not provide quantification of the microbial challenge. Therefore, sterilization with a liquid chemical sterilant may not convey the same sterility assurance as other sterilization methods.

One of the differences between thermal and liquid chemical processes for sterilization of devices is the accessibility of microorganisms to the sterilant. Heat can penetrate barriers, such as biofilm, tissue, and blood, to attain organism kill, whereas liquids cannot adequately penetrate these barriers. In addition, the viscosity of some liquid chemical sterilants impedes their access to organisms in the narrow lumens and mated surfaces of devices⁹²². Another limitation to sterilization of devices with liquid chemical germicides is the post-processing environment of the device. Devices cannot be wrapped or adequately contained during processing in a liquid chemical sterilant to maintain sterility following processing and during storage. Furthermore, devices may require rinsing following exposure to the liquid chemical sterilant with water that typically is not sterile. Therefore, due to the inherent limitations of using liquid chemical sterilants, their use should be restricted to reprocessing critical devices that are heat-sensitive and incompatible with other sterilization methods.

Several published studies compare the sporicidal effect of liquid chemical germicides against spores of *Bacillus* and *Clostridium*^{78, 659, 660, 715}.

Performic Acid. Performic acid is a fast-acting sporicide that was incorporated into an automated endoscope reprocessing system⁴⁰⁰. Systems using performic acid are not currently FDA cleared.

Filtration. Although filtration is not a lethality-based process and is not an FDA-cleared sterilization method, this technology is used to remove bacteria from thermolabile pharmaceutical fluids that cannot be purified by any other means. In order to remove bacteria, the membrane pore size (e.g., 0.22 μ m)

must be smaller than the bacteria and uniform throughout⁹²³. Some investigators have appropriately questioned whether the removal of microorganisms by filtration really is a sterilization method because of slight bacterial passage through filters, viral passage through filters, and transference of the sterile filtrate into the final container under aseptic conditions entail a risk of contamination⁹²⁴.

Microwave. Microwaves are used in medicine for disinfection of soft contact lenses, dental instruments, dentures, milk, and urinary catheters for intermittent self-catheterization⁹²⁵⁻⁹³¹. However, microwaves must only be used with products that are compatible (e.g., do not melt) 931. Microwaves are radio-frequency waves, which are usually used at a frequency of 2450 MHz. The microwaves produce friction of water molecules in an alternating electrical field. The intermolecular friction derived from the vibrations generates heat and some authors believe that the effect of microwaves depends on the heat produced while others postulate a nonthermal lethal effect⁹³²⁻⁹³⁴. The initial reports showed microwaves to be an effective microbicide. The microwaves produced by a "home-type" microwave oven (2.45 GHz) completely inactivate bacterial cultures, mycobacteria, viruses, and G. stearothermophilus spores within 60 seconds to 5 minutes depending on the challenge organism^{933, 935-937}. Another study confirmed these results but also found that higher power microwaves in the presence of water may be needed for sterilization⁹³². Complete destruction of Mycobacterium bovis was obtained with 4 minutes of microwave exposure (600W, 2450 MHz)⁹³⁷. The effectiveness of microwave ovens for different sterilization and disinfection purposes should be tested and demonstrated as test conditions affect the results (e.g., presence of water, microwave power). Sterilization of metal instruments can be accomplished but requires certain precautions.⁹²⁶. Of concern is that home-type microwave ovens may not have even distribution of microwave energy over the entire dry device (there may be hot and cold spots on solid medical devices); hence there may be areas that are not sterilized or disinfected. The use of microwave ovens to disinfect intermittent-use catheters also has been suggested. Researchers found that test bacteria (e.g., E. coli, Klebsiella pneumoniae, Candida albicans) were eliminated from red rubber catheters within 5 minutes ⁹³¹. Microwaves used for sterilization of medical devices have not been FDA cleared.

Glass Bead "Sterilizer". Glass bead "sterilization" uses small glass beads (1.2-1.5 mm diameter) and high temperature (217 °C -232°C) for brief exposure times (e.g., 45 seconds) to inactivate microorganisms. These devices have been used for several years in the dental profession⁹³⁸⁻⁹⁴⁰. FDA believes there is a risk of infection with this device because of potential failure to sterilize dental instruments and their use should be discontinued until the device has received FDA clearance.

Vaporized Hydrogen Peroxide (VHP®). Hydrogen peroxide solutions have been used as chemical sterilants for many years. However, the VHP® was not developed for the sterilization of medical equipment until the mid-1980s. One method for delivering VHP to the reaction site uses a deep vacuum to pull liquid hydrogen peroxide (30-35% concentration) from a disposable cartridge through a heated vaporizer and then, following vaporization, into the sterilization chamber. A second approach to VHP delivery is the flow-through approach in which the VHP is carried into the sterilization chamber by a carrier gas such as air using either a slight negative pressure (vacuum) or slight positive pressure. Applications of this technology include vacuum systems for industrial sterilization of medical devices and atmospheric systems for decontaminating for large and small areas⁸⁵³. VHP offers several appealing features that include rapid cycle time (e.g., 30-45 minutes); low temperature; environmentally safe by-products (H₂O, oxygen [O₂]); good material compatibility; and ease of operation, installation and monitoring. VHP has limitations including that cellulose cannot be processed; nylon becomes brittle; and VHP penetration capabilities are less than those of ETO. VHP has not been cleared by FDA for sterilization of medical devices in healthcare facilities.

The feasibility of utilizing vapor-phase hydrogen peroxide as a surface decontaminant and sterilizer was evaluated in a centrifuge decontamination application. In this study, vapor-phase hydrogen peroxide was shown to possess significant sporicidal activity ⁹⁴¹. In preliminary studies, hydrogen peroxide vapor decontamination has been found to be a highly effective method of eradicating MRSA, *Serratia*

marcescens, Clostridium botulinum spores and *Clostridium difficile* from rooms, furniture, surfaces and/or equipment; however, further investigation of this method to demonstrate both safety and effectiveness in reducing infection rates are required⁹⁴²⁻⁹⁴⁵.

Ozone. Ozone has been used for years as a drinking water disinfectant. Ozone is produced when O_2 is energized and split into two monatomic (O_1) molecules. The monatomic oxygen molecules then collide with O_2 molecules to form ozone, which is O_3 . Thus, ozone consists of O_2 with a loosely bonded third oxygen atom that is readily available to attach to, and oxidize, other molecules. This additional oxygen atom makes ozone a powerful oxidant that destroys microorganisms but is highly unstable (i.e., half-life of 22 minutes at room temperature).

A new sterilization process, which uses ozone as the sterilant, was cleared by FDA in August 2003 for processing reusable medical devices. The sterilizer creates its own sterilant internally from USP grade oxygen, steam-quality water and electricity; the sterilant is converted back to oxygen and water vapor at the end of the cycle by a passing through a catalyst before being exhausted into the room. The duration of the sterilization cycle is about 4 h and 15 m, and it occurs at 30-35°C. Microbial efficacy has been demonstrated by achieving a SAL of 10⁻⁶ with a variety of microorganisms to include the most resistant microorganism, *Geobacillus stearothermophilus*.

The ozone process is compatible with a wide range of commonly used materials including stainless steel, titanium, anodized aluminum, ceramic, glass, silica, PVC, Teflon, silicone, polypropylene, polyethylene and acrylic. In addition, rigid lumen devices of the following diameter and length can be processed: internal diameter (ID): > 2 mm, length \leq 25 cm; ID > 3 mm, length \leq 47 cm; and ID > 4 mm, length \leq 60 cm.

The process should be safe for use by the operator because there is no handling of the sterilant, no toxic emissions, no residue to aerate, and low operating temperature means there is no danger of an accidental burn. The cycle is monitored using a self-contained biological indicator and a chemical indicator. The sterilization chamber is small, about 4 ft³ (Written communication, S Dufresne, July 2004).

A gaseous ozone generator was investigated for decontamination of rooms used to house patients colonized with MRSA. The results demonstrated that the device tested would be inadequate for the decontamination of a hospital room⁹⁴⁶.

Formaldehyde Steam. Low-temperature steam with formaldehyde is used as a low-temperature sterilization method in many countries, particularly in Scandinavia, Germany, and the United Kingdom. The process involves the use of formalin, which is vaporized into a formaldehyde gas that is admitted into the sterilization chamber. A formaldehyde concentration of 8-16 mg/l is generated at an operating temperature of 70-75°C. The sterilization cycle consists of a series of stages that include an initial vacuum to remove air from the chamber and load, followed by steam admission to the chamber with the vacuum pump running to purge the chamber of air and to heat the load, followed by a series of pulses of formaldehyde gas, followed by steam. Formaldehyde is removed from the sterilizer and load by repeated alternate evacuations and flushing with steam and air. This system has some advantages, e.g., the cycle time for formaldehyde gas is faster than that for ETO and the cost per cycle is relatively low. However, ETO is more penetrating and operates at lower temperatures than do steam/formaldehyde sterilizers. Low-temperature steam formaldehyde sterilization has been found effective against vegetative bacteria, mycobacteria, *B. atrophaeus* and *G. stearothermophilus* spores and *Candida albicans*⁹⁴⁷⁻⁹⁴⁹.

Formaldehyde vapor cabinets also may be used in healthcare facilities to sterilize heat-sensitive medical equipment⁹⁵⁰. Commonly, there is no circulation of formaldehyde and no temperature and humidity controls. The release of gas from paraformaldehyde tablets (placed on the lower tray) is slow and produces a low partial pressure of gas. The microbicidal quality of this procedure is unknown⁹⁵¹.

Reliable sterilization using formaldehyde is achieved when performed with a high concentration of

gas, at a temperature between 60° and 80°C and with a relative humidity of 75 to 100%.

Studies indicate that formaldehyde is a mutagen and a potential human carcinogen, and OSHA regulates formaldehyde. The permissible exposure limit for formaldehyde in work areas is 0.75 ppm measured as a 8-hour TWA. The OSHA standard includes a 2 ppm STEL (i.e., maximum exposure allowed during a 15-minute period). As with the ETO standard, the formaldehyde standard requires that the employer conduct initial monitoring to identify employees who are exposed to formaldehyde at or above the action level or STEL. If this exposure level is maintained, employers may discontinue exposure monitoring until there is a change that could affect exposure levels or an employee reports formaldehyde-related signs and symptoms^{269, 578}. The formaldehyde steam sterilization system has not been FDA cleared for use in healthcare facilities.

Gaseous chlorine dioxide. A gaseous chlorine dioxide system for sterilization of healthcare products was developed in the late 1980s^{853, 952, 953}. Chlorine dioxide is not mutagenic or carcinogenic in humans. As the chlorine dioxide concentration increases, the time required to achieve sterilization becomes progressively shorter. For example, only 30 minutes were required at 40 mg/l to sterilize the 10⁶ *B. atrophaeus* spores at 30° to 32°C⁹⁵⁴. Currently, no gaseous chlorine dioxide system is FDA cleared.

Vaporized Peracetic Acid. The sporicidal activity of peracetic acid vapor at 20, 40, 60, and 80% relative humidity and 25°C was determined on *Bacillus atrophaeus* spores on paper and glass surfaces. Appreciable activity occurred within 10 minutes of exposure to 1 mg of peracetic acid per liter at 40% or higher relative humidity⁹⁵⁵. No vaporized peracetic acid system is FDA cleared.

Infrared radiation. An infrared radiation prototype sterilizer was investigated and found to destroy *B. atrophaeus* spores. Some of the possible advantages of infrared technology include short cycle time, low energy consumption, no cycle residuals, and no toxicologic or environmental effects. This may provide an alternative technology for sterilization of selected heat-resistant instruments but there are no FDA-cleared systems for use in healthcare facilities ⁹⁵⁶.

The other sterilization technologies mentioned above may be used for sterilization of critical medical items if cleared by the FDA and ideally, the microbicidal effectiveness of the technology has been published in the scientific literature. The selection and use of disinfectants, chemical sterilants and sterilization processes in the healthcare field is dynamic, and products may become available that are not in existence when this guideline was written. As newer disinfectants and sterilization processes become available, persons or committees responsible for selecting disinfectants and sterilization processes should be guided by products cleared by FDA and EPA as well as information in the scientific literature.

Sterilizing Practices

Overview. The delivery of sterile products for use in patient care depends not only on the effectiveness of the sterilization process but also on the unit design, decontamination, disassembling and packaging of the device, loading the sterilizer, monitoring, sterilant quality and quantity, and the appropriateness of the cycle for the load contents, and other aspects of device reprocessing. Healthcare personnel should perform most cleaning, disinfecting, and sterilizing of patient-care supplies in a central processing department in order to more easily control quality. The aim of central processing is the orderly processing of medical and surgical instruments to protect patients from infections while minimizing risks to staff and preserving the value of the items being reprocessed⁹⁵⁷. Healthcare facilities should promote the same level of efficiency and safety in the preparation of supplies in other areas (e.g., operating room, respiratory therapy) as is practiced in central processing.

Ensuring consistency of sterilization practices requires a comprehensive program that ensures operator competence and proper methods of cleaning and wrapping instruments, loading the sterilizer, operating the sterilizer, and monitoring of the entire process. Furthermore, care must be consistent from

an infection prevention standpoint in all patient-care settings, such as hospital and outpatient facilities.

Sterilization Cycle Verification. A sterilization process should be verified before it is put into use in healthcare settings. All steam, ETO, and other low-temperature sterilizers are tested with biological and chemical indicators upon installation, when the sterilizer is relocated, redesigned, after major repair and after a sterilization failure has occurred to ensure they are functioning prior to placing them into routine use. Three consecutive empty steam cycles are run with a biological and chemical indicator in an appropriate test package or tray. Each type of steam cycle used for sterilization (e.g., vacuum-assisted, gravity) is tested separately. In a prevacuum steam sterilizer three consecutive empty cycles are also run with a Bowie-Dick test. The sterilizer is not put back into use until all biological indicators are negative and chemical indicators show a correct end-point response^{811-814, 819, 958}.

Biological and chemical indicator testing is also done for ongoing quality assurance testing of representative samples of actual products being sterilized and product testing when major changes are made in packaging, wraps, or load configuration. Biological and chemical indicators are placed in products, which are processed in a full load. When three consecutive cycles show negative biological indicators and chemical indicators with a correct end point response, you can put the change made into routine use^{811-814, 958}. Items processed during the three evaluation cycles should be quarantined until the test results are negative.

Physical Facilities. The central processing area(s) ideally should be divided into at least three areas: decontamination, packaging, and sterilization and storage. Physical barriers should separate the decontamination area from the other sections to contain contamination on used items. In the decontamination area reusable contaminated supplies (and possibly disposable items that are reused) are received, sorted, and decontaminated. The recommended airflow pattern should contain contaminates within the decontamination area and minimize the flow of contaminates to the clean areas. The American Institute of Architects ⁹⁵⁹ recommends negative pressure and no fewer than six air exchanges per hour in the decontamination area (AAMI recommends 10 air changes per hour) and 10 air changes per hour with positive pressure in the sterilizer equipment room. The packaging area is for inspecting, assembling, and packaging clean, but not sterile, material. The sterile storage area should be a limited access area with a controlled temperature (may be as high as 75°F) and relative humidity (30-60% in all works areas except sterile storage, where the relative humidity should not exceed 70%)⁸¹⁹. The floors and walls should be constructed of materials capable of withstanding chemical agents used for cleaning or disinfecting. Ceilings and wall surfaces should be constructed of non-shedding materials. Physical arrangements of processing areas are presented schematically in four references^{811, 819, 920, 957}.

Cleaning. As repeatedly mentioned, items must be cleaned using water with detergents or enzymatic cleaners ^{465, 466, 468} before processing. Cleaning reduces the bioburden and removes foreign material (i.e., organic residue and inorganic salts) that interferes with the sterilization process by acting as a barrier to the sterilization agent^{179, 426, 457, 911, 912}. Surgical instruments are generally presoaked or prerinsed to prevent drying of blood and tissue. Precleaning in patient-care areas may be needed on items that are heavily soiled with feces, sputum, blood, or other material. Items sent to central processing without removing gross soil may be difficult to clean because of dried secretions and excretions. Cleaning and decontamination should be done as soon as possible after items have been used.

Several types of mechanical cleaning machines (e.g., utensil washer-sanitizer, ultrasonic cleaner, washer-sterilizer, dishwasher, washer-disinfector) may facilitate cleaning and decontamination of most items. This equipment often is automated and may increase productivity, improve cleaning effectiveness, and decrease worker exposure to blood and body fluids. Delicate and intricate objects and heat- or moisture-sensitive articles may require careful cleaning by hand. All used items sent to the central processing area should be considered contaminated (unless decontaminated in the area of origin), handled with gloves (forceps or tongs are sometimes needed to avoid exposure to sharps), and decontaminated by one of the aforementioned methods to render them safer to handle. Items composed of more than one removable part should be disassembled. Care should be taken to ensure that all parts

are kept together, so that reassembly can be accomplished efficiently⁸¹¹.

Investigators have described the degree of cleanliness by visual and microscopic examination. One study found 91% of the instruments to be clean visually but, when examined microscopically, 84% of the instruments had residual debris. Sites that contained residual debris included junctions between insulating sheaths and activating mechanisms of laparoscopic instruments and articulations and grooves of forceps. More research is needed to understand the clinical significance of these findings ⁹⁶⁰ and how to ensure proper cleaning.

Personnel working in the decontamination area should wear household-cleaning-type rubber or plastic gloves when handling or cleaning contaminated instruments and devices. Face masks, eye protection such as goggles or full-length faceshields, and appropriate gowns should be worn when exposure to blood and contaminated fluids may occur (e.g., when manually cleaning contaminated devices)⁹⁶¹. Contaminated instruments are a source of microorganisms that could inoculate personnel through nonintact skin on the hands or through contact with the mucous membranes of eyes, nose, or mouth^{214, 811, 813}. Reusable sharps that have been in contact with blood present a special hazard. Employees must not reach with their gloved hands into trays or containers that hold these sharps to retrieve them²¹⁴. Rather, employees should use engineering controls (e.g., forceps) to retrieve these devices.

Packaging. Once items are cleaned, dried, and inspected, those requiring sterilization must be wrapped or placed in rigid containers and should be arranged in instrument trays/baskets according to the guidelines provided by the AAMI and other professional organizations^{454, 811-814, 819, 836, 962}. These guidelines state that hinged instruments should be opened; items with removable parts should be disassembled unless the device manufacturer or researchers provide specific instructions or test data to the contrary¹⁸¹; complex instruments should be prepared and sterilized according to device manufacturer's instructions and test data; devices with concave surfaces should be positioned to facilitate drainage of water; heavy items should be positioned not to damage delicate items; and the weight of the instrument set should be based on the design and density of the instruments and the distribution of metal mass^{811, 962}. While there is no longer a specified sterilization weight limit for surgical sets, heavy metal mass is a cause of wet packs (i.e., moisture inside the case and tray after completion of the sterilization cycle)⁹⁶³. Other parameters that may influence drying are the density of the wraps and the design of the set⁹⁶⁴.

There are several choices in methods to maintain sterility of surgical instruments, including rigid containers, peel-open pouches (e.g., self-sealed or heat-sealed plastic and paper pouches), roll stock or reels (i.e., paper-plastic combinations of tubing designed to allow the user to cut and seal the ends to form a pouch) ⁴⁵⁴ and sterilization wraps (woven and nonwoven). Healthcare facilities may use all of these packaging options. The packaging material must allow penetration of the sterilant, provide protection against contact contamination during handling, provide an effective barrier to microbial penetration, and maintain the sterility of the processed item after sterilization ⁹⁶⁵. An ideal sterilization wrap would successfully address barrier effectiveness, penetrability (i.e., allows sterilant to penetrate), aeration (e.g., allows ETO to dissipate), ease of use, drapeability, flexibility, puncture resistance, tear strength, toxicity, odor, waste disposal, linting, cost, and transparency⁹⁶⁶. Unacceptable packaging for use with ETO (e.g., foil, polyvinylchloride, and polyvinylidene chlorine [kitchen-type transparent wrap]) ⁸¹⁴ or hydrogen peroxide gas plasma (e.g., linens and paper) should not be used to wrap medical items.

In central processing, double wrapping can be done sequentially or nonsequentially (i.e., simultaneous wrapping). Wrapping should be done in such a manner to avoid tenting and gapping. The sequential wrap uses two sheets of the standard sterilization wrap, one wrapped after the other. This procedure creates a package within a package. The nonsequential process uses two sheets wrapped at the same time so that the wrapping needs to be performed only once. This latter method provides multiple layers of protection of surgical instruments from contamination and saves time since wrapping is

done only once. Multiple layers are still common practice due to the rigors of handling within the facility even though the barrier efficacy of a single sheet of wrap has improved over the years⁹⁶⁶. Written and illustrated procedures for preparation of items to be packaged should be readily available and used by personnel when packaging procedures are performed⁴⁵⁴.

Loading. All items to be sterilized should be arranged so all surfaces will be directly exposed to the sterilizing agent. Thus, loading procedures must allow for free circulation of steam (or another sterilant) around each item. Historically, it was recommended that muslin fabric packs should not exceed the maximal dimensions, weight, and density of 12 inches wide × 12 inches high × 20 inches long, 12 lbs, and 7.2 lbs per cubic foot, respectively. Due to the variety of textiles and metal/plastic containers on the market, the textile and metal/plastic container manufacturer and the sterilizer manufacturers should be consulted for instructions on pack preparation and density parameters⁸¹⁹.

There are several important basic principles for loading a sterilizer: allow for proper sterilant circulation; perforated trays should be placed so the tray is parallel to the shelf; nonperforated containers should be placed on their edge (e.g., basins); small items should be loosely placed in wire baskets; and peel packs should be placed on edge in perforated or mesh bottom racks or baskets^{454, 811, 836}.

Storage. Studies in the early 1970s suggested that wrapped surgical trays remained sterile for varying periods depending on the type of material used to wrap the trays. Safe storage times for sterile packs vary with the porosity of the wrapper and storage conditions (e.g., open versus closed cabinets). Heat-sealed, plastic peel-down pouches and wrapped packs sealed in 3-mil (3/1000 inch) polyethylene overwrap have been reported to be sterile for as long as 9 months after sterilization. The 3-mil polyethylene is applied after sterilization to extend the shelf life for infrequently used items⁹⁶⁷. Supplies wrapped in double-thickness muslin comprising four layers, or equivalent, remain sterile for at least 30 days. Any item that has been sterilized should not be used after the expiration date has been exceeded or if the sterilized package is wet, torn, or punctured.

Although some hospitals continue to date every sterilized product and use the time-related shelf-life practice, many hospitals have switched to an event-related shelf-life practice. This latter practice recognizes that the product should remain sterile until some event causes the item to become contaminated (e.g., tear in packaging, packaging becomes wet, seal is broken)⁹⁶⁸. Event-related factors that contribute to the contamination of a product include bioburden (i.e., the amount of contamination in the environment), air movement, traffic, location, humidity, insects, vermin, flooding, storage area space, open/closed shelving, temperature, and the properties of the wrap material^{966, 969}. There are data that support the event-related shelf-life practice⁹⁷⁰⁻⁹⁷². One study examined the effect of time on the sterile integrity of paper envelopes, peel pouches, and nylon sleeves. The most important finding was the absence of a trend toward an increased rate of contamination over time for any pack when placed in covered storage⁹⁷¹. Another evaluated the effectiveness of event-related outdating by microbiologically testing sterilized items. During the 2-year study period, all of the items tested were sterile⁹⁷². Thus, contamination of a sterile item is event-related and the probability of contamination increases with increased handling⁹⁷³.

Following the sterilization process, medical and surgical devices must be handled using aseptic technique in order to prevent contamination. Sterile supplies should be stored far enough from the floor (8 to 10 inches), the ceiling (5 inches unless near a sprinkler head [18 inches from sprinkler head]), and the outside walls (2 inches) to allow for adequate air circulation, ease of cleaning, and compliance with local fire codes (e.g., supplies must be at least 18 inches from sprinkler heads). Medical and surgical supplies should not be stored under sinks or in other locations where they can become wet. Sterile items that become wet are considered contaminated because moisture brings with it microorganisms from the air and surfaces. Closed or covered cabinets are ideal but open shelving may be used for storage. Any package that has fallen or been dropped on the floor must be inspected for damage to the packaging and contents (if the items are breakable). If the package is heat-sealed in impervious plastic and the seal is

still intact, the package should be considered not contaminated. If undamaged, items packaged in plastic need not be reprocessed.

Monitoring. The sterilization procedure should be monitored routinely by using a combination of mechanical, chemical, and biological indicators to evaluate the sterilizing conditions and indirectly the microbiologic status of the processed items. The mechanical monitors for steam sterilization include the daily assessment of cycle time and temperature by examining the temperature record chart (or computer printout) and an assessment of pressure via the pressure gauge. The mechanical monitors for ETO include time, temperature, and pressure recorders that provide data via computer printouts, gauges, and/or displays⁸¹⁴. Generally, two essential elements for ETO sterilization (i.e., the gas concentration and humidity) cannot be monitored in healthcare ETO sterilizers.

Chemical indicators are convenient, are inexpensive, and indicate that the item has been exposed to the sterilization process. In one study, chemical indicators were more likely than biological indicators to inaccurately indicate sterilization at marginal sterilization times (e.g., 2 minutes)⁸⁴⁷. Chemical indicators should be used in conjunction with biological indicators, but based on current studies should not replace them because they indicate sterilization at marginal sterilization time and because only a biological indicator consisting of resistant spores can measure the microbial killing power of the sterilization process.^{847, 974}. Chemical indicators are affixed on the outside of each pack to show that the package has been processed through a sterilization cycle, but these indicators do not prove sterilization has been achieved. Preferably, a chemical indicator also should be placed on the inside of each pack to verify sterilant penetration. Chemical indicators usually are either heat-or chemical-sensitive inks that change color when one or more sterilization parameters (e.g., steam-time, temperature, and/or saturated steam; ETO-time, temperature, relative humidity and/or ETO concentration) are present. Chemical indicators have been grouped into five classes based on their ability to monitor one or multiple sterilization parameters^{813, 819}. If the internal and/or external indicator suggests inadequate processing, the item should not be used⁸¹⁵. An air-removal test (Bowie-Dick Test) must be performed daily in an empty dynamic-air-removal sterilizer (e.g., prevacuum steam sterilizer) to ensure air removal.

Biological indicators are recognized by most authorities as being closest to the ideal monitors of the sterilization process ^{974, 975} because they measure the sterilization process directly by using the most resistant microorganisms (i.e., *Bacillus* spores), and not by merely testing the physical and chemical conditions necessary for sterilization. Since the *Bacillus* spores used in biological indicators are more resistant and present in greater numbers than are the common microbial contaminants found on patient-care equipment, the demonstration that the biological indicator has been inactivated strongly implies that other potential pathogens in the load have been killed⁸⁴⁴.

An ideal biological monitor of the sterilization process should be easy to use, be inexpensive, not be subject to exogenous contamination, provide positive results as soon as possible after the cycle so that corrective action may be accomplished, and provide positive results only when the sterilization parameters (e.g., steam-time, temperature, and/or saturated steam; ETO-time, temperature, relative humidity and/or ETO concentration) are inadequate to kill microbial contaminates⁸⁴⁷.

Biological indicators are the only process indicators that directly monitor the lethality of a given sterilization process. Spores used to monitor a sterilization process have demonstrated resistance to the sterilizing agent and are more resistant than the bioburden found on medical devices^{179, 911, 912}. *B. atrophaeus* spores (10⁶) are used to monitor ETO and dry heat, and *G. stearothermophilus* spores (10⁵) are used to monitor ETO and dry heat, and *G. stearothermophilus* spores (10⁵) are used to monitor steam sterilization, hydrogen peroxide gas plasma, and liquid peracetic acid sterilizers. *G. stearothermophilus* is incubated at 55-60°C, and *B. atrophaeus* is incubated at 35-37°C. Steam and low temperature sterilizers (e.g., hydrogen peroxide gas plasma, peracetic acid) should be monitored at least weekly with the appropriate commercial preparation of spores. If a sterilizer is used frequently (e.g., several loads per day), daily use of biological indicators allows earlier discovery of equipment malfunctions or procedural errors and thus minimizes the extent of patient surveillance and

product recall needed in the event of a positive biological indicator⁸¹¹. Each load should be monitored if it contains implantable objects. If feasible, implantable items should not be used until the results of spore tests are known to be negative.

Originally, spore-strip biological indicators required up to 7 days of incubation to detect viable spores from marginal cycles (i.e., when few spores remained viable). The next generation of biological indicator was self-contained in plastic vials containing a spore-coated paper strip and a growth media in a crushable glass ampoule. This indicator had a maximum incubation of 48 hours but significant failures could be detected in \leq 24 hours. A rapid-readout biological indicator that detects the presence of enzymes of *G. stearothermophilus* by reading a fluorescent product produced by the enzymatic breakdown of a nonfluorescent substrate has been marketed for the more than 10 years. Studies demonstrate that the sensitivity of rapid-readout tests for steam sterilization (1 hour for 132°C gravity sterilizers, 3 hrs for 121°C gravity and 132°C vacuum sterilizers) parallels that of the conventional sterilization-specific biological indicators ^{846, 847, 976, 977} and the fluorescent rapid readout results reliably predict 24- and 48-hour and 7- day growth⁹⁷⁸. The rapid-readout biological indicator is a dual indicator system as it also detects acid metabolites produced during growth of the *G. stearothermophilus* spores. This system is different from the indicator system consisting of an enzyme system of bacterial origin without spores. Independent comparative data using suboptimal sterilization cycles (e.g., reduced time or temperature) with the enzyme-based indicator system have not been published⁹⁷⁹.

A new rapid-readout ETO biological indicator has been designed for rapid and reliable monitoring of ETO sterilization processes. The indicator has been cleared by the FDA for use in the United States⁴⁰⁰. The rapid-readout ETO biological indicator detects the presence of *B. atrophaeus* by detecting a fluorescent signal indicating the activity of an enzyme present within the *B. atrophaeus* organism, beta-glucosidase. The fluorescence indicates the presence of an active spore-associated enzyme and a sterilization process failure. This indicator also detects acid metabolites produced during growth of the *B. atrophaeus* spore. Per manufacturer's data, the enzyme always was detected whenever viable spores were present. This was expected because the enzyme is relatively ETO resistant and is inactivated at a slightly longer exposure time than the spore. The rapid-readout ETO biological indicator can be used to monitor 100% ETO, and ETO-HCFC mixture sterilization cycles. It has not been tested in ETO-CO₂ mixture sterilization cycles.

The standard biological indicator used for monitoring full-cycle steam sterilizers does not provide reliable monitoring flash sterilizers⁹⁸⁰. Biological indicators specifically designed for monitoring flash sterilization are now available, and studies comparing them have been published^{846, 847, 981}.

Since sterilization failure can occur (about 1% for steam)⁹⁸², a procedure to follow in the event of positive spore tests with steam sterilization has been provided by CDC and the Association of periOperative Registered Nurses (AORN). The 1981 CDC recommendation is that "objects, other than implantable objects, do not need to be recalled because of a single positive spore test unless the steam sterilizer or the sterilization procedure is defective." The rationale for this recommendation is that single positive spore tests in sterilizers occur sporadically. They may occur for reasons such as slight variation in the resistance of the spores⁹⁸³, improper use of the sterilizer, and laboratory contamination during culture (uncommon with self-contained spore tests). If the mechanical (e.g., time, temperature, pressure in the steam sterilizer) and chemical (internal and/or external) indicators suggest that the sterilizer was functioning properly, a single positive spore test probably does not indicate sterilizer malfunction but the spore test should be repeated immediately 983. If the spore tests remain positive, use of the sterilizer should be discontinued until it is serviced¹. Similarly, AORN states that a single positive spore test does not necessarily indicate a sterilizer failure. If the test is positive, the sterilizer should immediately be rechallenged for proper use and function. Items, other than implantable ones, do not necessarily need to be recalled unless a sterilizer malfunction is found. If a sterilizer malfunction is discovered, the items must be considered nonsterile, and the items from the suspect load(s) should be recalled, insofar as possible, and reprocessed ⁹⁸⁴. A suggested protocol for management of positive biological indicators is shown in

Table 12⁸³⁹. A more conservative approach also has been recommended ⁸¹³ in which any positive spore test is assumed to represent sterilizer malfunction and requires that all materials processed in that sterilizer, dating from the sterilization cycle having the last negative biologic indicator to the next cycle showing satisfactory biologic indicator challenge results, must be considered nonsterile and retrieved, if possible, and reprocessed. This more conservative approach should be used for sterilization methods other than steam (e.g., ETO, hydrogen peroxide gas plasma). However, no action is necessary if there is strong evidence for the biological indicator being defective ⁹⁸³ or the growth medium contained a *Bacillus* contaminant⁹⁸⁵.

If patient-care items were used before retrieval, the infection control professional should assess the risk of infection in collaboration with central processing, surgical services, and risk management staff. The factors that should be considered include the chemical indicator result (e.g., nonreactive chemical indicator may indicate temperature not achieved); the results of other biological indicators that followed the positive biological indicator (e.g., positive on Tuesday, negative on Wednesday); the parameters of the sterilizer associated with the positive biological indicator (e.g., reduced time at correct temperature); the time-temperature chart (or printout); and the microbial load associated with decontaminated surgical instruments (e.g., 85% of decontaminated surgical instruments have less than 100 CFU). The margin of safety in steam sterilization is sufficiently large that there is minimal infection risk associated with items in a load that show spore growth, especially if the item was properly cleaned and the temperature was achieved (e.g., as shown by acceptable chemical indicator or temperature chart). There are no published studies that document disease transmission via a nonretrieved surgical instrument following a sterilization cycle with a positive biological indicator.

False-positive biological indicators may occur from improper testing or faulty indicators. The latter may occur from improper storage, processing, product contamination, material failure, or variation in resistance of spores. Gram stain and subculture of a positive biological indicator may determine if a contaminant has created a false-positive result^{839, 986}. However, in one incident, the broth used as growth medium contained a contaminant, *B. coagulans,* which resulted in broth turbidity at 55°C⁹⁸⁵. Testing of paired biological indicators from different manufacturers can assist in assessing a product defect⁸³⁹. False-positive biological indicators due to extrinsic contamination when using self-contained biological indicator until a thorough analysis of the entire sterilization process shows this to be likely.

The size and composition of the biological indicator test pack should be standardized to create a significant challenge to air removal and sterilant penetration and to obtain interpretable results. There is a standard 16-towel pack recommended by AAMI for steam sterilization ^{813, 819, 987} consisting of 16 clean, preconditioned, reusable huck or absorbent surgical towels each of which is approximately 16 inches by 26 inches. Each towel is folded lengthwise into thirds and then folded widthwise in the middle. One or more biological indicators are placed between the eight and ninth towels in the approximate geometric center of the pack. When the towels are folded and placed one on top of another, to form a stack (approximately 6 inch height) it should weigh approximately 3 pounds and should have a density of approximately 11.3 pounds per cubic foot⁸¹³. This test pack has not gained universal use as a standard pack that simulates the actual in-use conditions of steam sterilizers. Commercially available disposable test packs that have been shown to be equivalent to the AAMI 16 towel test pack also may be used. The test pack should be placed flat in an otherwise fully loaded sterilizer chamber, in the area least favorable to sterilization (i.e., the area representing the greatest challenge to the biological indicator). This area is normally in the front, bottom section of the sterilizer, near the drain^{811, 813}. A control biological indicator from the lot used for testing should be left unexposed to the sterilant, and then incubated to verify the presterilization viability of the test spores and proper incubation. The most conservative approach would be to use a control for each run; however, less frequent use may be adequate (e.g., weekly). There also is a routine test pack for ETO where a biological indicator is placed in a plastic syringe with plunger, then placed in the folds of a clean surgical towel, and wrapped. Alternatively, commercially available disposal test packs that have been shown to be equivalent to the AAMI test pack may be used. The test pack is placed in the center of the sterilizer load⁸¹⁴. Sterilization records (mechanical, chemical, and biological)

should be retained for a time period in compliance with standards (e.g., Joint Commission for the Accreditation of Healthcare Facilities requests 3 years) and state and federal regulations.

In Europe, biological monitors are not used routinely to monitor the sterilization process. Instead, release of sterilizer items is based on monitoring the physical conditions of the sterilization process that is termed "parametric release." Parametric release requires that there is a defined quality system in place at the facility performing the sterilization and that the sterilization process be validated for the items being sterilized. At present in Europe, parametric release is accepted for steam, dry heat, and ionizing radiation processes, as the physical conditions are understood and can be monitored directly⁹⁸⁸. For example, with steam sterilizers the load could be monitored with probes that would yield data on temperature, time, and humidity at representative locations in the chamber and compared to the specifications developed during the validation process.

Periodic infection control rounds to areas using sterilizers to standardize the sterilizer's use may identify correctable variances in operator competence; documentation of sterilization records, including chemical and biological indicator test results; sterilizer maintenance and wrapping; and load numbering of packs. These rounds also may identify improvement activities to ensure that operators are adhering to established standards⁹⁸⁹.

Reuse of Single-Use Medical Devices

The reuse of single-use medical devices began in the late 1970s. Before this time most devices were considered reusable. Reuse of single-use devices increased as a cost-saving measure. Approximately 20 to 30% of U.S. hospitals reported that they reuse at least one type of single-use device. Reuse of singleuse devices involves regulatory, ethical, medical, legal and economic issues and has been extremely controversial for more than two decades⁹⁹⁰. The U.S. public has expressed increasing concern regarding the risk of infection and injury when reusing medical devices intended and labeled for single use. Although some investigators have demonstrated it is safe to reuse disposable medical devices such as cardiac electrode catheters, 991-993 additional studies are needed to define the risks 994 and document the benefits. In August 2000, FDA released a guidance document on single-use devices reprocessed by third parties or hospitals⁹⁹⁵. In this guidance document, FDA states that hospitals or third-party reprocessors will be considered "manufacturers" and regulated in the same manner. A reused single-use device will have to comply with the same regulatory requirements of the device when it was originally manufactured. This document presents FDA's intent to enforce premarket submission requirements within 6 months (February 2001) for class III devices (e.g., cardiovascular intra-aortic balloon pump, transluminal coronary angioplasty catheter); 12 months (August 2001) for class II devices (e.g., blood pressure cuff, bronchoscope biopsy forceps); and 18 months (February 2002) for class I devices (e.g., disposable medical scissors, ophthalmic knife). FDA uses two types of premarket requirements for nonexempt class I and II devices, a 510(k) submission that may have to show that the device is as safe and effective as the same device when new, and a premarket approval application. The 510(k) submission must provide scientific evidence that the device is safe and effective for its intended use. FDA allowed hospitals a year to comply with the nonpremarket requirements (registration and listing, reporting adverse events associated with medical devices, quality system regulations, and proper labeling). The options for hospitals are to stop reprocessing single-use devices, comply with the rule, or outsource to a third-party reprocessor. FDA guidance document does not apply to permanently implantable pacemakers, hemodialyzers, opened but unused single-use devices, or healthcare settings other than acute-care hospitals. The reuse of single use medical devices continues to be an evolving area of regulations. For this reason, healthcare workers should refer to FDA (http://www.fda.gov/) for the latest guidance 996.

Conclusion

When properly used, disinfection and sterilization can ensure the safe use of invasive and noninvasive medical devices. However, current disinfection and sterilization guidelines must be strictly followed.

Web-Based Disinfection and Sterilization Resources

Additional information about disinfection and sterilization is available at the following dedicated websites:

Food and Drug Administration, Rockville, Maryland [This link is no longer active: http://www.fda.gov/dcrh/ode/germlab.html.]

Environmental Protection Agency, Washington, D.C <u>Selected EPA-registered Disinfectants</u> (<u>https://www.epa.gov/pesticide-registration/selected-epa-registered-disinfectants</u>)

Centers for Disease Control and Prevention, Atlanta, Georgia <u>Healthcare-associated Infections (HAI)</u> (https://www.cdc.gov/HAI/index.html)

University of North Carolina, Chapel Hill, North Carolina Disinfection and Sterilization (http://www.disinfectionandsterilization.org)

Recommendations for Disinfection and Sterilization in Healthcare Facilities

A. Rationale

The ultimate goal of the Recommendations for Disinfection and Sterilization in Health-Care Facilities, 2008, is to reduce rates of health-care associated infections through appropriate use of both disinfection and sterilization. Each recommendation is categorized according to scientific evidence, theoretical rationale, applicability, and federal regulations. Examples are included in some recommendations to aid the reader; however, these examples are not intended to define the only method of implementing the recommendation. The CDC system for categorizing recommendations is defined in the following (Rankings) section.

B. Rankings

Category IA. Strongly recommended for implementation and strongly supported by well-designed experimental, clinical, or epidemiologic studies.

Category IB. Strongly recommended for implementation and supported by some experimental, clinical, or epidemiologic studies, and by a strong theoretical rationale.

Category IC. Required by state or federal regulations. Because of state differences, readers should not assume that the absence of an *IC* recommendation implies the absence of state regulations.

Category II. Suggested for implementation and supported by suggestive clinical or epidemiologic studies or by a theoretical rationale.

No recommendation. Unresolved issue. These include practices for which insufficient evidence or no consensus exists regarding efficacy.

C. Recommendations

1. Occupational Health and Exposure

- a. Inform each worker of the possible health effects of his or her exposure to infectious agents (e.g., hepatitis B virus [HBV], hepatitis C virus, human immunodeficiency virus [HIV]), and/or chemicals (e.g., EtO, formaldehyde). The information should be consistent with Occupational Safety and Health Administration (OSHA) requirements and identify the areas and tasks in which potential exists for exposure. *Category II, IC*^{214, 320, 959, 997, 998}
- b. Educate health-care workers in the selection and proper use of personal protective equipment (PPE). *Category II, IC*
- c. Ensure that workers wear appropriate PPE to preclude exposure to infectious agents or chemicals through the respiratory system, skin, or mucous membranes of the eyes, nose, or mouth. PPE can include gloves, gowns, masks, and eye protection. The exact type of PPE depends on the infectious or chemical agent and the anticipated duration of exposure. The employer is responsible for making such equipment and training available. *Category II, IC.* ^{214, 997-999}
- d. Establish a program for monitoring occupational exposure to regulated chemicals (e.g., formaldehyde, EtO) that adheres to state and federal regulations. *Category II, IC*. ^{997, 1000, 1001}
- e. Exclude healthcare workers with weeping dermatitis of hands from direct contact with patientcare equipment. *Category IB*. ^{1002, 1003}

2. Cleaning of Patient-Care Devices

- a. In hospitals, perform most cleaning, disinfection, and sterilization of patient-care devices in a central processing department in order to more easily control quality. *Category II*. ^{454, 836, 959}
- b. Meticulously clean patient-care items with water and detergent, or with water and enzymatic cleaners before high-level disinfection or sterilization procedures. *Category IB*. ^{6, 83, 101, 104-106, 124,}

179, 424-426, 436, 465, 471, 911-913, 1004

- *i.* Remove visible organic residue (e.g., residue of blood and tissue) and inorganic salts with cleaning. Use cleaning agents that are capable of removing visible organic and inorganic residues. *Category IB*. ^{424-426, 466, 468, 469, 471, 908, 910}
- *ii.* Clean medical devices as soon as practical after use (e.g., at the point of use) because soiled materials become dried onto the instruments. Dried or baked materials on the instrument make the removal process more difficult and the disinfection or sterilization process less effective or ineffective. *Category IB.* ^{55, 56, 59, 291, 465, 1005, 1006}
- c. Perform either manual cleaning (i.e., using friction) or mechanical cleaning (e.g., with ultrasonic cleaners, washer-disinfector, washer-sterilizers). *Category IB*. ^{426, 456, 471, 999}
- d. If using an automatic washer/disinfector, ensure that the unit is used in accordance with the manufacturer's recommendations. *Category IB*. ^{7, 133, 155, 725}
- e. Ensure that the detergents or enzymatic cleaners selected are compatible with the metals and other materials used in medical instruments. Ensure that the rinse step is adequate for removing cleaning residues to levels that will not interfere with subsequent disinfection/sterilization processes. *Category II.* ^{836, 1004}
- f. Inspect equipment surfaces for breaks in integrity that would impair either cleaning or disinfection/sterilization. Discard or repair equipment that no longer functions as intended or cannot be properly cleaned, and disinfected or sterilized. *Category II.* ⁸⁸⁸

3. Indications for Sterilization, High-Level Disinfection, and Low-Level Disinfection

- a. Before use on each patient, sterilize critical medical and surgical devices and instruments that enter normally sterile tissue or the vascular system or through which a sterile body fluid flows (e.g., blood). See recommendation 7g for exceptions. *Category IA*. ^{179, 497, 821, 822, 907, 911, 912}
- b. Provide, at a minimum, high-level disinfection for semicritical patient-care equipment (e.g., gastrointestinal endoscopes, endotracheal tubes, anesthesia breathing circuits, and respiratory therapy equipment) that touches either mucous membranes or nonintact skin. *Category IA*. ^{6-8, 17, 20, 99, 101, 108, 113-115, 129, 138, 139, 147, 152-154, 471, 1007}
- c. Perform low-level disinfection for noncritical patient-care surfaces (e.g., bedrails, over-the-bed table) and equipment (e.g., blood pressure cuff) that touch intact skin (see Recommendation 5g). *Category II*. ^{17, 46-48, 50-52, 67, 68, 372, 373, 378, 382, 401}

4. Selection and Use of Low-Level Disinfectants for Noncritical Patient-Care Devices

- a. Process noncritical patient-care devices using a disinfectant and the concentration of germicide listed in Table 1. *Category IB*. ^{17, 46-48, 50-52, 67, 68, 378, 382, 401}
- b. Disinfect noncritical medical devices (e.g., blood pressure cuff) with an EPA-registered hospital disinfectant using the label's safety precautions and use directions. Most EPA-registered hospital disinfectants have a label contact time of 10 minutes. However, multiple scientific studies have demonstrated the efficacy of hospital disinfectants against pathogens with a contact time of at least 1 minute. By law, all applicable label instructions on EPA-registered products must be followed. If the user selects exposure conditions that differ from those on the EPA-registered product label, the user assumes liability from any injuries resulting from off-label use and is potentially subject to enforcement action under FIFRA. *Category IB*. ^{17, 47, 48, 50, 51, 53-57, 59, 60, 62-64, 355, 378, 382}
- c. Ensure that, at a minimum, noncritical patient-care devices are disinfected when visibly soiled and on a regular basis (such as after use on each patient or once daily or once weekly). *Category II.* ^{378, 380, 1008}
- d. If dedicated, disposable devices are not available, disinfect noncritical patient-care equipment after using it on a patient who is on contact precautions before using this equipment on another patient. *Category IB*. ^{47, 67, 391, 1009}

5. Cleaning and Disinfecting Environmental Surfaces in Healthcare Facilities

- Clean housekeeping surfaces (e.g., floors, tabletops) on a regular basis, when spills occur, and a. when these surfaces are visibly soiled. Category II. 23, 378, 380, 382, 1008, 1010
- Disinfect (or clean) environmental surfaces on a regular basis (e.g., daily, three times per week) b. and when surfaces are visibly soiled. Category II. 378, 380, 402, 1008
- c. Follow manufacturers' instructions for proper use of disinfecting (or detergent) products --- such as recommended use-dilution, material compatibility, storage, shelf-life, and safe use and disposal. Category II. 327, 365, 404
- Clean walls, blinds, and window curtains in patient-care areas when these surfaces are visibly d. contaminated or soiled. Category II. 1011
- Prepare disinfecting (or detergent) solutions as needed and replace these with fresh solution e. frequently (e.g., replace floor mopping solution every three patient rooms, change no less often than at 60-minute intervals), according to the facility's policy. Category IB. 68, 379
- f. Decontaminate mop heads and cleaning cloths regularly to prevent contamination (e.g., launder and dry at least daily). Category II. 68, 402, 403
- Use a one-step process and an EPA-registered hospital disinfectant designed for housekeeping g. purposes in patient care areas where
 - 1. uncertainty exists about the nature of the soil on the surfaces (e.g., blood or body fluid contamination versus routine dust or dirt); or
 - 2. uncertainty exists about the presence of multidrug resistant organisms on such surfaces. See 5n for recommendations requiring cleaning and disinfecting blood-contaminated surfaces.

Category II. 23, 47, 48, 51, 214, 378, 379, 382, 416, 1012

- Detergent and water are adequate for cleaning surfaces in nonpatient-care areas (e.g., h. administrative offices). Category II. 23
- Do not use high-level disinfectants/liquid chemical sterilants for disinfection of non-critical i. surfaces. Category IB. 23, 69, 318
- j. Wet-dust horizontal surfaces regularly (e.g., daily, three times per week) using clean cloths moistened with an EPA-registered hospital disinfectant (or detergent). Prepare the disinfectant (or detergent) as recommended by the manufacturer. Category II. 68, 378, 380, 402, 403, 1008
- Disinfect noncritical surfaces with an EPA-registered hospital disinfectant according to the label's k. safety precautions and use directions. Most EPA-registered hospital disinfectants have a label contact time of 10 minutes. However, many scientific studies have demonstrated the efficacy of hospital disinfectants against pathogens with a contact time of at least 1 minute. By law, the user must follow all applicable label instructions on EPA-registered products. If the user selects exposure conditions that differ from those on the EPA-registered product label, the user assumes liability for any injuries resulting from off-label use and is potentially subject to enforcement action under FIFRA. Category II, IC. 17, 47, 48, 50, 51, 53-57, 59, 60, 62-64, 355, 378, 382
- I. Do not use disinfectants to clean infant bassinets and incubators while these items are occupied. If disinfectants (e.g., phenolics) are used for the terminal cleaning of infant bassinets and incubators, thoroughly rinse the surfaces of these items with water and dry them before these items are reused. Category IB. 17, 739, 740
- Promptly clean and decontaminate spills of blood and other potentially infectious materials. m. Discard blood-contaminated items in compliance with federal regulations. Category IB. IC. 214
- For site decontamination of spills of blood or other potentially infectious materials (OPIM), n. implement the following procedures. Use protective gloves and other PPE (e.g., when sharps are involved use forceps to pick up sharps, and discard these items in a puncture-resistant container) appropriate for this task. Disinfect areas contaminated with blood spills using an EPAregistered tuberculocidal agent, a registered germicide on the EPA Lists D and E (i.e., products with specific label claims for HIV or HBV or freshly diluted hypochlorite solution. Category II, IC. 214, 215, 557, 1013

- * If sodium hypochlorite solutions are selected use a 1:100 dilution (e.g., 1:100 dilution of a 5.25-6.15% sodium hypochlorite provides 525-615 ppm available chlorine) to decontaminate nonporous surfaces after a small spill (e.g., <10 mL) of either blood or OPIM. If a spill involves large amounts (e.g., >10 mL) of blood or OPIM, or involves a culture spill in the laboratory, use a 1:10 dilution for the first application of hypochlorite solution before cleaning in order to reduce the risk of infection during the cleaning process in the event of a sharp injury. Follow this decontamination process with a terminal disinfection, using a 1:100 dilution of sodium hypochlorite. *Category IB, IC.* ^{63, 215, 557}
- o. If the spill contains large amounts of blood or body fluids, clean the visible matter with disposable absorbent material, and discard the contaminated materials in appropriate, labeled containment. *Category II, IC.*^{44, 214}
- p. Use protective gloves and other PPE appropriate for this task. Category II, IC. 44, 214
- q. In units with high rates of endemic *Clostridium difficile* infection or in an outbreak setting, use dilute solutions of 5.25%–6.15% sodium hypochlorite (e.g., 1:10 dilution of household bleach) for routine environmental disinfection. Currently, no products are EPA-registered specifically for inactivating *C. difficile* spores. *Category II.* ²⁵⁷⁻²⁵⁹
- r. If chlorine solution is not prepared fresh daily, it can be stored at room temperature for up to 30 days in a capped, opaque plastic bottle with a 50% reduction in chlorine concentration after 30 days of storage (e.g., 1000 ppm chlorine [approximately a 1:50 dilution] at day 0 decreases to 500 ppm chlorine by day 30). *Category IB*. ^{327, 1014}
- s. An EPA-registered sodium hypochlorite product is preferred, but if such products are not available, generic versions of sodium hypochlorite solutions (e.g., household chlorine bleach) can be used. *Category II.*⁴⁴

6. Disinfectant Fogging

a. Do not perform disinfectant fogging for routine purposes in patient-care areas. Category II. ^{23, 228}

Environmental Fogging Clarification Statement [December 2009]: CDC and HICPAC have recommendations in the 2008 Guideline for Disinfection and Sterilization in Healthcare Facilities that state that the CDC does not support disinfectant fogging.

These recommendations refer to the spraying or fogging of chemicals (e.g., formaldehyde, phenolbased agents, or quaternary ammonium compounds) as a way to decontaminate environmental surfaces or disinfect the air in patient rooms. The recommendation against fogging was based on studies in the 1970's that reported a lack of microbicidal efficacy (e.g., use of quaternary ammonium compounds in mist applications) but also adverse effects on healthcare workers and others in facilities where these methods were utilized. Furthermore, some of these chemicals are not EPA-registered for use in fogging-type applications.

These recommendations do not apply to newer technologies involving fogging for room decontamination (e.g., ozone mists, vaporized hydrogen peroxide) that have become available since the 2008 recommendations were made. These newer technologies were assessed by CDC and HICPAC in the 2011 Guideline for the Prevention and Control of Norovirus Gastroenteritis Outbreaks in Healthcare Settings, which makes the recommendation:

"More research is required to clarify the effectiveness and reliability of fogging, UV irradiation, and ozone mists to reduce norovirus environmental contamination. (No recommendation/ unresolved issue)"

The 2008 recommendations still apply; however, CDC does not yet make a recommendation regarding these newer technologies. This issue will be revisited as additional evidence becomes available.

7. High-Level Disinfection of Endoscopes

- a. To detect damaged endoscopes, test each flexible endoscope for leaks as part of each reprocessing cycle. Remove from clinical use any instrument that fails the leak test, and repair this instrument. *Category II*. ^{113, 115, 116}
- b. Immediately after use, meticulously clean the endoscope with an enzymatic cleaner that is

compatible with the endoscope. Cleaning is necessary before both automated and manual disinfection. *Category IA.* ^{83, 101, 104-106, 113, 115, 116, 124, 126, 456, 465, 466, 471, 1015}

- c. Disconnect and disassemble endoscopic components (e.g., suction valves) as completely as possible and completely immerse all components in the enzymatic cleaner. Steam sterilize these components if they are heat stable. *Category IB*. ^{115, 116, 139, 465, 466}
- d. Flush and brush all accessible channels to remove all organic (e.g., blood, tissue) and other residue. Clean the external surfaces and accessories of the devices by using a soft cloth or sponge or brushes. Continue brushing until no debris appears on the brush. *Category IA* ^{6, 17, 108, 113, 115, 116, 137, 145, 147, 725, 856, 903}.
- e. Use cleaning brushes appropriate for the size of the endoscope channel or port (e.g., bristles should contact surfaces). Cleaning items (e.g., brushes, cloth) should be disposable or, if they are not disposable, they should be thoroughly cleaned and either high-level disinfected or sterilized after each use. *Category II.* ^{113, 115, 116, 1016}
- f. Discard enzymatic cleaners (or detergents) after each use because they are not microbicidal and, therefore, will not retard microbial growth. *Category IB*. ^{38, 113, 115, 116, 466}
- g. Process endoscopes (e.g., arthroscopes, cystoscope, laparoscopes) that pass through normally sterile tissues using a sterilization procedure before each use; if this is not feasible, provide at least high-level disinfection. High-level disinfection of arthroscopes, laparoscopes, and cystoscopes should be followed by a sterile water rinse. *Category IB*. ^{1, 17, 31, 32, 35, 89, 90, 113, 554}
- h. Phase out endoscopes that are critical items (e.g., arthroscopes, laparoscopes) but cannot be steam sterilized. Replace these endoscopes with steam sterilizable instruments when feasible. *Category II.*
- i. Mechanically clean reusable accessories inserted into endoscopes (e.g., biopsy forceps or other cutting instruments) that break the mucosal barrier (e.g., ultrasonically clean biopsy forceps) and then sterilize these items between each patient. *Category IA*. ^{1, 6, 8, 17, 108, 113, 115, 116, 138, 145, 147, 153, 278}
- j. Use ultrasonic cleaning of reusable endoscopic accessories to remove soil and organic material from hard-to-clean areas. *Category II*. ^{116, 145, 148}
- k. Process endoscopes and accessories that contact mucous membranes as semicritical items, and use at least high-level disinfection after use on each patient. *Category IA*. ^{1, 6, 8, 17, 108, 113, 115, 116, 129, 138, 145-148, 152-154, 278}
- I. Use an FDA-cleared sterilant or high-level disinfectant for sterilization or high-level disinfection (Table 1). *Category IA*. ^{1, 6-8, 17, 85, 108, 113, 115, 116, 147}
- m. After cleaning, use formulations containing glutaraldehyde, glutaraldehyde with phenol/phenate, ortho-phthalaldehyde, hydrogen peroxide, and both hydrogen peroxide and peracetic acid to achieve high-level disinfection followed by rinsing and drying (see Table 1 for recommended concentrations). *Category IB*. ^{1, 6-8, 17, 38, 85, 108, 113, 145-148}
- n. Extend exposure times beyond the minimum effective time for disinfecting semicritical patient-care equipment cautiously and conservatively because extended exposure to a high-level disinfectant is more likely to damage delicate and intricate instruments such as flexible endoscopes. The exposure times vary among the Food and Drug Administration (FDA)-cleared high-level disinfectants (Table 2). *Category IB*. ^{17, 69, 73, 76, 78, 83}
- o. Federal regulations are to follow the FDA-cleared label claim for high-level disinfectants. The FDA-cleared labels for high-level disinfection with >2% glutaraldehyde at 25°C range from 20-90 minutes, depending upon the product based on three tier testing which includes AOAC sporicidal tests, simulated use testing with mycobacteriał and in-use testing. *Category IC*.
- p. Several scientific studies and professional organizations support the efficacy of >2% glutaraldehyde for 20 minutes at 20°C; that efficacy assumes adequate cleaning prior to disinfection, whereas the FDA-cleared label claim incorporates an added margin of safety to accommodate possible lapses in cleaning practices. Facilities that have chosen to apply the 20 minute duration at 20°C have done so based on the IA recommendation in the July 2003 SHEA position paper, "Multi-society Guideline for Reprocessing Flexible Gastrointestinal Endoscopes." ^{12.}

17, 19, 26, 27, 49, 55, 57, 58, 60, 73, 76, 79-81, 83-85, 93, 94, 104-106, 110, 111, 115-121, 124, 125, 233, 235, 236, 243, 265, 266, 609

Update [June 2011]: <u>Multisociety guideline on reprocessing flexible gastrointestinal</u> <u>endoscopes: 2011</u>

(<u>http://www.asge.org/uploadedFiles/Publications_and_Products/Practice_Guidelines/Multisociety%</u> 20guideline%20on%20reprocessing%20flexible%20gastrointestinal.pdf).

- q. When using FDA-cleared high-level disinfectants, use manufacturers' recommended exposure conditions. Certain products may require a shorter exposure time (e.g., 0.55% ortho-phthalaldehyde for 12 minutes at 20°C, 7.35% hydrogen peroxide plus 0.23% peracetic acid for 15 minutes at 20°C) than glutaraldehyde at room temperature because of their rapid inactivation of mycobacteria or reduced exposure time because of increased mycobactericidal activity at elevated temperature (e.g., 2.5% glutaraldehyde at 5 minutes at 35°C). *Category IB.* ^{83, 100, 689, 694, 700}
- r. Select a disinfectant or chemical sterilant that is compatible with the device that is being reprocessed. Avoid using reprocessing chemicals on an endoscope if the endoscope manufacturer warns against using these chemicals because of functional damage (with or without cosmetic damage). *Category IB*. ^{69, 113, 116}
- s. Completely immerse the endoscope in the high-level disinfectant, and ensure all channels are perfused. As soon as is feasible, phase out nonimmersible endoscopes. *Category IB*. ^{108, 113-116, 137, 725, 856, 882}
- t. After high-level disinfection, rinse endoscopes and flush channels with sterile water, filtered water, or tapwater to prevent adverse effects on patients associated with disinfectant retained in the endoscope (e.g., disinfectant induced colitis). Follow this water rinse with a rinse with 70% 90% ethyl or isopropyl alcohol. *Category IB.* ^{17, 31-35, 38, 39, 108, 113, 115, 116, 134, 145-148, 620-622, 624-630, 1017}
- u. After flushing all channels with alcohol, purge the channels using forced air to reduce the likelihood of contamination of the endoscope by waterborne pathogens and to facilitate drying. *Category IB*. 39, 113, 115, 116, 145, 147
- v. Hang endoscopes in a vertical position to facilitate drying. Category II. 17, 108, 113, 115, 116, 145, 815
- w. Store endoscopes in a manner that will protect them from damage or contamination. *Category II*. 17, 108, 113, 115, 116, 145
- x. Sterilize or high-level disinfect both the water bottle used to provide intraprocedural flush solution and its connecting tube at least once daily. After sterilizing or high-level disinfecting the water bottle, fill it with sterile water. *Category IB.* ^{10, 31-35, 113, 116, 1017}
- y. Maintain a log for each procedure and record the following: patient's name and medical record number (if available), procedure, date, endoscopist, system used to reprocess the endoscope (if more than one system could be used in the reprocessing area), and serial number or other identifier of the endoscope used. *Category II*. ^{108, 113, 115, 116}
- z. Design facilities where endoscopes are used and disinfected to provide a safe environment for healthcare workers and patients. Use air-exchange equipment (e.g., the ventilation system, outexhaust ducts) to minimize exposure of all persons to potentially toxic vapors (e.g., glutaraldehyde vapor). Do not exceed the allowable limits of the vapor concentration of the chemical sterilant or high-level disinfectant (e.g., those of ACGIH and OSHA). Category IB, IC. ^{116, 145, 318, 322, 577, 652}
- aa. Routinely test the liquid sterilant/high-level disinfectant to ensure minimal effective concentration of the active ingredient. Check the solution each day of use (or more frequently) using the appropriate chemical indicator (e.g., glutaraldehyde chemical indicator to test minimal effective concentration of glutaraldehyde) and document the results of this testing. Discard the solution if the chemical indicator shows the concentration is less than the minimum effective concentration. Do not use the liquid sterilant/high-level disinfectant beyond the reuse-life recommended by the manufacturer (e.g., 14 days for *ortho*-phthalaldehyde). *Category IA*. ^{76, 108, 113, 115, 116, 608, 609}
- ab. * Provide personnel assigned to reprocess endoscopes with device-specific reprocessing instructions to ensure proper cleaning and high-level disinfection or sterilization. Require competency testing on a regular basis (e.g., beginning of employment, annually) of all personnel who reprocess endoscopes. *Category IA*. ^{6-8, 108, 113, 115, 116, 145, 145, 145, 155}

- ac. * Educate all personnel who use chemicals about the possible biologic, chemical, and environmental hazards of performing procedures that require disinfectants. Category IB, IC. ^{116, 997, 998, 1018, 1019}
- ad. * Make PPE (e.g., gloves, gowns, eyewear, face mask or shields, respiratory protection devices) available and use these items appropriately to protect workers from exposure to both chemicals and microorganisms (e.g., HBV). *Category IB, IC*. ^{115, 116, 214, 961, 997, 998, 1020, 1021}
- ae. * If using an automated endoscope reprocessor (AER), place the endoscope in the reprocessor and attach all channel connectors according to the AER manufacturer's instructions to ensure exposure of all internal surfaces to the high-level disinfectant/chemical sterilant. *Category IB*. ^{7, 8,} 115, 116, 155, 725, 903
- af. * If using an AER, ensure the endoscope can be effectively reprocessed in the AER. Also, ensure any required manual cleaning/disinfecting steps are performed (e.g., elevator wire channel of duodenoscopes might not be effectively disinfected by most AERs). *Category IB*. ^{7, 8, 115, 116, 155, 725}
- ag. * Review the FDA advisories and the scientific literature for reports of deficiencies that can lead to infection because design flaws and improper operation and practices have compromised the effectiveness of AERs. *Category II.* ^{7, 98, 133, 134, 155, 725}
- ah. * Develop protocols to ensure that users can readily identify an endoscope that has been properly processed and is ready for patient use. *Category II*.
- ai. * Do not use the carrying case designed to transport clean and reprocessed endoscopes outside of the healthcare environment to store an endoscope or to transport the instrument within the healthcare environment. *Category II.*
- aj. * No recommendation is made about routinely performing microbiologic testing of either endoscopes or rinse water for quality assurance purposes. *Unresolved Issue*. ^{116, 164}
- ak. * If environmental microbiologic testing is conducted, use standard microbiologic techniques. *Category II.*^{23, 116, 157, 161, 167}
- al. * If a cluster of endoscopy-related infections occurs, investigate potential routes of transmission (e.g., person-to-person, common source) and reservoirs. *Category IA*. ^{8, 1022}
- am.* Report outbreaks of endoscope-related infections to persons responsible for institutional infection control and risk management and to FDA. *Category IB*. ^{6, 7, 113, 116, 1023}
 - 1. * Notify the local and the state health departments, CDC, and the manufacturer(s). Category II.
- an. * No recommendation is made regarding the reprocessing of an endoscope again immediately before use if that endoscope has been processed after use according to the recommendations in this guideline. *Unresolved issue*. ¹⁵⁷
- ao. * Compare the reprocessing instructions provided by both the endoscope's and the AER's manufacturer's instructions and resolve any conflicting recommendations. *Category IB*. ^{116, 155}

8. Management of Equipment and Surfaces in Dentistry

- a. Dental instruments that penetrate soft tissue or bone (e.g., extraction forceps, scalpel blades, bone chisels, periodontal scalers, and surgical burs) are classified as critical and should be sterilized after each use or discarded. In addition, after each use, sterilize dental instruments that are not intended to penetrate oral soft tissue or bone (e.g., amalgam condensers, air-water syringes) but that might contact oral tissues and are heat-tolerant, although classified as semicritical. Clean and, at a minimum, high-level disinfect heat-sensitive semicritical items. *Category IA*. ^{43, 209-211}
- b. Noncritical clinical contact surfaces, such as uncovered operatory surfaces (e.g., countertops, switches, light handles), should be barrier-protected or disinfected between patients with an intermediate-disinfectant (i.e., EPA-registered hospital disinfectant with a tuberculocidal claim) or low-level disinfectant (i.e., EPA-registered hospital disinfectant with HIV and HBV claim). *Category IB.* ^{43, 209-211}
- *c.* Barrier protective coverings can be used for noncritical clinical contact surfaces that are touched frequently with gloved hands during the delivery of patient care, that are likely to become

contaminated with blood or body substances, or that are difficult to clean. Change these coverings when they are visibly soiled, when they become damaged, and on a routine basis (e.g., between patients). Disinfect protected surfaces at the end of the day or if visibly soiled. *Category II.* ^{43, 210}

- 9. Processing Patient-Care Equipment Contaminated with Bloodborne Pathogens (HBV, Hepatitis C Virus, HIV), Antibiotic-Resistant Bacteria (e.g., Vancomycin-Resistant Enterococci, Methicillin-Resistant Staphylococcus aureus, Multidrug Resistant Tuberculosis), or Emerging Pathogens (e.g., Cryptosporidium, Helicobacter pylori, Escherichia coli O157:H7, Clostridium difficile, Mycobacterium tuberculosis, Severe Acute Respiratory Syndrome Coronavirus), or Bioterrorist Agents
 - a. Use standard sterilization and disinfection procedures for patient-care equipment (as recommended in this guideline), because these procedures are adequate to sterilize or disinfect instruments or devices contaminated with blood or other body fluids from persons infected with bloodborne pathogens or emerging pathogens, with the exception of prions. No changes in these procedures for cleaning, disinfecting, or sterilizing are necessary for removing bloodborne and emerging pathogens other than prions. *Category IA*. ^{22, 53, 60-62, 73, 79-81, 105, 118-121, 125, 126, 221, 224-234, 236, 244, 265, 266, 271-273, 279, 282, 283, 354-357, 666}

10. Disinfection Strategies for Other Semicritical Devices

- a. Even if probe covers have been used, clean and high-level disinfect other semicritical devices such as rectal probes, vaginal probes, and cryosurgical probes with a product that is not toxic to staff, patients, probes, and retrieved germ cells (if applicable). Use a high-level disinfectant at the FDA-cleared exposure time. (See Recommendation 7p for exceptions.) *Category IB*. ^{6-8, 17, 69}
- When probe covers are available, use a probe cover or condom to reduce the level of microbial contamination. *Category II.* ¹⁹⁷⁻²⁰¹
 Do not use a lower category of disinfection or cease to follow the appropriate disinfectant recommendations when using probe covers because these sheaths and condoms can fail. *Category IB* ¹⁹⁷⁻²⁰¹
- c. After high-level disinfection, rinse all items. Use sterile water, filtered water or tapwater followed by an alcohol rinse for semicritical equipment that will have contact with mucous membranes of the upper respiratory tract (e.g., nose, pharynx, esophagus). *Category II*. ^{10, 31-35, 1017}
- d. There is no recommendation to use sterile or filtered water rather than tapwater for rinsing semicritical equipment that contact the mucous membranes of the rectum (e.g., rectal probes, anoscope) or vagina (e.g., vaginal probes). *Unresolved issue.*¹¹
- e. Wipe clean tonometer tips and then disinfect them by immersing for 5-10 minutes in either 5000 ppm chlorine or 70% ethyl alcohol. None of these listed disinfectant products are FDA-cleared high-level disinfectants. *Category II.* ^{49, 95, 185, 188, 293}

11. Disinfection by Healthcare Personnel in Ambulatory Care and Home Care

- a. Follow the same classification scheme described above (i.e., that critical devices require sterilization, semicritical devices require high-level disinfection, and noncritical equipment requires low-level disinfection) in the ambulatory-care (outpatient medical/surgical facilities) setting because risk for infection in this setting is similar to that in the hospital setting (see Table 1). *Category IB.* ^{6-8, 17, 330}
- When performing care in the home, clean and disinfect reusable objects that touch mucous membranes (e.g., tracheostomy tubes) by immersing these objects in a 1:50 dilution of 5.25%-6.15% sodium hypochlorite (household bleach) (3 minutes), 70% isopropyl alcohol (5 minutes), or 3% hydrogen peroxide (30 minutes) because the home environment is, in most instances, safer than either hospital or ambulatory care settings because person-to-person transmission is less likely. *Category II*. ^{327, 328, 330, 331}
- c. Clean noncritical items that would not be shared between patients (e.g., crutches, blood

pressure cuffs) in the home setting with a detergent or commercial household disinfectant. *Category II.* ^{53, 330}

12. Microbial Contamination of Disinfectants

- a. Institute the following control measures to reduce the occurrence of contaminated disinfectants:
 - 1. prepare the disinfectant correctly to achieve the manufacturer's recommended use-dilution; and
 - prevent common sources of extrinsic contamination of germicides (e.g., container contamination or surface contamination of the healthcare environment where the germicide are prepared and/or used).
 Category IB. ^{404, 406, 1024}

13. Flash Sterilization

- a. Do not flash sterilize implanted surgical devices unless doing so is unavoidable. Category IB. ^{849,} ⁸⁵⁰
- b. Do not use flash sterilization for convenience, as an alternative to purchasing additional instrument sets, or to save time. Category II. ^{817, 962}
- c. When using flash sterilization, make sure the following parameters are met:
 - 1. clean the item before placing it in the sterilizing container (that are FDA cleared for use with flash sterilization) or tray;
 - 2. prevent exogenous contamination of the item during transport from the sterilizer to the patient; and 3) monitor sterilizer function with mechanical, chemical, and biologic monitors. *Category IB.* ^{812, 819, 846, 847, 962}
- d. Do not use packaging materials and containers in flash sterilization cycles unless the sterilizer and the packaging material/container are designed for this use. *Category IB*. ^{812, 819, 1025}
- e. When necessary, use flash sterilization for patient-care items that will be used immediately (e.g., to reprocess an inadvertently dropped instrument). *Category IB*. ^{812, 817, 819, 845}
- f. When necessary, use flash sterilization for processing patient-care items that cannot be packaged, sterilized, and stored before use. *Category IB*. ^{812, 819}

14. Methods of Sterilization

- a. Steam is the preferred method for sterilizing critical medical and surgical instruments that are not damaged by heat, steam, pressure, or moisture. *Category IA*. ^{181, 271, 425, 426, 827, 841, 1026, 1027}
- b. Cool steam- or heat-sterilized items before they are handled or used in the operative setting. *Category IB*. ⁸⁵⁰
- c. Follow the sterilization times, temperatures, and other operating parameters (e.g., gas concentration, humidity) recommended by the manufacturers of the instruments, the sterilizer, and the container or wrap used, and that are consistent with guidelines published by government agencies and professional organizations. *Category IB.* ^{811-814, 819, 825, 827, 841, 1026-1028}
- d. Use low-temperature sterilization technologies (e.g., EtO, hydrogen peroxide gas plasma) for reprocessing critical patient-care equipment that is heat or moisture sensitive. *Category IA*^{469,}721, 825, 856, 858, 878, 879, 881, 882, 890, 891, 1027.
- e. Completely aerate surgical and medical items that have been sterilized in the EtO sterilizer (e.g., polyvinylchloride tubing requires 12 hours at 50°C, 8 hours at 60°C) before using these items in patient care. *Category IB*. ⁸¹⁴
- f. Sterilization using the peracetic acid immersion system can be used to sterilize heat-sensitive immersible medical and surgical items. *Category IB*. ^{90, 717-719, 721-724}
- g. Critical items that have been sterilized by the peracetic acid immersion process must be used immediately (i.e., items are not completely protected from contamination, making long-term storage unacceptable). *Category II.*^{817, 825}
- h. Dry-heat sterilization (e.g., 340°F for 60 minutes) can be used to sterilize items (e.g., powders, oils) that can sustain high temperatures. *Category IB*. ^{815, 827}
- i. Comply with the sterilizer manufacturer's instructions regarding the sterilizer cycle parameters

(e.g., time, temperature, concentration). Category IB. 155, 725, 811-814, 819

j. Because narrow-lumen devices provide a challenge to all low-temperature sterilization technologies and direct contact is necessary for the sterilant to be effective, ensure that the sterilant has direct contact with contaminated surfaces (e.g., scopes processed in peracetic acid must be connected to channel irrigators). *Category IB.* ^{137, 725, 825, 856, 890, 891, 1029}

15. Packaging

- a. Ensure that packaging materials are compatible with the sterilization process and have received FDA 510[k] clearance. *Category IB*. ^{811-814, 819, 966}
- b. Ensure that packaging is sufficiently strong to resist punctures and tears to provide a barrier to microorganisms and moisture. *Category IB*. ^{454, 811-814, 819, 966}

16. Monitoring of Sterilizers

- a. Use mechanical, chemical, and biologic monitors to ensure the effectiveness of the sterilization process. *Category IB*. 811-815, 819, 846, 847, 975-977
- b. Monitor each load with mechanical (e.g., time, temperature, pressure) and chemical (internal and external) indicators. If the internal chemical indicator is visible, an external indicator is not needed. *Category II*. ^{811-815, 819, 846, 847, 975-977, 980}
- c. Do not use processed items if the mechanical (e.g., time, temperature, pressure) or chemical (internal and/or external) indicators suggest inadequate processing. *Category IB*^{811-814, 819}.
- d. Use biologic indicators to monitor the effectiveness of sterilizers at least weekly with an FDAcleared commercial preparation of spores (e.g., *Geobacillus stearothermophilus* for steam) intended specifically for the type and cycle parameters of the sterilizer. *Category IB*. ^{1, 811, 813-815, 819, 846, 847, 976, 977}
- e. After a single positive biologic indicator used with a method other than steam sterilization, treat as nonsterile all items that have been processed in that sterilizer, dating from the sterilization cycle having the last negative biologic indicator to the next cycle showing satisfactory biologic indicator results. These nonsterile items should be retrieved if possible and reprocessed. *Category II.*¹
- f. After a positive biologic indicator with steam sterilization, objects other than implantable objects do not need to be recalled because of a single positive spore test unless the sterilizer or the sterilization procedure is defective as determined by maintenance personnel or inappropriate cycle settings. If additional spore tests remain positive, consider the items nonsterile and recall and reprocess the items from the implicated load(s). *Category II.* ¹
- g. Use biologic indicators for every load containing implantable items and quarantine items, whenever possible, until the biologic indicator is negative. *Category IB*. ^{811-814, 819}

17. Load Configuration.

a. Place items correctly and loosely into the basket, shelf, or cart of the sterilizer so as not to impede the penetration of the sterilant. *Category IB*. ^{445, 454, 811, 813, 819, 836}

18. Storage of Sterile Items

- a. Ensure the sterile storage area is a well-ventilated area that provides protection against dust, moisture, insects, and temperature and humidity extremes. *Category II.* ^{454, 819, 836, 969}
- b. Store sterile items so the packaging is not compromised (e.g., punctured, bent). *Category II.* ^{454,} 816, 819, 968, 969, 1030
- c. Label sterilized items with a load number that indicates the sterilizer used, the cycle or load number, the date of sterilization, and, if applicable, the expiration date. *Category IB*. ^{811, 812, 814, 816, 819}
- d. The shelf life of a packaged sterile item depends on the quality of the wrapper, the storage conditions, the conditions during transport, the amount of handling, and other events (moisture) that compromise the integrity of the package. If event-related storage of sterile items is used, then packaged sterile items can be used indefinitely unless the packaging is compromised (see f

and g below). Category IB. 816, 819, 836, 968, 973, 1030, 1031

- e. Evaluate packages before use for loss of integrity (e.g., torn, wet, punctured). The pack can be used unless the integrity of the packaging is compromised. *Category II*. ^{819, 968}
- f. If the integrity of the packaging is compromised (e.g., torn, wet, or punctured), repack and reprocess the pack before use. *Category II*. ^{819, 1032}
- g. If time-related storage of sterile items is used, label the pack at the time of sterilization with an expiration date. Once this date expires, reprocess the pack. *Category II.*^{819, 968}

19. Quality Control

- a. Provide comprehensive and intensive training for all staff assigned to reprocess semicritical and critical medical/surgical instruments to ensure they understand the importance of reprocessing these instruments. To achieve and maintain competency, train each member of the staff that reprocesses semicritical and/or critical instruments as follows:
 - 1. provide hands-on training according to the institutional policy for reprocessing critical and semicritical devices;
 - 2. supervise all work until competency is documented for each reprocessing task;
 - 3. conduct competency testing at beginning of employment and regularly thereafter (e.g., annually); and
 - 4. review the written reprocessing instructions regularly to ensure they comply with the scientific literature and the manufacturers' instructions. *Category IB.* ^{6-8, 108, 114, 129, 155, 725, 813, 819}
- b. Compare the reprocessing instructions (e.g., for the appropriate use of endoscope connectors, the capping/noncapping of specific lumens) provided by the instrument manufacturer and the sterilizer manufacturer and resolve any conflicting recommendations by communicating with both manufacturers. *Category IB*. ^{155, 725}
- c. Conduct infection control rounds periodically (e.g., annually) in high-risk reprocessing areas (e.g., the Gastroenterology Clinic, Central Processing); ensure reprocessing instructions are current and accurate and are correctly implemented. Document all deviations from policy. All stakeholders should identify what corrective actions will be implemented. *Category IB*. ^{6-8, 129}
- d. Include the following in a quality control program for sterilized items: a sterilizer maintenance contract with records of service; a system of process monitoring; air-removal testing for prevacuum steam sterilizers; visual inspection of packaging materials; and traceability of load contents. *Category II*^{811-814, 819}.
- e. For each sterilization cycle, record the type of sterilizer and cycle used; the load identification number; the load contents; the exposure parameters (e.g., time and temperature); the operator's name or initials; and the results of mechanical, chemical, and biological monitoring. *Category II* 811-814, 819
- f. Retain sterilization records (mechanical, chemical, and biological) for a time period that complies with standards (e.g., 3 years), statutes of limitations, and state and federal regulations. *Category II, IC*. ¹⁰³³
- g. Prepare and package items to be sterilized so that sterility can be achieved and maintained to the point of use. Consult the Association for the Advancement of Medical Instrumentation or the manufacturers of surgical instruments, sterilizers, and container systems for guidelines for the density of wrapped packages. *Category II.* ^{811-814, 819}
- h. Periodically review policies and procedures for sterilization. Category II. ¹⁰³³
- i. Perform preventive maintenance on sterilizers by qualified personnel who are guided by the manufacturer's instruction. *Category II*. ^{811-814, 819}

20. Reuse of Single-Use Medical Devices

a. Adhere to the FDA enforcement document for single-use devices reprocessed by hospitals. FDA considers the hospital that reprocesses a single-use device as the manufacturer of the device and regulates the hospital using the same standards by which it regulates the original equipment manufacturer. *Category II, IC.* ⁹⁹⁵

Performance Indicators

- 1. Monitor adherence to high-level disinfection and/or sterilization guidelines for endoscopes on a regular basis. This monitoring should include ensuring the proper training of persons performing reprocessing and their adherence to all endoscope reprocessing steps, as demonstrated by competency testing at commencement of employment and annually.
- 2. Develop a mechanism for the occupational health service to report all adverse health events potentially resulting from exposure to disinfectants and sterilants; review such exposures; and implement engineering, work practice, and PPE to prevent future exposures.
- 3. Monitor possible sterilization failures that resulted in instrument recall. Assess whether additional training of personnel or equipment maintenance is required.

Acknowledgements

The authors gratefully acknowledge Eva P. Clontz, M.S., for her assistance in referencing this guideline. The Healthcare Infection Control Practices Advisory Committee thanks the following experts for reviewing a draft of this guideline: Martin S. Favero, Ph.D., Syed A. Sattar, Ph.D., A. Denver Russell, D.Sc., and Martin Exner, M.D. The opinions of the reviewers might not be reflected in all the recommendations contained in this document.

Glossary

Action level: concentration of a regulated substance (e.g., ethylene oxide, formaldehyde) within the employee breathing zone, above which OSHA requirements apply.

Activation of a sterilant: process of mixing the contents of a chemical sterilant that come in two containers (small vial with the activator solution; container of the chemical). Keeping the two chemicals separate until use extends the shelf life of the chemicals.

Aeration: method by which ethylene oxide (EtO) is removed from EtO-sterilized items by warm air circulation in an enclosed cabinet specifically designed for this purpose.

Antimicrobial agent: any agent that kills or suppresses the growth of microorganisms.

Antiseptic: substance that prevents or arrests the growth or action of microorganisms by inhibiting their activity or by destroying them. The term is used especially for preparations applied topically to living tissue.

Asepsis: prevention of contact with microorganisms.

Autoclave: device that sterilizes instruments or other objects using steam under pressure. The length of time required for sterilization depends on temperature, vacuum, and pressure.

Bacterial count: method of estimating the number of bacteria per unit sample. The term also refers to the estimated number of bacteria per unit sample, usually expressed as number of colony-forming units.

Bactericide: agent that kills bacteria.

Bioburden: number and types of viable microorganisms with which an item is contaminated; also called bioload or microbial load.

Biofilm: accumulated mass of bacteria and extracellular material that is tightly adhered to a surface and cannot be easily removed.

Biologic indicator: device for monitoring the sterilization process. The device consists of a standardized, viable population of microorganisms (usually bacterial spores) known to be resistant to the sterilization process being monitored. Biologic indicators are intended to demonstrate whether conditions were adequate to achieve sterilization. A negative biologic indicator does not prove that all items in the load are sterile or that they were all exposed to adequate sterilization conditions.

Bleach: Household bleach (that includes 5.25% or 6.00%-6.15% sodium hypochlorite depending on manufacturer) is usually diluted in water at 1:10 or 1:100. Approximate dilutions are 1.5 cups of bleach in a gallon of water for a 1:10 dilution (~6,000 ppm) and 0.25 cup of bleach in a gallon of water for a 1:100 dilution (~600 ppm). Sodium hypochlorite products that make pesticidal claims, such as sanitization or disinfection, must be registered by EPA and be labeled with an EPA Registration Number.

Expected Chlorine Concentrations by Various Dulitions of Househod Bleach (5.25-6.15% sodium hypochlorite)

Dilution	Chlorine (ppm)
None	52,500-61,500
1:10	5,250-6,150
1:100	525-615
1:1000	53-62



Format Change [February 2017]: The format of this section was changed to improve readability and accessibility. The content is unchanged.

Bowie-Dick test: diagnostic test of a sterilizer's ability to remove air from the chamber of a prevacuum steam sterilizer. The air-removal or Bowie-Dick test is not a test for sterilization.

Ceiling limit: concentration of an airborne chemical contaminant that should not be exceeded during any part of the workday. If instantaneous monitoring is not feasible, the ceiling must be assessed as a 15-minute time-weighted average exposure.

Centigrade or **Celsius**: a temperature scale (0°C = freezing point of water; 100°C = boiling point of water at sea level). Equivalents mentioned in the guideline are as follows: 20°C = 68°F; 25°C = 77°F; 121°C = 250°F; 132°C = 270°F; 134°C = 273°F. For other temperatures the formula is: F° = (C° × 9/5) + 32 or C° = (F° - 32) × 5/9.

Central processing or **Central service department**: the department within a health-care facility that processes, issues, and controls professional supplies and equipment, both sterile and nonsterile, for some or all patient-care areas of the facility.

Challenge test pack: pack used in installation, qualification, and ongoing quality assurance testing of health-care facility sterilizers.

Chemical indicator: device for monitoring a sterilization process. The device is designed to respond with a characteristic chemical or physical change to one or more of the physical conditions within the sterilizing chamber. Chemical indicators are intended to detect potential sterilization failures that could result from incorrect packaging, incorrect loading of the sterilizer, or malfunctions of the sterilizer. The "pass" response of a chemical indicator does not prove the item accompanied by the indicator is necessarily sterile. The Association for the Advancement of Medical Instrumentation has defined five classes of chemical indicators: Class 1 (process indicator); Class 2 (Bowie-Dick test indicator); Class 3 (single-parameter indicator); Class 5 (integrating indicator).

Contact time: time a disinfectant is in direct contact with the surface or item to be disinfected. For surface disinfection, this period is framed by the application to the surface until complete drying has occurred.

Container system, rigid container: sterilization containment device designed to hold medical devices for sterilization, storage, transportation, and aseptic presentation of contents.

Contaminated: state of having actual or potential contact with microorganisms. As used in health care, the term generally refers to the presence of microorganisms that could produce disease or infection.

Control, positive: biologic indicator, from the same lot as a test biologic indicator, that is left unexposed to the sterilization cycle and then incubated to verify the viability of the test biologic indicator.

Cleaning: removal, usually with detergent and water or enzyme cleaner and water, of adherent visible soil, blood, protein substances, microorganisms and other debris from the surfaces, crevices, serrations, joints, and lumens of instruments, devices, and equipment by a manual or mechanical process that prepares the items for safe handling and/or further decontamination.

Culture: growth of microorganisms in or on a nutrient medium; to grow microorganisms in or on such a medium.

Culture medium: substance or preparation used to grow and cultivate microorganisms.

Cup: 8 fluid ounces.

Decontamination: according to OSHA, "the use of physical or chemical means to remove, inactivate, or destroy bloodborne pathogens on a surface or item to the point where they are no longer capable of transmitting infectious particles and the surface or item is rendered safe for handling, use, or disposal" [29 CFR 1910.1030]. In health-care facilities, the term generally refers to all pathogenic organisms.

Decontamination area: area of a health-care facility designated for collection, retention, and cleaning of soiled and/or contaminated items.

Detergent: cleaning agent that makes no antimicrobial claims on the label. They comprise a hydrophilic component and a lipohilic component and can be divided into four types: anionic, cationic, amphoteric, and non-ionic detergents.

Disinfectant: usually a chemical agent (but sometimes a physical agent) that destroys disease-causing pathogens or other harmful microorganisms but might not kill bacterial spores. It refers to substances applied to inanimate objects. EPA groups disinfectants by product label claims of "limited," "general," or "hospital" disinfection.

Disinfection: thermal or chemical destruction of pathogenic and other types of microorganisms. Disinfection is less lethal than sterilization because it destroys most recognized pathogenic microorganisms but not necessarily all microbial forms (e.g., bacterial spores).

D value: time or radiation dose required to inactivate 90% of a population of the test microorganism under stated exposure conditions.

Endoscope: an instrument that allows examination and treatment of the interior of the body canals and hollow organs.

Enzyme cleaner: a solution used before disinfecting instruments to improve removal of organic material (e.g., proteases to assist in removing protein).

EPA Registration Number or **EPA Reg. No.**: a hyphenated, two- or three-part number assigned by EPA to identify each germicidal product registered within the United States. The first number is the company identification number, the second is the specific product number, and the third (when present) is the company identification number for a supplemental registrant.

Exposure time: period in a sterilization process during which items are exposed to the sterilant at the specified sterilization parameters. For example, in a steam sterilization process, exposure time is the period during which items are exposed to saturated steam at the specified temperature.

Flash sterilization: process designed for the steam sterilization of unwrapped patient-care items for immediate use (or placed in a specially designed, covered, rigid container to allow for rapid penetration of steam).

Fungicide: agent that destroys fungi (including yeasts) and/or fungal spores pathogenic to humans or other animals in the inanimate environment.

General disinfectant: EPA-registered disinfectant labeled for use against both gram-negative and grampositive bacteria. Efficacy is demonstrated against both *Salmonella choleraesuis* and *Staphylococcus aureus*. Also called *broad-spectrum disinfectant*.

Germicide: agent that destroys microorganisms, especially pathogenic organisms.

Germicidal detergent: detergent that also is EPA-registered as a disinfectant.

High-level disinfectant: agent capable of killing bacterial spores when used in sufficient concentration under suitable conditions. It therefore is expected to kill all other microorganisms.

Hospital disinfectant: disinfectant registered for use in hospitals, clinics, dental offices, and any other medical-related facility. Efficacy is demonstrated against *Salmonella choleraesuis*, *Staphylococcus aureus*, and *Pseudomonas aeruginosa*. EPA has registered approximately 1,200 hospital disinfectants.

Huck towel: all-cotton surgical towel with a honey-comb weave; both warp and fill yarns are tightly twisted. Huck towels can be used to prepare biologic indicator challenge test packs.

Implantable device: according to FDA, "device that is placed into a surgically or naturally formed cavity of the human body if it is intended to remain there for a period of 30 days or more" [21 CFR 812.3(d)].

Inanimate surface: nonliving surface (e.g., floors, walls, furniture).

Incubator: apparatus for maintaining a constant and suitable temperature for the growth and cultivation of microorganisms.

Infectious microorganisms: microorganisms capable of producing disease in appropriate hosts.

Inorganic and organic load: naturally occurring or artificially placed inorganic (e.g., metal salts) or organic (e.g., proteins) contaminants on a medical device before exposure to a microbicidal process.

Intermediate-level disinfectant: agent that destroys all vegetative bacteria, including tubercle bacilli, lipid and some nonlipid viruses, and fungi, but not bacterial spores.

Limited disinfectant: disinfectant registered for use against a specific major group of organisms (gramnegative or gram-positive bacteria). Efficacy has been demonstrated in laboratory tests against either *Salmonella choleraesuis* or *Staphylococcus aureus* bacteria.

Lipid virus: virus surrounded by an envelope of lipoprotein in addition to the usual core of nucleic acid surrounded by a coat of protein. This type of virus (e.g., HIV) is generally easily inactivated by many types of disinfectants. Also called *enveloped* or *lipophilic virus*.

Low-level disinfectant: agent that destroys all vegetative bacteria (except tubercle bacilli), lipid viruses, some nonlipid viruses, and some fungi, but not bacterial spores.

Mechanical indicator: devices that monitor the sterilization process (e.g., graphs, gauges, printouts).

Medical device: instrument, apparatus, material, or other article, whether used alone or in combination, including software necessary for its application, intended by the manufacturer to be used for human beings for

- diagnosis, prevention, monitoring treatment, or alleviation of disease;
- diagnosis, monitoring, treatment, or alleviation of or compensation for an injury or handicap;
- investigation, replacement, or modification of the anatomy or of a physiologic process; or
- control of conception

and that does not achieve its primary intended action in or on the human body by pharmacologic, immunologic, or metabolic means but might be assisted in its function by such means.

Microbicide: any substance or mixture of substances that effectively kills microorganisms.

Microorganisms: animals or plants of microscopic size. As used in health care, generally refers to bacteria, fungi, viruses, and bacterial spores.

Minimum effective concentration (MEC): the minimum concentration of a liquid chemical germicide needed to achieve the claimed microbicidal activity as determined by dose-response testing. Sometimes used interchangeably with *minimum recommended concentration*.

Muslin: loosely woven (by convention, 140 threads per square inch), 100% cotton cloth. Formerly used as a wrap for sterile packs or a surgical drape. Fabric wraps used currently consist of a cotton-polyester blend.

Mycobacteria: bacteria with a thick, waxy coat that makes them more resistant to chemical germicides than other types of vegetative bacteria.

Nonlipid viruses: generally considered more resistant to inactivation than lipid viruses. Also called nonenveloped or hydrophilic viruses.

One-step disinfection process: simultaneous cleaning and disinfection of a noncritical surface or item.

Pasteurization: process developed by Louis Pasteur of heating milk, wine, or other liquids to 65–77°C (or the equivalent) for approximately 30 minutes to kill or markedly reduce the number of pathogenic and spoilage organisms other than bacterial spores.

Parametric release: declaration that a product is sterile on the basis of physical and/or chemical process data rather than on sample testing or biologic indicator results.

Penicylinder: carriers inoculated with the test bacteria for in vitro tests of germicides. Can be constructed of stainless steel, porcelain, glass, or other materials and are approximately 8 x 10 mm in diameter.

Permissible exposure limit (PEL): time-weighted average maximum concentration of an air contaminant to which a worker can be exposed, according to OSHA standards. Usually calculated over 8 hours, with exposure considered over a 40-hour work week.

Personal protective equipment (PPE): specialized clothing or equipment worn by an employee for protection against a hazard. General work clothes (e.g., uniforms, pants, shirts) not intended to function as protection against a hazard are not considered to be PPE.

Parts per million (ppm): common measurement for concentrations by volume of trace contaminant gases in the air (or chemicals in a liquid); 1 volume of contaminated gas per 1 million volumes of contaminated air or 1¢ in \$10,000 both equal 1 ppm. Parts per million = $\mu g/mL$ or mg/L.

Prions: transmissible pathogenic agents that cause a variety of neurodegenerative diseases of humans and animals, including sheep and goats, bovine spongiform encephalopathy in cattle, and Creutzfeldt-Jakob disease in humans. They are unlike any other infectious pathogens because they are composed of an abnormal conformational isoform of a normal cellular protein, the prion protein (PrP). Prions are extremely resistant to inactivation by sterilization processes and disinfecting agents.

Process challenge device (PCD): item designed to simulate product to be sterilized and to constitute a defined challenge to the sterilization process and used to assess the effective performance of the process. A PCD is a challenge test pack or test tray that contains a biologic indicator, a Class 5 integrating indicator, or an enzyme-only indicator.

QUAT: abbreviation for *quaternary ammonium compound*, a surface-active, water-soluble disinfecting substance that has four carbon atoms linked to a nitrogen atom through covalent bonds.

Recommended exposure limit (REL): occupational exposure limit recommended by NIOSH as being protective of worker health and safety over a working lifetime. Frequently expressed as a 40-hour time-weighted-average exposure for up to 10 hours per day during a 40-work week.

Reprocess: method to ensure proper disinfection or sterilization; can include: cleaning, inspection, wrapping, sterilizing, and storing.

Sanitizer: agent that reduces the number of bacterial contaminants to safe levels as judged by public health requirements. Commonly used with substances applied to inanimate objects. According to the protocol for the official sanitizer test, a sanitizer is a chemical that kills 99.999% of the specific test bacteria in 30 seconds under the conditions of the test.

Shelf life: length of time an undiluted or use dilution of a product can remain active and effective. Also refers to the length of time a sterilized product (e.g., sterile instrument set) is expected to remain sterile.

Spaulding classification: strategy for reprocessing contaminated medical devices. The system classifies a medical device as critical, semicritical, or noncritical on the basis of risk to patient safety from contamination on a device. The system also established three levels of germicidal activity (sterilization, high-level disinfection, and low-level disinfection) for strategies with the three classes of medical devices (critical, semicritical, and noncritical).

Spore: relatively water-poor round or elliptical resting cell consisting of condensed cytoplasm and nucleus surrounded by an impervious cell wall or coat. Spores are relatively resistant to disinfectant and sterilant activity and drying conditions (specifically in the genera *Bacillus* and *Clostridium*).

Spore strip: paper strip impregnated with a known population of spores that meets the definition of biological indicators.

Steam quality: steam characteristic reflecting the dryness fraction (weight of dry steam in a mixture of dry saturated steam and entrained water) and the level of noncondensable gas (air or other gas that will not condense under the conditions of temperature and pressure used during the sterilization process). The dryness fraction (i.e., the proportion of completely dry steam in the steam being considered) should not fall below 97%.

Steam sterilization: sterilization process that uses saturated steam under pressure for a specified exposure time and at a specified temperature, as the sterilizing agent.

Steam sterilization, dynamic air removal type: one of two types of sterilization cycles in which air is removed from the chamber and the load by a series of pressure and vacuum excursions (prevacuum cycle) or by a series of steam flushes and pressure pulses above atmospheric pressure (steam-flush-pressure-pulse cycle).

Sterile or **Sterility**: state of being free from all living microorganisms. In practice, usually described as a probability function, e.g., as the probability of a microorganism surviving sterilization being one in one million.

Sterility assurance level (SAL): probability of a viable microorganism being present on a product unit after sterilization. Usually expressed as 10^{-6} ; a SAL of 10^{-6} means $\leq 1/1$ million chance that a single viable microorganism is present on a sterilized item. A SAL of 10^{-6} generally is accepted as appropriate for items intended to contact compromised tissue (i.e., tissue that has lost the integrity of the natural body barriers). The sterilizer manufacturer is responsible for ensuring the sterilizer can achieve the desired SAL. The user is responsible for monitoring the performance of the sterilizer to ensure it is operating in conformance to the manufacturer's recommendations.

Sterilization: validated process used to render a product free of all forms of viable microorganisms. In a sterilization process, the presence of microorganisms on any individual item can be expressed in terms of probability. Although this probability can be reduced to a very low number, it can never be reduced to zero.

Sterilization area: area of a health-care facility designed to house sterilization equipment, such as steam ethylene oxide, hydrogen peroxide gas plasma, or ozone sterilizers.

Sterilizer: apparatus used to sterilize medical devices, equipment, or supplies by direct exposure to the sterilizing agent.

Sterilizer, gravity-displacement type: type of steam sterilizer in which incoming steam displaces residual air through a port or drain in or near the bottom (usually) of the sterilizer chamber. Typical operating temperatures are 121–123°C (250–254°F) and 132–135°C (270–275°F).

Sterilizer, prevacuum type: type of steam sterilizer that depends on one or more pressure and vacuum excursions at the beginning of the cycle to remove air. This method of operation results in shorter cycle times for wrapped items because of the rapid removal of air from the chamber and the load by the vacuum system and because of the usually higher operating temperature (132–135°C [270–275°F]; 141–144°C [285–291°F]). This type of sterilizer generally provides for shorter exposure time and accelerated drying of fabric loads by pulling a further vacuum at the end of the sterilizing cycle.

Sterilizer, steam-flush pressure-pulse type: type of sterilizer in which a repeated sequence consisting of a steam flush and a pressure pulse removes air from the sterilizing chamber and processed materials using steam at above atmospheric pressure (no vacuum is required). Like a prevacuum sterilizer, a steam-flush pressure-pulse sterilizer rapidly removes air from the sterilizing chamber and wrapped items; however, the system is not susceptible to air leaks because air is removed with the sterilizing chamber pressure at above atmospheric pressure. Typical operating temperatures are 121–123°C (250–254°F), 132–135°C (270–275°F), and 141–144°C (285–291°F).

Surfactant: agent that reduces the surface tension of water or the tension at the interface between water and another liquid; a wetting agent found in many sterilants and disinfectants.

Tabletop steam sterilizer: a compact gravity-displacement steam sterilizer that has a chamber volume of not more than 2 cubic feet and that generates its own steam when distilled or deionized water is added.

Time-weighted average (TWA): an average of all the concentrations of a chemical to which a worker has been exposed during a specific sampling time, reported as an average over the sampling time. For example, the permissible exposure limit for ethylene oxide is 1 ppm as an 8-hour TWA. Exposures above the ppm limit are permitted if they are compensated for by equal or longer exposures below the limit during the 8-hour workday as long as they do not exceed the ceiling limit; short-term exposure limit; or, in the case of ethylene oxide, excursion limit of 5 ppm averaged over a 15-minute sampling period.

Tuberculocide: an EPA-classified hospital disinfectant that also kills *Mycobacterium tuberculosis* (tubercle bacilli). EPA has registered approximately 200 tuberculocides. Such agents also are called *mycobactericides*.

Use-life: the length of time a diluted product can remain active and effective. The stability of the chemical and the storage conditions (e.g., temperature and presence of air, light, organic matter, or metals) determine the use-life of antimicrobial products.

Vegetative bacteria: bacteria that are devoid of spores and usually can be readily inactivated by many types of germicides.

Virucide: an agent that kills viruses to make them noninfective.

Adapted from Association for the Advancement of Medical Instrumentation; ^{811-814, 819} Association of periOperating Registered Nurses (AORN), ⁸¹⁵ American Hospital Association, ³¹⁹ and Block. ^{16, 1034}

Tables and Figure

Table 1. Methods of sterilization and disinfection

Sterilization

• Critical items (will enter tissue or vascular system or blood will flow through them)

Disinfection

- High-level (semicritical items; [except dental] will come in contact with mucous membrane or nonintact skin)
- Intermediate-level (some semicritical items1 and noncritical items)
- Low-level (noncritical items; will come in contact with intact skin)

Format Change [February 2017]: The format of this section was changed to improve readability and accessibility. The content is unchanged.

Table 1A	Smooth	Hard	Surface ^{1,4}
	SHIUUUII,	riaiu	Juliace

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Intermediate Disinfection (exposure time ≥1 m) ⁹	Low-level Disinfection (exposure time ≥1 m) ⁹
Heat sterilization, including steam or hot air (see manufacturer's recommendations, steam sterilization processing time from 3-30 minutes)	Manufacturer's recommend- ations	Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	Ethyl or isopropyl alcohol (70–90%)	Ethyl or isopropyl alcohol (70– 90%)
Ethylene oxide gas (see manufacturer's recommendations, generally 1–6 hours processing time plus aeration time of 8–12 hours at 50–60°C)	Manufacturer's recommend- ations	Ortho-phthalaldehyde (OPA) 0.55%	Sodium hypochlorite (5.25–6.15% household bleach diluted 1:500 provides >100 ppm available chlorine) ⁵	Sodium hypochlorite (5.25–6.15% household bleach diluted 1:500 provides >100 ppm available chlorine)
Hydrogen peroxide gas plasma (see manufacturer's recommendations for internal diameter and length restrictions, processing time between 45–72 minutes).	Manufacturer's recommend- ations	Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	Phenolic germicidal detergent solution (follow product label for use- dilution)	Phenolic germicidal detergent solution (follow product label for use-dilution)
Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	10 h at 20– 25°C	Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	lodophor germicidal detergent solution (follow product label for use- dilution)	lodophor germicidal detergent solution (follow product label for use-dilution)

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Intermediate Disinfection (exposure time ≥1 m) ⁹	Low-level Disinfection (exposure time ≥1 m) ⁹
Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	6 h	Wet pasteurization at 70°C for 30 minutes with detergent cleaning ⁶	n/a	Quaternary ammonium germicidal detergent solution (follow product label for use-dilution)
Peracetic acid, concentration variable but 0.2% or greater is sporicidal. Peracetic acid immersion system operates at 50–56°C	12 m at 50–56°C	Hypochlorite, single use chlorine generated on-site by electrolyzing saline containing >650–675 active free chlorine; (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	3–8 h	n/a	n/a	n/a

Table 1B. Rubber tubing and catheters^{3,4}

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Intermediate Disinfection (exposure time ≥1 m) ⁹	Low-level Disinfection (exposure time ≥1 m) ⁹
Heat sterilization, including steam or hot air (see manufacturer's recommendations, steam sterilization processing time from 3-30 minutes)	Manufacturer's recommend- ations	Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	n/a	n/a
Ethylene oxide gas (see manufacturer's recommendations, generally 1–6 hours processing time plus aeration time of 8–12 hours at 50–60°C)	Manufacturer's recommend- ations	Ortho-phthalaldehyde (OPA) 0.55%	n/a	n/a
Hydrogen peroxide gas plasma (see manufacturer's recommendations for internal diameter and length restrictions, processing time between 45–72 minutes).	Manufacturer's recommend- ations	Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	n/a	n/a

	Sterilization Exposure	High-level Disinfection (exposure time 12-30 m at	Intermediate Disinfection (exposure time	Low-level Disinfection (exposure time
Sterilization Procedure	Time	≥20°C) ^{2,3}	≥1 m) ⁹	≥1 m) ⁹
Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	10 h at 20–25°C	Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	6 h	Wet pasteurization at 70°C for 30 minutes with detergent cleaning ⁶	n/a	n/a
Peracetic acid, concentration variable but 0.2% or greater is sporicidal. Peracetic acid immersion system operates at 50–56°C ^G	12 m at 50–56°C	Hypochlorite, single use chlorine generated on-site by electrolyzing saline containing >650–675 active free chlorine; (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	3–8 h	n/a	n/a	n/a

Table 1C. Polyethylene tubing and catheters^{3,4,7}

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Intermediate Disinfection (exposure time ≥1 m) ⁹	Low-level Disinfection (exposure time ≥1 m) ⁹
Heat sterilization, including steam or hot air (see manufacturer's recommendations, steam sterilization processing time from 3–30 minutes)	recommend- ations	Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	n/a	n/a
Ethylene oxide gas (see manufacturer's recommendations, generally 1–6 hours processing time plus aeration time of 8-12 hours at 50–60°C)	Manufacturer's recommend- ations	Ortho-phthalaldehyde (OPA) 0.55%	n/a	n/a
Hydrogen peroxide gas plasma (see manufacturer's recommendations for internal diameter and length restrictions, processing time between 45–72 minutes).		Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	n/a	n/a

	Sterilization Exposure	High-level Disinfection (exposure time 12-30 m at	Intermediate Disinfection (exposure time	Low-level Disinfection (exposure time
Sterilization Procedure	Time	≥20°C) ^{2,3}	≥1 m) ⁹	≥1 m) ⁹
Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	10 h at 20– 25℃	Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	6 h	Wet pasteurization at 70°C for 30 minutes with detergent cleaning ⁶	n/a	n/a
Peracetic acid, concentration variable but 0.2% or greater is sporicidal. Peracetic acid immersion system operates at 50–56°C	12 m at 50– 56°C	Hypochlorite, single use chlorine generated on-site by electrolyzing saline containing >650–675 active free chlorine; (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	3–8 h	n/a	n/a	n/a

Table 1D. Lensed instruments⁴

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Intermediate Disinfection (exposure time ≥1 m) ⁹	Low-level Disinfection (exposure time ≥1 m) ⁹
Heat sterilization, including steam or hot air (see manufacturer's recommendations, steam sterilization processing time from 3–30 minutes)	Manufacturer's recommend- ations	Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	n/a	n/a
Ethylene oxide gas (see manufacturer's recommendations, generally 1–6 hours processing time plus aeration time of 8–12 hours at 50–60°C)	recommend- ations	Ortho-phthalaldehyde (OPA) 0.55%	n/a	n/a
Hydrogen peroxide gas plasma (see manufacturer's recommendations for internal diameter and length restrictions, processing time between 45–72 minutes).		Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	n/a	n/a

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Intermediate Disinfection (exposure time ≥1 m) ⁹	Low-level Disinfection (exposure time ≥1 m) ⁹
Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	10 h at 20– 25°C	Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	6 h	Hypochlorite, single use chlorine generated on-site by electrolyzing saline containing >650–675 active free chlorine; (will corrode metal instruments)	n/a	n/a
Peracetic acid, concentration variable but 0.2% or greater is sporicidal. Peracetic acid immersion system operates at 50–56°C	12 m at 50–56°C	n/a	n/a	n/a
Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	3–8 h	n/a	n/a	n/a

Table 1E. Thermometers (Oral and Rectal)⁸

Sterilization Procedure	Sterilization Exposure Time	High-level Disinfection (exposure time 12-30 m at ≥20°C) ^{2,3}	Low-level Disinfection (exposure time ≥1 m) ⁹
n/a	n/a	n/a	Ethyl or isopropyl alcohol (70-90%) ⁸

Table 1F. Hinged Instruments⁴

Sterilization Presedure	Sterilization Exposure Time	High-level Disinfection (expective time 12.20 m et ≥20°C)23	Intermediate Disinfection (exposure time	Low-level Disinfection (exposure time
Sterilization Procedure Heat sterilization, including steam or hot air (see manufacturer's recommendations, steam sterilization processing time from 3–30 minutes)		(exposure time 12-30 m at \geq 20°C) ^{2,3} Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	<u>≥1 m) ⁹</u> n/a	<u>≥1 m) ⁹</u> n/a
Ethylene oxide gas (see manufacturer's recommendations, generally 1–6 hours processing time plus aeration time of 8–12 hours at 50–60°C)	Manufacturer's recommend- ations		n/a	n/a
Hydrogen peroxide gas plasma (see manufacturer's recommendations for internal diameter and length restrictions, processing time between 45–72 minutes).	recommend- ations	Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	n/a	n/a
Glutaraldehyde-based formulations (>2% glutaraldehyde, caution should be exercised with all glutaraldehyde formulations when further in-use dilution is anticipated); glutaraldehyde (1.12%) and 1.93% phenol/phenate. One glutaraldehyde-based product has a high-level disinfection claim of 5 minutes at 35°C.	10 h at 20– 25°C	Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide 7.5% (will corrode copper, zinc, and brass)	6 h	Wet pasteurization at 70°C for 30 minutes with detergent cleaning ⁶	n/a	n/a
Peracetic acid, concentration variable but 0.2% or greater is sporicidal. Peracetic acid immersion system operates at 50–56°C	12 m at 50- 56°C	Hypochlorite, single use chlorine generated on-site by electrolyzing saline containing >650–675 active free chlorine; (will corrode metal instruments)	n/a	n/a
Hydrogen peroxide (7.35%) and 0.23% peracetic acid; hydrogen peroxide 1% and peracetic acid 0.08% (will corrode metal instruments)	3-8 h	n/a	n/a	n/a

Modified from Rutala and Simmons. ^{15, 17, 18, 421} The selection and use of disinfectants in the healthcare field is dynamic, and products may become available that are not in existence when this guideline was written. As newer disinfectants become available, persons or committees responsible for selecting disinfectants and sterilization processes should be guided by products cleared by the FDA and the EPA as well as information in the scientific literature.

- ¹ See text for discussion of hydrotherapy.
- ² The longer the exposure to a disinfectant, the more likely it is that all microorganisms will be eliminated. Follow the

FDA-cleared high-level disinfection claim. Ten-minute exposure is not adequate to disinfect many objects, especially those that are difficult to clean because they have narrow channels or other areas that can harbor organic material and bacteria. Twenty-minute exposure at 20°C is the minimum time needed to reliably kill *M. tuberculosis* and nontuberculous mycobacteria with a 2% glutaraldehyde. Some high-level disinfectants have a reduced exposure time (e.g., ortho-phthalaldehyde at 12 minutes at 20°C) because of their rapid activity against mycobacteria or reduced exposure time due to increased mycobactericidal activity at elevated temperature (e.g., 2.5% glutaraldehyde at 5 minutes at 35°C, 0.55% OPA at 5 min at 25°C in automated endoscope reprocessor).

- ³ Tubing must be completely filled for high-level disinfection and liquid chemical sterilization; care must be taken to avoid entrapment of air bubbles during immersion.
- ⁴ Material compatibility should be investigated when appropriate.
- ⁵ A concentration of 1000 ppm available chlorine should be considered where cultures or concentrated preparations of microorganisms have spilled (5.25% to 6.15% household bleach diluted 1:50 provides > 1000 ppm available chlorine). This solution may corrode some surfaces.
- ⁶ Pasteurization (washer-disinfector) of respiratory therapy or anesthesia equipment is a recognized alternative to highlevel disinfection. Some data challenge the efficacy of some pasteurization units.
- ⁷ Thermostability should be investigated when appropriate.
- ⁸ Do not mix rectal and oral thermometers at any stage of handling or processing.
- ⁹ By law, all applicable label instructions on EPA-registered products must be followed. If the user selects exposure conditions that differ from those on the EPA-registered products label, the user assumes liability from any injuries resulting from off-label use and is potentially subject to enforcement action under FIFRA.

Table 2. Properties of an ideal disinfectant

- Broad spectrum: should have a wide antimicrobial spectrum
- Fast acting: should produce a rapid kill
- Not affected by environmental factors: should be active in the presence of organic matter (e.g., blood, sputum, feces) and compatible with soaps, detergents, and other chemicals encountered in use
- Nontoxic: should not be harmful to the user or patient
- Surface compatibility: should not corrode instruments and metallic surfaces and should not cause the deterioration of cloth, rubber, plastics, and other materials
- Residual effect on treated surfaces: should leave an antimicrobial film on the treated surface
- Easy to use with clear label directions
- Odorless: should have a pleasant odor or no odor to facilitate its routine use
- Economical: should not be prohibitively high in cost
- Solubility: should be soluble in water
- Stability: should be stable in concentrate and use-dilution
- Cleaner: should have good cleaning properties
- Environmentally friendly: should not damage the environment on disposal

Modified from Molinari¹⁰³⁵.

Table 3. Epidemiologic evidence associated with the use of surface disinfectants or detergents on noncritical environmental surfaces

Justification for Use of Disinfectants for Noncritical Environmental Surfaces

- Surfaces may contribute to transmission of epidemiologically important microbes (e.g., vancomycin-resistant Enterococci, methicillin-resistant S. aureus, viruses)
- Disinfectants are needed for surfaces contaminated by blood and other potentially infective material
- Disinfectants are more effective than detergents in reducing microbial load on floors
- Detergents become contaminated and result in seeding the patient's environment with bacteria
- Disinfection of noncritical equipment and surfaces is recommended for patients on isolation precautions by the Centers for Disease Control and Prevention.
- Advantage of using a single product for decontamination of noncritical surfaces, both floors and equipment
- Some newer disinfectants have persistent antimicrobial activity

Justification for Using a Detergent on Noncritical Environmental Surfaces

- Noncritical surfaces contribute minimally to endemic healthcare-associated infections
- No difference in healthcare-associated infection rates when floors are cleaned with detergent versus disinfectant
- No environmental impact (aquatic or terrestrial) issues with disposal
- No occupational health exposure issues
- Lower costs
- Use of antiseptics/disinfectants selects for antibiotic-resistant bacteria (?)
- More aesthetically pleasing floor

Modified from Rutala³⁷⁸.

Figure 1. Decreasing order of resistance of microorganisms to disinfection and sterilization (with the disinfection levels indicated)

Format Change [February 2017]: The format of this section was changed to improve readability and accessibility. The content is unchanged.

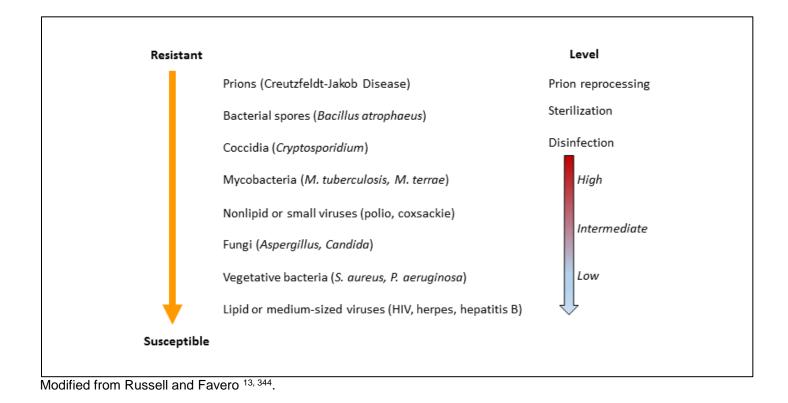


Table 4. Comparison of the characteristics of selected chemicals used as highlevel disinfectants or chemical Sterilants



Format Change [February 2017]: The format of this section was changed to improve readability and accessibility. The content is unchanged.

Chemical Charasterics	Hydrogen Peroxide (7.5%)	Peracetic Acid (0.2%)	Glutaraldehyde (≥2.0%)	OPA (0.55%)	Hydrogen Peroxide / Peracetic Acid (7.35%/0.23%)
High-level disinfectant	30 minutes @	Not Applicable	20-90 minutes @	12 minutes @ 20°C,	15 minutes @
claim	20°C		20°-25°C	5 minutes @ 25°C in AER	20°C
Sterilization Claim	6 hours @ 20°C	12 minutes @ 50-56°C	10 hours @ 20º-25°C	None	3 h @ 20°C
Activation	No	No	Yes (alkaline glutaraldehyde)	No	No
Reuse life (number of days a product can be reused as determined by re-use protocol)	21 days	Single use	14-30 days	14 days	14 days
Shelf life stability (time a product can remain in storage (unused))	2 years	6 months	2 years	2 years	2 years
Disposal Restrictions	None	None	Local (no U.S. EPA regulations exist but some states and local authorities have disposal restrictions)	Local (no U.S. EPA regulations exist but some states and local authorities have disposal restrictions)	None
Materials Compatibility	Good	Good	Excellent	Excellent	No data
Monitor MEC of solution	Yes (6%)	No	Yes (1.5% or higher)	Yes (0.3% OPA)	No
Safety	Serious eye irritant (safety glasses)	Serious eye and skin irritant (concentrated solution)	Respiratory irritant	Eye irritant, stains skin	Eye irritant
Processing	Manual or automated	Automated	Manual or automated	Manual or automated	Manual
Organic material resistance	Yes	Yes	Yes	Yes	Yes
OSHA exposure limit	1 ppm TWA	None	None (The ceiling limit recommended by the American Conference of Governmental Industrial Hygienists is 0.05 ppm.)	None	Hydrogen Peroxide - 1 ppm (time- weighted average for a conventional 8- hour workday.)
Cost profile (per	+ (manual)	+++++	+ (manual)	++ (manual)	++ (manual)
cycle) ¹ Modified from Butele ⁶⁹	++ (automated)	(automated)	++ (automated)		

Modified from Rutala 69.

Abbreviations and Footnotes:

ortho-phthalaldehyde (FDA cleared as a high-level disinfectant, included for comparison to other chemical agents OPA used for high-level disinfection)

Automated Endoscope Reprocessor AER

minimum effective concentration is the lowest concentration of active ingredients at which the product is still effective MEC per cycle cost profile considers cost of the processing solution (suggested list price to healthcare facilities in August 1 2001) and assumes maximum use life (e.g., 21 days for hydrogen peroxide, 14 days for glutaraldehyde), 5 reprocessing cycles per day, 1-gallon basin for manual processing, and 4-gallon tank for automated processing.

least expensive;

+++ most expensive

Sterilization Method	Advantages	Disadvantages
Peracetic Acid/Hydrogen Peroxide	 No activation required Odor or irritation not significant 	 Materials compatibility concerns (lead, brass, copper, zinc) both cosmetic and functional Limited clinical experience Potential for eye and skin damage
Glutaraldehyde	 Numerous use studies published Relatively inexpensive Excellent materials compatibility 	 Respiratory irritation from glutaraldehyde vapor Pungent and irritating odor Relatively slow mycobactericidal activity Coagulates blood and fixes tissue to surfaces Allergic contact dermatitis Glutaraldehyde vapor monitoring recommended
Hydrogen Peroxide	 No activation required May enhance removal of organic matter and organisms No disposal issues No odor or irritation issues Does not coagulate blood or fix tissues to surfaces Inactivates <i>Cryptosporidium</i> Use studies published 	 Material compatibility concerns (brass, zinc, copper, and nickel/silver plating) both cosmetic and functional Serious eye damage with contact
Ortho- phthalaldehyde	 Fast acting high-level disinfectant No activation required Odor not significant Excellent materials compatibility claimed Does not coagulate blood or fix tissues to surfaces claimed 	 Stains skin, mucous membranes, clothing, and environmental surfaces Repeated exposure may result in hypersensitivity in some patients with bladder cancer More expensive than glutaraldehyde Eye irritation with contact Slow sporicidal activity
Peracetic Acid	 Rapid sterilization cycle time (30-45 minutes) Low temperature (50-55°C) liquid immersion sterilization Environmental friendly by-products (acetic acid, O₂, H₂0) Fully automated Single-use system eliminates need for concentration testing Standardized cycle May enhance removal of organic material and endotoxin No adverse health effects to operators under normal operating conditions Compatible with many materials and instruments Does not coagulate blood or fix tissues to surfaces Sterilant flows through scope facilitating salt, protein, and microbe removal Rapidly sporicidal Provides procedure standardization (constant dilution, perfusion of channel, temperatures, exposure) 	 Potential material incompatibility (e.g., aluminum anodized coating becomes dull) Used for immersible instruments only Biological indicator may not be suitable for routine monitoring One scope or a small number of instruments can be processed in a cycle More expensive (endoscope repairs, operating costs, purchase costs) than high-level disinfection Serious eye and skin damage (concentrated solution) with contact Point-of-use system, no sterile storage

Table 5. Summary of advantages and disadvantages of chemical agents used as chemical sterilants¹ or as high-level disinfectants

Modified from Rutala⁶⁹.

¹All products effective in presence of organic soil, relatively easy to use, and have a broad spectrum of antimicrobial activity (bacteria, fungi, viruses, bacterial spores, and mycobacteria). The above characteristics are documented in the literature; contact the manufacturer of the instrument and sterilant for additional information. All products listed above are FDA-cleared as chemical sterilants except OPA, which is an FDA-cleared high-level disinfectant.

Table 6. Summary of advantages and disadvantages of commonly used sterilization technologies

Sterilization Method	Advantages	Disadvantages
Steam	 Nontoxic to patient, staff, environment Cycle easy to control and monitor Rapidly microbicidal Least affected by organic/inorganic soils among sterilization processes listed Rapid cycle time Penetrates medical packing, device lumens 	 Deleterious for heat-sensitive instruments Microsurgical instruments damaged by repeated exposure May leave instruments wet, causing them to rust Caution: Potential for burns
Hydrogen Peroxide Gas Plasma	 Safe for the environment Leaves no toxic residuals Cycle time is 28-75 minutes (varies with model type) and no aeration necessary Used for heat- and moisture-sensitive items since process temperature <50°C Simple to operate, install (208 V outlet), and monitor Compatible with most medical devices Only requires electrical outlet 	 Cellulose (paper), linens and liquids cannot be processed Sterilization chamber size from 1.8-9.4 ft³ total volume (varies with model type) Some endoscopes or medical devices with long or narrow lumens cannot be processed at this time in the United States (see manufacturer's recommendations for internal diameter and length restrictions) Requires synthetic packaging (polypropylene wraps, polyolefin pouches) and special container tray Caution: Hydrogen peroxide may be toxic at levels greater than 1 ppmTWA
100% Ethylene Oxide (ETO)	 Penetrates packaging materials, device lumens Single-dose cartridge and negative- pressure chamber minimizes the potential for gas leak and ETO exposure Simple to operate and monitor Compatible with most medical materials 	 Requires aeration time to remove ETO residue Sterilization chamber size from 4.0-7.9 ft³ total volume (varies with model type) ETO emission regulated by states but catalytic cell removes 99.9% of ETO and converts it to CO₂ and H₂O ETO cartridges should be stored in flammable liquid storage cabinet Lengthy cycle/aeration time Caution: ETO is toxic, a carcinogen, and flammable
ETO Mixtures 8.6% ETO/91.4% HCFC 10% ETO/90% HCFC 8.5% ETO/91.5% CO ₂	 Penetrates medical packaging and many plastics Compatible with most medical materials Cycle easy to control and monitor 	 Some states (e.g., CA, NY, MI) require ETO emission reduction of 90-99.9% Chlorofluorocarbon (CFC) (inert gas that eliminates explosion hazard) banned in 1995 Potential hazards to staff and patients Lengthy cycle/aeration time Caution: ETO is toxic, a carcinogen, and flammable
Peracetic Acid	 Rapid cycle time (30-45 minutes) Low temperature (50-55°C) liquid immersion sterilization Environmental friendly by-products Sterilant flows through endoscope which facilitates salt, protein and microbe removal 	 Point-of-use system, no sterile storage Biological indicator may not be suitable for routine monitoring Used for immersible instruments only Some material incompatibility (e.g., aluminum anodized coating becomes dull) One scope or a small number of instruments processed in a cycle Caution: Potential for serious eye and skin damage (concentrated solution) with contact

Modified from Rutala. 825

Abbreviations:

HCFC=hydrochlorofluorocarbon.

Table 7. Minimum cycle times for steam sterilization cycles

Type of Sterilizer	ltem	Exposure time at 250°F (121°C)	Exposure time at 270°F (132°C)	Drying time
Gravity displacement	Wrapped instruments	30 min	15 min	15-30 min
Gravity displacement	Textile packs	30 min	25 min	15 min
Gravity displacement	Wrapped utensils	30 min	15 min	15-30 min
Dynamic-air-removal (e.g., prevacuum)	Wrapped instruments	Not Applicable	4 min	20-30 min
Dynamic-air-removal (e.g., prevacuum)	Textile packs	Not Applicable	4 min	5-20 min
Dynamic-air-removal (e.g., prevacuum)	Wrapped utensils	Not Applicable	4 min	20 min

Modified from Association for the Advancement of Medical Instrumentation. 813, 819

Type of Sterilizer	Load configuration	Temperature	Time
Gravity displacement	Nonporous items only (i.e., routine metal instruments, no lumens)	132°C (270°F)	3 minutes
Gravity displacement	Nonporous and porous items (e.g., rubber or plastic items, items with lumens) sterilized together		10 minutes
Dynamic-air-removal (e.g., prevacuum)	Nonporous items only (i.e., routine metal instruments, no lumens)	132°C (270°F)	3 minutes
Dynamic-air-removal (e.g., prevacuum)	Nonporous and porous items (e.g., rubber or plastic items, items with lumens) sterilized together	132°C (270°F)	4 minutes
Steam-flush pressure-pulse	Nonporous or mixed nonporous/porous items	132°C (270°F) Manufacturers' instruction	4 minutes

Table 8. Examples of flash steam sterilization parameters

Modified from Association for the Advancement of Medical Instrumentation. 812, 819

Table 9. Characteristics of an ideal low-temperature sterilization process

- High efficacy: the agent should be virucidal, bactericidal, tuberculocidal, fungicidal and sporicidal
- Rapid activity: ability to quickly achieve sterilization
- Strong penetrability: ability to penetrate common medical-device packaging materials and penetrate into the interior of device lumens
- Material compatibility: produces only negligible changes in the appearance or the function of processed items and packaging materials even after repeated cycling
- Nontoxic: presents no toxic health risk to the operator or the patient and poses no hazard to the environment
- Organic material resistance: withstands reasonable organic material challenge without loss of efficacy
- Adaptability: suitable for large or small (point of use) installations
- Monitoring capability: monitored easily and accurately with physical, chemical, and biological process monitors
- **Cost effectiveness:** reasonable cost for installation and for routine operation

Modified from Schneider. 851

Factors	Effect
Cleaning ¹	Failure to adequately clean instrument results in higher bioburden, protein load, and salt
	concentration. These will decrease sterilization efficacy.
Bioburden ¹	The natural bioburden of used surgical devices is 10 ⁰ to 10 ³ organisms (primarily vegetative
	bacteria), which is substantially below the 10 ⁵ -10 ⁶ spores used with biological indicators.
Pathogen type	Spore-forming organisms are most resistant to sterilization and are the test organisms
	required for FDA clearance. However, the contaminating microflora on used surgical
	instruments consists mainly of vegetative bacteria.
Protein ¹	Residual protein decreases efficacy of sterilization. However, cleaning appears to rapidly
	remove protein load.
Salt ¹	Residual salt decreases efficacy of sterilization more than does protein load. However,
	cleaning appears to rapidly remove salt load.
Biofilm accumulation ¹	Biofilm accumulation reduces efficacy of sterilization by impairing exposure of the sterilant
	to the microbial cell.
Lumen length	Increasing lumen length impairs sterilant penetration. May require forced flow through
	lumen to achieve sterilization.
Lumen diameter	Decreasing lumen diameter impairs sterilant penetration. May require forced flow through
	lumen to achieve sterilization.
Restricted flow	Sterilant must come into contact with microorganisms. Device designs that prevent or
	inhibit this contact (e.g., sharp bends, blind lumens) will decrease sterilization efficacy.
Device design and	Materials used in construction may affect compatibility with different sterilization processes
construction	and affect sterilization efficacy. Design issues (e.g., screws, hinges) will also affect
	sterilization efficacy.
Madifiad frame Alfa and D	

Modified from Alfa and Rutala. 470, 825

¹ Factor only relevant for reused surgical/medical devices

Table 11. Comparative evaluation of the microbicidal activity of low-temperature sterilization technology (based on tests of inoculated carriers in lab setting)

Challenge	Ethylene Oxide 12/88	100% Ethylene Oxide	HCFC - Ethylene Oxide	Hydrogen Peroxide Gas Plasma 100	Hydrogen Peroxide Gas Plasma 100S	Peracetic Acid	Reference
No salt or serum ¹	100%	100%	96%	100%	No Data	No Data	Alfa 721
10% serum and 0.65% salt ²	97%	60%	95%	37%	No Data	No Data	Alfa ⁷²¹
Lumen (125 cm long × 3 mm wide) without serum or salt ¹	No Data	96%	96%	No Data	No Data	No Data	Alfa ⁷²¹
Lumen (125 cm long × 3 mm wide) with 10% serum and 0.65% salt ²	44%	40%	49%	35%	No Data	100% ¹	Alfa ⁷²¹
Lumen (40 cm long × 3 mm wide) ³	No Data	No Data	100%	95%	100%	8%	Rutala ⁸⁵⁶
Lumen (40 cm long × 2 mm wide) ³	No Data	No Data	100%	93%	100%	No Data	Rutala 856
Lumen (40 cm long × 1 mm wide) ³	No Data	No Data	100%	26%	100%	No Data	Rutala ⁸⁵⁶
Lumen (40 cm long × 3 mm wide) ⁴	No Data	No Data	100%	100%	100%	No Data	Rutala ⁸⁵⁶

Modified from Rutala. 825

HCFC=hydrochlorofluorocarbon

- ¹ Test organisms included Enterococcus faecalis, Mycobacterium chelonae, and Bacillus atrophaeus spores.
- ² Test organisms included *E. faecalis*, *P. aeruginosa*, *E. coli*, *M. chelonae*, *B. atrophaeus* spores, *G. stearothermophilus* spores, and *B. circulans* spores.
- ³ Test organism was *G. stearothermophilus* spores . The lumen test units had a removable 5 cm center piece (1.2 cm diameter) of stainless steel sealed to the narrower steel tubing by hard rubber septums.
- ⁴ Test organism was *G. stearothermophilus* spores. The lumen test unit was a straight stainless steel tube.

Table 12. Suggested protocol for management of positive biological indicator in a steam sterilizer

- Format Change [February 2017]: The format of this section was changed to improve readability and accessibility. The content is unchanged.
- 1. Take the sterilizer out of service. Notify area supervisor and infection control department.
- 2. Objects, other than implantable objects, do not need to be recalled because of a single positive spore test unless the sterilizer or the sterilization procedure is defective. As soon as possible, repeat biological indicator test in three consecutive sterilizer cycles. If additional spore tests remain positive, the items should be considered nonsterile, and supplies processed since the last acceptable (negative) biological indicator should be recalled. The items from the suspect load(s) should be recalled and reprocessed.
- 3. Check to ensure the sterilizer was used correctly (e.g., verify correct time and temperature setting). If not, repeat using appropriate settings and recall and reprocess all inadequately processed items.
- 4. Check with hospital maintenance for irregularities (e.g., electrical) or changes in the hospital steam supply (i.e., from standard ≥97% steam, <3% moisture). Any abnormalities should be reported to the person who performs sterilizer maintenance (e.g., medical engineering, sterilizer manufacturer).
- 5. Check to ensure the correct biological indicator was used and appropriately interpreted. If not, repeat using appropriate settings.
- If steps 1 through 5 resolve the problem
 If all three repeat biological indicators from three consecutive sterilizer cycles (step 2 above) are negative, put the
 sterilizer back in service.
- 7. If one or both biological indicators are positive, do one or more of the following until problem is resolved.
 - A. Request an inspection of the equipment by sterilizer maintenance personnel.
 - B. Have hospital maintenance inspect the steam supply lines.
 - C. Discuss the abnormalities with the sterilizer manufacturer.
 - D. Repeat the biological indicator using a different manufacturer's indicator.

If step 7 does not resolve the problem

Close sterilizer down until the manufacturer can assure that it is operating properly. Retest at that time with biological indicators in three consecutive sterilizer cycles.

Modified from Bryce. 839

Disclosure of Financial Interests and Relationships (2000- July 2004)

William A. Rutala: Honoraria from Advanced Sterilization Products, Kimberly-Clark; consultation with Advanced Sterilization Products, Aesculap, Clorox, 3M, SC Johnson, Intelligent Biocides, Metrex; and an educational grant from Consumer Specialty Products Association, Kimberly-Clark.

David J. Weber: Honoraria from Consumer Specialty Products Association; consultation with Clorox; and educational grant from Consumer Specialty Products Association.

References

- 1. Garner JS, Favero MS. CDC Guideline for handwashing and hospital environmental control, 1985. Infect. Control 1986;7:231-43.
- 2. Centers for Disease Control. Ambulatory and inpatient procedures in the United States, 1996. Atlanta, GA, 1998:1-39.
- 3. Uttley AH, Simpson RA. Audit of bronchoscope disinfection: a survey of procedures in England and Wales and incidents of mycobacterial contamination. J. Hosp. Infect. 1994;26:301-8.
- 4. Zaidi M, Angulo M, Sifuentes-Osornio J. Disinfection and sterilization practices in Mexico. J. Hosp. Infect. 1995;31:25-32.
- 5. McCarthy GM, Koval JJ, John MA, MacDonald JK. Infection control practices across Canada: do dentists follow the recommendations? J. Can. Dent. Assoc. 1999;65:506-11.
- 6. Spach DH, Silverstein FE, Stamm WE. Transmission of infection by gastrointestinal endoscopy and bronchoscopy. Ann. Intern. Med. 1993;118:117-28.
- 7. Weber DJ, Rutala WA. Lessons from outbreaks associated with bronchoscopy. Infect. Control Hosp. Epidemiol. 2001;22:403-8.
- 8. Weber DJ, Rutala WA, DiMarino AJ, Jr. The prevention of infection following gastrointestinal endoscopy: the importance of prophylaxis and reprocessing. In: DiMarino AJ, Jr, Benjamin SB, eds. Gastrointestinal diseases: an endoscopic approach. Thorofare, NJ: Slack Inc., 2002:87-106.
- 9. Meyers H, Brown-Elliott BA, Moore D, et al. An outbreak of *Mycobacterium chelonae* infection following liposuction. Clin. Infect. Dis. 2002;34:1500-7.
- 10. Lowry PW, Jarvis WR, Oberle AD, et al. *Mycobacterium chelonae* causing otitis media in an ear-nose-and-throat practice. N. Engl. J. Med. 1988;319:978-82.
- 11. Centers for Disease Control and Prevention. *Pseudomonas aeruginosa* infections associated with transrectal ultrasound-guided prostate biopsies--Georgia, 2005. MMWR CDC Surveill. Summ. 2006;55:776-7.
- 12. Mehta AC, Prakash UBS, Garland R, et al. Prevention of flexible bronchoscopy-associated infection. Chest 2006;128:1742-55.
- 13. Favero MS, Bond WW. Chemical disinfection of medical and surgical materials. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:881-917.
- 14. Spaulding EH. Chemical disinfection of medical and surgical materials. In: Lawrence C, Block SS, eds. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1968:517-31.
- 15. Simmons BP. CDC guidelines for the prevention and control of nosocomial infections. Guideline for hospital environmental control. Am. J. Infect. Control 1983;11:97-120.
- 16. Block SS. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001.
- 17. Rutala WA, 1994, 1995, and 1996 APIC Guidelines Committee. APIC guideline for selection and use of disinfectants. Association for Professionals in Infection Control and Epidemiology, Inc. Am. J. Infect. Control 1996;24:313-42.
- 18. Rutala WA. Disinfection, sterilization and waste disposal. In: Wenzel RP, ed. Prevention and control of nosocomial infections. Baltimore: Williams and Wilkins, 1997:539-93.
- 19. Rutala WA. APIC guideline for selection and use of disinfectants. Am. J. Infect. Control 1990;18:99-117.
- 20. Association of peri-Operative Registered Nurses. Recommended practices for high-level disinfection. AORN J. 2005;81:402-12.
- 21. Garner JS, Favero MS. CDC guidelines for the prevention and control of nosocomial infections. Guideline for handwashing and hospital environmental control, 1985. Supersedes guideline for hospital environmental control published in 1981. Am. J. Infect. Control 1986;14:110-29.
- 22. Centers for Disease Control. Guidelines for prevention of transmission of human immunodeficiency virus and hepatitis B virus to health-care and public-safety workers. MMWR 1989;38:1-37.
- 23. Centers for Disease Control. Guidelines for Environmental Infection Control in Health-Care Facilities, 2003. MMWR 2003;52 (No. RR-10):1-44.
- 24. Bhattachatyya M, Kepnes LJ. The effectiveness of immersion disinfection for flexible fiberoptic laryngoscopes. Otolaryngol Head Neck 2004;130:681-5.
- 25. Hamasuna R, Nose K, Sueyoshi T, Nagano M, Hasui Y, Osada Y. High-level disinfection of cystoscopic equipment with ortho-phthalaldehyde solution. J. Hosp. Infect. 2004;57:346-8.
- 26. Foliente RL KB, Aprecio RM, Bains HJ, Kettering JD, Chen YK. Efficacy of high-level disinfectants for

reprocessing gastrointestinal endoscopes in simulated-use testing. Gastrointest. Endosc. 2001;53:456-62.

- 27. Kovacs BJ, Chen YK, Kettering JD, Aprecio RM, Roy I. High-level disinfection of gastrointestinal endoscopes: are current guidelines adequate? Am. J. Gastroenterol. 1999;94:1546-50.
- 28. Rutala WA, Clontz EP, Weber DJ, Hoffmann KK. Disinfection practices for endoscopes and other semicritical items. Infect. Control Hosp. Epidemiol. 1991;12:282-8.
- 29. Phillips J, Hulka B, Hulka J, Keith D, Keith L. Laparoscopic procedures: The American Association of Gynecologic Laparoscopists' Membership Survey for 1975. J. Reprod. Med. 1977;18:227-32.
- 30. Muscarella LF. Current instrument reprocessing practices: Results of a national survey. Gastrointestinal Nursing 2001;24:253-60.
- 31. Wright EP, Collins CH, Yates MD. *Mycobacterium xenopi* and *Mycobacterium kansasii* in a hospital water supply. J. Hosp. Infect. 1985;6:175-8.
- 32. Wallace RJ, Jr., Brown BA, Driffith DE. Nosocomial outbreaks/pseudo-outbreaks caused by nontuberculous mycobacteria. Annu. Rev. Microbiol. 1998;52:453-90.
- 33. Mitchell DH, Hicks LJ, Chiew R, Montanaro JC, Chen SC. Pseudoepidemic of *Legionella pneumophila* serogroup 6 associated with contaminated bronchoscopes. J. Hosp. Infect. 1997;37:19-23.
- 34. Meenhorst PL, Reingold AL, Groothuis DG, et al. Water-related nosocomial pneumonia caused by *Legionella pneumophila* serogroups 1 and 10. J. Infect. Dis. 1985;152:356-64.
- 35. Atlas RM. *Legionella*: from environmental habitats to disease pathology, detection and control. Environ. Microbiol. 1999;1:283-93.
- 36. Rutala WA, Weber DJ. Water as a reservoir of nosocomial pathogens. Infect. Control Hosp. Epidemiol. 1997;18:609-16.
- 37. Weber DJ, Rutala WA. Environmental issues and nosocomial infections. In: Wenzel RP, ed. Prevention and control of nosocomial infections. Baltimore: Williams and Wilkins, 1997:491-514.
- 38. Society of Gastroenterology Nurses and Associates. Standards for infection control and reprocessing of flexible gastrointestinal endoscopes. Gastroenterol. Nurs. 2000;23:172-9.
- 39. Gerding DN, Peterson LR, Vennes JA. Cleaning and disinfection of fiberoptic endoscopes: evaluation of glutaraldehyde exposure time and forced-air drying. Gastroenterology 1982;83:613-8.
- 40. Society of Gastroenterology Nurses and Associates. Guideline for the use of high-level disnfectants and sterilants in reprocessing of flexible gastrointestinal endoscopes.
- 41. Turner AG, Higgins MM, Craddock JG. Disinfection of immersion tanks (Hubbard) in a hospital burn unit. Arch. Environ. Health 1974;28:101-4.
- 42. Rutala DR, Rutala WA, Weber DJ, Thomann CA. Infection risks associated with spirometry. Infect. Control Hosp. Epidemiol. 1991;12:89-92.
- 43. Kohn WG, Collins AS, Cleveland JL, Harte JA, Eklund KJ, Malvitz DM. Guidelines for infection control in dental health-care settings-2003. MMWR 2003;52 (no. RR-17):1-67.
- 44. Sehulster L, Chinn RYW, Healthcare Infection Control Practices Advisory Committee. Guidelines for environmental infection control in health-care facilities. MMWR 2003;52:1-44.
- 45. Rutala WA, White MS, Gergen MF, Weber DJ. Bacterial contamination of keyboards: Efficacy and functional impact of disinfectants. Infect Control Hosp Epidemiol 2006;27:372-7.
- Sattar SA, Lloyd-Evans N, Springthorpe VS, Nair RC. Institutional outbreaks of rotavirus diarrhoea: potential role of fomites and environmental surfaces as vehicles for virus transmission. J. Hyg. (Lond). 1986;96:277-89.
- 47. Weber DJ, Rutala WA. Role of environmental contamination in the transmission of vancomycin-resistant enterococci. Infect. Control Hosp. Epidemiol. 1997;18:306-9.
- 48. Ward RL, Bernstein DI, Knowlton DR, et al. Prevention of surface-to-human transmission of rotaviruses by treatment with disinfectant spray. J. Clin. Microbiol. 1991;29:1991-6.
- 49. Tyler R, Ayliffe GA, Bradley C. Virucidal activity of disinfectants: studies with the poliovirus. J. Hosp. Infect. 1990;15:339-45.
- 50. Gwaltney JM, Jr., Hendley JO. Transmission of experimental rhinovirus infection by contaminated surfaces. Am. J. Epidemiol. 1982;116:828-33.
- 51. Sattar SA, Jacobsen H, Springthorpe VS, Cusack TM, Rubino JR. Chemical disinfection to interrupt transfer of rhinovirus type 14 from environmental surfaces to hands. Appl. Environ. Microbiol. 1993;59:1579-85.
- 52. Sattar SA, Jacobsen H, Rahman H, Cusack TM, Rubino JR. Interruption of rotavirus spread through chemical disinfection. Infect. Control Hosp. Epidemiol. 1994;15:751-6.
- 53. Rutala WA, Barbee SL, Aguiar NC, Sobsey MD, Weber DJ. Antimicrobial activity of home disinfectants and

natural products against potential human pathogens. Infect. Control Hosp. Epidemiol. 2000;21:33-8.

- 54. Silverman J, Vazquez JA, Sobel JD, Zervos MJ. Comparative in vitro activity of antiseptics and disinfectants versus clinical isolates of *Candida* species. Infect. Control Hosp. Epidemiol. 1999;20:676-84.
- 55. Best M, Sattar SA, Springthorpe VS, Kennedy ME. Efficacies of selected disinfectants against *Mycobacterium tuberculosis.* J. Clin. Microbiol. 1990;28:2234-9.
- 56. Best M, Kennedy ME, Coates F. Efficacy of a variety of disinfectants against *Listeria* spp. Appl. Environ. Microbiol. 1990;56:377-80.
- 57. Best M, Springthorpe VS, Sattar SA. Feasibility of a combined carrier test for disinfectants: studies with a mixture of five types of microorganisms. Am. J. Infect. Control 1994;22:152-62.
- 58. Mbithi JN, Springthorpe VS, Sattar SA. Chemical disinfection of hepatitis A virus on environmental surfaces. Appl. Environ. Microbiol. 1990;56:3601-4.
- 59. Springthorpe VS, Grenier JL, Lloyd-Evans N, Sattar SA. Chemical disinfection of human rotaviruses: efficacy of commercially-available products in suspension tests. J. Hyg. (Lond). 1986;97:139-61.
- 60. Akamatsu T, Tabata K, Hironga M, Kawakami H, Uyeda M. Transmission of *Helicobacter pylori* infection via flexible fiberoptic endoscopy. Am. J. Infect. Control 1996;24:396-401.
- 61. Sattar SA, Springthorpe VS. Survival and disinfectant inactivation of the human immunodeficiency virus: a critical review. Rev. Infect. Dis. 1991;13:430-47.
- 62. Resnick L, Veren K, Salahuddin SZ, Tondreau S, Markham PD. Stability and inactivation of HTLV-III/LAV under clinical and laboratory environments. JAMA 1986;255:1887-91.
- 63. Weber DJ, Barbee SL, Sobsey MD, Rutala WA. The effect of blood on the antiviral activity of sodium hypochlorite, a phenolic, and a quaternary ammonium compound. Infect. Control Hosp. Epidemiol. 1999;20:821-7.
- 64. Rice EW, Clark RM, Johnson CH. Chlorine inactivation of *Escherichia coli* O157:H7. Emerg. Infect. Dis. 1999;5:461-3.
- 65. Pentella MA, Fisher T, Chandler S, Britt-Ohrmund T, Kwa BH, Yangco BG. Are disinfectants accurately prepared for use in hospital patient care areas? Infect. Control Hosp. Epidemiol. 2000;21:103.
- 66. Bhalla A, Pultz NJ, Gries DM, et al. Acquisition of nosocomial pathogens on hands after contact with environmental surfaces near hospitalized patients. Infect Control Hosp Epidemiol 2004;25:164-7.
- 67. Ray AJ, Hoyen CK, Taub TF, Eckstein EC, Donskey CJ. Nosocomial transmission of vancomycin-resistant enterococci from surfaces. JAMA 2002;287:1400-1.
- 68. Westwood JC, Mitchell MA, Legace S. Hospital sanitation: the massive bacterial contamination of the wet mop. Appl. Microbiol. 1971;21:693-7.
- 69. Rutala WA, Weber DJ. Disinfection of endoscopes: review of new chemical sterilants used for high-level disinfection. Infect. Control Hosp. Epidemiol. 1999;20:69-76.
- 70. Russell AD. Bacterial spores and chemical sporicidal agents. Clin. Microbiol. Rev. 1990;3:99-119.
- 71. Terleckyj B, Axler DA. Quantitative neutralization assay of fungicidal activity of disinfectants. Antimicrob. Agents Chemother. 1987;31:794-8.
- 72. Klein M, DeForest A. The inactivation of viruses by germicides. Chem. Specialists Manuf. Assoc. Proc. 1963;49:116-8.
- 73. Rutala WA, Cole EC, Wannamaker NS, Weber DJ. Inactivation of *Mycobacterium tuberculosis* and *Mycobacterium bovis* by 14 hospital disinfectants. Am. J. Med. 1991;91:267S-271S.
- 74. Robison RA, Bodily HL, Robinson DF, Christensen RP. A suspension method to determine reuse life of chemical disinfectants during clinical use. Appl. Environ. Microbiol. 1988;54:158-64.
- 75. Isenberg HD, Giugliano ER, France K, Alperstein P. Evaluation of three disinfectants after in-use stress. J. Hosp. Infect. 1988;11:278-85.
- 76. Cole EC, Rutala WA, Nessen L, Wannamaker NS, Weber DJ. Effect of methodology, dilution, and exposure time on the tuberculocidal activity of glutaraldehyde-based disinfectants. Appl. Environ. Microbiol. 1990;56:1813-7.
- 77. Power EG, Russell AD. Sporicidal action of alkaline glutaraldehyde: factors influencing activity and a comparison with other aldehydes. J. Appl. Bacteriol. 1990;69:261-8.
- 78. Rutala WA, Gergen MF, Weber DJ. Sporicidal activity of chemical sterilants used in hospitals. Infect. Control Hosp. Epidemiol. 1993;14:713-8.
- 79. Rutala WA, Gergen MF, Weber DJ. Inactivation of *Clostridium difficile* spores by disinfectants. Infect. Control Hosp. Epidemiol. 1993;14:36-9.
- 80. Ascenzi JM, Ezzell RJ, Wendt TM. A more accurate method for measurement of tuberculocidal activity of

disinfectants. Appl. Environ. Microbiol. 1987;53:2189-92.

- 81. Collins FM. Use of membrane filters for measurement of mycobactericidal activity of alkaline glutaraldehyde solution. Appl. Environ. Microbiol. 1987;53:737-9.
- 82. Rubbo SD, Gardner JF, Webb RL. Biocidal activities of glutaraldehyde and related compounds. J. Appl. Bacteriol. 1967;30:78-87.
- 83. Rutala WA, Weber DJ. FDA labeling requirements for disinfection of endoscopes: a counterpoint. Infect. Control Hosp. Epidemiol. 1995;16:231-5.
- 84. Collins FM. Kinetics of the tuberculocidal response by alkaline glutaraldehyde in solution and on an inert surface. J. Appl. Bacteriol. 1986;61:87-93.
- 85. Food and Drug Administration. 2005. FDA-cleared sterilants and high-level disinfectants with general claims for processing reusable medical and dental devices, May 13, 2005. [This link is no longer active: www.fda.gov/cdrh/ode/germlab.html].
- 86. Crow S, Metcalf RW, Beck WC, Birnbaum D. Disinfection or sterilization? Four views on arthroscopes. AORN J. 1983;37:854-9, 862-8.
- 87. Loffer FD. Disinfection vs. sterilization of gynecologic laparoscopy equipment. The experience of the Phoenix Surgicenter. J. Reprod. Med. 1980;25:263-6.
- 88. Johnson LL, Shneider DA, Austin MD, Goodman FG, Bullock JM, DeBruin JA. Two per cent glutaraldehyde: a disinfectant in arthroscopy and arthroscopic surgery. J. Bone Joint Surg. 1982;64:237-9.
- 89. Burns S, Edwards M, Jennings J, et al. Impact of variation in reprocessing invasive fiberoptic scopes on patient outcomes. Infect. Control Hosp. Epidemiol. 1996;17(suppl):P42.
- 90. Fuselier HA, Jr., Mason C. Liquid sterilization versus high level disinfection in the urologic office. Urology 1997;50:337-40.
- 91. Muscarella LF. High-level disinfection or "sterilization" of endoscopes? Infect. Control Hosp. Epidemiol. 1996;17:183-7.
- 92. Miles RS. What standards should we use for the disinfection of large equipment? J. Hosp. Infect. 1991;18:264-73.
- 93. Lee RM, Kozarek RA, Sumida SE, Raltz SL. Risk of contamination of sterile biopsy forceps in disinfected endoscopes. Gastrointest. Endosc. 1998;47:377-81.
- 94. Kinney TP, Kozarek RA, Raltz S, Attia F. Contamination of single-use biopsy forceps: A prospective in vitro analysis. Gastrointest. Endosc. 2002;56:209-12.
- 95. Centers for Disease Control. Recommendations for preventing possible transmission of human Tlymphotropic virus type III/lymphadenopathy-associated virus from tears. MMWR 1985;34:533-4.
- 96. Lettau LA, Bond WW, McDougal JS. Hepatitis and diaphragm fitting. JAMA 1985;254:752.
- 97. Schembre DB. Infectious complications associated with gastrointestinal endoscopy. Gastrointest. Endosc. Clin. N. Am. 2000;10:215-32.
- 98. Nelson DB. Infectious disease complications of GI endoscopy: Part II, exogenous infections. Gastrointest. Endosc. 2003;57:695-711.
- 99. Chu NS, Favero M. The microbial flora of the gastrointestinal tract and the cleaning of flexible endoscopes. Gastrointest. Endosc. Clin. N. Am. 2000;10:233-44.
- 100. Alfa MJ, Sitter DL. In-hospital evaluation of orthophthalaldehyde as a high level disinfectant for flexible endoscopes. J. Hosp. Infect. 1994;26:15-26.
- 101. Vesley D, Melson J, Stanley P. Microbial bioburden in endoscope reprocessing and an in-use evaluation of the high-level disinfection capabilities of Cidex PA. Gastroenterol. Nurs. 1999;22:63-8.
- 102. Chu NS, McAlister D, Antonoplos PA. Natural bioburden levels detected on flexible gastrointestinal endoscopes after clinical use and manual cleaning. Gastrointest. Endosc. 1998;48:137-42.
- Rutala WA, Weber DJ. Reprocessing endoscopes: United States perspective. J. Hosp. Infect. 2004;56:S27-S39.
- 104. Hanson PJ, Gor D, Clarke JR, et al. Contamination of endoscopes used in AIDS patients. Lancet 1989;2:86-8.
- 105. Hanson PJ, Gor D, Clarke JR, et al. Recovery of the human immunodeficiency virus from fibreoptic bronchoscopes. Thorax 1991;46:410-2.
- 106. Chaufour X, Deva AK, Vickery K, et al. Evaluation of disinfection and sterilization of reusable angioscopes with the duck hepatitis B model. J. Vasc. Surg. 1999;30:277-82.
- 107. Cheung RJ, Ortiz D, DiMarino AJ, Jr. GI endoscopic reprocessing practices in the United States. Gastrointest. Endosc. 1999;50:362-8.

- 108. American Society for Gastrointestinal Endoscopy. Position statement: reprocessing of flexible gastrointestinal endoscopes. Gastrointest. Endosc. 1996;43:541-6.
- 109. Food and Drug Administration. Content and format of premarket notification [510(k] submissions for liquid chemical sterilants/high level disinfectants. [This link is no longer active: www.fda.gov/cdrh/ode/397] 2000.
- 110. Urayama S, Kozarek RA, Sumida S, Raltz S, Merriam L, Pethigal P. Mycobacteria and glutaraldehyde: is high-level disinfection of endoscopes possible? Gastrointest. Endosc. 1996;43:451-6.
- 111. Jackson J, Leggett JE, Wilson DA, Gilbert DN. *Mycobacterium gordonae* in fiberoptic bronchoscopes. Am. J. Infect. Control 1996;24:19-23.
- 112. Martiny H, Floss H, Zuhlsdorf B. The importance of cleaning for the overall results of processing endoscopes. J. Hosp. Infect. 2004;56:S16-S22.
- 113. Alvarado CJ, Reichelderfer M. APIC guideline for infection prevention and control in flexible endoscopy. Association for Professionals in Infection Control. Am. J. Infect. Control 2000;28:138-55.
- 114. Society of Gastroenterology Nurses and Associates. Guideline for the use of high-level disinfectants and sterilants for reprocessing of flexible gastrointestinal endoscopes. Gastroenterol. Nurs. 2000;23:180-7.
- 115. Society of Gastroenterology Nurses and Associates. Standards of infection control in reprocessing of flexible gastrointestinal endoscopes. Gastroenterol. Nurs. 2006;29:142-8.
- 116. Nelson DB, Jarvis WR, Rutala WA, et al. Multi-society guideline for reprocessing flexible gastrointestinal endoscopes. Infect Control Hosp Epidemiol 2003;24:532-537.
- 117. Martin MA, Reichelderfer M, 1991, and 1993 APIC Guidelines Committee. APIC guidelines for infection prevention and control in flexible endoscopy. Am. J. Infect. Control 1994;22:19-38.
- 118. Rey JF, Halfon P, Feryn JM, Khiri H, Masseyeff MF, Ouzan D. Risk of transmission of hepatitis C virus by digestive endoscopy. Gastroenterol. Clin. Biol. 1995;19:346-9.
- 119. Cronmiller JR, Nelson DK, Jackson DK, Kim CH. Efficacy of conventional endoscopic disinfection and sterilization methods against *Helicobacter pylori* contamination. Helicobacter 1999;4:198-203.
- 120. Sartor C, Charrel RN, de Lamballerie X, Sambuc R, De Micco P, Boubli L. Evaluation of a disinfection procedure for hysteroscopes contaminated by hepatitis C virus. Infect. Control Hosp. Epidemiol. 1999;20:434-6.
- Hanson PJ, Chadwick MV, Gaya H, Collins JV. A study of glutaraldehyde disinfection of fibreoptic bronchoscopes experimentally contaminated with *Mycobacterium tuberculosis*. J. Hosp. Infect. 1992;22:137-42.
- 122. Merighi A, Contato E, Scagliarini R, et al. Quality improvement in gastrointestinal endoscopy: microbiologic surveillance of disinfection. Gastrointest. Endosc. 1996;43:457-62.
- 123. Bond WW. Endoscope reprocessing: Problems and solutions. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:151-163.
- Deva AK, Vickery K, Zou J, West RH, Harris JP, Cossart YE. Establishment of an in-use testing method for evaluating disinfection of surgical instruments using the duck hepatitis B model. J. Hosp. Infect. 1996;33:119-30.
- 125. Hanson PJ, Gor D, Jeffries DJ, Collins JV. Elimination of high titre HIV from fibreoptic endoscopes. Gut 1990;31:657-9.
- 126. Wu MS, Wang JT, Yang JC, et al. Effective reduction of *Helicobacter pylori* infection after upper gastrointestinal endoscopy by mechanical washing of the endoscope. Hepatogastroenterology. 1996;43:1660-4.
- 127. Kirschke DL, Jones TF, Craig AS, et al. *Pseudomonas aeruginosa* and *Serratia marcescens* contamination associated with a manufacturing defect in bronchoscopes. N. Engl. J. Med. 2003;348:214-20.
- 128. Srinivasan A, Wolfenden LL, Song X, et al. An outbreak of *Pseudomonas aeruginosa* infections associated with flexible bronchoscopes. N. Engl. J. Med. 2003;348:221-7.
- 129. Kaczmarek RG, Moore RM, Jr., McCrohan J, et al. Multi-state investigation of the actual disinfection/sterilization of endoscopes in health care facilities. Am. J. Med. 1992;92:257-61.
- 130. Bradley CR, Babb JR. Endoscope decontamination: automated vs. manual. J. Hosp. Infect. 1995;30:537-42.
- 131. Muscarella LF. Advantages and limitations of automatic flexible endoscope reprocessors. Am. J. Infect. Control 1996;24:304-9.
- Muscarella LF. Automatic flexible endoscope reprocessors. Gastrointest. Endosc. Clin. N. Am. 2000;10:245-57.
- 133. Alvarado CJ, Stolz SM, Maki DG. Nosocomial infections from contaminated endoscopes: a flawed automated endoscope washer. An investigation using molecular epidemiology. Am. J. Med. 1991;91:272S-280S.

- 134. Fraser VJ, Jones M, Murray PR, Medoff G, Zhang Y, Wallace RJ, Jr. Contamination of flexible fiberoptic bronchoscopes with *Mycobacterium chelonae* linked to an automated bronchoscope disinfection machine. Am. Rev. Respir. Dis. 1992;145:853-5.
- 135. Cooke RP, Whymant-Morris A, Umasankar RS, Goddard SV. Bacteria-free water for automatic washerdisinfectors: an impossible dream? J. Hosp. Infect. 1998;39:63-5.
- 136. Muscarella LF. Deja Vu...All over again? The importance of instrument drying. Infect. Control Hosp. Epidemiol. 2000;21:628-9.
- 137. Rutala WA, Weber DJ. Importance of lumen flow in liquid chemical sterilization. Am. J. Infect. Control 1999;20:458-9.
- 138. Dwyer DM, Klein EG, Istre GR, Robinson MG, Neumann DA, McCoy GA. *Salmonella newport* infections transmitted by fiberoptic colonoscopy. Gastrointest. Endosc. 1987;33:84-7.
- 139. Wheeler PW, Lancaster D, Kaiser AB. Bronchopulmonary cross-colonization and infection related to mycobacterial contamination of suction valves of bronchoscopes. J. Infect. Dis. 1989;159:954-8.
- 140. Bond WW. Virus transmission via fiberoptic endoscope: recommended disinfection. JAMA 1987;257:843-4.
- 141. Lynch DA, Porter C, Murphy L, Axon AT. Evaluation of four commercial automatic endoscope washing machines. Endoscopy 1992;24:766-70.
- 142. Bond WW. Disinfection and endoscopy: microbial considerations. J. Gastroenterol. Hepatol. 1991;6:31-6.
- 143. Nelson D. Newer technologies for endoscope disinfection: electrolyzed acid water and disposablecomponent endoscope systems. Gastrointest. Endosc. Clin. N. Am. 2000;10:319-28.
- 144. Silberman HD. Non-inflatable sterile sheath for introduction of the flexible nasopharyngolaryngoscope. Ann Otol, Rhinol, Laryngol 2001;110:385-7.
- 145. Kruse A, Rey JF. Guidelines on cleaning and disinfection in GI endoscopy. Update 1999. The European Society of Gastrointestinal Endoscopy. Endoscopy 2000;32:77-80.
- 146. British Thoracic Society. British Thoracic Society guidelines on diagnostic flexible bronchoscopy. Thorax 2001;56:1-21.
- 147. Association of Operating Room Nurses. Recommended practices for use and care of endoscopes. 2000 standards, recommended practices, and guidelines. Denver, CO: AORN, 2000:243-7.
- 148. British Society of Gastroenterology. Cleaning and disinfection of equipment for gastrointestinal endoscopy. Report of a working party of the British Society of Gastroenterology Endoscope Committee. Gut 1998;42:585-93.
- 149. Jackson FW, Ball MD. Correction of deficiencies in flexible fiberoptic sigmoidoscope cleaning and disinfection technique in family practice and internal medicine offices. Arch. Fam. Med. 1997;6:578-82.
- 150. Orsi GB, Filocamo A, Di Stefano L, Tittobello A. Italian National Survey of Digestive Endoscopy Disinfection Procedures. Endoscopy 1997;29:732-8; quiz 739-40.
- 151. Honeybourne D, Neumann CS. An audit of bronchoscopy practice in the United Kingdom: a survey of adherence to national guidelines. Thorax 1997;52:709-13.
- 152. Michele TM, Cronin WA, Graham NM, et al. Transmission of *Mycobacterium tuberculosis* by a fiberoptic bronchoscope. Identification by DNA fingerprinting. JAMA 1997;278:1093-5.
- 153. Bronowicki JP, Venard V, Botte C, et al. Patient-to-patient transmission of hepatitis C virus during colonoscopy. N. Engl. J. Med. 1997;337:237-40.
- 154. Agerton T, Valway S, Gore B, et al. Transmission of a highly drug-resistant strain (strain W1) of *Mycobacterium tuberculosis*. Community outbreak and nosocomial transmission via a contaminated bronchoscope. JAMA 1997;278:1073-7.
- 155. Food and Drug Administration, Centers for Disease Control and Prevention. FDA and CDC public health advisory: Infections from endoscopes inadequately reprocessed by an automated endoscope reprocessing system, Food and Drug Administration, Rockville, MD. 1999.
- 156. Nelson DB, Muscarella LF. Current issues in endoscope reprocessing and infection control during gastrointestinal endoscopy. World J Gastroenterol 2006;12:3953-64.
- 157. Riley R, Beanland C, Bos H. Establishing the shelf life of flexible colonoscopes. Gastroenterol. Nurs. 2002;25:114-9.
- 158. Rejchrt S, Cermak P, Pavlatova L, Mickova E, Bures J. Bacteriologic testing of endoscopes after high-level disinfection. Gastrointest. Endosc. 2004;60:76-8.
- 159. Willis C. Bacteria-free endoscopy rinse water a realistic aim? Epidemiol. Infect. 2005;134:279-84.
- 160. Humphreys H, McGrath H, McCormick PA, Walsh C. Quality of final rinse water used in washer-disinfectors for endoscopes. J. Hosp. Infect. 2002;51:151-3.

- 161. Pang J, Perry P, Ross A, Forbes GM. Bacteria-free rinse water for endoscope disinfection. Gastrointest. Endosc. 2002;56:402-6.
- 162. Leung J, Vallero R, Wilson R. Surveillance cultures to monitor quality of gastrointestinal endoscope reprocessing. Am. J. Gastroenterol. 2003;98.
- 163. Moses FM, Lee J. Surveillance cultures to monitor quality of gastrointestinal endoscope reprocessing. Am. J. Gastroenterol. 2003;98:77-81.
- 164. Tunuguntla A, Sullivan MJ. Monitoring quality of flexible endoscopic disinfection by microbiologic surveillance cultures. Tennessee Med 2004;October:453-6.
- 165. Muscarella LF. Application of environmental sampling to flexible endoscope reprocessing: The importance of monitoring the rinse water. Infect . Control Hosp. Epidemiol. 2002;23:285-9.
- 166. Fraser TG, Reiner S, Malcznski M, Yarnold PR, Warren J, Noskin GA. Multidrug-resistant *Pseudomonas* aeruginosa cholangiopancreatography: Failure of routine endoscope cultures to prevent an outbreak. Infect Control Hosp Epidemiol 2004;25:856-9.
- 167. Bond WW, Hedrick ER. Microbiological culturing of environmental and medical-device surfaces. In: Isenberg HD, and M.J.R. Gilchrist, ed. Clinical Microbiology Procedures Handbook, Section 11, Epidemiologic and Infection Control Microbiology. Washington, DC: American Society for Microbiology, 1992:11.10.1-11.10.9.
- Murray PR, Baron EJ, Pfaller MA, Jorgensen JH, Yolken RH. Manual of Clinical Microbiology. In: Murray PR, Baron EJ, Pfaller MA, Jorgensen JH, Yolken RH, eds. Washington, D.C.: American Society for Microbiology Press, 2003.
- 169. Blob R, Kampf G. Test models to determine cleaning efficacy with different types of bioburden and its clinical correlation. J. Hosp. Infect. 2004;56 (suppl):S44-S48.
- 170. Obee PC, Griffith CJ, Cooper RA, Cooke RP, Bennion NE, Lewis M. Real-time monitoring in managing the decontamination of flexible gastrointestinal endoscopes. Am. J. Infect. Control 2005;33:202-6.
- 171. Sciortino CV, Xia EL, Mozee A. Assessment of a novel approach to evaluate the outcome of endoscope reprocessing. Infect Control Hosp Epidemiol 2004;25:284-90.
- 172. Murphy C. Inactivated glutaraldehyde: Lessons for infection control. Am. J. Infect. Control 1998;26:159-60.
- 173. Carsauw H, Debacker N. Recall of patients after use of inactive batch of Cidex disinfection solution in Belgian hospitals, Fifth International Conference of the Hospital Infection Society, Edinburgh, September 15-18, 2002. Hospital Infections Society.
- 174. Ad hoc Committee on Infection Control in the Handling of Endoscopic Equipment. Guidelines for preparation of laparoscopic instrumentation. AORN J. 1980;32:65-6, 70, 74, 76.
- 175. Taylor EW, Mehtar S, Cowan RE, Feneley RC. Endoscopy: disinfectants and health. Report of a meeting held at the Royal College of Surgeons of England, February 1993. J. Hosp. Infect. 1994;28:5-14.
- 176. Hulka JF, Wisler MG, Bruch C. A discussion: laparoscopic instrument sterilization. Med. Instrum. 1977;11:122-3.
- 177. Corson SL, Block S, Mintz C, Dole M, Wainwright A. Sterilization of laparoscopes. Is soaking sufficient? J. Reprod. Med. 1979;23:49-56.
- 178. Corson SL, Dole M, Kraus R, Richards L, Logan B. Studies in sterilization of the laparoscope: II. J. Reprod. Med. 1979;23:57-9.
- 179. Chan-Myers H, McAlister D, Antonoplos P. Natural bioburden levels detected on rigid lumened medical devices before and after cleaning. Am. J. Infect. Control 1997;25:471-6.
- 180. Rodrigues C, Mehta AC, Jha U, Bharucha M, Dastur FD, Udwadia TE. Nosocomial *Mycobacterium chelonae* infection in laparoscopic surgery. Infect . Control Hosp. Epidemiol. 2001;22:474-5.
- 181. Marshburn PB, Rutala WA, Wannamaker NS, Hulka JF. Gas and steam sterilization of assembled versus disassembled laparoscopic equipment. Microbiologic studies. J. Reprod. Med. 1991;36:483-7.
- 182. Bernhang AM. Clostridium pyoarthrosis following arthroscopy. Arthroscopy 1987;3:56-8.
- 183. D'Angelo GL, Ogilvie-Harris DJ. Septic arthritis following arthroscopy, with cost/benefit analysis of antibiotic prophylaxis. Arthroscopy 1988;4:10-4.
- 184. Weber DJ, Rutala WA. Nosocomial ocular infections. In: Mayhall CG, ed. Infect. Control and Hosp. Epidemiol. Philadelphia: Lippincott Williams & Wilkins, 1999:287-99.
- 185. Rutala WA, Peacock JE, Gergen MF, Sobsey MD, Weber DJ. Efficacy of hospital germicides against adenovirus 8, a common cause of epidemic keratoconjunctivitis in health care facilities. Antimicrob. Agents Chemother. 2006;50:1419-24.
- 186. Sattar SA, Springthorpe VS, Karim Y, Loro P. Chemical disinfection of non-porous inanimate surfaces experimentally contaminated with four human pathogenic viruses. Epidemiol. Infect. 1989;102:493-505.

- 187. Chronister CL. Structural damage to Schiotz tonometers after disinfection with solutions. Optom. Vis. Sci. 1997;74:164-6.
- 188. Nagington J, Sutehall GM, Whipp P. Tonometer disinfection and viruses. Br. J. Ophthalmol. 1983;67:674-6.
- 189. Craven ER, Butler SL, McCulley JP, Luby JP. Applanation tonometer tip sterilization for adenovirus type 8. Ophthalmology 1987;94:1538-40.
- 190. American Academy of Ophthalmology. Updated recommendations for ophthalmic practice in relation to the human immunodeficiency virus. American Academy of Ophthalmology, San Francisco, CA, 1988.
- 191. Pepose JS, Linette G, Lee SF, MacRae S. Disinfection of Goldmann tonometers against human immunodeficiency virus type 1. Arch. Ophthalmol. 1989;107:983-5.
- 192. Ventura LM, Dix RD. Viability of herpes simplex virus type 1 on the applanation tonometer. Am. J. Ophthalmol. 1987;103:48-52.
- 193. Koo D, Bouvier B, Wesley M, Courtright P, Reingold A. Epidemic keratoconjunctivitis in a university medical center ophthalmology clinic; need for re-evaluation of the design and disinfection of instruments. Infect. Control Hosp. Epidemiol. 1989;10:547-52.
- 194. Jernigan JA, Lowry BS, Hayden FG, et al. Adenovirus type 8 epidemic keratoconjunctivitis in an eye clinic: risk factors and control. J. Infect. Dis. 1993;167:1307-13.
- 195. Fritz S, Hust MH, Ochs C, Gratwohl I, Staiger M, Braun B. Use of a latex cover sheath for transesophageal echocardiography (TEE) instead of regular disinfection of the echoscope? Clin. Cardiol. 1993;16:737-40.
- 196. Lawrentschuk N, Chamberlain M. Sterile disposable sheath sytsem for flexible cytoscopes. Urology 2005;66:1310-3.
- 197. Milki AA, Fisch JD. Vaginal ultrasound probe cover leakage: implications for patient care. Fertil. Steril. 1998;69:409-11.
- 198. Storment JM, Monga M, Blanco JD. Ineffectiveness of latex condoms in preventing contamination of the transvaginal ultrasound transducer head. South. Med. J. 1997;90:206-8.
- 199. Hignett M, Claman P. High rates of perforation are found in endovaginal ultrasound probe covers before and after oocyte retrieval for *in vitro* fertilization-embryo transfer. J. Assist. Reprod. Genet. 1995;12:606-9.
- 200. Amis S, Ruddy M, Kibbler CC, Economides DL, MacLean AB. Assessment of condoms as probe covers for transvaginal sonography. J. Clin. Ultrasound 2000;28:295-8.
- 201. Rooks VJ, Yancey MK, Elg SA, Brueske L. Comparison of probe sheaths for endovaginal sonography. Obstet. Gynecol. 1996;87:27-9.
- 202. Odwin CS, Fleischer AC, Kepple DM, Chiang DT. Probe covers and disinfectants for transvaginal transducers. J. Diagnostic Med. Sonography 1990;6:130-5.
- 203. Benson WG. Exposure to glutaraldehyde. J. Soc. Occup. Med. 1984;34:63-4.
- 204. Garland SM, de Crespigny L. Prevention of infection in obstetrical and gynaecological ultrasound practice. Aust. N. Z. J. Obstet Gynaecol. 1996;36:392-5.
- 205. Fowler C, McCracken D. US probes: risk of cross infection and ways to reduce it--comparison of cleaning methods. Radiology 1999;213:299-300.
- 206. Muradali D, Gold WL, Phillips A, Wilson S. Can ultrasound probes and coupling gel be a source of nosocomial infection in patients undergoing sonography? An in vivo and in vitro study. AJR. Am. J. Roentgenol. 1995;164:1521-4.
- 207. Lewis DL, Arens M, Appleton SS, et al. Cross-contamination potential with dental equipment. Lancet 1992;340:1252-4.
- 208. Lewis DL, Boe RK. Cross-infection risks associated with current procedures for using high-speed dental handpieces. J. Clin. Microbiol. 1992;30:401-6.
- 209. American Dental Association. Infection control recommendations for the dental office and the dental laboratory. JADA 1996;127:672-80.
- 210. Centers for Disease Control. Recommended Infection-Control Practices for Dentistry, 1993. MMWR 1993;41:1-12.
- 211. Department of Health and Human Services. Food and Drug Administration. Dental handpiece sterilization, Food and Drug Administration, Rockville, MD, 1992.
- 212. Silverstone SE, Hill DE. Evaluation of sterilization of dental handpieces by heating in synthetic compressor lubricant. Gen. Dent. 1999;47:158-60.
- 213. Goodman HS, Carpenter RD, Cox MR. Sterilization of dental instruments and devices: an update. Am. J. Infect. Control 1994;22:90-4.
- 214. Occupational Safety and Health Administration. Occupational exposure to bloodborne pathogens; final rule.

Fed. Regist. 1991;56:64003-182.

- 215. Occupational Safety and Health Administration. OSHA Memorandum from Stephen Mallinger. EPAregistered disinfectants for HIV/HBV. Washington, DC, 1997.
- 216. Gurevich I, Dubin R, Cunha BA. Dental instrument and device sterilization and disinfection practices. J. Hosp. Infect. 1996;32:295-304.
- 217. Smith A, Dickson M, Aitken J, Bagg J. Contaminated dental instruments. J. Hosp. Infect. 2002;51:233-5.
- 218. Hastreiter RJ, Molinari JA, Falken MC, Roesch MH, Gleason MJ, Merchant VA. Effectiveness of dental office instrument sterilization procedures. J. Am. Dent. Assoc. 1991;122:51-6.
- 219. Andres MT, Tejerina JM, Fierro JF. Reliability of biologic indicators in a mail-return sterilization-monitoring service: a review of 3 years. Quintessence Int. 1995;26:865-70.
- 220. Miller CH, Sheldrake MA. The ability of biological indicators to detect sterilization failures. Am. J. Dent. 1994;7:95-7.
- 221. Sarin PS, Scheer DI, Kross RD. Inactivation of human T-cell lymphotropic retrovirus (HTLV-III) by LD. N. Engl. J. Med. 1985;313:1416.
- 222. Sarin PS, Scheer DI, Kross RD. Inactivation of human T-cell lymphotropic retrovirus. Environ Microbiol 1990;56:1423-8.
- 223. Ascenzi JM. Standardization of tuberculocidal testing of disinfectants. J. Hosp. Infect. 1991;18:256-63.
- 224. Bond WW, Favero MS, Petersen NJ, Ebert JW. Inactivation of hepatitis B virus by intermediate-to-high-level disinfectant chemicals. J. Clin. Microbiol. 1983;18:535-8.
- 225. Kobayashi H, Tsuzuki M. The effect of disinfectants and heat on hepatitis B virus. J. Hosp. Infect. 1984;5:93 4.
- 226. Spire B, Barre-Sinoussi F, Montagnier L, Chermann JC. Inactivation of lymphadenopathy associated virus by chemical disinfectants. Lancet 1984;2:899-901.
- 227. Martin LS, McDougal JS, Loskoski SL. Disinfection and inactivation of the human T lymphotropic virus type III/Lymphadenopathy-associated virus. J. Infect. Dis. 1985;152:400-3.
- 228. Centers for Disease Control. Recommendations for prevention of HIV transmission in health-care settings. MMWR 1987;36:S3-S18.
- 229. Prince DL, Prince HN, Thraenhart O, Muchmore E, Bonder E, Pugh J. Methodological approaches to disinfection of human hepatitis B virus. J. Clin. Microbiol. 1993;31:3296-304.
- 230. Prince DL, Prince RN, Prince HN. Inactivation of human immunodeficiency virus type 1 and herpes simplex virus type 2 by commercial hospital disinfectants. Chemical Times and Trends 1990;13:13-16.
- 231. Sattar SA, Springthorpe VS, Conway B, Xu Y. Inactivation of the human immunodeficiency virus: an update. Rev. Med. Microbiol. 1994;5:139-150.
- 232. Kaplan JC, Crawford DC, Durno AG, Schooley RT. Inactivation of human immunodeficiency virus by Betadine. Infect. Control 1987;8:412-4.
- 233. Hanson PJ, Gor D, Jeffries DJ, Collins JV. Chemical inactivation of HIV on surfaces. Br. Med. J. 1989;298:862-4.
- 234. Hanson PJ, Jeffries DJ, Collins JV. Viral transmission and fibreoptic endoscopy. J. Hosp. Infect. 1991;18:136-40.
- 235. Payan C, Cottin J, Lemarie C, Ramont C. Inactivation of hepatitis B virus in plasma by hospital in-use chemical disinfectants assessed by a modified HepG2 cell culture. J. Hosp. Infect. 2001;47:282-87.
- 236. Chanzy B, Duc-Bin DL, Rousset B, et al. Effectiveness of a manual disinfection procedure in eliminating hepatitis C virus from experimentally contaminated endoscopes. Gastrointest. Endosc. 1999;50:147-51.
- Druce JD, Russell JS, Birch CJ, Yates LA, Harper RW, Smolich JJ. A decontamination and sterilization protocol employed during reuse of cardiac electrophysiology catheters inactivates human immunodeficiency virus. Infect. Control Hosp. Epidemiol. 2003;24:184-90.
- 238. Payan C, Pivert A, Kampf G, Ramont C, Cottin J, Lemarie C. Assessment of new chemical disinfectants for HBV virucidal activity in a cell culture model. J. Hosp. Infect. 2004;56 (suppl):S58-S63.
- 239. Reynolds CD, Rhinehart E, Dreyer P, Goldmann DA. Variability in reprocessing policies and procedures for flexible fiberoptic endoscopes in Massachusetts hospitals. Am. J. Infect. Control 1992;20:283-90.
- 240. Handsfield HH, Cummings MJ, Swenson PD. Prevalence of antibody to human immunodeficiency virus and hepatitis B surface antigen in blood samples submitted to a hospital laboratory. Implications for handling specimens. JAMA 1987;258:3395-7.
- 241. Baker JL, Kelen GD, Sivertson KT, Quinn TC. Unsuspected human immunodeficiency virus in critically ill emergency patients. JAMA 1987;257:2609-11.

- 242. Kelen GD, Fritz S, Qaqish B, et al. Unrecognized human immunodeficiency virus infection in emergency department patients. N. Engl. J. Med. 1988;318:1645-50.
- 243. Ishino Y, Ido K, Sugano K. Contamination with hepatitis B virus DNA in gastrointestinal endoscope channels: Risk of infection on reuse after on-site cleaning. Endoscopy 2005;37:548-51.
- 244. Agolini G, Russo A, Clementi M. Effect of phenolic and chlorine disinfectants on hepatitis C virus binding and infectivity. Am. J. Infect. Control 1999;27:236-9.
- 245. Alter MJ, Tokars JI, Arduino MJ, Favero MS. Nosocomial infections with hemodialysis. In: Mayhall CG, ed. Infect. Control and Hosp. Epidemiol. Philadelphia: Lippincott Williams & Wilkins, 2004:1139-60.
- 246. Centers for Disease Control. Recommendations for preventing transmission of infections among chronic hemodialysis patients. MMWR. 2001;50:1-43.
- 247. Velandia M, Fridkin SK, Cardenas V, et al. Transmission of HIV in dialysis centre. Lancet 1995;345:1417-22.
- 248. Guinto CH, Bottone EJ, Raffalli JT, Montecalvo MA, Wormser GP. Evaluation of dedicated stethoscopes as a potential source of nosocomial pathogens. Am. J. Infect. Control 2002;30:499-502.
- 249. Tokars JI, Miller ER, Alter MJ, Arduino MJ. National surveillance of dialysis-associated diseases in the United States, 1997. Semin. Dialysis 2000;13:75-85.
- 250. Amato RL, Curtis JM. The practical application of ozone in dialysis. Nephrol. News Issues 2002;September 27-9.
- 251. Smeets E, Koonman J, van der Sande F, et al. Prevention of biofilm formation in dialysis water treatment systems. Kidney Int. 2003;63:1574-6.
- 252. Finelli L, Miller JT, Tokars JI, Alter MJ, Arduino MJ. National surveillance of dialysis-associated diseases in the United States, 2002. Semin Dialysis 2005;18:52-61.
- 253. Association for the Advancement of Medical Instrumentation. Reuse of hemodialyzers: Association for the Advancement of Medical Instrumentation, Arlington VA, 2002/2003:ANSI/AAMI RD47:2002 & RD47:2002/A1:2003; 1-32.
- 254. Kim KH, Fekety R, Batts DH, et al. Isolation of *Clostridium difficile* from the environment and contacts of patients with antibiotic-associated colitis. J. Infect. Dis. 1981;143:42-50.
- 255. Skoutelis AT, Westenfelder GO, Beckerdite M, Phair JP. Hospital carpeting and epidemiology of *Clostridium difficile*. Am. J. Infect. Control 1994;22:212-7.
- 256. Wilcox MH, Fawley WN. Hospital disinfectants and spore formation by *Clostridium difficile*. Lancet 2000;356:1324.
- 257. Kaatz GW, Gitlin SD, Schaberg DR, et al. Acquisition of *Clostridium difficile* from the hospital environment. Am. J. Epidemiol. 1988;127:1289-94.
- 258. Wilcox MH, Fawley WN, Wigglesworth N, Parnell P, Verity P, Freeman J. Comparison of the effect of detergent versus hypochlorite cleaning on environmental contamination and incidence of *Clostridium difficile* infection. J. Hosp. Infect. 2003;54:109-14.
- 259. Mayfield JL, Leet T, Miller J, Mundy LM. Environmental control to reduce transmission of *Clostridium difficile*. Clin. Infect. Dis. 2000;31:995-1000.
- 260. Marler LM, Siders JA, Wolters LC, Pettigrew Y, Skitt BL, Allen SD. Comparison of five cultural procedures for isolation of *Clostridium difficile* from stools. J. Clin. Microbiol. 1992;30:514-6.
- 261. Brazier JG. The diagnosis of *Clostridium difficile*-associated disease. J. Antimicrob. Chemother. 1998;41 (suppl):29-40.
- 262. Perez J, Springthorpe S, Sattar SA. Activity of selected oxidizing microbicides against spores of *Clostridium difficile*: Relevance to environmental control. Am. J. Infect. Control 2005;33:320-5.
- 263. McFarland LV, Mulligan ME, Kwok RY, Stamm WE. Nosocomial acquisition of *Clostridium difficile* infection. N. Engl. J. Med. 1989;320:204-10.
- 264. Jernigan JA, Siegman-Igra Y, Guerrant RC, Farr BM. A randomized crossover study of disposable thermometers for prevention of *Clostridium difficile* and other nosocomial infections. Infect Control Hosp Epidemiol 1998;19:494-9.
- 265. Hughes CE, Gebhard RL, Peterson LR, Gerding DN. Efficacy of routine fiberoptic endoscope cleaning and disinfection for killing *Clostridium difficile*. Gastrointest. Endosc. 1986;32:7-9.
- 266. Dyas A, Das BC. The activity of glutaraldehyde against *Clostridium difficile*. J. Hosp. Infect. 1985;6:41-5.
- 267. Wullt M, Odenholt I, Walder M. Activity of three disinfectants and acidified nitrite against *Clostridium difficile* spores. Infect Control Hosp Epidemiol 2003;24:765-8.
- 268. Block C. The effect of Perasafe and sodium dichloroisocyanurate (NaDCC) against spores of *Clostridium difficile* and *Bacillus atrophaeus* on stainless steel and polyvinyl chloride surfaces. J. Hosp. Infect.

2004;57:144-8.

- 269. Occupational Safety and Health Administration. OSHA instruction CPL 2-2.44C. Office of Health Compliance Assistance. Washington, DC, 1992.
- 270. Rutala WA, Weber DJ. Infection control: the role of disinfection and sterilization. J. Hosp. Infect. 1999;43:S43-55.
- 271. Barbee SL, Weber DJ, Sobsey MD, Rutala WA. Inactivation of *Cryptosporidium parvum* oocyst infectivity by disinfection and sterilization processes. Gastrointest. Endosc. 1999;49:605-11.
- 272. Wilson JA, Margolin AB. The efficacy of three common hospital liquid germicides to inactivate *Cryptosporidium parvum* oocysts. J. Hosp. Infect. 1999;42:231-7.
- 273. Fayer R, Graczyk TK, Cranfield MR, Trout JM. Gaseous disinfection of *Cryptosporidium parvum* oocysts. Appl. Environ. Microbiol. 1996;62:3908-9.
- Venkitanarayanan KS, Ezeike GO, Hung YC, Doyle MP. Inactivation of *Escherichia coli* O157:H7 and *Listeria monocytogenes* on plastic kitchen cutting boards by electrolyzed oxidizing water. J. Food Prot. 1999;62:857-60.
- 275. Taormina PJ, Beuchat LR. Behavior of enterohemorrhagic *Escherichia coli* O157:H7 on alfalfa sprouts during the sprouting process as influenced by treatments with various chemicals. J. Food Prot. 1999;62:850-6.
- 276. Taormina PJ, Beuchat LR. Comparison of chemical treatments to eliminate enterohemorrhagic *Escherichia coli* O157:H7 on alfalfa seeds. J. Food Prot. 1999;62:318-24.
- 277. Castillo A, Lucia LM, Kemp GK, Acuff GR. Reduction of *Escherichia coli* O157:H7 and *Salmonella typhimurium* on beef carcass surfaces using acidified sodium chlorite. J. Food Prot. 1999;62:580-4.
- 278. Graham DY, Osato MS. Disinfection of biopsy forceps and culture of *Helicobacter pylori* from gastric mucosal biopsies. Am. J. Gastroenterol. 1999;94:1422-3.
- 279. Kaneko H, Mitsuma T, Kotera H, Uchida K, Furusawa A, Morise K. Are routine cleaning methods sufficient to remove *Helicobacter pylori* from endoscopic equipment? Endoscopy 1993;25:435.
- 280. Langenberg W, Rauws EA, Oudbier JH, Tytgat GN. Patient-to-patient transmission of *Campylobacter pylori* infection by fiberoptic gastroduodenoscopy and biopsy. J. Infect. Dis. 1990;161:507-11.
- 281. Miyaji H, Kohli Y, Azuma T, et al. Endoscopic cross-infection with *Helicobacter pylori*. Lancet 1995;345:464.
- 282. Fantry GT, Zheng QX, James SP. Conventional cleaning and disinfection techniques eliminate the risk of endoscopic transmission of *Helicobacter pylori*. Am. J. Gastroenterol. 1995;90:227-32.
- 283. Shimada T, Terano A, Ota S, Takikawa H, Sumino S. Risk of iatrogenic transmission of *Helicobacter pylori* by gastroscopes. Lancet 1996;347:1342-3.
- 284. Roosendaal R, Kuipers EJ, van den Brule AJ, et al. Detection of *Helicobacter pylori* DNA by PCR in gastrointestinal equipment. Lancet 1993;341:900.
- 285. Johnson CH, Rice EW, Reasoner DJ. Inactivation of *Helicobacter pylori* by chlorination. Appl. Environ. Microbiol. 1997;63:4969-70.
- 286. Chapin M, Yatabe J, Cherry JD. An outbreak of rotavirus gastroenteritis on a pediatric unit. Am. J. Infect. Control 1983;11:88-91.
- 287. Keswick BH, Pickering LK, DuPont HL, Woodward WE. Survival and detection of rotaviruses on environmental surfaces in day care centers. Appl. Environ. Microbiol. 1983;46:813-16.
- 288. Ansari SA, Spingthorpe S, Sattar SA. Survival and vehicular spread of human rotaviruses: Possible relation to seasonality of outbreaks. Rev. Infect. Dis. 1991;13:448-61.
- 289. Ansari SA, Sattar SA, Springthorpe VS, Wells GA, Tostowaryk W. Rotavirus survival on human hands and transfer of infectious virus to animate and nonporous inanimate surfaces. J. Clin. Microbiol. 1988;26:1513-8.
- 290. Sattar SA, Raphael RA, Lochnan H, Springthorpe VS. Rotavirus inactivation by chemical disinfectants and antiseptics used in hospitals. Can. J. Microbiol. 1983;29:1464-9.
- 291. Lloyd-Evans N, Springthorpe VS, Sattar SA. Chemical disinfection of human rotavirus-contaminated inanimate surfaces. J. Hyg. (Lond). 1986;97:163-73.
- 292. Tan JA, Schnagl RD. Inactivation of a rotavirus by disinfectants. Med. J. Aust. 1981;1:19-23.
- 293. Sattar SA, Springthorpe VS, Karim Y, Loro P. Chemical disinfection of non-porous inanimate surfaces experimentally contaminated with four human pathogenic viruses. Epidemiol Infect 1989;102:493-505.
- 294. Green J, Wright PA, Gallimore CI, Mitchell O, Morgan-Capner P, Brown DWG. The role of environmental contamination with small round structured viruses in a hospital outbreak investigated by reverse-transcriptase polymerase chain reaction assay. J. Hosp. Infect. 1998;39:39-45.
- 295. Evans MR, Meldrum R, Lane W, et al. An outbreak of viral gastroenteritis following environmental contamination at a concert hall. Epidemiol & Infect 2002;129:355-360.

- 296. Marks PJ, Vipond IB, Regan FM, Wedgwood K, Fey RE, Caul EO. A school outbreak of Norwalk-like virus: Evidence for airborne transmission. Epidemiol Infect 2003;131:727-36.
- 297. Doultree JC, Druce JD, Birch CJ, Bowden DS, Marshall JA. Inactivation of feline calicivirus, a Norwalk virus surrogate. J. Hosp. Infect. 1999;41:51-7.
- 298. Sattar SA. Microbicides and the environmental control of nosocomial viral infections. J. Hosp. Infect. 2004;56 (suppl):S64-S69.
- 299. Jimenez L, Chiang M. Virucidal activity of a quaternary ammonium compound disinfectant against feline calicivirus: A surrogate for norovirus. Am. J. Infect. Control 2006;34:269-73.
- Gehrke C, Steinmann J, Goroncy-Bermes P. Inactivation of feline calicivirus, a surrogate of norovirus (formerly Norwalk-like viruses), by different types of alcohol in vitro and in vivo. J. Hosp. Infect. 2004;56:49-55.
- 301. Centers for Disease Control and Prevention. Update: Severe acute respiratory syndrome United States, May 14, 2003. MMWR 2003;52:436-8.
- 302. Saknimit M, Inatsuki I, Sugiyama Y, Yagami K. Virucidal efficacy of physico-chemical treatments against coronaviruses and parvoviruses of laboratory animals. Jikken Dobutsu. 1988;37:341-5.
- 303. Sizun J, Yu MW, Talbot PJ. Survival of human coronaviruses 229E and OC43 in suspension and after drying onsurfaces: a possible source of hospital-acquired infections. J. Hosp. Infect. 2000;46:55-60.
- 304. Kariwa H, Fujii N, Takashima I. Inactivation of SARS coronavirus by means of povidone-iodine, physical conditions, and chemical reagents. Jpn. J. Vet. Res. 2004;52:105-12.
- 305. Greub G, Raoult D. Biocides currently used for bronchoscope decontamination are poorly effective against free-living amoebae. Infect Control Hosp Epidemiol 2003;24:784-6.
- 306. Leggiadro RJ. The threat of biological terrorism: A public health and infection control reality. Infect . Control Hosp. Epidemiol. 2000;21:53-6.
- 307. Henderson DA. The looming threat of bioterrorism. Science 1999;283:1279-82.
- 308. Centers for Disease Control. Biological and chemical terrorism: strategic plan for preparedness and response. MMWR 2000;49 (no. RR-4):1-14.
- 309. Weber DJ, Rutala WA. Disinfection and sterilization of potential bioterrorism agents. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices and new research. Washington DC: Association for Professionals in Infection Control and Epidemiology, 2004:86-103.
- 310. Ferrier A, Garin D, Crance JM. Rapid inactivation of vaccinia virus in suspension and dried on surfaces. J. Hosp. Infect. 2004;57:73-9.
- 311. Butcher W, Ulaeto D. Contact inactivation of orthopoxviruses by household disinfectants. J. Appl. Microbiol. 2005;99:279-84.
- 312. Brazis AR, Leslie JE, PW K, RL W. The inactivation of spores of *Bacillus globigii* and *Bacillus anthracis* by free available chlorine. Appl. Microbiol. 1958;6:338-342.
- 313. Sattar SA, Springthorpe VS, Adegbunrin O. Is *Bacillus subtilis* (ATCC 19659) a suitable surrogate for evaluating microbicides against Bacillus anthracis., Association for Official Analytical Chemists International Annual Meeting, St. Louis, Missouri, 2004.
- 314. Whitney EAS, Beatty ME, Taylor TH Jr, et al. Inactivation of *Bacillus anthracis* spores. Emerg. Infect. Dis. 2003;9:623-7.
- 315. Weber DJ, Rutala WA. Risks and prevention of nosocomial transmission of rare zoonotic diseases. Clin. Infect. Dis. 2001;32:446-456.
- 316. Chataigner D, Garnier R, Sans S, Efthymiou ML. [Acute accidental poisoning with hospital disinfectant. 45 cases of which 13 with fatal outcome]. Presse Med. 1991;20:741-3.
- 317. Hess JA, Molinari JA, Gleason MJ, Radecki C. Epidermal toxicity of disinfectants. Am. J. Dent. 1991;4:51-6.
- 318. Weber DJ, Rutala WA. Occupational risks associated with the use of selected disinfectants and sterilants. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:211-26.
- 319. Cokendolpher JC, Haukos JF. The Practical Application of Disinfection and Sterilization in Health Care Facilities. Chicago: American Hospital Association, 1996.
- 320. Rideout K, Teschke K, Dimich-Ward H, Kennedy SM. Considering risks to healthcare workers from glutaraldehyde alternatives in high-level dinfection. J. Hosp. Infect. 2005;59:4-11.
- 321. Oie S, Kamiya A. Assessment of and intervention for the misuse of aldehyde disinfectants in Japan. Infect . Control Hosp. Epidemiol. 2002;23:98-9.
- 322. American Conference of Governmental Industrial Hygienists (ACGIH). Threshold Limit Values for Chemical

Substances and Physical Agents and Biological Exposure Indices. Cincinnati: ACGIH, 2001.

- 323. Jordan SLP, Russo MR, Blessing RL, Grab LA. Glutaraldehyde safety: inactivation and disposal. Abstract. Am. J. Infect. Control 1997;25:154-55.
- 324. Jordan SL. The correct use of glutaraldehyde in the healthcare environment. Gastroenterol. Nurs. 1995;18:143-5.
- 325. Cheung HY, Brown MR. Evaluation of glycine as an inactivator of glutaraldehyde. J Pharmacy Pharmacol 1982;34:211-4.
- 326. Daschner F. The hospital and pollution: Role of the hospital epidemiologist in protecting the environment. In: Wenzel RP, ed. Prevention and control of nosocomial infections. Baltimore: Williams and Wilkins, 1997:595-605.
- 327. Rutala WA, Cole EC, Thomann CA, Weber DJ. Stability and bactericidal activity of chlorine solutions. Infect. Control Hosp. Epidemiol. 1998;19:323-7.
- 328. Rutala WA, Weber DJ. Uses of inorganic hypochlorite (bleach) in health-care facilities. Clin. Microbiol. Rev. 1997;10:597-610.
- 329. Dychdala GR. Chlorine and chlorine compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:135-157.
- 330. Rutala WA, Weber DJ. Principles of disinfecting patient-care items. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:133-49.
- 331. Luebbert P. Home care. In: Pfeiffer JA, ed. APIC text of infection control and epidemiology. Vol. 1. Washington: Association for Professionals in Infection control and epidemiology, 2000:44-7.
- 332. Parnes CA. Efficacy of sodium hypochlorite bleach and "alternative" products in preventing transfer of bacteria to and from inanimate surfaces. Environ. Health 1997;59:14-20.
- 333. Karapinar M, Gonul SA. Effects of sodium bicarbonate, vinegar, acetic and citric acids on growth and survival of *Yersinia enterocolitica*. Int. J. Food Microbiol. 1992;16:343-7.
- 334. McMurry LM, Oethinger M, Levy SB. Triclosan targets lipid synthesis. Nature 1998;394:531-2.
- 335. Moken MC, McMurry LM, Levy SB. Selection of multiple-antibiotic-resistant (mar) mutants of *Escherichia coli* by using the disinfectant pine oil: roles of the *mar* and *acr*AB loci. Antimicrob. Agents Chemother. 1997;41:2770-2.
- 336. Scott E, Bloomfield SF, Barlow CG. An investigation of microbial contamination in the home. J. Hyg. (Lond). 1982;89:279-93.
- 337. Rusin P, Orosz-Coughlin P, Gerba C. Reduction of faecal coliform, coliform and heterotrophic plate count bacteria in the household kitchen and bathroom by disinfection with hypochlorite cleaners. J. Appl. Microbiol. 1998;85:819-28.
- 338. Gilbert P, McBain AJ. Potential impact of increased use of biocides in consumer products on prevalence of antibiotic resistance. Clin Microbiol Reviews 2003;16:189-208.
- 339. Bueumer R, Bloomfield SF, Exner M, Fara G, Scott EA. The need for a home hygiene policy and guidelines on home hygiene. Ann. Ig. 1999;11:11-26.
- 340. <u>International Scientific Forum on Home Hygiene</u>. Accessed August 2016. [Current version of this document may differ from original.].
- Russell AD, Russell NJ. Biocides: activity, action and resistance. In: Hunter PA, Darby GK, Russell NJ, eds. Fifty years of antimicrobials: past perspectives and future trends. England: Cambridge University Press, 1995:327-65.
- 342. Russell AD. Bacterial resistance to disinfectants: Present knowledge and future problems. J. Hosp. Infect. 1998;43:S57-S68.
- 343. Russell AD. Plasmids and bacterial resistance to biocides. J. Appl. Microbiol. 1997;83:155-65.
- 344. Russell AD. Bacterial resistance to disinfectants: present knowledge and future problems. J. Hosp. Infect. 1998;43:S57-68.
- 345. Russell AD. Principles of antimicrobial activity and resistance. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:31-55.
- 346. McDonnell G, Russell AD. Antiseptics and disinfectants: activity, action, and resistance. Clin. Microbiol. Rev. 1999;12:147-79.
- 347. Gerba CP, Rusin P. Relationship between the use of antiseptics/disinfectants and the development of antimicrobial resistance. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:187-94.

- 348. Townsend DE, Ashdown N, Greed LC, Grubb WB. Transposition of gentamicin resistance to staphylococcal plasmids encoding resistance to cationic agents. J. Antimicrob. Chemother. 1984;14:115-24.
- 349. Brumfitt W, Dixson S, Hamilton-Miller JM. Resistance to antiseptics in methicillin and gentamicin resistant *Staphylococcus aureus*. Lancet 1985;1:1442-3.
- 350. Al-Masaudi SB, Day MJ, Russell AD. Sensitivity of methicillin-resistant *Staphylococcus aureus* strains to some antibiotics, antiseptics and disinfectants. J. Appl. Bacteriol. 1988;65:329-37.
- 351. Tennent JM, Lyon BR, Midgley M, Jones IG, Purewal AS, Skurray RA. Physical and biochemical characterization of the qacA gene encoding antiseptic and disinfectant resistance in *Staphylococcus aureus*. J. Gen. Microbiol. 1989;135:1-10.
- 352. Kaulfers PM, Laufs R. [Transmissible formaldehyde resistance in *Serratia marcescens*]. Zentralblatt fur Bakteriologie, Mikrobiologie und Hygiene 1 Abt Originale B, Hygiene 1985;181:309-19.
- 353. Tennent JM, Lyon BR, Gillespie MT, May JW, Skurray RA. Cloning and expression of *Staphylococcus aureus* plasmid-mediated quaternary ammonium resistance in *Escherichia coli*. Antimicrob. Agents Chemother. 1985;27:79-83.
- 354. Rutala WA, Stiegel MM, Sarubbi FA, Weber DJ. Susceptibility of antibiotic-susceptible and antibiotic-resistant hospital bacteria to disinfectants. Infect. Control Hosp. Epidemiol. 1997;18:417-21.
- 355. Anderson RL, Carr JH, Bond WW, Favero MS. Susceptibility of vancomycin-resistant enterococci to environmental disinfectants. Infect. Control Hosp. Epidemiol. 1997;18:195-9.
- 356. Sakagami Y, Kajimura K. Bactericidal activities of disinfectants against vancomycin-resistant enterococci. J. Hosp. Infect. 2002;50:140-4.
- 357. Sehulster LM, Anderson RL. Susceptibility of glycopeptide-intermediate resistant *Staphylococcus aureus* (GISA) to surface disinfectants, hand washing chemicals, and a skin antiseptic. Abstract Y-3. 98th General Meeting of American Society for Microbiology, May, 1998:547.
- 358. Rutala WA, Weber DJ, Gergen MF. Studies on the disinfection of VRE-contaminated surfaces. Infect. Control Hosp. Epidemiol. 2000;21:548.
- 359. Byers KE, Durbin LJ, Simonton BM, Anglim AM, Adal KA, Farr BM. Disinfection of hospital rooms contaminated with vancomycin-resistant *Enterococcus faecium*. Infect. Control Hosp. Epidemiol. 1998;19:261-4.
- 360. Carling PC, Briggs JL, Perkins J, Highlander D. Improved cleaning of patient rooms using a new targeting method. Clin Inf Dis 2006;42:385-8.
- 361. Russell AD, Suller MT, Maillard JY. Do antiseptics and disinfectants select for antibiotic resistance? J. Med. Microbiol. 1999;48:613-5.
- 362. Russell AD. Bacterial adaptation and resistance to antiseptics, disinfectants and preservatives is not a new phenomenon. J. Hosp. Infect. 2004;57:97-104.
- 363. Levy SB. The challenge of antibiotic resistance. Scientific Am. 1998;278:46-53.
- 364. Jones RD, Jampani HB, Newman JL, Lee AS. Triclosan: a review of effectiveness and safety in health care settings. Am. J. Infect. Control 2000;28:184-96.
- 365. Russell AD, McDonnell G. Concentration: a major factor in studying biocidal action. J. Hosp. Infect. 2000;44:1-3.
- 366. Russell AD, Maillard JY. Reaction and response-relationship between antibiotic resistance and resistance to antiseptics and disinfectants. Am. J. Infect. Control 2000;28:204-6.
- 367. Murtough SM, Hiom SJ, Palmer M, Russell AD. Biocide rotation in the healthcare setting: is there a case for policy implementation? J. Hosp. Infect. 2001;48:1-6.
- 368. Murtough SM, Hiom SJ, Palmer M, Russell AD. A survey of rotational use of biocides in hospital pharmacy aseptic units. J. Hosp. Infect. 2002;50:228-31.
- 369. Gebel J, Sonntag H-G, Werner H-P, Vavata V, Exner M, Kistemann T. The higher disinfectant resistance of nosocomial isolates of *Klebsiella oxytoca*: How reliable are indicator organisms in disinfectant testing? J. Hosp. Infect. 2002;50:309-11.
- 370. Ruden H, Daschner F. Should we routinely disinfect floors? J. Hosp. Infect. 2002;51:309.
- 371. Rutala WA, DJ W. Should we routinely disinfect floors? Reply to Professor F. Daschner. J. Hosp. Infect. 2002;51:309-11.
- 372. Rutala WA, Weber DJ. The benefits of surface disinfection. Am. J. Infect. Control 2004;32:226-31.
- Dettenkofer M, Wenzler S, Amthor S, Antes G, Motschall E, Daschner FD. Does disinfection of environmental surfaces influence nosocomial infection rates? A systematic review. Am. J. Infect. Control 2004;32:84-9.

- 374. Daschner F, Schuster A. Disinfection and the prevention of infectious disease-no adverse effects? Am. J. Infect. Control 2004;32:224-5.
- 375. Cozad A, Jones RD. Disinfection and the prevention of infectious disease. Am. J. Infect. Control 2003;31:243-54.
- 376. Favero MS, Bond WW. Chemical disinfection of medical and surgical materials. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1991:617-41.
- 377. Rheinbaben FV, Schunemann S, Grob T, Wolff MH. Transmission of viruses via contact in a household setting: experiments using bacteriophage OX174 as a model virus. J. Hosp. Infect. 2000;46:61-66.
- 378. Rutala WA, Weber DJ. Surface disinfection: should we do it? J. Hosp. Infect. 2001;48 (supplement A):S64-S68.
- 379. Ayliffe GAJ, Collins DM, Lowbury EJL. Cleaning and disinfection of hospital floors. Brit. Med. J. 1966;2:442-5.
- 380. Ayliffe GA, Collins BJ, Lowbury EJ, Babb JR, Lilly HA. Ward floors and other surfaces as reservoirs of hospital infection. J. Hyg. (Lond). 1967;65:515-36.
- 381. Exner M, Vacata V, Hornei B, Dietlein E, Gebel J. Household cleaning and surface disinfection: New insights and strategies. J. Hosp. Infect. 2004;56 (suppl):S70-S75.
- 382. Dharan S, Mourouga P, Copin P, Bessmer G, Tschanz B, Pittet D. Routine disinfection of patients' environmental surfaces. Myth or reality? J. Hosp. Infect. 1999;42:113-7.
- 383. Engelhart S KL, Glasmacher A, Fischnaller E, Marklein G, Exner M. *Pseudomonas aeruginosa* outbreak in a haematology-oncology unit associated with contaminated surface cleaning equipment. J. Hosp. Infect. 2002;52:93-98.
- 384. Denton M, Wilcox MH, Parnell P, et al. Role of environmental cleaning in controlling an outbreak of *Acinetobacter baumanni* on a neurosurgical intensive care unit. J. Hosp. Infect. 2004;56:106-10.
- 385. Barker J, Vipond IB, Bloomfield SF. Effects of cleaning and disinfection in reducing the spread of Norovirus contamination via environmental surfaces. J. Hosp. Infect. 2004;58:42-9.
- 386. Garner JS. Guideline for isolation precautions in hospitals. The Hospital Infection Control Practices Advisory Committee. Infect. Control Hosp. Epidemiol. 1996;17:53-80.
- 387. Maki DG, Alvarado CJ, Hassemer CA, Zilz MA. Relation of the inanimate hospital environment to endemic nosocomial infection. N. Engl. J. Med. 1982;307:1562-6.
- 388. Daschner F, Rabbenstein G, Langmaack H. [Surface decontamination in the control of hospital infections: comparison of different methods (author's transl]. Dtsch. Med. Wochenschr. 1980;105:325-9.
- 389. Danforth D, Nicolle LE, Hume K, Alfieri N, Sims H. Nosocomial infections on nursing units with floors cleaned with a disinfectant compared with detergent. J. Hosp. Infect. 1987;10:229-35.
- 390. Smith TL, Iwen PC, Olson SB, Rupp ME. Environmental contamination with vancomycin-resistant enterococci in an outpatient setting. Infect. Control Hosp. Epidemiol. 1998;19:515-8.
- Boyce JM, Potter-Bynoe G, Chenevert C, King T. Environmental contamination due to methicillin-resistant Staphylococcus aureus: possible infection control implications. Infect. Control Hosp. Epidemiol. 1997;18:622-7.
- 392. Bonten MJM, Hayden MJ, Nathan C, et al. Epidemiology of colonisation of patients and environment with vancomycin-resistant enterococci. Lancet 1996;348:1615-9.
- 393. Hardy KJ, Oppenheim BA, Gossain S, Gao F, Hawkey PM. A study of the relationship between environmental contamination with methicillin-resistant *Staphylococcus aureus* (MRSA) and patients' acquisition of MRSA. Infect Control Hosp Epidemiol 2006;27:127-32.
- 394. Hota B. Contamination, disinfection, and cross-contamination: Are hospital surfaces reservoirs for nosocomial infection? Clin. Infect. Dis. 2004;39:1182-9.
- 395. Neely AN, Maley MP. Survival of enterococci and staphylococci on hospital fabrics and plastic. J. Clin. Microbiol. 2000;38:724-6.
- 396. Wendt C, Wiensenthal B, Dietz E, Ruden H. Survival of enterococci on dry surfaces. J. Clin. Microbiol. 1998;36:3734-6.
- 397. Neely AN, Maley MP. The 1999 Lindberg award. 3% hydrogen peroxide for the gram-positive disinfection of fabrics. J. Burn Care Rehabil. 1999;20:471-7.
- 398. Griffith CJ, Cooper RA, Gilmore J, Davies C, Lewis M. An evaluation of hospital cleaning regimes and standards. J. Hosp. Infect. 2000;45:19-28.
- 399. Tiller JC, Liao CJ, Lewis K, Klibanov AM. Designing surfaces that kill bacteria on contact. Proc. Natl. Acad. Sci. 2001;98:5981-5.
- 400. Rutala WA, Weber DJ. New disinfection and sterilization methods. Emerg. Inf. Dis. 2001;7:348-53.

- 401. Whitby JL, Rampling A. *Pseudomonas aeruginosa* contamination in domestic and hospital environments. Lancet 1972;1:15-7.
- 402. Scott E, Bloomfield SF. The survival and transfer of microbial contamination via cloths, hand and utensils. J. Appl. Bacteriol. 1990;68:271-8.
- 403. Scott E, Bloomfield SF. Investigations of the effectiveness of detergent washing, drying and chemical disinfection on contamination of cleaning cloths. J. Appl. Bacteriol. 1990;68:279-83.
- 404. Rutala WA, Cole EC. Antiseptics and disinfectants--safe and effective? Infect. Control 1984;5:215-8.
- 405. Oie S, Huang Y, Kamiya A, Konishi H, Nakazawa T. Efficacy of disinfectants against biofilm cells of methicillin-resistant *Staphylococcus aureus*. Microbios 1996;85:223-30.
- 406. Sartor C, Jacomo V, Duvivier C, Tissot-Dupont H, Sambuc R, Drancourt M. Nosocomial *Serratia marcescens* infections associated with extrinsic contamination of a liquid nonmedicated soap. Infect. Control Hosp. Epidemiol. 2000;21:196-9.
- 407. Reiss I, Borkhardt A, Fussle R, Sziegoleit A, Gortner L. Disinfectant contaminated with *Klebsiella oxytoca* as a source of sepsis in babies. Lancet 2000;356:310.
- 408. O'Rourke E, Runyan D, O'Leary J, Stern J. Contaminated iodophor in the operating room. Am. J. Infect. Control 2003;31:255-6.
- 409. Chuanchuen R, Karkhoff-Schweizer RR, Schweizer HP. High-level triclosan resistance in *Pseudomonas aeruginosa* is solely a result of efflux. Am. J. Infect. Control 2003;31:124-7.
- 410. Newman KA, Tenney JH, Oken HA, Moody MR, Wharton R, Schimpff SC. Persistent isolation of an unusual *Pseudomonas* species from a phenolic disinfectant system. Infect. Control 1984;5:219-22.
- 411. Bean HS. Types and characteristics of disinfectants. J. Appl. Bacteriol. 1967;30:6-16.
- 412. Russell AD, Hugo WB, Ayliffe GAJ. Principles and Practice of Disinfection, Preservation and Sterilization. Oxford, England: Blackwell Scientific Publications, 1999.
- 413. Russell AD. Factors influencing the efficacy of germicides. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington DC: Association for Professionals in Infection Control and Epidemiology, 2004:162-70.
- 414. Gillis RJ, Schmidt WC. Scanning electron microscopy of spores on inoculated product surfaces. MD 1983:46-9.
- 415. Favero MS, Petersen NJ, Carson LA, Bond WW, Hindman SH. Gram-negative water bacteria in hemodialysis systems. Health Lab. Sci. 1975;12:321-34.
- 416. Rutala WA, Cole EC. Ineffectiveness of hospital disinfectants against bacteria: a collaborative study. Infect. Control 1987;8:501-6.
- 417. Favero MS. Naturally occurring microrganisms and their resistance to physical and chemical agents. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington DC: Association for Professionals in Infection Control and Epidemiology, 2004:1-14.
- 418. Lee DH, Miles RJ, Perry BF. The mycoplasmacidal properties of sodium hypochlorite. J. Hyg. (Lond). 1985;95:243-53.
- 419. Scott GH, Williams JC. Susceptibility of *Coxiella burnetii* to chemical disinfectants. Ann. N. Y. Acad. Sci. 1990;590:291-6.
- 420. Russell AD. Factors influencing the efficacy of antimicrobial agents. In: Russell AD, Hugo WB, Ayliffe GAJ, eds. Principles and practice of disinfection, preservation and sterilization. Oxford: Blackwell Science, 1999:95-123.
- 421. Rutala WA. Selection and use of disinfectants in healthcare. In: Mayhall CG, ed. Infect. Control and Hosp. Epidemiol. Philadelphia: Lippincott Williams & Wilkins, 1999:1161-87.
- 422. Lewis DL, Arens M. Resistance of microorganisms to disinfection in dental and medical devices. Nat. Med. 1995;1:956-8.
- 423. Muscarella LF. Sterilizing dental equipment. Nat. Med. 1995;1:1223-5.
- 424. Abbott CF, Cockton J, Jones W. Resistance of crystalline substances to gas sterilization. J. Pharm. Pharmacol. 1956;8:709-20.
- 425. Doyle JE, Ernst RR. Resistance of *Bacillus subtilis* var. *niger* spores occluded in water-insoluble crystals to three sterilization agents. Appl. Microbiol. 1967;15:726-30.
- 426. Jacobs P. Cleaning: Principles, methods and benefits. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:165-81.
- 427. Gorham RA, Jacobs P, Roberts CG. Laboratory artifacts due to protein and salt crystals on the inactivation of *Bacillus stearothermophilus*. J. Hosp. Infect. 1998;40:abstract P.9.2.2.

- 428. Cole EC, Rutala WA, Carson JL, Alfano EM. *Pseudomonas* pellicle in disinfectant testing: electron microscopy, pellicle removal, and effect on test results. Appl. Environ. Microbiol. 1989;55:511-3.
- 429. Anderson RL, Holland BW, Carr JK, Bond WW, Favero MS. Effect of disinfectants on pseudomonads colonized on the interior surface of PVC pipes. Am. J. Public Health 1990;80:17-21.
- 430. Anderson RL, Vess RW, Carr JH, Bond WW, Panlilio AL, Favero MS. Investigations of intrinsic *Pseudomonas cepacia* contamination in commercially manufactured povidone-iodine. Infect. Control Hosp. Epidemiol. 1991;12:297-302.
- 431. LeChevallier MW, Cawthon CD, Lee RG. Inactivation of biofilm bacteria. Appl. Environ. Microbiol. 1988;54:2492-9.
- 432. LeChevallier MW, Cawthon CD, Lee RG. Factors promoting survival of bacteria in chlorinated water supplies. Appl. Environ. Microbiol. 1988;54:649-54.
- 433. Costerton JS, Steward PS, Greenberg EP. Bacterial biofilms: a common cause of persistent infections. Science 1999;284:1318-22.
- 434. Donlan RM, Costerton JW. Biofilms: Survival mechanisms of clinically relevant mirocorganisms. Clin. Microbiol. Rev. 2002;15:167-93.
- 435. Dunne WM. Bacterial adhesion: Seen any good biofilms lately? Clin. Microbiol. Rev. 2002;15:155-66.
- 436. Vickery K, Pajkos A, Cossart Y. Removal of biofilm from endoscopes: Evaluation of detergent efficiency. Am. J. Infect. Control 2004;32:170-6.
- 437. Marion K, Freney J, James G, Bergeron E, Renaud FNR, Costerton JW. Using an efficient biofilm detaching agent: An essential step for the improvement of endoscope reprocessing protocols. J. Hosp. Infect. 2006;In press.
- 438. Marion-Ferey K, Pasmore M, Stoodley P, Wilson S, Husson GP, Costerton JW. Biofilm removal from silicone tubing: an assessment of the efficacy of dialysis machine decontamination procedures using an in vitro model. J. Hosp. Infect. 2003;53:64-71.
- 439. Brown ML, Aldrich HC, Gauthier JJ. Relationship between glycocalyx and povidone-iodine resistance in *Pseudomonas aeruginosa* (ATCC 27853) biofilms. Appl. Environ. Microbiol. 1995;61:187-93.
- 440. Price D, Ahearn DG. Incidence and persistence of *Pseudomonas aeruginosa* in whirlpools. J. Clin. Microbiol. 1988;26:1650-4.
- 441. Anonymous. Dental Unit Waterlines: Approaching the Year 2000. ADA Council on Scientific Affairs. JADA 1999;130:1653-64.
- 442. Donlan RM. Biofilms: a source of infection? In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:219-26.
- 443. Loukili NH, Zink E, Grandadam S, Bientz M, Meunier O. Effectiveness of detergent-disinfecting agents on *Escherichia coli* 54127 biofilm. J. Hosp. Infect. 2004;57:175-8.
- 444. Johansen C, Falholt P, Gram L. Enzymatic removal and disinfection of bacterial biofilms. Appl. Environ. Microbiol. 1997;63:3724-8.
- 445. Reichert M. Preparation of supplies for terminal sterilization. In: Reichert M, Young JH, eds. Sterilization technology for the health care facility. Gaithersburg, MD: Aspen Publication, 1997:36-50.
- 446. Miller CH, Riggen SD, Sheldrake MA, Neeb JM. Presence of microorganisms in used ultrasonic cleaning solutions. Am. J. Dent. 1993;6:27-31.
- 447. Jatzwauk L, Schone H, Pietsch H. How to improve instrument disinfection by ultrasound. J. Hosp. Infect. 2001;48 (Supple):S80-S83.
- 448. Richburg FA, Reidy JJ, Apple DJ, Olson RJ. Sterile hypopyon secondary to ultrasonic cleaning solution. J. Cataract Refract. Surg. 1986;12:248-51.
- 449. Schultz JK. Decontamination alternative. Infect. Control Hosp. Epidemiol. 1990;11:8-9.
- 450. Rutala WA, Shafer KM. General information on cleaning, disinfection, and sterilization. In: Pfeiffer JA, ed. APIC infection control and applied epidemiology: principles and practice,. St. Louis: Mosby, 1996:15.1-15.17.
- 451. Leonard DL, Mills SE. Comparison of automated instrument cleaning: preliminary results. Infect. Control Steril. Technol. 1997:20-23, 26-28.
- 452. Ransjo U, Engstrom L, Hakansson P, et al. A test for cleaning and disinfection processes in a washerdisinfector. APMIS 2001;109:299-304.
- 453. American Society for Hospital Central Service Personnel. Training manual for central service technicians. Chicago: American Hospital Association, 2001:1-271.
- 454. Ninemeier JD. Central service technical manual. Chicago: International Association of Healthcare Central

Service Materiel Management, 1998.

- 455. Reichert M, Young JH. Sterilization technology for the health care facility. Gaithersburg: Aspen Publication, 1997:307.
- 456. Vesley D, Norlien KG, Nelson B, Ott B, Streifel AJ. Significant factors in the disinfection and sterilization of flexible endoscopes. Am. J. Infect. Control 1992;20:291-300.
- 457. Roberts CG. Studies on the bioburden on medical devices and the importance of cleaning. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:63-9.
- 458. Baxter RL, Baxter HC, Campbell GA, et al. Quantitaive analysis of residual protein contamination on reprocessed surgical instruments. J. Hosp. Infect. 2006;63:439-44.
- 459. Murdoch H, Taylor D, Dickinson J, et al. Surface decontamination of surgical instruments: An ongoing dilemma. J. Hosp. Infect. 2006;63:432-8.
- 460. Alfa MJ, Nemes R. Manual versus automated methods for cleaning reusable accessory devices used for minimally invasine surgical procedures. J. Hosp. Infect. 2004;58:50-8.
- 461. Alfa MJ, Nemes R, Olson N, Mulaire A. Manual methods are suboptimal compared with automated methods for cleaning of single-use biopsy forceps. Infect Control Hosp Epidemiol 2006;27:841-6.
- 462. Lee CH, Cheng SM, Humar A, et al. Acute febrile reactions with hypotension temporally associated with the introduction of a concentrated bioenzyme preparation in the cleaning and sterilization process of endomyocardial bioptones. Infect. Control Hosp. Epidemiol. 2000;21:102.
- 463. Hutchisson B, LeBlanc C. The truth and consequences of enzymatic detergents. Gastroenterol. Nurs. 2005;28:372-6.
- 464. Zuhlsdorf B EM, Floss H, Martiny H,. Cleaning efficacy of nine different cleaners in a washer-disinfector designed for flexible endoscopes. J. Hosp. Infect. 2002;52:206-11.
- 465. Merritt K, Hitchins VM, Brown SA. Safety and cleaning of medical materials and devices. J. Biomed. Mater. Res. 2000;53:131-6.
- 466. Babb JR, Bradley CR. Endoscope decontamination: where do we go from here? J. Hosp. Infect. 1995;30:543-51.
- 467. Zuhlsdorf B, Floss H, Martiny H. Efficacy of 10 different cleaning processes in a washer-disinfector for flexible endoscopes. J. Hosp. Infect. 2004;56:305-11.
- 468. Alfa MJ, Jackson M. A new hydrogen peroxide-based medical-device detergent with germicidal properties: Comparison with enzymatic cleaners. Am. J. Infect. Control 2001;29:168-77.
- 469. Alfa MJ, DeGagne P, Olson N, Puchalski T. Comparison of ion plasma, vaporized hydrogen peroxide and 100% ethylene oxide sterilizers to the 12/88 ethylene oxide gas sterilizer. Infect. Control Hosp. Epidemiol. 1996;17:92-100.
- 470. Alfa MJ. Flexible endoscope reprocessing. Infect. Control Steril. Technol. 1997;3:26-36.
- 471. Alfa MJ, Degagne P, Olson N. Worst-case soiling levels for patient-used flexible endoscopes before and after cleaning. Am. J. Infect. Control 1999;27:392-401.
- 472. Rutala WA, Weber DJ. Low-temperature sterilization technology: Do we need to redefine sterilization? Infect. Control Hosp. Epidemiol. 1996;17:89-91.
- 473. Dancer SJ. How do we assess hospital cleaning? A proposal for microbiological standards for surface hygiene in hospitals. J. Hosp. Infect. 2004;56:10-5.
- 474. Pfeifer M. Standardised test soil blood 1:Composition, preparation, application. Zentr. Steril. 1998;6:381-5.
- 475. Pfeifer M. Blood as a soil on surgical instruments: Chemical profile, cleaning, detection. Zentr. Steril. 1998;6:304-10.
- 476. Fengier TW, Pahike H, Bisson S, Michels W. Are processed surgical instruments free of protein? Zentr. Steril. 2001;9:20-32.
- 477. Takashina M. Application of a bioluminescent method for checking cleaning results. Zentr. Steril. 2001;9:248-58.
- 478. Lipscomb IP, Sihota AK, Botham M, Harris KL, Keevil CW. Rapid method for the sensitive detection of protein contamination on surgical instruments. J. Hosp. Infect. 2006;62:141-8.
- 479. Malik RE, Cooper RA, Griffith CJ. Use of audit tools to evaluate the efficacy of cleaning systems in hospitals. Am. J. Infect. Control 2003;31:181-7.
- 480. Hansen KS. Occupational dermatoses in hospital cleaning women. Contact Dermatitis 1983;9:343-51.
- 481. Melli MC, Giorgini S, Sertoli A. Sensitization from contact with ethyl alcohol. Contact Dermatitis 1986;14:315.
- 482. Spaulding EH. Alcohol as a surgical disinfectant. AORN J. 1964;2:67-71.

- 483. Morton HE. The relationship of concentration and germicidal efficiency of ethyl alcohol. Ann N.Y. Acad. Sci. 1950;53:191-96.
- 484. Ali Y, Dolan MJ, Fendler EJ, Larson EL. Alcohols. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:229-54.
- 485. Morton HE. Alcohols. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1983:225-239.
- 486. Sykes G. The influence of germicides on the dehydrogenases of *Bact. coli*. Part I. The succinic acid dehydrogenase of *Bact. coli*. J. Hyg. (Camb) 1939;39:463-69.
- 487. Dagley S, Dawes EA, Morrison GA. Inhibition of growth of *Aerobacter aerogenes*: the mode of action of phenols, alcohols, acetone and ethyl acetate. J. Bacteriol. 1950;60:369-78.
- 488. Tilley FW, Schaffer JM. Relation between the chemical constitution and germicidal activity of the monohydric alcohols and phenols. J. Bacteriol. 1926;12:303-9.
- 489. Coulthard CE, Sykes G. The germicidal effect of alcohol with special reference to its action on bacterial spores. Pharmaceutical J. 1936;137:79-81.
- 490. Tyler R, Ayliffe GA. A surface test for virucidal activity of disinfectants: preliminary study with herpes virus. J. Hosp. Infect. 1987;9:22-9.
- 491. Kurtz JB, Lee TW, Parsons AJ. The action of alcohols on rotavirus, astrovirus and enterovirus. J. Hosp. Infect. 1980;1:321-5.
- 492. Smith CR. Alcohol as a disinfectant against the tubercle bacillus. Public Health Rep. 1947;62:1285-95.
- 493. Kruse RH, Green TD, Chambers RC, Jones MW. Disinfection of aerosolized pathogenic fungi on laboratory surfaces. 1. Tissue phase. Appl. Microbiol. 1963;11:436-45.
- 494. Kruse RH, Green TD, Chambers RC, Jones MW. Disinfection of aerosolized pathogenic fungi on laboratory surfaces. II Culture phase. Appl. Microbiol. 1964;12:155-60.
- 495. Connor CG, Hopkins SL, Salisbury RD. Effectivity of contact lens disinfection systems against *Acanthamoeba culbertsoni*. Optom. Vis. Sci. 1991;68:138-41.
- 496. Turner NA, Russell AD, Furr JR, Lloyd D. *Acanthamoeba* spp., antimicrobial agents and contact lenses. Sci. Prog. 1999;82:1-8.
- 497. Nye RN, Mallory TB. A note on the fallacy of using alcohol for the sterilization of surgical instruments. Boston Med. Surg. J. 1923;189:561-3.
- 498. Frobisher M, Sommermeyer L, Blackwell MJ. Studies on disinfection of clinical thermometers. I. Oral thermometers. Appl. Microbiol. 1973;1:187-94.
- 499. Sommermeyer L, Frobisher M. Laboratory studies on disinfection of rectal thermometers. Nurs. Res. 1973;2:85-9.
- 500. Singh D, Kaur H, Gardner WG, Treen LB. Bacterial contamination of hospital pagers. Infect . Control Hosp. Epidemiol. 2002;23:274-6.
- 501. Embil JM, Zhanel GG, Plourde J, Hoban D. Scissors: A potential source of nosocomial infection. Infect . Control Hosp. Epidemiol. 2002;23:147-51.
- 502. Zachary KC, Bayne PS, Morrison VJ, Ford DS, Silver LC, Hooper DC. Contamination of gowns, gloves, and stethoscopes with vancomycin-resistant enterococci. Infect . Control Hosp. Epidemiol. 2001;22:560-4.
- 503. Babb JR, Bradley CR, Deverill CE, Ayliffe GA, Melikian V. Recent advances in the cleaning and disinfection of fibrescopes. J. Hosp. Infect. 1981;2:329-40.
- 504. Garcia de Cabo A, Martinez Larriba PL, Checa Pinilla J, Guerra Sanz F. A new method of disinfection of the flexible fibrebronchoscope. Thorax 1978;33:270-2.
- 505. Elson CO, Hattori K, Blackstone MO. Polymicrobial sepsis following endoscopic retrograde cholangiopancreatography. Gastroenterology 1975;69:507-10.
- 506. Weber DJ, Wilson MB, Rutala WA, Thomann CA. Manual ventilation bags as a source for bacterial colonization of intubated patients. Am. Rev. Respir. Dis. 1990;142:892-4.
- 507. Cavagnolo RZ. Inactivation of herpesvirus on CPR manikins utilizing a currently recommended disinfecting procedure. Infect. Control 1985;6:456-8.
- 508. Ohara T, Itoh Y, Itoh K. Ultrasound instruments as possible vectors of staphylococcal infection. J. Hosp. Infect. 1998;40:73-7.
- 509. Talbot GH, Skros M, Provencher M. 70% alcohol disinfection of transducer heads: experimental trials. Infect. Control 1985;6:237-9.
- 510. Platt R, Lehr JL, Marino S, Munoz A, Nash B, Raemer DB. Safe and cost-effective cleaning of pressuremonitoring transducers. Infect. Control Hosp. Epidemiol. 1988;9:409-16.

- 511. Beck-Sague CM, Jarvis WR. Epidemic bloodstream infections associated with pressure transducers: a persistent problem. Infect. Control Hosp. Epidemiol. 1989;10:54-9.
- 512. Chronister CL, Russo P. Effects of disinfecting solutions on tonometer tips. Optom. Vis. Sci. 1990;67:818-21.
- 513. Lingel NJ, Coffey B. Effects of disinfecting solutions recommended by the Centers for Disease Control on Goldmann tonometer biprisms. J. Am. Optom. Assoc. 1992;63:43-8.
- 514. Soukiasian SH, Asdourian GK, Weiss JS, Kachadoorian HA. A complication from alcohol-swabbed tonometer tips. Am. J. Ophthalmol. 1988;105:424-5.
- 515. Jakobsson SW, Rajs J, Jonsson JA, Persson H. Poisoning with sodium hypochlorite solution. Report of a fatal case, supplemented with an experimental and clinico-epidemiological study. Am. J. Forensic Med. Pathol. 1991;12:320-7.
- 516. Heidemann SM, Goetting MG. Treatment of acute hypoxemic respiratory failure caused by chlorine exposure. Pediatr. Emerg. Care 1991;7:87-8.
- 517. Hoy RH. Accidental systemic exposure to sodium hypochlorite (Clorox) during hemodialysis. Am. J. Hosp. Pharm. 1981;38:1512-4.
- 518. Landau GD, Saunders WH. The effect of chlorine bleach on the esophagus. Arch. Otolaryngol. 1964;80:174-6.
- 519. French RJ, Tabb HG, Rutledge LJ. Esophageal stenosis produced by ingestion of bleach: report of two cases. South. Med. J. 1970;63:1140-4.
- 520. Ward MJ, Routledge PA. Hypernatraemia and hyperchloraemic acidosis after bleach ingestion. Hum. Toxicol. 1988;7:37-8.
- 521. Ingram TA. Response of the human eye to accidental exposure to sodium hypochlorite. J Endodontics 1990;16:235-8.
- 522. Haag JR, Gieser RG. Effects of swimming pool water on the cornea. JAMA 1983;249:2507-8.
- 523. Mrvos R, Dean BS, Krenzelok EP. Home exposures to chlorine/chloramine gas: review of 216 cases. South. Med. J. 1993;86:654-7.
- 524. Reisz GR, Gammon RS. Toxic pneumonitis from mixing household cleaners. Chest 1986;89:49-52.
- 525. Gapany-Gapanavicius M, Yellin A, Almog S, Tirosh M. Pneumomediastinum. A complication of chlorine exposure from mixing household cleaning agents. JAMA 1982;248:349-50.
- 526. Hoffman PN, Death JE, Coates D. The stability of sodium hypochlorite solutions. In: Collins CH, Allwood MC, Bloomfield SF, Fox A, eds. Disinfectants: their use and evaluation of effectiveness. London: Academic Press, 1981:77-83.
- 527. Gamble MR. Hazard: formaldehyde and hypochlorites. Lab. Anim. 1977;11:61.
- 528. Helms C, Massanari R, Wenzel R, et al. Control of epidemic nosocomial legionellosis: a 5 year progress report on continuous hyperchlorination of a water distribution system. Abstracts of 27th Interscience Conference of Antimicrobial Agents and Chemotherapy, 1987:349, p.158.
- 529. Environmental Protection Agency. R.E.D. Facts sodium and calcium hypochlorite salts. 1991.
- 530. Coates D. Comparison of sodium hypochlorite and sodium dichloroisocyanurate disinfectants: neutralization by serum. J. Hosp. Infect. 1988;11:60-7.
- 531. Coates D. A comparison of sodium hypochlorite and sodium dichloroisocyanurate products. J. Hosp. Infect. 1985;6:31-40.
- 532. Coates D, Wilson M. Use of sodium dichloroisocyanurate granules for spills of body fluids. J. Hosp. Infect. 1989;13:241-51.
- 533. Bloomfield SF, Uso EE. The antibacterial properties of sodium hypochlorite and sodium dichloroisocyanurate as hospital disinfectants. J. Hosp. Infect. 1985;6:20-30.
- 534. Coates D. An evaluation of the use of chlorine dioxide (Tristel One-Shot) in an automated washer/disinfector (Medivator) fitted with a chlorine dioxide generator for decontamination of flexible endoscopes. J. Hosp. Infect. 2001;48:55-65.
- 535. Isomoto H, Urata M, Kawazoe K, et al. Endoscope disinfection using chlorine dioxide in an automated washer-disinfector. J. Hosp. Infect. 2006;63:298-305.
- 536. Sampson MN MA. Not all super-oxidized waters are the same. J. Hosp. Infect. 2002;52:227-8.
- 537. Selkon JB, Babb JR, Morris R. Evaluation of the antimicrobial activity of a new super-oxidized water, Sterilox®, for the disinfection of endoscopes. J. Hosp. Infect. 1999;41:59-70.
- 538. Fraise AP. Choosing disinfectants. J. Hosp. Infect. 1999;43:255-64.
- 539. Tanaka H, Hirakata Y, Kaku M, et al. Antimicrobial activity of superoxidized water. J. Hosp. Infect. 1996;34:43-9.

- 540. Tanaka N, Fujisawa T, Daimon T, Fujiwara K, Yamamoto M, Abe T. The use of electrolyzed solutions for the cleaning and disinfecting of dialyzers. Artif. Organs 2000;24:921-8.
- 541. Williams ND, Russell AD. The effects of some halogen-containing compounds on *Bacillus subtilis* endospores. J. Appl. Bacteriol. 1991;70:427-36.
- 542. Babb JR, Bradley CR, Ayliffe GAJ. Sporicidal activity of glutaraldehydes and hypochlorites and other factors influencing their selection for the treatment of medical equipment. J. Hosp. Infect. 1980;1:63-75.
- 543. Brown DG, Skylis TP, Fekety FR. Comparison of chemical sterilant/disinfectant solutions against spores of *Clostridium difficile*. Abstracts of the American Society for Microbiology, 1983:Q39,267.
- 544. Grant D, Venneman M, Burns RM. Mycobactericidal activity of Alcide an experimental liquid sterilant. Abstracts of the Annual Meeting of the American Society of Microbiology, 1982: Q101,226.
- 545. Korich DG, Mead JR, Madore MS, Sinclair NA, Sterling CR. Effects of ozone, chlorine dioxide, chlorine, and monochloramine on *Cryptosporidium parvum* oocyst viability. Appl. Environ. Microbiol. 1990;56:1423-8.
- 546. Griffiths PA, Babb JR, Fraise AP. Mycobactericidal activity of selected disinfectants using a quantitative suspension test. J. Hosp. Infect. 1999;41:111-21.
- 547. Centers for Disease Control. Bacteremia associated with reuse of disposable hollow-fiber hemodialyzers. MMWR 1986;35:417-8.
- 548. Bloomfield SF, Miller EA. A comparison of hypochlorite and phenolic disinfectants for disinfection of clean and soiled surfaces and blood spillages. J. Hosp. Infect. 1989;13:231-9.
- 549. Shetty N, Srinivasan S, Holton J, Ridgway GL. Evaluation of microbicidal activity of a new disinfectant: Sterilox® 2500 against *Clostridium difficile* spores, *Helicobacter pylori*, vancomycin resistant *Enterococcus* species, *Candida albicans* and several *Mycobacterium* species. J. Hosp. Infect. 1999;41:101-5.
- 550. Urata M, Isomot H, Murase K, et al. Comparison of the microbicidal activities of superoxidized and ozonated water in the disinfection of endoscopes. J Intern Med Res 2003;31:299-306.
- 551. Tsuji S, Kawano S, Oshita M, et al. Endoscope disinfection using acidic electrolytic water. Endoscopy 1999;31:528-35.
- 552. Lee JH, Rhee PL, Kim JH, et al. Efficacy of electolyzed acid water in reprocessing patient-used flexible upper endoscopes: Comparison with 2% alkaline glutaraldehyde. J. Gastroenterol. Hepatol. 2004;19:897-904.
- 553. Centers for Disease Control. Acquired immune deficiency syndrome (AIDS): precautions for clinical and laboratory staffs. MMWR 1982;31:577-80.
- 554. Garner JS, Simmons BP. Guideline for isolation precautions in hospitals. Infect. Control 1983;4:245-325.
- 555. Van Bueren J, Simpson RA, Salman H, Farrelly HD, Cookson BD. Inactivation of HIV-1 by chemical disinfectants: sodium hypochlorite. Epidemiol. Infect. 1995;115:567-79.
- 556. Coates D. Disinfection of spills of body fluids: how effective is a level of 10,000 ppm available chlorine? J. Hosp. Infect. 1991;18:319-22.
- 557. Chitnis V, Chitnis S, Patil S, Chitnis D. Practical limitations of disinfection of body fluid spills with 10,000 ppm sodium hypochlorite (NaOCI). Am. J. Infect. Control 2004;32:306-8.
- 558. Anonymous. Recommendations for decontaminating manikins used in cardiopulmonary resuscitation training, 1983 update. Infect . Control 1984;5:399-401.
- 559. Centers for Disease Control. Use of bleach for disinfection of drug injection equipment. MMWR 1993;42:418-9.
- 560. Shapshak P, McCoy CB, Rivers JE, et al. Inactivation of human immunodeficiency virus-1 at short time intervals using undiluted bleach. J. Acquir. Immune Defic. Syndr. 1993;6:218-9.
- 561. Shapshak P, McCoy CB, Shah SM, et al. Preliminary laboratory studies of inactivation of HIV-1 in needles and syringes containing infected blood using undiluted household bleach. J. Acquir. Immune Defic. Syndr. 1994;7:754-9.
- 562. Brystrom A, Sundqvist G. Bacteriologic evaluation of the effect of 0.5 percent sodium hypochlorite in endodontic therapy. Oral Surg. Oral Med. Oral Pathol. 1983;55:307-12.
- 563. Favero MJ, Tokars JI, Arduino MJ, Alter MJ. Nosocomial infections associated with hemodialysis. In: Mayhall CG, ed. Infect. Control and Hosp. Epidemiol. Philadelphia: Lippincott Williams & Wilkins, 1999:897-917.
- 564. Helms CM, Massanari RM, Zeitler R, et al. Legionnaires' disease associated with a hospital water system: a cluster of 24 nosocomial cases. Ann. Intern. Med. 1983;99:172-8.
- 565. Heffelfinger JD, Kool JL, Fridkin S, et al. Risk of hospital-acquired legionnaires' disease in cities using monchloramine versus other water disinfectants. Infect Control Hosp Epidemiol 2003;24:569-74.
- 566. Moore MR, Pryor M, Fields B, Lucas C, Phelan M, Besser RE. Introduction of monochloramine into a municipal water system: Impact on colonization of buildings by *Legionella spp*. Appl. Environ. Microbiol.

2006;72:378-83.

- 567. Srinivasan A, Bova G, Ross T, et al. A 17-month evaluation of a chlorine dioxide water treatment system to control *Legionella* species in a hospital water supply. Infect Control Hosp Epidemiol 2003;24:575-9.
- 568. Steve L, Goodhart P, Alexander J. Hydrotherapy burn treatment: use of chloramine-T against resistant microorganisms. Arch. Phys. Med. Rehabilit. 1979;60:301-3.
- 569. Coates D, Wilson M. Powders, composed of chlorine-releasing agent acrylic resin mixtures or based on peroxygen compounds, for spills of body fluids. J. Hosp. Infect. 1992;21:241-52.
- 570. Tulis JJ. Formaldehyde as a gas. In: Phillips GB, Miller WS, eds. Industrial Sterilization. Durham: Duke University Press, 1972:209-38.
- 571. Emmons CW. Fungicidal action of some common dsinfectants on two dermatophytes. Arch. Dermatol. Syphil. 1933;28:15-21.
- 572. McCulloch EC, Costigan S. A comparison of the efficiency of phenol, liquor cresolis, formaldehyde, sodium hypochlorite and sodium hydroxide against *Eberthella typhi* at various temperatures. J. Infect. Dis. 1936;59:281-4.
- 573. Sagripanti JL, Eklund CA, Trost PA, et al. Comparative sensitivity of 13 species of pathogenic bacteria to seven chemical germicides. Am. J. Infect. Control 1997;25:335-9.
- 574. NIOSH. Formaldehyde: evidence of carcinogenicity. NIOSH Current Intelligence Bulletin 34. DHEW (NIOSH) Publication No. 81-111. 1981.
- 575. Occupational Safety and Health Administration. OSHA amends formaldehyde standard. Occupational Safety and Health News 1991:1.
- 576. Occupational Health and Safety Administration. OSHA Fact Sheet: Formaldehyde: Occupational Safety and Health Administration, U.S. Department of Labor, 2002.
- 577. Occupational Safety and Health Administration. Air Contaminants Final Rule. Fed. Regist. 1993;58:35338-51.
- 578. Occupational Safety and Health Administration. Formaldehyde: OSHA Fact Sheet: Occupational Safety and Health Administration, 2002.
- 579. Centers for Disease Control. Occupational exposures to formaldehyde in dialysis units. MMWR 1986;35:399-01.
- 580. Centers for Disease Control. Formaldehyde exposures in a gross anatomy laboratory Colorado. MMWR 1983;52:698-700.
- 581. Tokars JI, Miller ER, Alter MJ, Arduino MJ. National surveillance of dialysis associated diseases in the United States, 1995. ASAIO J. 1998;44:98-107.
- 582. Favero MS, Alter MJ, Tokars JI, Bland LA. Dialysis-associated disease and their control. In: Bennett JV, Brachman PS, eds. Hospital Infections. Boston: Little, Brown and Company, 1998:357-80.
- 583. Bland LA, Favero MS. Microbial contamination control strategies for hemodialysis system. Plant, Technology & Safety Management Series: infection control issues in PTSM 1990. Oakbrook Terrace, Illinois.
- 584. Boucher RM. Potentiated acid 1,5 pentanedial solution--a new chemical sterilizing and disinfecting agent. Am. J. Hosp. Pharm. 1974;31:546-57.
- 585. Miner NA, McDowell JW, Willcockson GW, Bruckner NI, Stark RL, Whitmore EJ. Antimicrobial and other properties of a new stabilized alkaline glutaraldehyde disinfectant/sterilizer. Am. J. Hosp. Pharm. 1977;34:376-82.
- 586. Pepper RE. Comparison of the activities and stabilities of alkaline glutaraldehyde sterilizing solutions. Infect. Control 1980;1:90-2.
- 587. Leach ED. A new synergized glutaraldehyde-phenate sterilizing solution and concentrated disinfectant. Infect. Control 1981;2:26-30.
- 588. Miner NA, Ross C. Clinical evaluation of ColdSpor, a glutaraldehyde-phenolic disinfectant. Respir. Care 1991;36:104-9.
- 589. Collins FM, Montalbine V. Mycobactericidal activity of glutaraldehyde solutions. J. Clin. Microbiol. 1976;4:408-12.
- 590. Masferrer R, Marquez R. Comparison of two activated glutaraldehyde solutions: Cidex Solution and Sonacide. Respir. Care 1977;22:257-62.
- 591. Jette LP, Ringuette L, Ishak M, Miller M, Saint-Antoine P. Evaluation of three glutaraldehyde-based disinfectants used in endoscopy. J. Hosp. Infect. 1995;30:295-303.
- 592. Scott EM, Gorman SP. Glutaraldehyde. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:361-81.
- 593. Scott EM, Gorman SP. Glutaraldehyde. In: Block SS, ed. Disinfection, Sterilization, and Preservation.

Philadelphia: Lea & Febiger, 1991:596-616.

- 594. Stonehill AA, Krop S, Borick PM. Buffered glutaraldehyde a new chemical sterilizing solution. Am. J. Hosp. Pharm. 1963;20:458-65.
- 595. Borick PM, Dondershine FH, Chandler VL. Alkalinized glutaraldehyde, a new antimicrobial agent. J. Pharm. Sci. 1964;53:1273-5.
- 596. Russell AD. Glutaraldehyde: current status and uses. Infect. Control Hosp. Epidemiol. 1994;15:724-33.
- 597. Hanson PJ, Bennett J, Jeffries DJ, Collins JV. Enteroviruses, endoscopy and infection control: an applied study. J. Hosp. Infect. 1994;27:61-7.
- 598. van Klingeren B, Pullen W. Glutaraldehyde resistant mycobacteria from endoscope washers. J. Hosp. Infect. 1993;25:147-9.
- 599. Griffiths PA, Babb JR, Bradley CR, Fraise AP. Glutaraldehyde-resistant *Mycobacterium chelonae* from endoscope washer disinfectors. J. Appl. Microbiol. 1997;82:519-26.
- 600. Dauendorffer JN, Laurain C, Weber M, Dailloux M. Evaluation of the bactericidal efficiency of a 2% alkaline glutaraldehyde solution on *Mycobacterium xenopi*. J. Hosp. Infect. 2000;46:73-6.
- 601. Nomura K, Ogawa M, Miyamoto H, Muratani T, Taniguchi H. Antibiotic susceptibility of glutaraldehydetolerant *Mycobacterium chelonae* from bronchoscope washing machines. J. Hosp. Infect. 2004;32:185-8.
- 602. Webster E, Ribner B, Streed LL, Hutton N. Microbial contamination of activated 2% glutaraldehyde used in high-level disinfection of endoscopes (abstract). Am. J. Infect. Control 1996;24:153.
- 603. Casemore DP, Blewett DA, Wright SE. Cleaning and disinfection of equipment for gastrointestinal flexible endoscopy: interim recommendations of a Working Party of the British Society of Gastroenterology. Gut 1989;30:1156-7.
- 604. Laskowski LF, Marr JJ, Spernoga JF, et al. Fastidious mycobacteria grown from porcine prosthetic-heartvalve cultures. N. Engl. J. Med. 1977;297:101-2.
- 605. Collins FM. Bactericidal activity of alkaline glutaraldehyde solution against a number of atypical mycobacterial species. J. Appl. Bacteriol. 1986;61:247-51.
- 606. Food and Drug Administration. Sterilants and high level disinfectants cleared by FDA in a 510(k) as of January 30, 2002 with general claims for processing reusable medical and dental devices, [This link is no longer active: www.fda.gov/cdrh/ode/germlab.html], 2001.
- 607. Hernandez A, Martro E, Pizo C, et al. In-use evaluation of Perasafe compared with Cidex in fibreoptic bronchoscope disinfection. J. Hosp. Infect. 2003;54:46-52.
- 608. Leong D, Dorsey G, Klapp M. Dilution of glutaraldehyde by automatic endoscope machine washers: the need for a quality control program. Abstracts of the 14th Annual Educational Conference of Association for Practitioners in Infection Control, 1987:108, p.130.
- 609. Mbithi JN, Springthorpe VS, Sattar SA, Pacquette M. Bactericidal, virucidal, and mycobactericidal activities of reused alkaline glutaraldehyde in an endoscopy unit. J. Clin. Microbiol. 1993;31:2988-95.
- 610. Kleier DJ, Averbach RE. Glutaraldehyde nonbiologic monitors. Infect. Control Hosp. Epidemiol. 1990;11:439-41.
- 611. Kleier DJ, Tucker JE, Averbach RE. Clinical evaluation of glutaraldehyde nonbiologic monitors. Quintessence Int. 1989;20:271-7.
- 612. Overton D, Burgess JO, Beck B, Matis B. Glutaraldehyde test kits: evaluation for accuracy and range. Gen. Dent. 1989;37:126, 128.
- 613. Cooke RPD, Goddard SV, Chatterley R, Whymant-Morris A, Cheale J. Monitoring glutaraldehyde dilution in automated washer/disinfectors. J. Hosp. Infect. 2001;48:242-6.
- 614. Ayliffe GA, Babb JR, Bradley CR. Disinfection of endoscopes. J. Hosp. Infect. 1986;7:296-9.
- 615. Centers for Disease Control. Federal regulatory action against sporicidin cold sterilizing solution. MMWR 1991;40:880-1.
- 616. Husni L, Kale E, Climer C, Bostwick B, Parker TF, 3rd. Evaluation of a new disinfectant for dialyzer reuse. Am. J. Kidney Dis. 1989;14:110-8.
- 617. Townsend TR, Wee SB, Koblin B. An efficacy evaluation of a synergized glutaraldehyde-phenate solution in disinfecting respiratory therapy equipment contaminated during patient use. Infect. Control 1982;3:240-4.
- 618. Petersen NJ, Carson LA, Doto IL, Aguero SM, Favero MS. Microbiologic evaluation of a new glutaraldehydebased disinfectant for hemodialysis systems. Trans. Am. Soc. Artif. Intern. Organs 1982;28:287-90.
- 619. Gundogdu H, Ocal K, Caglikulekci M, Karabiber N, Bayramoglu E, Karahan M. High-level disinfection with 2% alkalinized glutaraldehyde solution for reuse of laparoscopic disposable plastic trocars. J. Laparoendosc. Adv. Surg. Techniques. Part A 1998;8:47-52.

- 620. Castelli M, Qizilbash A, Seaton T. Post-colonoscopy proctitis. Am. J. Gastroenterol. 1986;81:887.
- 621. Jonas G, Mahoney A, Murray J, Gertler S. Chemical colitis due to endoscope cleaning solutions: a mimic of pseudomembranous colitis. Gastroenterology 1988;95:1403-8.
- 622. Levine DS. Proctitis following colonoscopy. Gastrointest. Endosc. 1988;34:269-72.
- 623. Riney S, Grimes M, Khalife K, Warbasse L, Massanari M. Diarrhea associated with disinfection of sigmoidoscopes. [abstract]. Am J Infect Control 1991;19:109.
- 624. Durante L, Zulty JC, Israel E, et al. Investigation of an outbreak of bloody diarrhea: association with endoscopic cleaning solution and demonstration of lesions in an animal model. Am. J. Med. 1992;92:476-80.
- 625. Burtin P, Ruget O, Petit R, Boyer J. Glutaraldehyde-induced proctitis after endorectal ultrasound examination: a higher risk of incidence than expected? Gastrointest. Endosc. 1993;39:859-60.
- 626. Babb RR, Paaso BT. Glutaraldehyde proctitis. West. J. Med. 1995;163:477-8.
- 627. Ryan CK, Potter GD. Disinfectant colitis. Rinse as well as you wash. J. Clin. Gastroenterol. 1995;21:6-9.
- 628. Rozen P, Somjen GJ, Baratz M, Kimel R, Arber N, Gilat T. Endoscope-induced colitis: description, probable cause by glutaraldehyde, and prevention. Gastrointest. Endosc. 1994;40:547-53.
- 629. West AB, Kuan SF, Bennick M, Lagarde S. Glutaraldehyde colitis following endoscopy: clinical and pathological features and investigation of an outbreak. Gastroenterology 1995;108:1250-5.
- 630. Dolce P, Gourdeau M, April N, Bernard PM. Outbreak of glutaraldehyde-induced proctocolitis. Am. J. Infect. Control 1995;23:34-9.
- 631. Farina A, Fievet MH, Plassart F, Menet MC, Thuillier A. Residual glutaraldehyde levels in fiberoptic endoscopes: measurement and implications for patient toxicity. J. Hosp. Infect. 1999;43:293-7.
- 632. Dailey JR, Parnes RE, Aminlari A. Glutaraldehyde keratopathy. Am. J. Ophthalmol. 1993;115:256-8.
- 633. Courtright P, Lewallen S, Holland SP, Wendt TM. Corneal decompensation after cataract surgery. An outbreak investigation in Asia. Ophthalmology 1995;102:1461-5.
- 634. Leinster P, Baum JM, Baxter PJ. An assessment of exposure to glutaraldehyde in hospitals: typical exposure levels and recommended control measures. Br. J. Ind. Med. 1993;50:107-11.
- 635. Beauchamp RO, St Clair MB, Fennell TR, Clarke DO, Morgan KT. A critical review of the toxicology of glutaraldehyde. Crit. Rev. Toxicol. 1992;22:143-74.
- 636. Corrado OJ, Osman J, Davies RJ. Asthma and rhinitis after exposure to glutaraldehyde in endoscopy units. Hum. Toxicol. 1986;5:325-8.
- 637. Norback D. Skin and respiratory symptoms from exposure to alkaline glutaraldehyde in medical services. Scand. J. Work, Environ. Health 1988;14:366-71.
- 638. Mwaniki DL, Guthua SW. Occupational exposure to glutaraldehyde in tropical climates. Lancet 1992;340:1476-7.
- 639. Centers for Disease Control. Symptoms of irritation associated with exposure to glutaraldehyde. MMWR 1987;36:190-1.
- 640. Wiggins P, McCurdy SA, Zeidenberg W. Epistaxis due to glutaraldehyde exposure. J. Occup. Med. 1989;31:854-6.
- 641. Di Prima T, De Pasquale R, Nigro M. Contact dermatitis from glutaraldehyde. Contact Dermatitis 1988;19:219-20.
- 642. Fowler JF, Jr. Allergic contact dermatitis from glutaraldehyde exposure. J. Occup. Med. 1989;31:852-3.
- 643. Fisher AA. Allergic contact dermatitis of the hands from Sporicidin (glutaraldehyde-phenate) used to disinfect endoscopes. Cutis 1990;45:227-8.
- 644. Nethercott JR, Holness DL, Page E. Occupational contact dermatitis due to glutaraldehyde in health care workers. Contact Dermatitis 1988;18:193-6.
- 645. Gannon PF, Bright P, Campbell M, O'Hickey SP, Burge PS. Occupational asthma due to glutaraldehyde and formaldehyde in endoscopy and x ray departments. Thorax 1995;50:156-9.
- 646. Chan-Yeung M, McMurren T, Catonio-Begley F, Lam S. Occupational asthma in a technologist exposed to glutaraldehyde. J. Allergy Clin. Immunol. 1993;91:974-8.
- 647. Schnuch A, Uter W, Geier J, Frosch PJ, Rustemeyer T. Contact allergies in healthcare workers. Results from the IVDK. Acta Derm. Venereol. 1998;78:358-63.
- 648. Wellons SL, Trawick EG, Stowers MF, Jordan SL, Wass TL. Laboratory and hospital evaluation of four personal monitoring methods for glutaraldehyde in ambient air. Am. Ind. Hyg. Assoc. J. 1998;59:96-103.
- 649. Newman MA, Kachuba JB. Glutaraldehyde: a potential health risk to nurses. Gastroenterol. Nurs. 1992;14:296-300, discussion 300-1.

- 650. Association for the Advancement of Medical Instrumentation. Safe use and handling of glutaraldehyde-based products in healthcare facilities. Arlington, VA: AAMI, 1995.
- 651. Anonymous. Glutaraldehyde. New York: Occupational Health Services, Inc., 1992.
- 652. Rutala WA, Hamory BH. Expanding role of hospital epidemiology: employee health--chemical exposure in the health care setting. Infect. Control Hosp. Epidemiol. 1989;10:261-6.
- 653. Turner FJ. Hydrogen peroxide and other oxidant disinfectants. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1983:240-50.
- 654. Block SS. Peroxygen compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:185-204.
- 655. Sattar SA, Springthorpe VS, Rochon M. A product based on accelerated and stabilized hydrogen peroxide: Evidence for broad-spectrum germicidal activity. Canadian J Infect Control 1998 (Winter):123-30.
- 656. Omidbakhsh N, Sattar SA. Broad-spectrum microbicidal activity, toxicologic assessment, and materials compatibility of a new generation of accelerated hydrogen peroxide-based environmental surface disinfectant. Am. J. Infect. Control 2006;34:251-7.
- 657. Schaeffer AJ, Jones JM, Amundsen SK. Bacterial effect of hydrogen peroxide on urinary tract pathogens. Appl. Environ. Microbiol. 1980;40:337-40.
- 658. Wardle MD, Renninger GM. Bactericidal effect of hydrogen peroxide on spacecraft isolates. Appl. Microbiol. 1975;30:710-1.
- 659. Sagripanti JL, Bonifacino A. Comparative sporicidal effect of liquid chemical germicides on three medical devices contaminated with spores of *Bacillus subtilis*. Am. J. Infect. Control 1996;24:364-71.
- 660. Sagripanti JL, Bonifacino A. Effects of salt and serum on the sporicidal activity of liquid disinfectants. J. AOAC Int. 1997;80:1198-207.
- 661. Saurina G, Landman D, Quale JM. Activity of disinfectants against vancomycin-resistant *Enterococcus faecium*. Infect. Control Hosp. Epidemiol. 1997;18:345-7.
- 662. Kilvington S. Moist-heat disinfection of *Acanthamoeba* cysts. Rev. Infect. Dis. 1991;13:S418.
- 663. Sattar SA, Adegbunrin O, Ramirez J. Combined application of simulated reuse and quantitative carrier test to assess high-level disinfection: Experiments with an accelerated hydrogen peroxide-based formulation. Am. J. Infect. Control 2002;30:449-57.
- 664. Leaper S. Influence of temperature on the synergistic sporicidal effect of peracetic acid plus hydrogen peroxide in *Bacillus subtilis* SA22(NCA 72-52). Food Microbiol. 1984;1:199-203.
- 665. Mentel R, Schmidt J. Investigations on rhinovirus inactivation by hydrogen peroxide. Acta Virol. 1973;17:3514.
- 666. Sattar SA. Effect of liquid chemical germicides on mycobacteria including multi-drug resistant isolates of *Mycobacteria tuberculosis*. Abstracts of the 37th Interscience Conference on Antimicrobial Agents of Chemotherapy; September 28-October 1, 1997; Toronto, Ontario, Canada; E166., 1997.
- 667. Reckitt & Colman. Sporox sterilant and high-level disinfectant technical report. Montvale, NJ: Reckitt & Colman, 1997:1-12.
- 668. Sattar SA, Taylor YE, Paquette M, Rubino J. In-hospital evaluation of 7.5% hydrogen peroxide as a disinfectant for flexible endoscopes. Can. J. Infect. Control 1996;11:51-4.
- 669. Hobson DW, Seal LA. Evaluation of a novel, rapid-acting, sterilizing solution at room temperature. Am. J. Infect. Control 2000;28:370-5.
- 670. Anonymous. Hydrogen peroxide, ACS reagent. Vol. 2001: Sigma Product Information Sheet, [This link is no longer active: http://www.sigma.sial.com/sigma/proddata/h0904.htm].
- 671. Silvany RE, Dougherty JM, McCulley JP, Wood TS, Bowman RW, Moore MB. The effect of currently available contact lens disinfection systems on *Acanthamoeba castellanii* and *Acanthamoeba polyphaga*. Ophthalmology 1990;97:286-90.
- 672. Moore MB. *Acanthamoeba keratitis* and contact lens wear: the patient is at fault. Cornea 1990;9:S33-5; discussion S39-40.
- 673. Judd PA, Tomlin PJ, Whitby JL, Inglis TC, Robinson JS. Disinfection of ventilators by ultrasonic nebulisation. Lancet 1968;2:1019-20.
- 674. Levenson JE. Corneal damage from improperly cleaned tonometer tips. Arch. Ophthalmol. 1989;107:1117.
- 675. Thompson RL, Haley CE, Searcy MA, et al. Catheter-associated bacteriuria. Failure to reduce attack rates using periodic instillations of a disinfectant into urinary drainage systems. JAMA 1984;251:747-51.
- 676. Bilotta JJ, Waye JD. Hydrogen peroxide enteritis: the "snow white" sign. Gastrointest. Endosc. 1989;35:428-30.

- 677. Gottardi W. Iodine and iodine compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1991:152-66.
- 678. Gottardi W. Iodine and iodine compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:159-84.
- 679. Craven DE, Moody B, Connolly MG, Kollisch NR, Stottmeier KD, McCabe WR. Pseudobacteremia caused by povidone-iodine solution contaminated with *Pseudomonas cepacia*. N. Engl. J. Med. 1981;305:621-3.
- 680. Berkelman RL, Lewin S, Allen JR, et al. Pseudobacteremia attributed to contamination of povidone-iodine with *Pseudomonas cepacia*. Ann. Intern. Med. 1981;95:32-6.
- 681. Parrott PL, Terry PM, Whitworth EN, et al. *Pseudomonas aeruginosa* peritonitis associated with contaminated poloxamer-iodine solution. Lancet 1982;2:683-5.
- 682. Favero MS. Iodine--champagne in a tin cup. Infect. Control 1982;3:30-2.
- 683. Berkelman RL, Holland BW, Anderson RL. Increased bactericidal activity of dilute preparations of povidoneiodine solutions. J. Clin. Microbiol. 1982;15:635-9.
- 684. Chang SL. Modern concept of disinfection. J. Sanit. Eng. Div. Proc. Am. Soc. Civ. Eng. 1971:689-705.
- 685. Wallbank AM, Drulak M, Poffenroth L, Barnes C, Kay C, Lebtag I. Wescodyne: lack of activity against poliovirus in the presence of organic matter. Health Lab. Sci. 1978;15:133-7.
- 686. Carson JA, Favero MS. Comparative resistance of nontuberculous mycobacteria to iodophor germicides. Abstracts of the Annual Meeting of the American Society for Microbiology, 1984:Q101, p.221.
- 687. Medcomp. Medcomp Frequently Asked Questions. [This link is no longer active: www.medcompnet.com/faq/faq/html], 2000.
- 688. Simons C, Walsh SE, Maillard JY, Russell AD. A note: ortho-phthalaldehyde: proposed mechanism of action of a new antimicrobial agent. Lett. Appl. Microbiol. 2000;31:299-302.
- 689. Walsh SE, Maillard JY, Simons C, Russell AD. Studies on the mechanisms of the antibacterial action of ortho-phthalaldehyde. J. Appl. Microbiol. 1999;87:702-10.
- 690. Fraud S, Hann AC, Maillard J-Y, Russell AD. Effects of ortho-phthalaldehyde, glutaraldehyde and chlorhexidine diacetate on *Mycobacterium chelonae* and *Mycobacterium abscessus* strains with modified permeability. J. Antimicrob. Chemother. 2003;51:575-84.
- 691. Cabrera-Martinez RM, Setlow B, Setlow P. Studies on the mechanisms of the sporicidal action of orthophthalaldehyde. J. Appl. Microbiol. 2002;92:675-80.
- 692. Gordon MD, Ezzell RJ, Bruckner NI, Ascenzi JM. Enhancement of mycobactericidal activity of glutaraldehyde with Unsaturated and aromatic aldehydes. J. Indust. Microbiol. 1994;13:77-82.
- 693. Gregory AW, Schaalje GB, Smart JD, Robison RA. The mycobactericidal efficacy of ortho-phthalaldehyde and the comparative resistances of *Mycobacterium bovis, Mycobacterium terrae*, and *Mycobacterium chelonae*. Infect. Control Hosp. Epidemiol. 1999;20:324-30.
- 694. Walsh SE, Maillard JY, Russell AD. Ortho-phthalaldehyde: a possible alternative to glutaraldehyde for high level disinfection. J. Appl. Microbiol. 1999;86:1039-46.
- 695. Roberts CG, Chan Myers H. Mycobactericidal activity of dilute ortho-phthalaldehyde solutions. In: Abstracts in Environmental and General Applied Microbiology, Q-265, ASM 98th General Meeting, Atlanta, Georgia, USA, 1998:464-5.
- 696. Chan-Myers H. Sporicidal activity of ortho-phthalaldehyde as a function of temperature (abstract). Infect. Control Hosp. Epidemiol. 2000;21:101.
- 697. Chan-Myers H, Roberts C. Effect of temperature and organic soil concentration on biocidal activity of orthophthalaldehyde solution (abstract). 2000 Education Meeting of the Association for Professional in Infection Control and Epidemiology, Minneapolis, MN, 2000:31.
- 698. Bruckner NI, Gordon MD, Howell RG. Odorless aromatic dialdehyde disinfecting and sterilizing composition. US Patent 4,851,449. July, 1989.
- 699. McDonnell G, Pretzer D. New and developing chemical antimicrobials. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:431-43.
- 700. Fraud S, Maillard J-Y, Russell AD. Comparison of the mycobactericidal activity of ortho-phthalaldehyde, glutaraldehyde, and other dialdehydes by a quantitative suspension test. J. Hosp. Infect. 2001;48:214-21.
- 701. Sattar SA, Springthorpe VS. New methods for efficacy testing of disinfectants and antiseptics. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:174-86.
- 702. Walsh SE, Maillard J-Y, Russell AD, Hann AC. Possible mechanisms for the relative efficiacies of orthophthalaldehyde and glutaraldehyde aganst glutaradehyde-resistant *Mycobacterium chelonae*. J. Appl.

Microbiol. 2001;91:80-92.

- 703. Herruzo-Cabrera R, Vizcaino-Alcaide MJ, Rodriquez J. Comparison of the microbicidal efficacy on germ carriers of several tertiary amine compounds with ortho-phthalaldehyde and Perasafe. J. Hosp. Infect. 2006;63:73-8.
- 704. Herruzo-Cabrera R, Vizcaino-Alcaide MJ, Fernandez-Acenero MJ. The influence of laboratory adaptation on test strains, such as *Pseudomonas aeruginosa*, in the evaluation of the antimicrobial efficacy of ortho-phthalaldehyde. J. Hosp. Infect. 2004;57:217-22.
- 705. Favero MS. Naturally occurring microorganisms and their resistance to physical and chemical agents. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington, DC: Association for Professionals in Infection Control and Epidemiology, 2004:1-15.
- 706. Cooke RPD, Goddard SV, Whymant-Morris A, Sherwood J, Chatterly R. An evaluation of Cidex OPA (0.55% ortho-phthalaldehyde) as an alternative to 2% glutaraldehyde for high-level disinfection of endoscopes. J. Hosp. Infect. 2003;54:226-31.
- 707. Streckenbach SC, Alston TA. Perioral stains after ortho-phthalaldehyde disinfectioon of echo probes. Anesthesiol 2003;99:1032.
- 708. Wardle E, Jones D. Determination of rinsing volumes following manual endoscope disinfection with orthophthalaldehyde (OPA). J Gastroenterol Nurses College Australia 2003;January:7-9.
- 709. Sokol WN. Nine episodes of anaphylaxis following cytoscopy caused by Cidex OPA (ortho-phthalaldehyde) high-level disinfectant in 4 patients after cystoscopy. J. Allergy Clin. Immunol. 2004;114:392-7.
- 710. Hession SM. Endoscopic disinfection by ortho-phthalaldehyde in a clinical setting: An evaluation of reprocessing time and costs compared with glutaraldehyde. Gastroenterol. Nurs. 2003;26:110-4.
- 711. Tucker RC, Lestini BJ, Marchant RE. Surface analysis of clinically used expanded PTFE endoscopic tubing treated by the STERIS PROCESS. ASAIO J. 1996;42:306-13.
- 712. Hernandez A, Martro E, Matas L, Ausina V. In-vitro evaluation of Pearsafe compared with 2% alkaline glutaraldehyde against *Mycobacterium* spp. J. Hosp. Infect. 2003;54:52-6.
- 713. Vizcaino-Alcaide MJ, Herruzo-Cabrera R, Fernandez-Acenero MJ. Comparison of the disinfectant efficacy of Persafe and 2% glutaraldehyde in in vitro tests. J. Hosp. Infect. 2003;53:124-8.
- 714. Lensing HH, Oei HL. Investigations on the sporicidal and fungicidal activity of disinfectants. Zentralblatt fur Bakteriologie, Mikrobiologie und Hygiene 1 Abt Originale B, Hygiene 1985;181:487-95.
- 715. Sagripanti JL, Bonifacino A. Comparative sporicidal effects of liquid chemical agents. Appl. Environ. Microbiol. 1996;62:545-51.
- 716. Crow S. Peracetic acid sterilization: a timely development for a busy healthcare industry. Infect. Control Hosp. Epidemiol. 1992;13:111-3.
- 717. Malchesky PS. Medical applications of peracetic acid. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:979-96.
- 718. Mannion PT. The use of peracetic acid for the reprocessing of flexible endoscopes and rigid cystoscopes and laparoscopes. J. Hosp. Infect. 1995;29:313-5.
- 719. Bradley CR, Babb JR, Ayliffe GA. Evaluation of the Steris System 1 Peracetic Acid Endoscope Processor. J. Hosp. Infect. 1995;29:143-51.
- 720. Duc DL, Ribiollet A, Dode X, Ducel G, Marchetti B, Calop J. Evaluation of the microbicidal efficacy of Steris System I for digestive endoscopes using GERMANDE and ASTM validation protocols. J. Hosp. Infect. 2001;48:135-41.
- 721. Alfa MJ, Olson N, Degagne P, Hizon R. New low temperature sterilization technologies: microbicidal activity and clinical efficacy. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:67-78.
- 722. Alfa MJ, DeGagne P, Olson N, Hizon R. Comparison of liquid chemical sterilization with peracetic acid and ethylene oxide sterilization for long narrow lumens. Am. J. Infect. Control 1998;26:469-77.
- 723. Seballos RJ, Walsh AL, Mehta AC. Clinical evaluation of a liquid chemical sterilization system for flexible bronchoscopes. J. Bronch. 1995;2:192-99.
- 724. Wallace CG, Agee PM, Demicco DD. Liquid chemical sterilization using peracetic acid. An alternative approach to endoscope processing. ASAIO J. 1995;41:151-4.
- 725. Centers for Disease Control and Prevention. Bronchoscopy-related infections and pseudoinfections New York, 1996 and 1998. MMWR 1999;48:557-60.
- 726. Middleton AM, Chadwick MV, Gaya H. Disinfection of bronchoscopes, contaminated in vitro with Mycobacterium tuberculosis, Mycobacterium avium-intracellulare and Mycobacterium chelonae in sputum,

using stabilized, buffered peracetic acid solution ('Nu-Cidex'). J. Hosp. Infect. 1997;37:137-43.

- 727. Holton J, Shetty N. In-use stability of Nu-Cidex. J. Hosp. Infect. 1997;35:245-8.
- 728. Alasri A, Roques C, Michel G, Cabassud C, Aptel P. Bactericidal properties of peracetic acid and hydrogen peroxide, alone and in combination, and chlorine and formaldehyde against bacterial water strains. Can. J. Microbiol. 1992;38:635-42.
- 729. Stanley P. Destruction of a glutaraldehyde-resistant mycobacterium by a per-oxygen disinfectant. (Abstract). Am. J. Infect. Control 1998;26:185.
- 730. Fleming SJ, Foreman K, Shanley K, Mihrshahi R, Siskind V. Dialyser reprocessing with Renalin. Am. J. Nephrol. 1991;11:27-31.
- 731. Kahn G. Depigmentation caused by phenolic detergent germicides. Arch. Dermatol. 1970;102:177-87.
- 732. Prindle RF. Phenolic compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1983:197-224.
- 733. Hegna IK. A comparative investigation of the bactericidal and fungicidal effects of three phenolic disinfectants. J. Appl. Bacteriol. 1977;43:177-81.
- 734. Hegna IK. An examination of the effect of three phenolic disinfectants on *Mycobacterium tuberculosis*. J. Appl. Bacteriol. 1977;43:183-7.
- 735. Bergan T, Lystad A. Antitubercular action of disinfectants. J. Appl. Bacteriol. 1971;34:751-6.
- 736. Narang HK, Codd AA. Action of commonly used disinfectants against enteroviruses. J. Hosp. Infect. 1983;4:209-12.
- 737. Cole EC, Rutala WA, Samsa GP. Disinfectant testing using a modified use-dilution method: collaborative study. J. Assoc. Off. Anal. Chem. 1988;71:1187-94.
- 738. Goddard PA, McCue KA. Phenolic compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:255-81.
- 739. Wysowski DK, Flynt JW, Jr., Goldfield M, Altman R, Davis AT. Epidemic neonatal hyperbilirubinemia and use of a phenolic disinfectant detergent. Pediatrics 1978;61:165-70.
- 740. Doan HM, Keith L, Shennan AT. Phenol and neonatal jaundice. Pediatrics 1979;64:324-5.
- 741. Shickman MD, Guze LB, Pearce ML. Bacteremia following cardiac catheterization. N. Engl. J. Med. 1959;260:1164-6.
- 742. Ehrenkranz NJ, Bolyard EA, Wiener M, Cleary TJ. Antibiotic-sensitive *Serratia marcescens* infections complicating cardiopulmonary operations: contaminated disinfectant as a reservoir. Lancet 1980;2:1289-92.
- 743. Shere L. Some comparisons of the disinfecting properties of hypochlorites and quaternary ammonium compounds. Milk Plant Monthly March 1948:66-9.
- 744. MacDougall KD, Morris C. Optimizing disinfectant application in healthcare facilities. Infect Control Today 2006;June:62-7.
- 745. Sykes G. Disinfection and sterilization. London: E & FN Spon Ltd, 1965.
- 746. Merianos JJ. Surface-active agents. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:283-320.
- 747. Purohit A, Kopferschmitt-Kubler MC, Moreau C, Popin E, Blaumeiser M, Pauli G. Quaternary ammonium compounds and occupational asthma. International Archives of Occupational & Environmental Health 2000;73:423-7.
- 748. Petrocci AN. Surface active agents: quaternary ammonium compounds. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1983:309-29.
- 749. Smith CR, Nishihara H, Golden F, Hoyt A, Guss CO, Kloetzel MC. The bactericidal effect of surface-active agents on tubercle bacilli. Public Health Rep. 1950;48:1588-1600.
- 750. Broadley SJ, Furr JR, Jenkins PA, Russell AD. Antimycobacterial activity of 'Virkon'. J. Hosp. Infect. 1993;23:189-97.
- 751. Angelillo IF, Bianco A, Nobile CG, Pavia M. Evaluation of the efficacy of glutaraldehyde and peroxygen for disinfection of dental instruments. Lett. Appl. Microbiol. 1998;27:292-6.
- 752. Coates D. Disinfectants and spills of body fluids. Nurs. RSA 1992;7:25-7.
- 753. Hamouda T, Hayes MM, Cao ZH, et al. A novel surfactant nanoemulsion with broad-spectrum sporicidal activity against *Bacillus* species. J. Infect. Dis. 1999;180:1939-49.
- 754. Hamouda T, Myc A, Donovan B, Shih AY, Reuter JD, Baker JR. A novel surfactant nanoemulsion with a unique non-irritant topical antimicrobial activity against bacteria, enveloped viruses and fungi. Microbiol. Res. 2001;156:1-7.
- 755. Hamouda T, Baker JR, Jr. Antimicrobial mechanism of action of surfactant lipid preparations in enteric Gram-

negative bacilli. J. Appl. Microbiol. 2000;89:397-403.

- 756. Widmer AF, Frei R. Antimicrobial activity of glucoprotamin: A clinical study of a new disinfectant for instruments. Infect Control Hosp Epidemiol 2003;24:762-4.
- 757. Wilson M. Light-activated antimicrobial coating for the continuous disinfection of surfaces. Infect Control Hosp Epidemiol 2003;24:782-4.
- 758. Schneider PM. New technologies for disinfection and sterilization. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington DC: Association for Professionals in Infection Control and Epidemiology, 2004:127-39.
- 759. Rotter ML. Handwashing, hand disinfection, and skin disinfection. In: Wenzel RP, ed. Prevention and control of nosocomial infections. Baltimore: Williams and Wilkins, 1997:691-709.
- 760. Mandel GL, Bennett JE, Dolin R. Principles and practices of infectious diseases. New York: Livingstone, 2000.
- 761. Weber DJ, Rutala WA. Use of metals and microbicides in the prevention of nosocomial infections. In: Rutala W, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1995:271-85.
- 762. Weber DJ, Rutala WA. Use of metals as microbicides in preventing infections in healthcare. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:415-30.
- 763. Brady MJ, Lisay CM, Yurkovetskiy AV, Sawan SP. Persistent silver disinfectant for the environmental control of pathogenic bacteria. Am. J. Infect. Control 2003;31:208-214.
- 764. Rusin P, Bright K, Gerba C. Rapid reduction of *Legionella pneumophila* on stainless steel with zeolite coatings containing silver and zinc ions. Lett. Appl. Microbiol. 2003;36:69-72.
- 765. Bright KR, Gerba CP, Rusin PA. Rapid reduction of *Staphylococcus aureus* populations on stainless stell surfaces by zeolite ceramic coatings containing silver and zinc ions. J. Hosp. Infect. 2002;52:307-9.
- 766. Landeen LK, Yahya MT, Gerba CP. Efficacy of copper and silver ions and reduced levels of free chlorine in inactivation of *Legionella pneumophila*. Appl. Environ. Microbiol. 1989;55:3045-50.
- 767. Pyle BH, Broadaway SC, McFeters GA. Efficacy of copper and silver ions with iodine in the inactivation of *Pseudomonas cepacia*. J. Appl. Bacteriol. 1992;72:71-9.
- 768. Yahya MT, Landeen LK, Messina MC, Kutz SM, Schulze R, Gerba CP. Disinfection of bacteria in water systems by using electrolytically generated copper:silver and reduced levels of free chlorine. Can. J. Microbiol. 1990;36:109-16.
- 769. Liu Z, Stout JE, Tedesco L, et al. Controlled evaluation of copper-silver ionization in eradicating *Legionella pneumophila* from a hospital water distribution system. J. Infect. Dis. 1994;169:919-22.
- 770. Noyce JO, Michels H, Keevil CW. Potential use of copper surfaces to reduce survival of epidemic methicillinresistant *Staphylococcus aureus* in the healthcare environment. J. Hosp. Infect. 2006;63:289-97.
- 771. Goetz A, Yu VL. Copper-silver ionization: cautious optimism for *Legionella* disinfection and implications for environmental culturing. Am. J. Infect. Control 1997;25:449-51.
- 772. Miuetzner S, Schwille RC, Farley A, et al. Efficacy of thermal treatment and copper-silver ionization for controlling *Legionella pneumophila* in high-volume hot water plumbing systems in hospitals. Am. J. Infect. Control 1997;25:452-7.
- 773. Stout JE, Lin YS, Goetz AM, Muder RR. Controlling *Legionella* in hospital water systems: experience with the superheat-and-flush method and copper-silver ionization. Infect. Control Hosp. Epidemiol. 1998;19:911-4.
- 774. Stout JE, Yu VL. Experience of the first 16 hospitals using copper-silver ionization for *Legionella* control: Implications for the evaluation of other disinfection modalities. Infect Control Hosp Epidemiol 2003;24:563-8.
- 775. Russell AD. Ultraviolet radiation. In: Russell AD, Hugo WB, Ayliffe GAJ, eds. Principles and practices of disinfection, preservation and sterilization. Oxford: Blackwell Science, 1999:688-702.
- 776. Hall KK, Giannetta ET, Getchell-White SI, Durbin LJ, Farr BM. Ultraviolet light disinfection of hospital water for preventing nosocomial *Legionella* infection: A 13-year follow-up. Infect Control Hosp Epidemiol 2003;24:580-3.
- 777. Singh S, Schaaf NG. Dynamic sterilization of titanium implants with ultraviolet light. Internat. J. Oral Maxillofac. Implants 1989;4:139-46.
- 778. Dolman PJ, Dobrogowski MJ. Contact lens disinfection by ultraviolet light. Am. J. Ophthalmol. 1989;108:665-9.
- 779. Shechmeister IL. Sterilization by ultraviolet irradiation. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1991:553-65.
- 780. National Research Council. Postoperative wound infections the influence of ultraviolet irradiation of the

operating room and of various other factors. Ann. Surg. 1964;160:1-125.

- 781. Sensakovic JW, Smith LG. Nosocomial ultraviolet keratoconjunctivitis. Infect. Control 1982;3:475-6.
- 782. Cefai C, Richards J, Gould FK, McPeake P. An outbreak of respiratory tract infection resulting from incomplete disinfection of ventilatory equipment. J. Hosp. Infect. 1990;15:177-82.
- 783. Gurevich I, Tafuro P, Ristuccia P, Herrmann J, Young AR, Cunha BA. Disinfection of respirator tubing: a comparison of chemical versus hot water machine-assisted processing. J. Hosp. Infect. 1983;4:199-208.
- 784. Rutala WA, Weber DJ, Gergen MF, Gratta AR. Efficacy of a washer-pasteurizer for disinfection of respiratory-care equipment. Infect. Control Hosp. Epidemiol. 2000;21:333-6.
- 785. Jette LP, Lambert NG. Evaluation of two hot water washer disinfectors for medical instruments. Infect. Control Hosp. Epidemiol. 1988;9:194-9.
- 786. Wang C-Y, Wu H-D, Lee L-N, et al. Pasteurization is effective against multidrug-resistant bacteria. Am. J. Infect. Control 2006;34:320-2.
- 787. Dempsey KM, Chiew RF, McKenzie JA, Mitchell DH. Evaluation of the cleaning and disinfection efficacy of the DEKO-190; award-based automated washer/disinfector. J. Hosp. Infect. 2000;46:50-4.
- 788. Kearns AM, Freeman R, Lightfoot NF. Nosocomial enterococci: resistance to heat and sodium hypochlorite. J. Hosp. Infect. 1995;30:193-9.
- 789. Bradley CR, Fraise AP. Heat and chemical resistance of enterococci. J. Hosp. Infect. 1996;34:191-6.
- 790. Chadwick PR, Oppenheim BA. Vancomycin-resistant enterococci and bedpan washer machines. Lancet 1994;344:685.
- 791. Nystrom B. New technology for sterilization and disinfection. Am. J. Med. 1991;91:264S-266S.
- 792. Sanders FT, Morrow MS. The EPA's role in the regulation of antimicrobial pesticides in the United States. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington, DC: Association for Professionals in Infection Control and Epidemiology, 2004:29-41.
- 793. Anonymous. Memorandum of understanding between the Food and Drug Administration, Public Health Service, and the Environmental Protection Agency, 1993.
- 794. Ulatowski TA. Current activities concerning the premarket evaluation of infection control devices at the Food and Drug Administration. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:1-7.
- 795. Cole EC, Rutala WA. Bacterial numbers on penicylinders used in disinfectant testing: use of 24 hour adjusted broth cultures. J. Assoc. Off. Anal. Chem. 1988;71:9-11.
- 796. Cole EC, Rutala WA, Alfano EM. Comparison of stainless steel penicylinders used in disinfectant testing. J. Assoc. Off. Anal. Chem. 1988;71:288-9.
- 797. Cole EC, Rutala WA, Carson JL. Evaluation of penicylinders used in disinfectant testing: bacterial attachment and surface texture. J. Assoc. Off. Anal. Chem. 1987;70:903-6.
- 798. Cole EC, Rutala WA, Samsa GP. Standardization of bacterial numbers of penicylinders used in disinfectant testing: interlaboratory study. J. Assoc. Off. Anal. Chem. 1987;70:635-7.
- 799. Alfano EM, Cole EC, Rutala WA. Quantitative evaluation of bacteria washed from stainless steel penicylinders during AOAC use-dilution method. J. Assoc. Off. Anal. Chem. 1988;71:868-71.
- 800. Favero MS, Groschel DHM. Chemical germicides in the health care field: current status and evaluation of efficacy and research needs. Washington, DC: American Society for Microbiology, 1987.
- 801. Sattar SA. Microbicidal testing of germicides: an update. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:227-40.
- 802. Best M. Development of a combined carrier test for disinfectant efficacy. Ottawa, Canada: University of Ottawa, 1994.
- 803. Sattar SA, Springthorpe VS. Recent developments in methods for testing the germicidal activity of disinfectants and antiseptics. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington, DC: Association for Professionals in Infection Control and Epidemiology, 2004:180-8.
- 804. Sanders FT. Environmental protection agency's role in the regulation of antimicrobial pesticides in the United States. In: Rutala WA, ed. Disinfection, Sterilization and Antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:28-40.
- 805. Groschel DHM. Caveat emptor: do your disinfectants work? Infect. Control 1983;4:144.
- 806. United States General Accounting Office. Disinfectants: EPA lacks assurance they work., 1990.
- 807. Johnston MD, Lambert RJW, Hanlon GW, Denyer SP. A rapid method for assessing the suitability of quenching agents for individual biocides as well as combinations. J. Appl. Microbiol. 2002;92:784-9.

- 808. Russell AD. Neutralization procedures in the evaluation of bactericidal activity. In: Collins CH, Allwood MC, Bloomfield SF, Fox A, eds. Disinfectants: their use and evaluation of effectiveness. London: Academic Press, 1981:45-59.
- 809. Russell AD, Ahonkhai I, Rogers DT. Microbiological applications of the inactivation of antibiotics and other antimicrobial agents. J. Appl. Bacteriol. 1979;46:207-45.
- 810. Engley FB, Jr, Dey BP. A universal neutralizing medium for antimicrobial chemicals. Chem. Specialists Manuf. Assoc. Proc. 1970:100-6.
- 811. Association for the Advancement of Medical Instrumentation. Good hospital practice: Steam sterilization and sterility assurance. AAMI. Arlington, VA, 1993.
- 812. Association for the Advancement of Medical Instrumentation. Flash sterilization: Steam sterilization of patient care items for immediate use. AAMI. Arlington, VA, 1996.
- 813. Association for the Advancement of Medical Instrumentation. Steam sterilization and sterility assurance in health care facilities. ANSI/AAMI ST46. Arlington, VA, 2002:ANSI/AAMI ST46:2002.
- 814. Association for the Advancement of Medical Instrumentation. Ethylene oxide sterilization in health care facilities: Safety and effectiveness. AAMI. Arlington, VA, 1999.
- 815. Association of Operating Room Nurses. Recommended practices for sterilization in perioperative practice settings. 2000 Standards, Recommended Practices, and Guidelines. Denver, CO: AORN, 2000:347-58.
- 816. Association for peri-Operative Registered Nurses. Recommended practices for cleaning and caring for surgical instruments and powered equipment. AORN J. 2002;75:727-41.
- Mangram AJ, Horan TC, Pearson ML, Silver LC, Jarvis WR. Guideline for prevention of surgical site infection, 1999. Hospital Infection Control Practices Advisory Committee. Infect. Control Hosp. Epidemiol. 1999;20:250-78.
- 818. Education Design. Best practices for the prevention of surgical site infection. Denver Colorado: Education Design, 1998.
- 819. Association for the Advancement of Medical Instrumentation. Comprehensive guide to steam sterilization and sterility assurance in health care facilities, ANSI/AAMI ST79. 2006.
- 820. Association for peri-Operative Registered Nurses. Recommended pracice for sterilization in the perioperative practice setting. AORN J. 2006;83:700-22.
- 821. Singh J, Bhatia R, Gandhi JC, et al. Outbreak of viral hepatitis B in a rural community in India linked to inadequately sterilized needles and syringes. Bull. World Health Organ. 1998;76:93-8.
- 822. Eickhoff TC. An outbreak of surgical wound infections due to *Clostridium perfringens*. Surg. Gynecol. Obstet. 1962;114:102-8.
- 823. Favero MS. Sterility assurance: Concepts for patient safety. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:110-9.
- 824. Oxborrow GS, Berube R. Sterility testing-validation of sterilization processes, and sporicide testing. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1991:1047-57.
- 825. Rutala WA, Weber DJ. Clinical effectiveness of low-temperature sterilization technologies. Infect. Control Hosp. Epidemiol. 1998;19:798-804.
- 826. Adler S, Scherrer M, Daschner FD. Costs of low-temperature plasma sterilization compared with other sterilization methods. J. Hosp. Infect. 1998;40:125-34.
- 827. Joslyn L. Sterilization by heat. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:695-728.
- 828. Bucx MJ, Veldman DJ, Beenhakker MM, Koster R. The effect of steam sterilization at 134 degrees C on light intensity provided by fibrelight Macintoch laryngoscopes. Anaesthesia 2000;55:185-6.
- 829. Gilbert JA, Phillips HO. The effect of steam sterilization on plaster casting material. Clinical Orthopaed Rel Res 1984:241-4.
- 830. Agalloco JP, Akers JE, Madsen RE. Moist heat sterilization--myths and realities. PDA J. Pharmaceutical Sci. Technol. 1998;52:346-50.
- 831. Rutala WA, Stiegel MM, Sarubbi FA, Jr. Decontamination of laboratory microbiological waste by steam sterilization. Appl. Environ. Microbiol. 1982;43:1311-6.
- 832. Lauer JL, Battles DR, Vesley D. Decontaminating infectious laboratory waste by autoclaving. Appl. Environ. Microbiol. 1982;44:690-4.
- 833. Rhodes P, Zelner L, Laufman H. A new disposable bowie-Dick-type test pack for prevacuum hightemperature sterilizers. Med. Instrum. 1982;16:117-20.

- 834. Association for the Advancement of Medical Instrumentation. Technical Information Report on process challenge devices/test packs for use in health care facilities, 2003.
- 835. Young M. Sterilization process monitoring. Managing Infect Control 2004; August: 70-6.
- 836. American Society for Healthcare Central Service Professionals. Training Manual for Health Care Central Service Technicians. In: Association AH, ed. Chicago: The Jossey-Bass/American Hospital Association Press Series, 2001:1-271.
- 837. Crow S. Steam sterilizers: an evolution in design. Infect. Control Hosp. Epidemiol. 1993;14:488-90.
- 838. Gurevich I, Jacobsen E, Cunha BA. Pseudoautoclave failure caused by differences in spore test steam sensitivities. Am. J. Infect. Control 1996;24:402-4.
- 839. Bryce EA, Roberts FJ, Clements B, MacLean S. When the biological indicator is positive: investigating autoclave failures. Infect. Control Hosp. Epidemiol. 1997;18:654-6.
- Barone MA, Faisel AJ, Andrews L, Ahmed J, Rashida B, Kristensen D. Adaptation and validation of a portable steam sterilizer for processing intrauterine device insertion instruments and supplies in low-resource settings. Am. J. Infect. Control 1997;25:350-6.
- 841. Young JH. Sterilization with steam under pressure. In: Morrissey RF, Phillips GB, eds. Sterilization technology: a practical guide for manufacturers and users of health care product. New York: Van Nostrand Reinhold, 1993:81-119.
- 842. Palenik CJ, Cumberlander ND. Effects of steam sterilization on the contents of sharps containers. Am. J. Infect. Control 1993;21:28-33.
- Rutala WA. Disinfection and flash sterilization in the operating room. J. Ophthal. Nurs. Technol. 1991;10:106-15.
- 844. Maki DG, Hassemer CA. Flash sterilization: carefully measured haste. Infect. Control 1987;8:307-10.
- 845. Barrett T. Flash sterilization: What are the risks? In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:70-6.
- 846. Vesley D, Langholz AC, Rohlfing SR, Foltz WE. Fluorimetric detection of a *Bacillus stearothermophilus* spore-bound enzyme, D-glucosidase, for rapid identification of flash sterilization failure. Appl. Environ. Microbiol. 1992;58:717-9.
- 847. Rutala WA, Gergen MF, Weber DJ. Evaluation of a rapid readout biological indicator for flash sterilization with three biological indicators and three chemical indicators. Infect. Control Hosp. Epidemiol. 1993;14:390-4.
- Strzelecki LR, Nelson JH. Evaluation of closed container flash sterilization system. Orthoped. Nurs. 1989;8:21-4.
- 849. Hood E, Stout N, Catto B. Flash sterilization and neurosurgical site infections: Guilt by association. Am. J. Infect. Control 1997;25:156.
- 850. Rutala WA, Weber DJ, Chappell KJ. Patient injury from flash-sterilized instruments. Infect. Control Hosp. Epidemiol. 1999;20:458.
- 851. Schneider PM. Low-temperature sterilization alternatives in the 1990s. Tappi J. 1994;77:115-9.
- 852. Environmental Protection Agency. Protection of stratospheric ozone; Proposed Rule. 40 CFR Part 82. Fed. Regist. 1993.
- 853. Schneider PM. Emerging low temperature sterilization technologies (non-FDA approved). In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:79-92.
- 854. Gross D. Ethylene oxide sterilization and alternative methods. Surg. Serv. Management 1995;1:16-7.
- 855. Holler C, Martiny H, Christiansen B, Ruden H, Gundermann KO. The efficacy of low temperature plasma (LTP) sterilization, a new sterilization technique. Zentralbl. Hyg. Umweltmed. 1993;194:380-91.
- 856. Rutala WA, Gergen MF, Weber DJ. Comparative evaluation of the sporicidal activity of new low-temperature sterilization technologies: ethylene oxide, 2 plasma sterilization systems, and liquid peracetic acid. Am. J. Infect. Control 1998;26:393-8.
- 857. Ernst RR, Doyle JE. Sterilization with gaseous ethylene oxide: a review of chemical and physical factors. Biotech. Bioeng. 1968;10.
- 858. Joslyn L. Gaseous chemical sterilization. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:337-60.
- 859. Fisher AA. Ethylene oxide dermatitis. Cutis 1984;34:20, 22, 24.
- 860. Jay WM, Swift TR, Hull DS. Possible relationship of ethylene oxide exposure to cataract formation. Am. J. Ophthalmol. 1982;93:727-32.

- 861. Salinas E, Sasich L, Hall DH, Kennedy RM, Morriss H. Acute ethylene oxide intoxication. Drug Intell. Clin. Pharm. 1981;15:384-6.
- 862. Marchand M, Delesvonx R, Claeys C. The toxicity of ethylene oxide and a report on three fatal cases of poisoning. Am. Arch. Indust. Health 1958;18:60.
- 863. Finelli PF, Morgan TF, Yaar I, Granger CV. Ethylene oxide-induced polyneuropathy. A clinical and electrophysiologic study. Arch. Neurol. 1983;40:419-21.
- 864. Estrin WJ, Becker CE. Evidence of neurologic dysfunction related to long-term ethylene oxide exposure. Arch. Neurol. 1987;44:1283-6.
- 865. Estrin WJ, Bowler RM, Lash A, Becker CE. Neurotoxicological evaluation of hospital sterilizer workers exposed to ethylene oxide. J. Toxicol. Clin. Toxicol. 1990;28:1-20.
- 866. Crystal HA, Schaumburg HH, Grober E, Fuld PA, Lipton RB. Cognitive impairment and sensory loss associated with chronic low-level ethylene oxide exposure. Neurology 1988;38:567-9.
- 867. Shaham J, Levi Z, Gurvich R, Shain R, Ribak J. Hematological changes in hospital workers due to chronic exposure to low levels of ethylene oxide. J. Occup. Environ. Med. 2000;42:843-50.
- 868. Lindbohm ML, Hemminki K, Bonhomme MG, et al. Effects of paternal occupational exposure on spontaneous abortions. Am. J. Public Health 1991;81:1029-33.
- 869. Hemminki K, Mutanen P, Saloniemi I, Niemei M-L, Vainio H. Spontaneous abortions in hospital staff engaged in sterilising instruments with chemical agents. Br. Med. J. 1982;285:1461-3.
- 870. Rowland AS, Baird DD, Shore DL, Darden B, Wilcox AJ. Ethylene oxide exposure may increase the risk of spontaneous abortion, preterm birth, and postterm birth. Epidemiology 1996;7:363-8.
- 871. <u>National Toxicology Program (https://ntp.niehs.nih.gov/</u>) [Current version of this document may differ from original.]
- 872. Anonymous. Ethylene oxide sterilization: How hospitals can adapt to the changes. Health Devices 1994;23:485-92.
- 873. Occupational Safety and Health Administration. Ethylene Oxide: OSHA Fact Sheet: Occupational Safety and Health Administration, 2002.
- 874. Cardenas-Camarena L. Ethylene oxide burns from improperly sterilized mammary implants. Ann. Plast. Surg. 1998;41:361-9.
- 875. Windebank AJ, Blexrud MD. Residual ethylene oxide in hollow fiber hemodialysis units is neurotoxic in vitro. Ann. Neurol. 1989;26:63-8.
- 876. Occupational Health and Safety Administration. Chemical sampling information-Ethylene chlorohydrin: Occupational Safety and Health Administration, 2002.
- 877. Parisi AN, Young WE. Sterilization with ethylene oxide and other gases. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lea & Febiger, 1991:580-95.
- 878. Ries MD, Weaver K, Beals N. Safety and efficacy of ethylene oxide sterilized polyethylene in total knee arthroplasty. Clin. Orthop. 1996:159-63.
- 879. Alfa MJ, DeGagne P, Olson N. Bacterial killing ability of 10% ethylene oxide plus 90% hydrochlorofluorocarbon sterilizing gas. Infect. Control Hosp. Epidemiol. 1997;18:641-5.
- 880. Parker HH, Johnson RB. Effectiveness of ethylene oxide for sterilization of dental handpieces. J. Dent. 1995;23:113-5.
- 881. Jacobs PT, Lin SM. Sterilization processes utilizing low-temperature plasma. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:747-63.
- 882. Rutala WA, Gergen MF, Weber DJ. Sporicidal activity of a new low-temperature sterilization technology: the Sterrad 50 sterilizer. Infect. Control Hosp. Epidemiol. 1999;20:514-6.
- 883. Kyi MS, Holton J, Ridgway GL. Assessment of the efficacy of a low temperature hydrogen peroxide gas plasma sterilization system. J. Hosp. Infect. 1995;31:275-84.
- 884. Jacobs PT, Smith D. The new Sterrad 100S sterilization system: Features and advantages. Zentr. Steril. 1998;6:86-94.
- 885. Rudolph H, Hilbert M. Practical testing of the new plasma sterilizer "Sterrad 100S" in the Diakonkrankenhaus Rotenburg. Zentr. Steril. 1997;5:207-15.
- 886. Bar W, Marquez de Bar G, Naumann A, Rusch-Gerdes S. Contamination of bronchoscopes with *Mycobacterium tuberculosis* and successful sterilization by low-temperature hydrogen peroxide plasma sterilization. Am. J. Infect. Control 2001;29:306-11.
- 887. Centers for Disease Control and Prevention. Corneal decompensation after intraocular ophthalmic surgery-Missouri, 1998. MMWR 1998;47:306-9.

- 888. Duffy RE, Brown SE, Caldwell KL, et al. An epidemic of corneal destruction caused by plasma gas sterilization. Arch. Ophthalmol. 2000;118:1167-76.
- 889. Jarvis WR. Hospital Infections Program, Centers for Disease Control and Prevention: On-site outbreak investigations, 1990-1999: How often are germicides or sterilants the source? In: Rutala WA, ed. Disinfection, sterilization and antisepsis: principles and practices in healthcare facilities. Washington, DC: Association for Professional in Infection Control and Epidemiology, 2001:41-8.
- 890. Borneff M, Ruppert J, Okpara J, et al. Efficacy testing of low-temperature plasma sterilization (LTP) with test object models simulating practice conditions. Zentr. Steril. 1995;3:361-71.
- 891. Borneff-Lipp M, Okpara J, Bodendorf M, Sonntag HG. Validation of low-temperature-plasma (LPT) sterilization systems: Comparison of two technical versions, the Sterrad 100, 1.8 and the 100S. Hygiene und Mikrobiologie 1997;3:21-8.
- 892. Roberts C, Antonoplos P. Inactivation of human immunodeficiency virus type 1, hepatitis A virus, respiratory syncytial virus, vaccinia virus, herpes simplex virus type 1, and poliovirus type 2 by hydrogen peroxide gas plasma sterilization. Am. J. Infect. Control 1998;26:94-101.
- Okpara-Hofmann J, Knoll M, Durr M, Schmitt B, Borneff-Lipp M. Comparison of low-temperature hydrogen peroxide gas plasma sterilization for endoscopes using various Sterrad models. J. Hosp. Infect. 2005;59:280-5.
- 894. Timm D, Gonzales D. Effect of sterilization on microstructure and function of microsurgical scissors. Surg. Serv. Management 1997;3:47-9.
- 895. Feldman LA, Hui HK. Compatibility of medical devices and materials with low-temperature hydrogen peroxide gas plasma. Med. Dev. Diag. Indust. 1997;19:57-62.
- 896. Muscarella LF. Leading a horse to water: Are crucial lessons in endoscopy and outbreak investigations being learned? Infect . Control Hosp. Epidemiol. 2002;23:358-60.
- 897. Gurevich I, Qadri SMH, Cunha BA. False-positive results of spore tests from improper clip use with the Steris chemical sterilant system. Infect. Control Hosp. Epidemiol. 1992;21:42-3.
- 898. Kralovic RC. Use of biological indicators designed for steam or ethylene oxide to monitor a liquid chemical sterilization process. Infect. Control Hosp. Epidemiol. 1993;14:313-9.
- 899. Bond WW. Biological indicators for a liquid chemical sterilizer. Infect. Control Hosp. Epidemiol. 1993;14:565.
- 900. Bond WW. Biological indicators for a liquid chemical sterilizer: a solution to the instrument reprocessing problem? Infect. Control Hosp. Epidemiol. 1993;14:309-12.
- 901. Malchesky PS. Biological indicators for a liquid chemical sterilizer. Infect. Control Hosp. Epidemiol. 1993;14:563-6.
- 902. Daschner F. STERIS SYSTEM 1 in Germany. Infect. Control Hosp. Epidemiol. 1994;15:294, 296.
- 903. Sorin M, Segal-Maurer S, Urban C, Combest A, Rahal JJ. Nosocomial transmission of imipenem-resistant *Pseudomonas aeruginosa* following bronchoscopy associated with improper connection to the STERIS System 1 Processor. Infect. Control Hosp. Epidemiol. 2001;22:409-13.
- 904. Food and Drug Administration, Division of General and Restorative Devices. Guidance on Premarket Notification [510(k] Submissions for Sterilizers Intended for Use in Health Care Facilities. Rockville, MD. 1993.
- 905. Vickery K, Deva AK, Zou J, Kumaradeva P, Bissett L, Cossart YE. Inactivation of duck hepatitis B virus by a hydrogen peroxide gas plasma sterilization system: laboratory and 'in use' testing. J. Hosp. Infect. 1999;41:317-22.
- 906. Vassal S, Favennec L, Ballet JJ, Brasseur P. Hydrogen peroxide gas plasma sterilization is effective against *Cryptosporidium parvum* oocysts. Am. J. Infect. Control 1998;26:136-8.
- 907. Penna TC, Ferraz CA, Cassola MA. The presterilization microbial load on used medical devices and the effectiveness of hydrogen peroxide gas plasma against *Bacillus subtilis* spores. Infect. Control Hosp. Epidemiol. 1999;20:465-72.
- 908. Bryce EA, Chia E, Logelin G, Smith JA. An evaluation of the AbTox Plazlyte Sterilization System. Infect. Control Hosp. Epidemiol. 1997;18:646-53.
- 909. Graham GS, Riley R. Sterilization manufacturers: Interactions with regulatory agencies. In: Rutala WA, ed. Disinfection, sterilization, and antisepsis in healthcare. Champlain, New York: Polyscience Publications, 1998:41-8.
- 910. Royce A, Bowler C. Ethylene oxide sterilisation-some experiences and some practical limitations. J. Pharm. Pharmacol. 1961;13:87t-94t.
- 911. Nystrom B. Disinfection of surgical instruments. J. Hosp. Infect. 1981;2:363-8.

- 912. Rutala WA, Gergen MF, Jones JF, Weber DJ. Levels of microbial contamination on surgical instruments. Am. J. Infect. Control 1998;26:143-5.
- 913. Alfa MJ, Nemes R. Inadequacy of manual cleaning for reprocessing single-use, triple-lumen sphinctertomes: Simulated-use testing comparing manual with automated cleaning methods. Am. J. Infect. Control 2003;31:193-207.
- 914. Alfa MJ, Nemes R. Reprocessing of lumened instruments. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington DC: Association for Professionals in Infection Control and Epidemiology, 2004:189-99.
- 915. Bargmann LS, Bargmann BC, Collier JP, Currier BH, Mayor MB. Current sterilization and packaging methods for polyethylene. Clin. Orthop. 1999:49-58.
- 916. Williams IR, Mayor MB, Collier JP. The impact of sterilization method on wear in knee arthroplasty. Clin. Orthop. 1998:170-80.
- 917. Russell AD. Ionizing radiation. In: Russell AD, Hugo WB, Ayliffe GAJ, eds. Principles and practices of disinfection, preservation and sterilization. Oxford: Blackwell Science, 1999:675-87.
- 918. Hansen JM, Shaffer HL. Sterilization and preservation by radiation sterilization. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:729-46.
- 919. Lagergren ER. Recent advances in sterilization. J Infect Control (Asia) 1998;1:11-3.
- 920. Perkins JJ. Principles and methods of sterilization in health sciences. Springfield, IL: Charles C Thomas, 1969.
- 921. Favero MS, Bond WW. The use of liquid chemical germicides. In: Morrissey RF, Phillips GB, eds. Sterilization technology: A practical guide for manufacturers and users of health care products. New York: Van Nostrand Reinhold, 1993:309-334.
- 922. Muscarella LF. Are all sterilization processes alike? AORN J. 1998;67:966-70, 973-6.
- 923. Levy RV. Sterile filtration of liquids and gases. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:795-822.
- 924. Wallhausser KH. Is the removal of microorganisms by filtration really a sterilization method? J. Parenter. Drug Assoc. 1979;33:156-70.
- 925. Webb BC, Thomas CJ, Harty DW, Willcox MD. Effectiveness of two methods of denture sterilization. J. Oral Rehabil. 1998;25:416-23.
- 926. Rohrer MD, Bulard RA. Microwave sterilization. J. Am. Dent. Assoc. 1985;110:194-8.
- 927. Rohrer MD, Terry MA, Bulard RA, Graves DC, Taylor EM. Microwave sterilization of hydrophilic contact lenses. Am. J. Ophthalmol. 1986;101:49-57.
- 928. Douglas C, Burke B, Kessler DL, Cicmanec JF, Bracken RB. Microwave: practical cost-effective method for sterilizing urinary catheters in the home. Urology 1990;35:219-22.
- 929. Kindle G, Busse A, Kampa D, Meyer-Konig U, Daschner FD. Killing activity of microwaves in milk. J. Hosp. Infect. 1996;33:273-8.
- 930. Harris MG, Rechberger J, Grant T, Holden BA. In-office microwave disinfection of soft contact lenses. Optom. Vis. Sci. 1990;67:129-32.
- 931. Mervine J, Temple R. Using a microwave oven to disinfect intermittent-use catheters. Rehabil. Nurs. 1997;22:318-20.
- 932. Najdovski L, Dragas AZ, Kotnik V. The killing activity of microwaves on some non-sporogenic and sporogenic medically important bacterial strains. J. Hosp. Infect. 1991;19:239-47.
- 933. Rosaspina S, Salvatorelli G, Anzanel D, Bovolenta R. Effect of microwave radiation on *Candida albicans*. Microbios 1994;78:55-9.
- 934. Welt BA, Tong CH, Rossen JL, Lund DB. Effect of microwave radiation on inactivation of *Clostridium sporogenes* (PA 3679) spores. Appl. Environ. Microbiol. 1994;60:482-8.
- 935. Latimer JM, Matsen JM. Microwave oven irradiation as a method for bacterial decontamination in a clinical microbiology laboratory. J. Clin. Microbiol. 1977;6:340-2.
- 936. Sanborn MR, Wan SK, Bulard R. Microwave sterilization of plastic tissue culture vessels for reuse. Appl. Environ. Microbiol. 1982;44:960-4.
- 937. Rosaspina S, Salvatorelli G, Anzanel D. The bactericidal effect of microwaves on *Mycobacterium bovis* dried on scalpel blades. J. Hosp. Infect. 1994;26:45-50.
- 938. Engelhardt JP, Grun L, Dahl HJ. Factors affecting sterilization in glass bead sterilizers. J. Endod. 1984;10:465-70.
- 939. Smith GE. Glass bead sterilization of orthodontic bands. Am. J. Orthod. Dentofacial Orthop. 1986;90:243-9.

- 940. Sisco V, Winters LL, Zange LL, Brennan PC. Efficacy of various methods of sterilization of acupuncture needles. J. Manip. Physiol. Therap. 1988;11:94-7.
- 941. Klapes NA, Vesley D. Vapor-phase hydrogen peroxide as a surface decontaminant and sterilant. Appl. Environ. Microbiol. 1990;56:503-6.
- 942. French GL, Otter JA, Shannon KP, Adams NMT, Watling D, Parks MJ. Tackling contamination of the hospital environment by methicillin-resistant *Staphylococcus aureus* (MRSA): A comparison between conventional terminal cleaning and hydrogen peroxide vapour decontamination. J. Hosp. Infect. 2004;57:31-7.
- 943. Jeanes A, Rao G, Osman M, Merrick P. Eradication of persistent environmental MRSA. J. Hosp. Infect. 2005;61:85-6.
- 944. Bates CJ, Pearse R. Use of hydrogen peroxide vapour for environmental control during a *Serratia* outbreak in a neonatal intensive care unit. J. Hosp. Infect. 2005;61:364-6.
- 945. Boyce JM, Havill NL, Otter JA, et al. Impact of hydrogen peroxide vapor room bio-decontamination on environmental contamination and nosocomial transmission of *Clostridium difficile*. The Society of Healthcare Epidemiology of America, 2006;Abstract 155:109.
- 946. Berrington AW, Pedler SJ. Investigation of gaseous ozone for MRSA decontamination of hospital side-rooms. J. Hosp. Infect. 1998;40:61-5.
- 947. Gaspar MC, Pelaez B, Fernandez C, Fereres J. Microbiological efficacy of Sterrad 100S and LTSF sterilisation systems compared to ethylene oxide. Zentr. Steril. 2002;10:91-9.
- 948. Kanemitsu K, Kunishima H, Imasaka T, et al. Evaluation of a low-temperature steam and formaldehyde sterilizer. J. Hosp. Infect. 2003;55:47-52.
- 949. Kanemitsu K, Imasaka T, Ishikawa S, et al. A comparative study of ethylene oxide gas, hydrogen peroxide gas plasma, and low-temperature steam formaldehyde sterilization. Infect Control Hosp Epidemiol 2005;26:486-9.
- 950. Roncoroni AJ, Casewell MW, Phillips I. The disinfection of clinically contaminated Matburn suction pumps and baby incubators in an 'Aseptor' formalin cabinet. J. Hosp. Infect. 1980;1:251-9.
- 951. Cumberland NS, Botting FG. Formaldehyde vapour cabinets. J. Hosp. Infect. 1991;19:67-70.
- 952. Jeng DK, Woodworth AG. Chlorine dioxide gas sterilization of oxygenators in an industrial scale sterilizer: a successful model. Artif. Organs 1990;14:361-8.
- 953. Knapp JE, Battisti DL. Chloride dioxide. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:215-27.
- 954. Kowalski JB. Sterilization of medical devices, pharmaceutical components, and barrier isolation systems with gaseous chlorine dioxide. In: Morrissey RF, Kowalski JB, eds. Sterilization of medical products,. Champlain, NY: Polyscience Publications, 1998:313-23.
- 955. Portner DM, Hoffman RK. Sporicidal effect of peracetic acid vapor. Appl. Microbiol. 1968;16:1782-5.
- 956. Mata-Portuguez VH, Perez LS, Acosta-Gio E. Sterilization of heat-resistant instruments with infrared radiation. Infect . Control Hosp. Epidemiol. 2002;23.
- 957. Frey R. The structural and functional prerequisites for a central sterile supply department (CSSD). Zentr. Steril. 2000;8:128-40.
- 958. Reich RR, Fleming W, Burgess DJ. Sterilization validation: it's not just for industry. Infect. Control Steril. Technol. 1996;2.
- 959. American Institute of Architects. Guidelines for design and construction of hospital and health care facilities. Washington, DC: The American Institute of Architects Press, 2001.
- 960. DesCoteaux JG, Poulin EC, Julien M, Guidoin R. Residual organic debris on processed surgical instruments. AORN J. 1995;62:23-30.
- 961. Rutala WA, Weber DJ. A review of the use of gowns and drapes (single use and reusable) in healthcare. Infect. Control Hosp. Epidemiol. 2001;22:248-57.
- 962. Association of peri-Operative Registered Nurses. Recommended practices for sterilization in the perioperative practice setting. AORN J. 2006;83:700-22.
- 963. Taurasi R. Comfortable PPE? Maximum tray weight? Healthcare Purchasing News 2004; July: 48.
- 964. Chobin N, Furr D, Nuyttens A. Wet packs and plastic accessory cases. Infect Control Today 2004;August:24, 28-30.
- 965. Dunkelberg H, Fleitmann-Glende F. Measurement of the microbial barrier effectiveness of sterilization containers in terms of the log reduction value for prevention of nosocomial infections. Am. J. Infect. Control 2006;34:285-9.
- 966. Rutala WA, Weber DJ. Choosing a sterilization wrap. Infect. Control Today 2000;4:64,70.

- 967. Maloney JM, Kohut RD. Infection control, barrier protection and the treatment environment. Dent. Hyg. (Chic). 1987;61:310-3.
- 968. Mayworm D. Sterile shelf life and expiration dating. J. Hosp. Supply, Process. Distri. 1984;2:32-5.
- 969. Cardo DM, Sehulster LM. Central sterile supply. In: Mayhall CG, ed. Infect. Control and Hosp. Epidemiol. Philadelphia: Lippincott Williams & Wilkins, 1999:1023-30.
- 970. Klapes NA, Greene VW, Langholz AC. Microbial contamination associated with routine aseptic practice. J. Hosp. Infect. 1987;10:299-304.
- 971. Butt WE, Bradley DV, Jr., Mayhew RB, Schwartz RS. Evaluation of the shelf life of sterile instrument packs. Oral Surg. Oral Med. Oral Pathol. 1991;72:650-4.
- 972. Webster J, Lloyd W, Ho P, Burridge C, George N. Rethinking sterilization practices: Evidence for eventrelated outdating. Infect Control Hosp Epidemiol 2003;24:622-4.
- 973. Widmer AF, Houston A, Bollinger E, Wenzel RP. A new standard for sterility testing for autoclaved surgical trays. J. Hosp. Infect. 1992;21:253-60.
- 974. Schneider PM, Reich RR, Kirckof SS, Foltz WG. Perfomance of various steam sterilization indicators under optimum and sub-optimum exposure conditions. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington DC: Association for Professionals in Infection Control and Epidemiology, 2004:200-23.
- 975. Greene VW. Control of sterilization process. In: Russell AD, Hugo WB, Ayliffe GAJ, eds. Principles and practice of disinfection, preservation and sterilization. Oxford, England: Blackwell Scientific Publications, 1992:605-24.
- 976. Vesley D, Nellis MA, Allwood PB. Evaluation of a rapid readout biological indicator for 121°C gravity and 132°C vacuum-assisted steam sterilization cycles. Infect. Control Hosp. Epidemiol. 1995;16:281-6.
- 977. Rutala WA, Jones SM, Weber DJ. Comparison of a rapid readout biological indicator for steam sterilization with four conventional biological indicators and five chemical indicators. Infect. Control Hosp. Epidemiol. 1996;17:423-8.
- 978. Alfa MJ, Olson N, DeGagne P, Jackson M. Evaluation of rapid readout biological indicators for 132°C gravity and 132°C vacuum-assisted steam sterilization cycles using a new automated fluorescent reader. Infect . Control Hosp. Epidemiol. 2002;23:388-92.
- 979. Koncur P, Janes JE, Ortiz PA. 20 second sterilization indicator tests equivalent to BIs. Infect. Control Steril. Technol. 1998:26-8, 30, 32-4.
- 980. Perkins RE, Bodman HA, Kundsin RB, Walter CW. Monitoring steam sterilization of surgical instruments: a dilemma. Appl. Environ. Microbiol. 1981;42:383-7.
- 981. Kotilainen HR, Gantz NM. An evaluation of three biological indicator systems in flash sterilization. Infect. Control 1987;8:311-6.
- 982. Kleinegger CL, Yeager DL, Huling JK, Drake DR. The effects of contamination on biological monitoring. Infect. Control Hosp. Epidemiol. 2001;22:391-2.
- 983. Centers for Disease Control. False-positive results of spore tests in ethylene oxide sterilizers Wisconsin. MMWR 1981;30:238-40.
- 984. Association of Operating Room Nurses. AORN standards and recommended practices for perioperative nursing. 1987:Section III:14.1-III:14.11, AORN, Denver, CO.
- 985. Gurevich I, Holmes JE, Cunha BA. Presumed autoclave failure due to false-positive spore strip tests. Infect. Control 1982;3:388-92.
- 986. Epstein BJ, Lattimer JM, Matsen JM, Garibaldi RA. False positive spore strip sterility tests with steam sterilization. Am. J. Infect. Control 1983;11:71-3.
- 987. Association for the Advancement of Medical Instrumentation. Good hospital practice: steam sterilization and sterility assurance. Arlington, VA: AAMI, 1988.
- 988. Baird RM. Sterility assurance: Concepts, methods and problems. In: Russell AD, Hugo WB, Ayliffe GAJ, eds. Principles and practice of disinfection, preservation and sterilization. Oxford, England: Blackwell Scientific Publications, 1999:787-99.
- 989. Coulter WA, Chew-Graham CA, Cheung SW, Burke FJT. Autoclave performance and operator knowledge of autoclave use in primary care: a survey of UK practices. J. Hosp. Infect. 2001;48:180-5.
- 990. Greene VW. Reuse of disposable devices. In: Mayhall CG, ed. Infect. Control and Hosp. Epidemiol. Philadelphia: Lippincott Williams & Wilkins, 1999:1201-8.
- 991. Avitall B, Khan M, Krum D, Jazayeri M, Hare J. Repeated use of ablation catheters: a prospective study. J. Am. Coll. Cardiol. 1993;22:1367-72.

- 992. Dunnigan A, Roberts C, McNamara M, Benson DW, Jr., Benditt DG. Success of re-use of cardiac electrode catheters. Am. J. Cardiol. 1987;60:807-10.
- 993. Aton EA, Murray P, Fraser V, Conaway L, Cain ME. Safety of reusing cardiac electrophysiology catheters. Am. J. Cardiol. 1994;74:1173-5.
- 994. Brown SA, Merritt K, Woods TO, McNamee SG, Hitchins VM. Effects of different disinfection and sterilization methods on tensile strength of materials used for single-use devices. Biomed. Instrum. Technol 2002;January/February:23-7.
- 995. Food and Drug Administration. Enforcement Priorities for Single-Use Devices Reprocessed by Third Parties and Hospitals, Rockville, MD., 2000.
- 996. Ulatowski TA. FDA: Reuse of Single-Use Devices. In: Rutala WA, ed. Disinfection, sterilization and antisepsis: Principles, practices, challenges, and new research. Washington, DC: Association for Professionals in Infection Control and Epidemiology, 2004:15-23.
- 997. Occupational Health and Safety Administration. Hazard Communication Standard.29 CFR 1910.1200, OSHA, Washington, DC.
- 998. Edens AL. Occupational Safety and Health Administration: Regulations affecting healthcare facilities. In: Rutala WA, ed. Disinfection, Sterilization and Antisepsis: Principles and practices in healthcare facilities. Washington, D.C,: Association for Professionals in Infection Control and Epidemiology, Inc., 2001:49-58.
- 999. Schultz JK. Decontamination: recommended practices. In: Reichert M, Young JH, eds. Sterilization technology for the health care facility. Gaithersburg, MD: Aspen Publication, 1997:10-20.
- 1000. Occupational Health and Safety Administration. Ethylene Oxide Standard. Vol. 29 CFR 1910.1047, OSHA, Washington, DC.
- 1001. Occupational Safety and Health Administration. Formaldehyde Standard. Vol. 29 CFR 1910.1048, Washington, DC.
- 1002. Buxton AE, Anderson RL, Werdegar D, Atlas E. Nosocomial respiratory tract infection and colonization with *Acinetobacter calcoaceticus*. Epidemiologic characteristics. Am. J. Med. 1978;65:507-13.
- 1003. Snydman DR. Hepatitis B infection from medical personnel. JAMA 1976;236:1009.
- 1004. Martiny H, Floss H. Residuals on medical devices following reprocessing. J. Hosp. Infect. 2001;48 (Supplement):S88-S92.
- 1005. Taylor DM. Inactivation of prions by physical and chemical means. J. Hosp. Infect. 1999;43 (supplement):S69-S76.
- 1006. Best M, Sattar SA, Springthorpe VS, Kennedy ME. Comparative mycobactericidal efficacy of chemical disinfectants in suspension and carrier tests. Appl. Environ. Microbiol. 1988;54:2856-8.
- 1007. Weber DJ, Rutala WA. Environmental issues and nosocomial infections. In: Wenzel RP, ed. Prevention and control of nosocomial infections. Baltimore: Williams and Wilkins, 1991:1042-64.
- 1008. Palmer PH, Yeoman DM. A study to assess the value of disinfectants when washing ward floors. Med. J. Aust. 1972;2:1237-9.
- 1009. Centers for Disease Control and Prevention. Preventing the spread of vancomycin resistance report from the Hospital Infection Control Practices Advisory Committee. Fed. Regist. 1994:25758-63.
- 1010. Ayliffe GA, Collins BJ, Lowbury EJ. Cleaning and disinfection of hospital floors. BMJ 1966;5511:442-5.
- 1011. Neely AN. A survey of gram-negative bacteria survival on hospital fabrics and plastics. J. Burn Care Rehabil. 2000;21:523-7.
- 1012. Rutala WA. Disinfection and sterilization of patient-care items. Infect. Control Hosp. Epidemiol. 1996;17:377-84.
- 1013. Environmental Protection Agency. Pesticides: Regulating Pesticides. [This link is no longer active: http://www.epa.gov/oppad001/chemregindex.htm], 2003.
- 1014. Hoffman PN, Layzell SK. Household bleach as disinfectant for use by injecting drug users. Lancet 1993;342:743.
- 1015. Chu NS, Chan-Myers H, Ghazanfari N, Antonoplos P. Levels of naturally occurring microorganisms on surgical instruments after clinical use and after washing. Am. J. Infect. Control 1999;27:315-9.
- 1016. Hoffmann KK, Weber DJ, Rutala WA. Pseudoepidemic of Rhodotorula rubra in patients undergoing fiberoptic bronchoscopy. Infect. Control Hosp. Epidemiol. 1989;10:511-4.
- 1017. Lowry PW, Jarvis WR. Use of tap water and disinfection practices in outpatient settings. A survey of otolaryngologists. Arch. Otolaryngol. Head Neck Surg. 1991;117:886-8.
- 1018. Fahey BJ, Koziol DE, Banks SM, Henderson DK. Frequency of nonparenteral occupational exposures to blood and body fluids before and after universal precautions training. Am. J. Med. 1991;90:145-53.

- 1019. Beekmann SE, Vlahov D, Koziol DE, McShalley ED, Schmitt JM, Henderson DK. Temporal association between implementation of universal precautions and a sustained, progressive decrease in percutaneous exposures to blood. Clin. Infect. Dis. 1994;18:562-9.
- 1020. Gerberding JL, Littell C, Tarkington A, Brown A, Schecter WP. Risk of exposure of surgical personnel to patients' blood during surgery at San Francisco General Hospital. N. Engl. J. Med. 1991;324:1788-93.
- 1021. Mast ST, Woolwine JD, Gerdberding JL. Efficacy of gloves in reducing blood volumes transferred during simulated needlestick injury. J. Infect. Dis. 1993;168:1589-92.
- 1022. Wendt C, Herwaldt LA. Epidemics: Identification and Management. In: Wenzel RP, ed. Prevention and Control of Nosocomial Infections. Baltimore: Williams & Wilkins, 1997:175-214.
- 1023. Feigal DW, Gardner SN, McClellan M. Ensuring safe and effective medical devices. N. Engl. J. Med. 2003;348:191-2.
- 1024. Oie S, Kamiya A. Microbial contamination of antiseptics and disinfectants. Am. J. Infect. Control 1996;24:389-95.
- 1025. Strzelecki LR, Nelson JH. Evaluation of closed container flash sterilization system. Orthop. Nurs. 1989;8:21-4.
- 1026. Burgess DJ, Reich RR. Industrial ethylene oxide sterilization. In: Morrissey RF, Phillips GB, eds. Sterilization technology: a practical guide for manufacturers and users of health care product. New York: Van Nostrand Reinhold, 1993:120-51.
- 1027. Conviser CA, CW. Ethylene oxide sterilization: sterilant alternatives. In: Reichert M, Young JH, eds. Sterilization technology for the health care facility. Gaithersburg, MD: Aspen Publication, 1997:189-99.
- 1028. Young JH. Steam sterilization: scientific principles. In: Reichert M, Young JH, eds. Sterilization technology for the health care facility. Gaithersburg, MD: Aspen Publication, 1997:123-144.
- 1029. Alfa MJ. Importance of lumen flow in liquid chemical sterilization. Am. J. Infect. Control 1999;27:373-5.
- 1030. Mallison GF, Standard PG. Safe storage times for sterile packs. Hospitals 1974;48:77-8, 80.
- 1031. Klapes NA, Greene VW, Langholz AC, Hunstiger C. Effect of long-term storage on sterile status of devices in surgical packs. Infect. Control 1987;8:289-93.
- 1032. Japp NF. Packaging: Shelf life. In: Reichert M, Young JH, eds. Sterilization Technology. Gaithersburg, Maryland: Aspen, 1997:99-102.
- 1033. Joint Commission for the Accreditation of Healthcare Organizations. Comprehensive accreditation manual for hospitals, JCAHO, Chicago, IL. 2003.
- 1034. Block SS. Definition of terms. In: Block SS, ed. Disinfection, sterilization, and preservation. Philadelphia: Lippincott Williams & Wilkins, 2001:19-28.
- 1035. Molinari JA, Gleason MJ, Cottone JA, Barrett ED. Comparison of dental surface disinfectants. Gen. Dent. 1987;35:171-5.