



# Softgel Polisher



# Softgel Polisher

CPP

## Features

- The capsule Polisher is a machine capable of polishing softgel capsules and can be operated in two modes: Manual or Automatic. Both modes can be combined if required at individuals times as semi-automatic.
- The basket can be operated in forward or reverse mode.
- Stack light and alarm buzzer to prompt the operator when the cycle is completed, including any other faults like basket motor.
- Reset button will start the cycle from the beginning once pressed.
- The basket lifting cover operates in safety conditions.
- Tumbler basket capacity: up to 120Kg or softgel capsules.
- Electrical panel incorporated within the machine.
- Sealed construction to prevent leakage and subsequent contamination.
- Aluminum frame construction prevents any possible rust and allows easy cleaning and maintenance.
- Outer covers made completely from stainless steel, designed according to GMP/FDA requirements.
- Ideal for both pharmaceutical and nutraceutical products.
- Easy access for set up and cleaning operations.
- No special electrical connections are required.
- Plug and play.
- Portable.


## Mechanism

- A variable-speed electrical motor rotates the basket using a combination of belt and plastic gear transmission.
- The capsules are evacuated by changing the direction of the rotation of the basket.
- A pneumatic cylinder is assisting the removal of the basket for cleaning.

# Softgel Polisher

CPP

## Specification

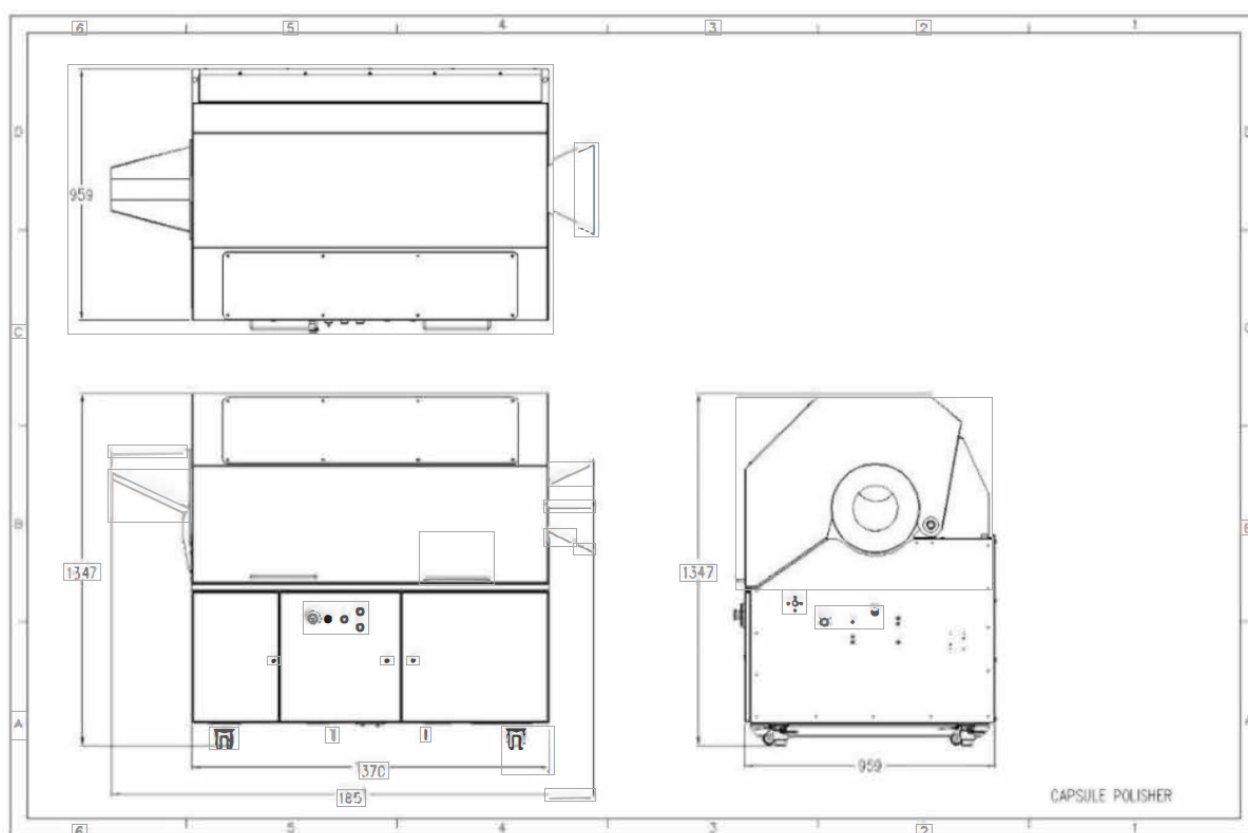

| CPP              |                           |
|------------------|---------------------------|
| Dimensions (mm)  | 1,851 W X 959 D X 1,347 H |
| Weight (kg)      | 500                       |
| Electrical Power | 305 KVA                   |

## TECHNOPHAR

Technophar Global Headquarters  
3293 St. Etienne Blvd.  
Windsor, Ontario, Canada N8W 5B1  
Phone: +1 519-946-0800  
[www.technophar.com](http://www.technophar.com)