



Softgel Tumble Dryer



TD24

Softgel Tumble Dryer

TD24

Features

- Allows longer drying time, thus extracting more water from the capsules and reducing the Drying Tunnels drying time, compare with the one row and less baskets Tumbler dryer model
- Anodized aluminum frame and stainless steel covers with number 4 finish for fast, thorough cleaning
- Maximum basket load - 12 kg.
- Built in heating available upon request, allowing for gentle heating of specific products
- Tumbler dryers are offered as stand alone units or in combination with Technophar machines
- Light-weight baskets with stainless steel perforated plate allow for quick changeovers, easy clean-up and inspection
- Lightweight, simple to install ,clear Lexan covers
- Easy lubrication from the side of the machine
- Discharge Conveyor allows transferring the capsules from the Tumbler Dryer to the Spreading Table.
- Blower for transferring the capsules from the Encapsulation Machine to the Tumbler Dryer available.


Mechanism

- Electrical panel and PLC are incorporated within the stand alone unit
- Two (2) variable speed direct drive electrical motors rotate the baskets
- Pneumatic actuated doors allow capsules to be transferred throughout the tumbler and be discharged at the end of the drying cycle
- Programmable drying and capsule transfer time

Softgel Tumble Dryer

TD24

Specification

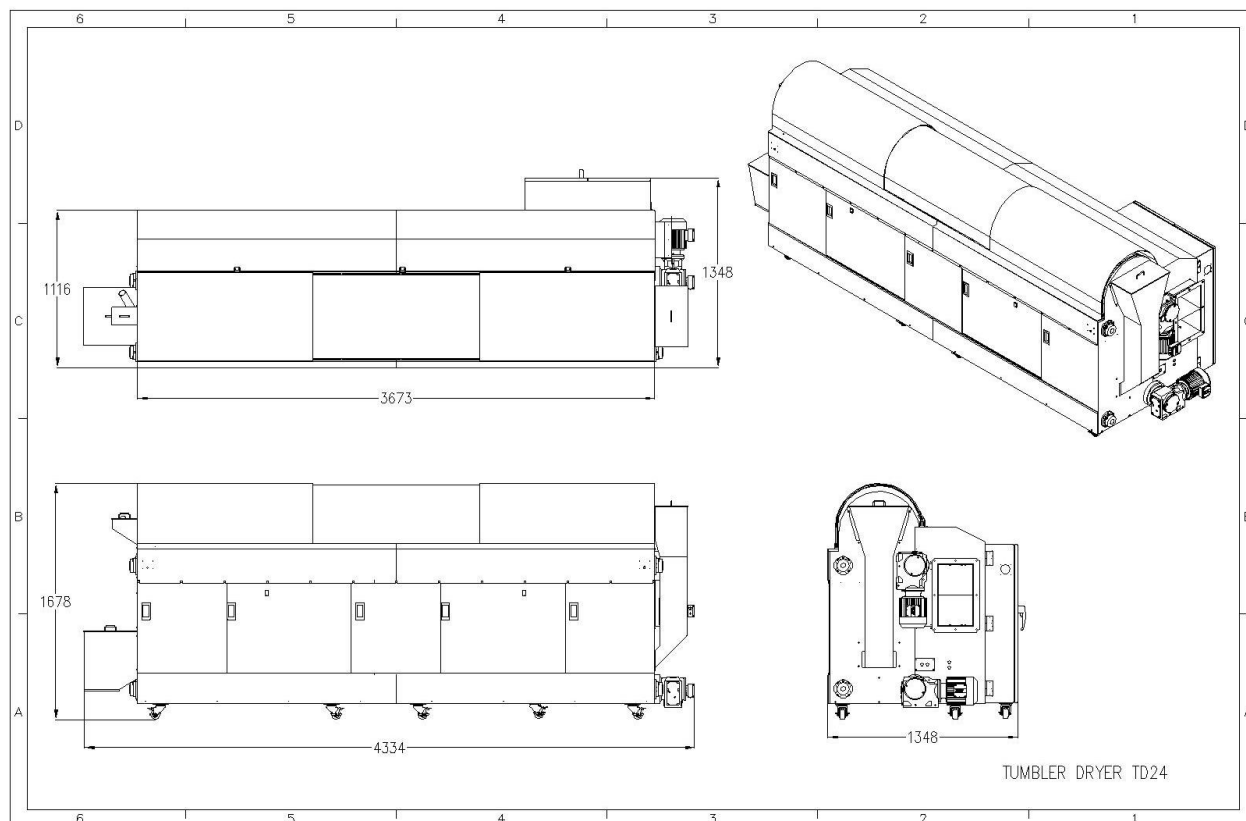

TD24	
Dimensions (mm)	4334 W X 1348 D X 1678 H
Weight (kg)	2500
TD Electrical Power	14 KVA without heaters
Compressed Air	8 Bar

TECHNOPHAR

Technophar Global Headquarters
3293 St. Etienne Blvd.
Windsor, Ontario, Canada N8W 5B1
Phone: +1 519-946-0800
www.technophar.com