



CMCcare Thumb Brace




The CMCcare Thumb Brace reduces thumb abduction and MCP hyperextension of a thumb “swan-neck” deformity. The embedded stay contours around the thenar eminence to realign the first metacarpal and stabilize the CMC joint. The low-profile design leaves the ulnar border of the hand free for unimpeded function and features a replaceable strap to secure the brace on the hand.

Features

- Easily customize fit by squeezing the built-in stay for increased stabilization at the CMC joint
- Low-profile design leaves the ulnar border free for unimpeded hand function
- Strap is adjustable at both ends allowing the user to determine how much force is needed to align the thumb in a functional position
- Durable TPE material resists stains and abrasions, can be worn in water and is ideal for all levels of activity
- Includes two closures tabs and a replacement strap; additional straps available
- Hand wash only; latex-free

Order Information

Measure around the hand just below the MP joints

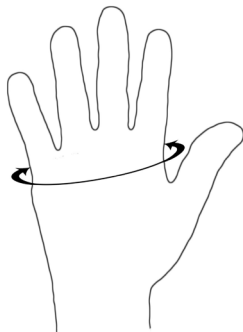

Small	6 1/4" to 7 5/8" (16.0 to 19.5CM)
Medium	7 5/8" to 8 1/8" (19.5 to 22.5CM)
Large	8 1/8" to 10 1/4" (22.5 to 26.0CM)

Size	Left	Right
Small	P2014-L2	P2014-R2
Medium	P2014-L3	P2014-R3
Large	P2014-L4	P2014-R4

Call 888-378-7763

Mon-Fri 8:30am - 5:00pm EST for ordering details

www.3pointproducts.com