



# CMCcare Thumb Brace




The CMCcare Thumb Brace reduces thumb abduction and MCP hyperextension of a thumb “swan-neck” deformity. The embedded stay contours around the thenar eminence to realign the first metacarpal and stabilize the CMC joint. The low-profile design leaves the ulnar border of the hand free for unimpeded function and features a replaceable strap to secure the brace on the hand.

## Features

- Easily customize fit by squeezing the built-in stay for increased stabilization at the CMC joint
- Low-profile design leaves the ulnar border free for unimpeded hand function
- Strap is adjustable at both ends allowing the user to determine how much force is needed to align the thumb in a functional position
- Durable TPU material resists stains and abrasions, can be worn in water and is ideal for all levels of activity
- Includes two closures tabs and a replacement strap; additional straps available
- Hand or machine wash only; latex-free


## Size Information

**A**= Measure across the MP joint (large knuckle)

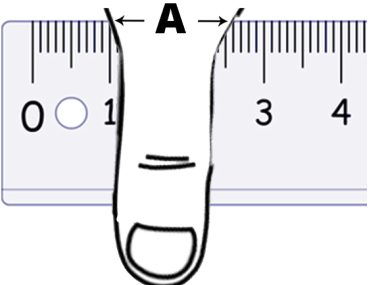

**Small** 5/8" - 1" (1.58 - 2.54 cm)  
**Medium** 7/8" - 1 1/8" (2.22 - 2.86 cm)  
**Large** 1" - 1 1/4" (2.54 - 3.18 cm)

**B**= Bend your thumb towards your palm. Measure from the top of the MP joint (large knuckle) to the wrist crease.

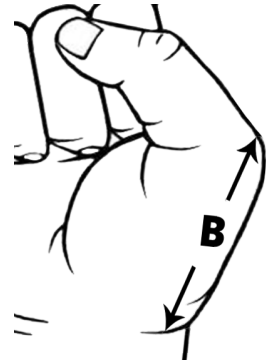

**Small** 1 1/2" - 2" (3.18 - 5.08 cm)  
**Medium** 1 7/8" - 2 1/4" (4.76 - 5.71 cm)  
**Large** 2 1/8" - 2 3/4" (5.08 - 6.98 cm)

The CMCcare sizes do overlap. If your measurements fall between sizes, follow these suggestions to choose the size most likely to fit:

- Choose the smaller size if your hands are thin and your thumb muscles are soft or flattened
- Choose the larger size if your hands are fuller and your thumb muscles are large and firm