


The CMCcare thumb brace reduces pain from CMC joint arthritis by realigning the first metacarpal and stabilizing the CMC joint. The embedded stay contours around the thumb for customized support. The molded brace is heat adjustable and the strap is adjustable and replaceable. The ideal support for comfortable, long-term reduction of CMC joint pain.

FEATURES

- Carefully contoured to support the CMC joint WITHOUT LIMITING WRIST motion or the ULNAR BORDER of the hand
- Fits below the MP crease to AVOID PRESSURE on the MP joint
- Contouring the EMBEDDED STAY allows for a CUSTOM FIT
- HEAT GUN ADJUSTABLE. PAD with self-adhesive SILICONE GEL
- The single STRAP IS ADJUSTABLE from both ends and it's REPLACEABLE
- TPU material is latex-free; hand or machine WASHABLE; stain and abrasion resistant; and WILL LAST FOR YEARS of regular use
- Includes REPLACEMENT STRAP and TWO CLOSURE STABS

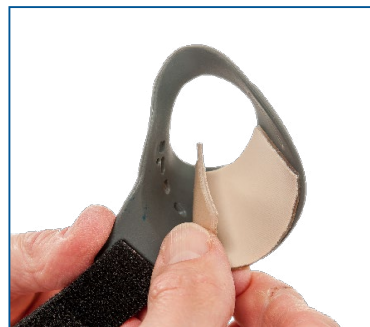

CMCcare Thumb Brace

Sizing Instructions



For the best fit and function you only need one measurement

Measure the length from your web space to your wrist crease

S	2" - 2 1/2" (5.1 - 6.4cm)
M	2 1/2" - 3" (6.4 - 7.6cm)
L	3" - 3 1/2" (7.6 - 8.9cm)

Sizing Tips:

If the measurement is near the end of a range, consider the following when choosing size:

Choose the smaller size if hands are thin and the thumb muscles are soft or flattened

Choose the larger size if hands are fuller and the thumb muscles are large and firm

It may be necessary to try two sizes on your patient for the best fit

TO ORDER

To order or request more information

888-378-7763

Mon-Fri 8:30 am-5:00pm EST

CMCcare thumb brace			
	Left		Right
Sm	P2014-L2	Sm	P2014-R2
Med	P2014-L3	Med	P2014-R3
Lrg	P2014-L4	Lrg	P2014-R4