



3pp[®] Side Step[™] Splint



NEW IMPROVED DESIGN!

Correct DIP or PIP joint deviation with the 3pp[®] Side Step[™] Splint. This adjustable tension splint applies a static progressive stretch to correct IP joint deviation caused by arthritis deformity or injury.

Features


- Cushioned slide adjusts to provide a stabilizing counterforce at the DIP or PIP.
- Aids in preventing the progression of Heberden's and Bouchard's nodes.
- Case hardened steel wire will not bend or break from tension.
- Apply static progressive force for short periods several times a day for successful correction.
- Soft foam-lined straps and counterforce pads are latex-free.

Order Information

Measure length from web space to fingertip and width at proximal phalanx. Fits right or left hand. If in between sizes, choose the larger size.

	Item Number	Web Space to Tip
Sm	P1201-2	1 1/2" - 2" (3.8-5cm)
Med	P1201-3	1 3/4" - 3 1/8" (4.4 - 8.3cm)
Lrg	P1201-4	3" - 3 1/2" (7.6-8.9cm)

Watch our **instructional video on the Side Step Splint** on our website: www.3pointproducts.com/product-videos-of-splints/3pp-side-step-video

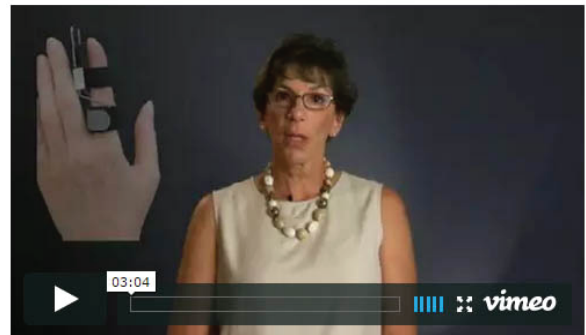

Call 888-378-7763

Mon-Fri 8:30am - 5:00pm EST for ordering details

www.3pointproducts.com