



3pp[®] Step Up[™] Splint



The 3pp[®] Step Up[™] Splint is designed to restore PIP and DIP extension using static progressive tension. Easily fit and adjustments can be made by the user for rapid reduction of contractures caused by jammed fingers, Boutonniere deformities or tendon injuries.

Features


- Accommodates contractures of 90 degrees or more and allows full correction to zero degrees.
- The angle of correction adjusts by repositioning the corrective strap along with the wire outrigger allowing the wearer to control the amount and angle of force.
- The palmar plate rotates to adjust the proper line of pull for any digit.
- Non-slip, foam-lined straps apply a secure and comfortable application of corrective force. Latex-free.

Order Information

Measure length from web space to fingertip and width at proximal phalanx. Fits right or left hand. If in between sizes, choose the larger size.

	Item Number	Web Space to Tip
Sm	P1200-2	1 1/2" - 2" (3.8-5cm)
Med	P1200-3	1 3/4" - 3 1/8" (4.4 - 8.3cm)
Lrg	P1200-4	3" - 3 1/2" (7.6-8.9cm)

Watch our **instructional video on the Step Up Splint** on our website: www.3pointproducts.com/product-videos-of-splints/3pp-step-up-video

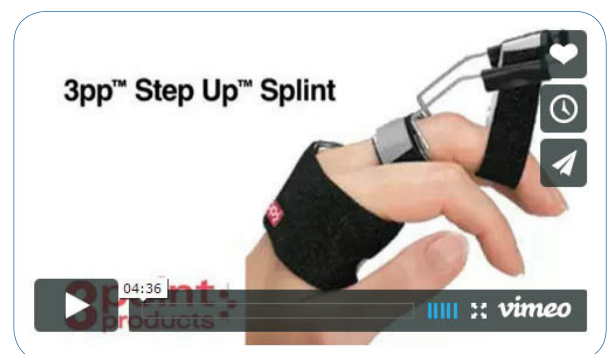

Call 888-378-7763

Mon-Fri 8:30am - 5:00pm EST for ordering details

www.3pointproducts.com