


DYNAWASH⁸²⁵ QUICK START

Before each use, inspect the lint filter and clean if necessary.

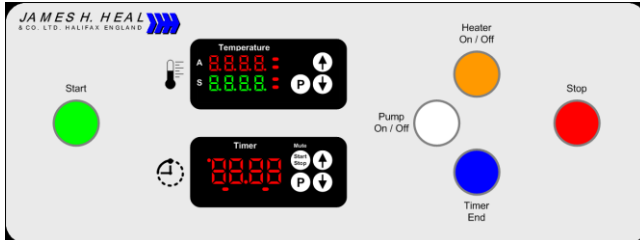


This is the **Dynawash** control panel.

The buttons and indicators are self-explanatory.

Ensure the bath is filled with water to the fill marker (approx. 40 litres).

Switch on the power and the display quickly goes through a self-test.

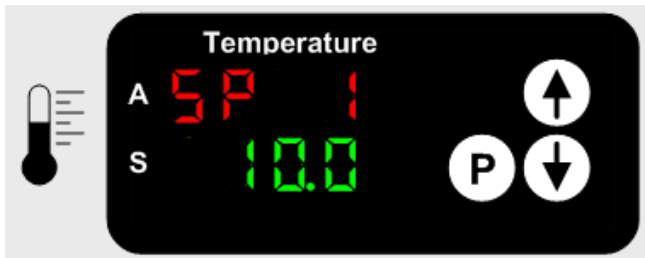


To set the temperature:

Press the **P** button on the Temperature controller.

SP 1 (set-point) is displayed.

Use the **↑** **↓** buttons to increase or decrease the temperature set-point.

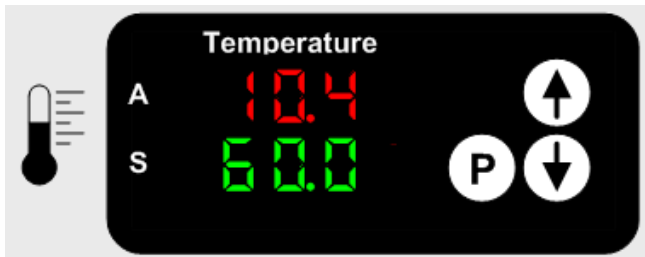

Press the **P** button again to Program the required temperature set-point.

A = Actual temperature of water

S = temperature **S**et-point





To set the timer:

Press the **P** button on the Timer controller. **t 1** is displayed.



Use the **↑** **↓** buttons to increase or decrease the time.

Press the **Start Stop** button to set the new Timer set-point.

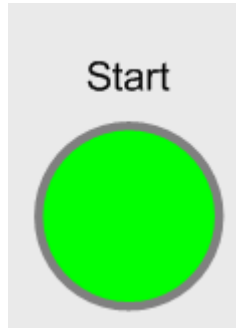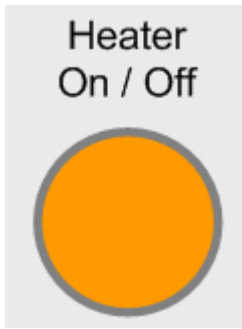


Press the **Heater** button to switch the electric heater on.

Press the **Start** button to start the impeller turning in the water bath.

Note: the impeller will not start if the lid is open.


Note: always run the impeller while heating.

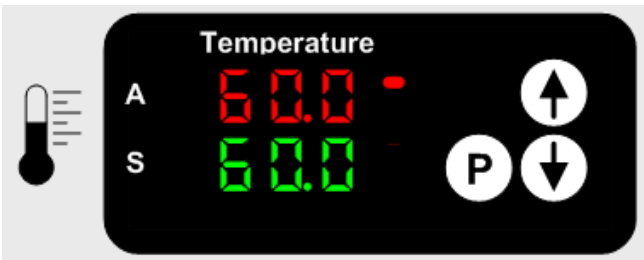

Prepare specimen, detergent and makeweights as required.

When set temperature is achieved, add detergent and run for approx. 1 minute.

If required, add makeweights and run for approx. 1 minute.

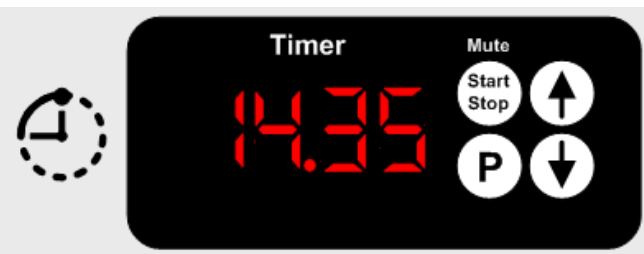


Add specimen, close lid and **Start** impeller.

Start the Timer by pressing the **Start Stop** button on the Timer.

The Timer will count down displaying *mm:ss*.

The Timer will count down to zero.

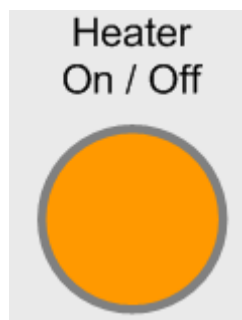

The **Timer End** indicator will illuminate and the buzzer will sound to alert the operator.


The impeller will **Stop** automatically.

To silence the buzzer press the **Mute (Start Stop)** button on the Timer.

Switch the electric heater off by pressing the **Heater** button.

Open the lid and remove the specimen first, then remove the makeweights.

Rinse and dry the specimen and makeweights, as required.

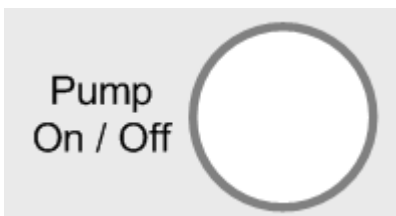

To empty the bath, ensure the heater is off.

Press the **Pump** button to drain the bath. This will take approx. 90s.

When the bath is empty, rinse with clean water while the pump is running.

When the bath is empty switch the **Pump** off.

Wipe the exterior of **Dynawash** with a soft dry cloth.