

Mirage Quattro™

FULL FACE MASK

FOURTH GENERATION TECHNOLOGY

Wide fit range: fourth generation full face mask from ResMed fits more than 95% of users.*

Mirage™ dual-wall cushion with spring air technology: unique cushion accommodates movement during sleep and provides an effective solution for jaw drop.

MicroFit dial: new forehead support adjustment offers personalised fitting and comfort.

Light, quiet and easy to use: these benefits are optimised with streamlined mask design, improved vent design and a range of user-friendly features.

MicroFit dial
provides better fit, seal and comfort for over 95% of users*

Set-and-forget headgear slots
save optimal fitting settings

Improved built-in vents
provide minimal noise output

Moving cushion membrane at nasal bridge
accommodates different types of nasal structures and prevents leak into the eyes


Cushion clip with tabs
allows easy removal of cushion from frame

Mirage dual-wall cushion with spring air technology
enhances seal where it's needed most

Moving cushion membrane at chin
maintains seal and accommodates jaw drop and movement

Set-and-forget headgear clips
allow convenient removal of mask

Simplified elbow with anti-asphyxia valve (AAV)
few elbow parts for easy assembly and disassembly; AAV enhances safety if the airflow to the mask is impeded

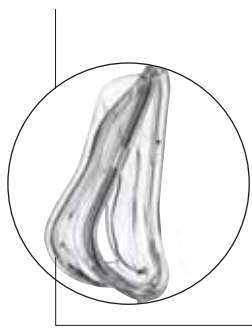

The membrane as seen on a side view of the mask cushion





TECHNOLOGY HIGHLIGHT

MicroFit dial

A choice of 24 positions at your fingertips—providing a comfortable, effective fit for over 95% of users*

- **Adjustable:** can be adjusted while on the user's face
- **Precise:** each turn of the dial adjusts the forehead support in small increments
- **Progressive:** moves in towards the face, allowing the user to sense when the fit is most comfortable for them
- **Works in combination with the cushion:** as the dial adjusts the angle of the mask to the face, the cushion adapts to enhance the seal where it's needed most—without the need to overtighten the headgear

USER BENEFITS

Personalised comfort

- MicroFit dial adjusts the mask to suit user's profile
- Dual-wall cushion with spring air technology adapts to facial contours and accommodates movement during sleep (including jaw drop)

Light and quiet

- Streamlined design feels light on the face
- Improved vent design reduces noise

Making it easier

- MicroFit dial and spring air cushion technology help achieve good mask fit—first time, every time
- Fewer parts for easy handling, cleaning, assembly and disassembly

PROFESSIONAL BENEFITS

Fit range increased to over 95% of users*

- MicroFit dial and spring air cushion technology adjust the mask to suit more facial structures
- New XS mask size

Customer satisfaction and compliance

- Comprehensive user support package with travel bag, mask care DVD, quick-fitting guide and user guide
- Increased therapy acceptance through first-time fit
- Personalised fit, comfort and seal encourage long-term compliance
- Simplified parts for quick and easy assembly and disassembly

Cost and time efficiencies

- Reduced inventory: 33% fewer product codes than Ultra Mirage™ Full Face Mask and fewer mask parts
- Reduced fitting and training time: improved cushion technology seals quickly
- Reduced callbacks: ease-of-use and fitting

PRODUCT CODES

Mirage Quattro Full Face Mask

Europe 1¹

X-Small	61204
Small	61205
Medium	61206
Large	61207

Europe 2²

X-Small	61208
Small	61209
Medium	61210
Large	61211

Europe 3³

X-Small	61228
Small	61229
Medium	61230
Large	61231

France

X-Small	61220
Small	61221
Medium	61222
Large	61223

UK

X-Small	61224
Small	61225
Medium	61226
Large	61227

Asia-Pacific

X-Small	61216
Small	61217
Medium	61218
Large	61219

ANZ

X-Small	61212
Small	61213
Medium	61214
Large	61215

¹ English, German, French, Spanish, Italian, Portuguese and Dutch

² English, Swedish, Norwegian, Danish, Finnish, Hungarian, Polish and Greek

³ English, Russian, Turkish, Czech, Estonian, Slovenian, Slovakian, Bulgarian, Romanian, Latvian and Lithuanian

* Source: internal ResMed anthropometric data