



caleva

TECHNICAL INFORMATION

LABORATORY

CALEVA

DATASTOR-2™

Easily save record operating parameters automatically from Caleva bench-top equipment to your PC or Laptop.
Protect your important trial results.

The Caleva Datastor-2 can be used* with the Caleva equipment shown below:-



Low pressure screen
Extruder 20



Variable Density Screw
Extruder (radial)



MBS Spheronizer 250



MBS Spheronizer 120



Variable Density Screw
Extruder (axial)

*The machines will require and options upgrade at the time of order. The user can use any modern PC or laptop. The software is supplied free of charge to users of compatible machines.

DATASTOR-2 - SAVE AND PROTECT YOUR LABORATORY TRIAL DATA AUTOMATICALLY


With the release of the new Datastor-2™ system Caleva laboratory extruders and spheronizers have (as an option) the ability to be connected by USB to a PC to record and store operational parameters.

The bench-top extruder or spheronizer is connected via USB to a suitably configured PC or laptop PC and the following data is recorded:-

- Power usage
- Operating speed in rpm
- Product temperature at a targeted location

The software can recognize and monitor up to two machines working independently at the same time from the same C.

Results can be opened by the Datastor software as a one page summary of relevant data. This can be printed as a permanent record. This enhances GLP and gives added quality assurance **automatically**.

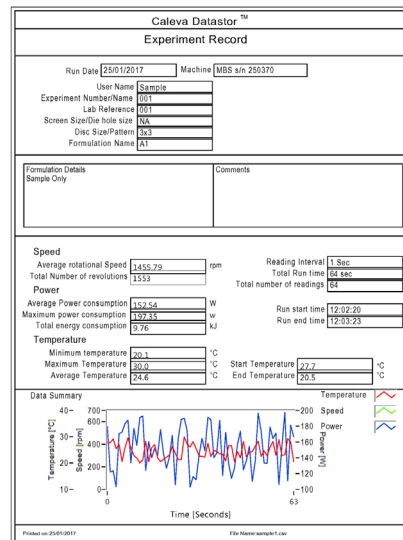

Results sheet produced automatically when the data files are opened.

CALEVA OFFER THE DATASTOR-2™

- A unique option that allows operators to connect an extruder and spheronizer directly to a PC or laptop
- Easy and intuitive operation
- Data readings automatically taken every second
- Permanent data storage on your PC without additional intervention
- Simultaneous operation with one PC on completely separate machines. This is normally one extruder and one spheronizer
- Results can be directly opened as a “read only” results sheet automatically
- It is not possible to retrofit the Datastor™ on extruders and spheronizers and this option must be chosen at time of purchase

TALK TO US

Please call us without obligation

+44 (0) 1258 471122

info@caleva.com



Cert No. 1503
ISO 9001