

PrimeSight™



Introducing PrimeSight™ Endoscopy

The PrimeSight name may be new but the results are not. For more than 20 years, PrimeSight endoscopy systems from Cogentix Medical have provided high-performance endoscopy products under the Vision-Sciences brands of EndoSheath®, Vision Systems® and through the 1000, 2000, 3000, 4000, 5000 and 7000 series of endoscopes.

Key Benefits of PrimeSight Endoscopy Systems

- » High-performance fiberoptic and video systems designed for optimal mobility and plug-and-play simplicity
- » Designed for use with the EndoSheath protective barrier to enable rapid, simplified reprocessing and maximum endoscopy efficiency
- » Supports full range of diagnostic and therapeutic procedures in urology, urogynecology, bariatric medicine, critical care, gastroenterology, otolaryngology and pulmonology

PrimeSight™

PrimeSight Flexible Endoscopes

Our 4000 Series of fiberoptic endoscopes are developed using state-of-the-art technology and offer optimal portability and utility across a variety of settings. In addition, our 5000 and 7000 series of video endoscopes offer CCD-based, high-resolution picture quality with LED based illumination for vibrant colors and tissue differentiation. The video series of endoscopes utilize our uniquely portable, all-in-one, easy-to-use PrimeSight Video System.



Used with the EndoSheath® protective barrier, PrimeSight flexible endoscopes can go anywhere, anytime for efficient and effective endoscopy.

Product	Specifications
Flexible Bronchoscopes	
BRS-4000 Fiberoptic Bronchoscope	Length 570mm, angulation 215° up / 140° down, field of view 95° Outer diameter varies based on choice of EndoSheath protective barrier
BRS-5100 Video Bronchoscope	Length 600mm, angulation 215° up / 140° down, field of view 110° Outer diameter varies based on choice of EndoSheath protective barrier
Flexible Cystoscopes	
CST-4000 Fiberoptic Cystoscope	Length: 370mm, angulation 215° up / 140° down, field of view 110° Outer diameter varies based on choice of EndoSheath protective barrier
CST-5000 Video Cystoscope	Length 370mm, angulation: 215° up / 140° down, field of view 110° Outer diameter varies based on choice of EndoSheath protective barrier
Flexible Esophagoscope	
TNE-5000 Video Esophagoscope	Length 650mm, angulation 140° up / 215° down, field of view 110° Outer diameter varies based on choice of EndoSheath protective barrier
Flexible Laryngoscope	
ENT-4000 Fiberoptic Laryngoscope	Length 300mm, angulation 140° up / down, field of view 90°, Distal tip diameter 3mm
ENT-4500 Fiberoptic Laryngoscope	Length 300mm, angulation 140° up / down, field of view 80° Distal tip diameter 2.4mm
ENT-5000 Video Laryngoscope	Length 300mm, angulation 140° up / down, field of view 90° Distal tip diameter 3mm
ENT-7000 Video Laryngoscope	Length 300mm, angulation 140° up / down, field of view 90° Distal tip diameter 2.4mm

PrimeSight Rigid Endoscopes

Cogentix Medical also offers a variety of HD otoscopes, sinusscopes, laryngoscopes and cystoscopes. Our high quality, stainless steel endoscopes have exceptional color reproduction, sharper image and greater depth of field compared to traditional endoscopes. PrimeSight rigid endoscopes are compatible with the PrimeSight video system.



PrimeSight Video System

The PrimeSight all-in-one video system is designed for efficient and portable high performance imaging with all PrimeSight video endoscopes and cameras. Featuring easy to use workflow simplicity features, audio/video recording and review, and a full range of connectivity options, the PrimeSight video system is ideal in any setting, from the physician's office to the large hospital networks.

PrimeSight Light Sources

Cogentix Medical offers a variety of light sources to support surgeon preference including halogen, xenon, and LED models. All light sources are designed with portability, durability and performance in mind.

- » Halogen Light Source
- » Xenon Light Source
- » LED Light Source
- » Battery Powered LED Light Source

PrimeSight Stroboscopy System

The PrimeSight stroboscopy system allows physicians an affordable way to integrate stroboscopy. Our compact unit replaces the large tower of equipment that was previously necessary for high quality, true flashlight stroboscopy. Using the PrimeSight stroboscopy unit with the PrimeSight video system, practitioners can easily move between standard visualization and stroboscopy with a touch of a single button, reducing the need to switch between flexible and rigid endoscopes.

PrimeSight Add-On Digital Camera

The PrimeSight digital camera is used with flexible fiberoptic or rigid endoscopes and plugs directly into a PrimeSight video system. The PrimeSight camera delivers optimal image quality by enabling the user to enhance brightness, adjust contrast, and control white balance and colors for sharp image representation and tissue differentiation.

Accessories

Cogentix Medical offers a variety of accessories to support your busy endoscopy practice from endoscope preparation through the procedures. These include:


- » Cytology brushes, aspirating needles, injection needles, forceps, graspers and baskets
- » Leak testers to verify the integrity of PrimeSight™ endoscopes before immersion into a chemical disinfection fluid or reprocessing machine
- » Light guide sleeves to ensure compatibility with virtually all manufacturers' light sources
- » EndoWipe® towelettes and sponges for quick and easy reprocessing of PrimeSight scopes in 10 minutes or less.

EndoSheath®

PROTECTIVE BARRIER

Designed for use with PrimeSight™ Endoscopy Systems

PrimeSight™ Scope

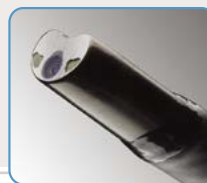

with EndoSheath®
Protective Barrier



The EndoSheath® protective barrier is a sterile, single-use product with an integrated working channel, that is used with PrimeSight™ flexible endoscopes. It is an always ready, always sterile solution.

Safety and efficiency? We've got you covered

The EndoSheath protective barrier covers the endoscope from end to end, isolating the endoscope from patient contact and protecting the endoscope controls from contamination during the procedure. As a result, there is no need for high level disinfection between cases and you can quickly reprocess a scope in less than 10 minutes.

Single-use working channel revolutionizes scope versatility

Since the disposable EndoSheath protective barrier includes a working/operative channel for suction, irrigation and tool passage, physicians can use PrimeSight endoscopes with the EndoSheath protective barrier for a full range of diagnostic and therapeutic procedures. The EndoSheath protective barrier is available in multiple models, each with a different size of working channel, allowing you to adapt a single endoscope to a wide array of procedures.

Product	Channel size	Compatible endoscope	Catalog number
For PrimeSight Bronchoscopes			
BF-2.1 EndoSheath Protective Barrier (10/box)	2.1mm	BRS-4000 Fiber Bronchoscope	04-4101
BF-1.5 EndoSheath Protective Barrier (10/box)	1.5mm	BRS-4000 Fiber Bronchoscope	04-4102
BF-0 EndoSheath Protective Barrier (10/box)	NA	BRS-4000 Fiber Bronchoscope	04-4105
BV-2.1 EndoSheath Protective Barrier (10/box)	2.1mm	BRS-5000 Video Bronchoscope	04-5101
BV-1.5 EndoSheath Protective Barrier (10/box)	1.5mm	BRS-5000 Video Bronchoscope	04-5102
BV-0 EndoSheath Protective Barrier (10/box)	NA	BRS-5000 Video Bronchoscope	04-5105
BVX-2.1 EndoSheath Protective Barrier (10/box)	2.1mm	BRS-5100 Video Bronchoscope	04-5111
BVX-1.5 EndoSheath Protective Barrier (10/box)	1.5mm	BRS-5100 Video Bronchoscope	04-5112
BVX-0 EndoSheath Protective Barrier (10/box)	NA	BRS-5100 Video Bronchoscope	04-5115
For PrimeSight Cystoscopes			
CF-2.1 EndoSheath Protective Barrier (10/box)	2.1mm	CST-4000 Fiber Cystoscope	08-4101
CF-1.5 EndoSheath Protective Barrier (10/box)	1.5mm	CST-4000 Fiber Cystoscope	08-4102
CV-2.1 EndoSheath Protective Barrier (10/box)	2.1mm	CST-5000 Video Cystoscope	08-5101
CV-1.5 EndoSheath Protective Barrier (10/box)	1.5mm	CST-5000 Video Cystoscope	08-5102
For PrimeSight Esophagoscopes			
TV-2.1 EndoSheath Protective Barrier (10/box)	2.1mm	TNE-5000 Video Esophagoscope	03-5101
TV-1.5 EndoSheath Protective Barrier (10/box)	1.5mm	TNE-5000 Video Esophagoscope	03-5102

Connect with Cogentix Medical



Rx Only. EndoSheath and EndoWipe are registered trademarks and PrimeSight is a trademark of Cogentix Medical. © 2015 Cogentix Medical. All rights reserved. 70020A 10/15

Cogentix
Medical

cogentixmedical.com

Cogentix Medical
5420 Feltl Road
Minnetonka, MN 55343

TEL 866.258.2182
FAX 866.255.4522